


LabOne Programming Manual

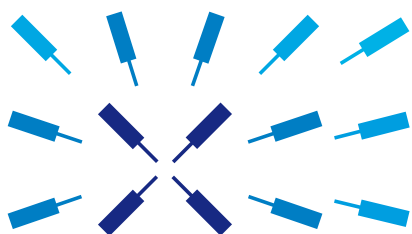

Zurich
Instruments

LabOne Programming Manual

Zurich Instruments AG

Publication date Revision 38200

Copyright © 2008-2016 Zurich Instruments AG

The contents of this document are provided by Zurich Instruments AG (ZI), “as is”. ZI makes no representations or warranties with respect to the accuracy or completeness of the contents of this publication and reserves the right to make changes to specifications and product descriptions at any time without notice.

LabVIEW is a registered trademark of National Instruments Inc. MATLAB is a registered trademark of The MathWorks, Inc. All other trademarks are the property of their respective owners.

Revision History

Revision 38200, 12-Jul-2016:

Update of the LabOne Programming Manual for LabOne Release 16.04.

- Update of Core Module parameter tables.

Revision 34390, 23-Dec-2015:

Update of the LabOne Programming Manual for LabOne Release 15.11.

- [Sweeper Module](#) parameter and description update; explanation of `sweep/inaccuracy`.
- LabOne Matlab, Python and C APIs: Documentation added describing the [APIs logging capabilities](#).
- LabOne C API Documentation added for [Core Module functionality](#).
- LabOne C API: Updates to the [Error Handling section](#) following the introduction of the `ziGetLastError` function.

Revision 31421, 8-Jul-2015:

Update of the LabOne Programming Manual for LabOne Release 15.05.

- The LabOne LabVIEW and C (ziAPI) APIs are now `ziCore`-based.
- Additions and modifications for MFLI support.

Revision 28870, 18-Mar-2015:

Update of the LabOne Programming Manual for LabOne Release 15.01.

- Added description of API Level 5.

Revision 26206, 01-Oct-2014:

Update of the LabOne Programming Manual for LabOne Release 14.08.

- Added PLL Advisor Module section to `ziCore` Modules.
- Consistency update of `ziCore` Module parameters.
- Improvements to plots in Software Trigger `ziCore` Module section.

Revision 23212, 23-Apr-2014:

- First release of the LabOne Programming Manual.
-

Table of Contents

I. LabOne Programming Concepts	4
1. Introduction	5
1.1. LabOne Software Architecture	6
1.2. Comparison of the LabOne APIs	8
1.3. Initializing a Connection to a Data Server	9
1.4. Configuring and Obtaining Data from an Instrument	13
1.5. Instrument-Specific Considerations	16
2. ziCore Programming Overview	17
2.1. An Introduction to ziCore-based APIs	18
2.2. Sweeper Module	21
2.3. zoomFFT Module	32
2.4. Software Trigger (Recorder) Module	34
2.5. Device Settings Module	40
2.6. PLL Advisor Module	41
2.7. Tips and Tricks	43
II. LabOne APIs	44
3. Matlab Programming	45
3.1. Installing the LabOne Matlab API	46
3.2. Getting Started with the LabOne Matlab API	49
3.3. LabOne Matlab API Tips and Tricks	54
3.4. Troubleshooting the LabOne Matlab API	56
3.5. LabOne Matlab API (ziDAQ) Command Reference	58
4. Python Programming	76
4.1. Installing the LabOne Python API	77
4.2. Getting Started with the LabOne Python API	81
4.3. LabOne Python API Tips and Tricks	85
4.4. LabOne Python API (ziPython) Command Reference	86
5. LabVIEW Programming	123
5.1. Installing the LabOne LabVIEW API	124
5.2. Getting Started with the LabOne LabVIEW API	126
5.3. LabVIEW Programming Tips and Tricks	131
6. C Programming	132
6.1. Getting Started	133
6.2. Module Documentation	135
6.3. Data Structure Documentation	248
6.4. File Documentation	303
Glossary	470
Index	476

Part I. LabOne Programming Concepts

This part of the manual provides an overview of LabOne programming and deals with generic concepts that apply to any of Zurich Instruments' APIs.

Refer to:

- [Chapter 1](#) for an overview of programming with Zurich Instruments LabOne.
- [Chapter 2](#) for an overview to working with `ziCore` Modules.

Chapter 1. Introduction

This chapter briefly describes the different possibilities to interface with a Zurich Instruments device, other than via the LabOne User Interface or ziControl (HF2 Series only). Zurich Instruments devices are designed with the concept that "the computer is the cockpit"; there are no controls on the front panel of the instrument, instead the user can configure their instrument from and stream data directly to their computer. The aim of this approach is to give the user the freedom to choose where they connect to, and how they control, their instrument.

As an example, the user can either work on a computer directly connected to the instrument via USB or remotely from a different computer on the network, away from their experimental setup. Then, on either computer, the user can configure and retrieve data from their instrument via a number of different software interfaces, i.e. via the web-based LabOne User Interface and/or their own custom programs. In this way the user can decide which connectivity setup and combination of interfaces best suits their experimental setup and data processing needs.

Refer to:

- [Section 1.1](#) for an overview of the [LabOne Software Architecture](#).
- [Section 1.2](#) for a [Comparison of the LabOne APIs](#).
- [Section 1.3](#) for help [Initializing a Connection to a Data Server](#).
- [Section 1.4](#) for help [Configuring and Obtaining Data from an Instrument](#).
- [Section 1.5](#) for [Instrument-Specific Considerations](#).

Note

New users could benefit by first familiarizing themselves with the instrument using the LabOne User Interface or ziControl; please refer to the appropriate user manual for your instrument for more details.

Note

The Real-time Option (RTK) for the HF2 Series is not a PC-based interface for controlling an instrument and is documented in the HF2 User Manual.

1.1. LabOne Software Architecture

Zurich Instruments devices uses a server-based connectivity methodology. Server-based means that all communication between the user and the instrument takes place via a computer program called a server, the Data Server. The Data Server recognizes available instruments and manages all communication between the instrument and the host computer on one side, and communication to all the connected clients on the other side. This allows for:

- A multi-client configuration: Multiple interfaces (even from multiple computers on the network) can access the settings and data on an instrument. Settings are synchronized between all interfaces by the single instance of the Data Server.
- A multi-device setup: Any of the Data Server's clients can access multiple devices simultaneously.

This software architecture is organized in layers, see [Figure 1.1](#) for a schematic of the software layers.

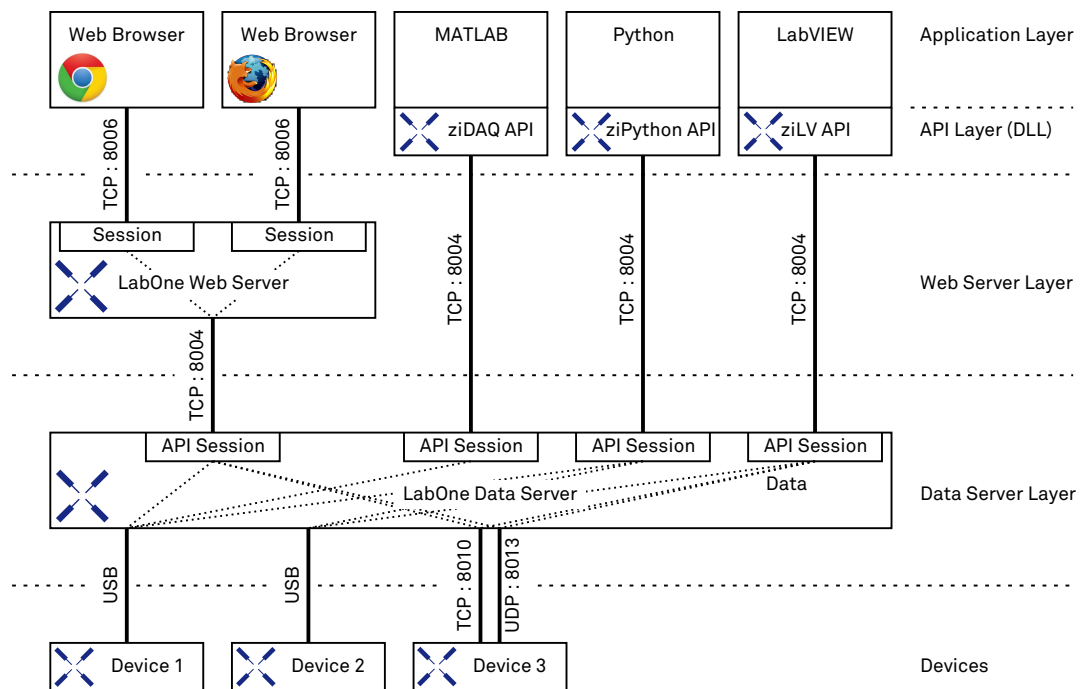


Figure 1.1. LabOne Software Architecture. The above diagram depicts the software architecture when using UHF and HF2 Instruments. In the case of MF Instruments the server runs on the device itself instead of on a PC; only one MF device can be accessed from the Data Server. Web Server and API usage for the MF is analogous to that of other instruments.

First, we briefly explain some terminology that is used throughout this manual.

- **Host computer:** The computer where the Data Server is running and that is directly connected to the instrument. Multiple remote computers on a local area network can access the instrument by creating an API connection to the Data Server running on the host computer.
- **Data Server:** A computer program that runs on the host computer and manages settings on, and data transfer to and from instruments by receiving commands from clients. It always has the most up-to-date configuration of the device and ensures that the configuration is synchronized between different clients.

- `ziServer.exe`: The Data Server that handles communication with HF2 Instruments.
- `ziDataServer.exe`: The Data Server that handles communication with UHFLI and MFLI Instruments. Note, in the case of MFLI Instruments the Data Server runs on the instrument itself.
- Remote computer: A computer, available on the same network as the host computer, that can communicate with an instrument via the Data Server program running on the host.
- Client: A computer program that communicates with an instrument via the Data Server. The client can be running either on the host or the remote computer.
- API (Application Programming Interface): a collection of functions and data structures which enable communication between software components. In our case, the various APIs (e.g., LabVIEW, Matlab®) provide functions to configure instruments and receive measured experimental data.
- Interface: Either a client or an API.
- GUI (Graphical User Interface): A computer program that the user can operate via images as opposed to text-based commands.
- LabOne User Interface: The browser-based user interface that connects to the Web Server.
- LabOne Web Server: The program that generates the browser-based LabOne User Interface.
- `ziControl`: The standard GUI shipped for use with HF2 Instruments (before software release 15.11). HF2 support was added to the LabOne User Interface for devices with the WEB Option installed in LabOne software release 15.11 .
- `ziCore`: The internal core C++ library upon which many APIs are based, see [Part II](#) of this document.
- Modules: `ziCore` software components that provide a unified interface to APIs to perform a specific high-level common task such as sweeping data.

1.2. Comparison of the LabOne APIs

The various software interfaces available in LabOne allow the user to pick a programming environment they are familiar with to achieve fast results. All other things being equal, here is a brief discussion of the merits of each interface.

- The [LabVIEW interface](#) allows for quick and efficient implementation of virtual instruments that run independently. These can easily be integrated in existing experiment control performed in LabVIEW. This interface requires a National Instruments LabVIEW license and LabVIEW 2009 (or higher).
- The [Matlab® interface](#) allows the user to directly obtain measurement data within the Matlab programming environment, where they can make use of the many built-in functions available. This interface requires a Mathworks Matlab license, but no additional Matlab Toolboxes.
- The [Python interface](#) allows the user to directly obtain measurement data within python. Python is available as free and open source software; no license is required to use it.
- The [C API](#), [ziAPI](#), is a very versatile interface that will run on most platforms. However, since C is a low-level programming language, the development cycle is slower than with the other programming environments.
- The text-based interface (HF2 Series only) allows the user to manually connect to the HF2 Data Server in a console via telnet. While this interface is a very useful tool for HF2 programmers to verify instrument configuration set by other interfaces, it is limited in terms of performance and maximum demodulator sample rate. See the HF2 User Manual for more details.

Note

From LabOne Release 15.05 onwards the high-level functionality provided by ziCore's Sweep and Software Trigger [Modules](#) is also available in the LabVIEW and C APIs, all modules are available from the Matlab and Python LabOne APIs.

1.3. Initializing a Connection to a Data Server

As described in [Section 1.1](#) an API client communicates with an instrument via a data server over a TCP/IP socket connection. As such, the first step towards communicating with an instrument is initializing an API session to the correct data server for the target device.

The choice of data server depends on the device class (and in general the user's network topology). Users must be aware that HF2 instruments operate via a different data server program than other instruments (UHF, MF) and users of MF instruments should be aware that the data server runs on the MF instrument itself and not on a separate PC. Regardless of which data server is used and where it is running the desired data server is specified by providing three parameters:

- the data server host's address (hostname),
- the data server port,
- the API level to use for the session.

1.3.1. Specifying the Data Server Hostname and Port

For users working with a single device, this section describes how to quickly connect to the correct data server by manually specifying the required data server's hostname and port and the required API Level. Each API has a connect function which takes these three parameters in order to initialize an API session, for example, in the LabOne Matlab API:

```
>>> ziDAQ('connect', serverHostname, serverPort, apiLevel);
```

Data Server Port

A LabOne API client connects to the correct Data Server for their instrument by specifying the appropriate port. By default, the data server programs for UHF and MF Instruments listen to port 8004 for API connections and the data server program for HF2 instruments listens to port 8005. The value of the port that the data server listens to can be changed using the `--port` command-line option when starting the data server.

Data Server Hostname (UHF, HF2 instruments)

In the simplest configuration for HF2 and UHF instruments, the instrument is attached to the same PC where both the data server and API client are running. Since the API client is running on the same PC as the data server, the `'localhost'` (equivalently, `'127.0.0.1'`) should be specified as the data server address, [Figure 1.2](#).

The API client may also connect to a data server running on a different PC from the client. In this case, the data server address should be the IP address (or hostname, if available) of the PC where the data server is running. Note, remote data server access is not enabled by default and the data server must be configured in order to listen to non-localhost connections by either enabling the `--open-override` command-line option when starting the data server or by setting the value of the server node `/zi/config/open` to 1 on a running data server (clearly only possible from a client running on the localhost). See [Section 1.4.1](#) for more information on nodes.

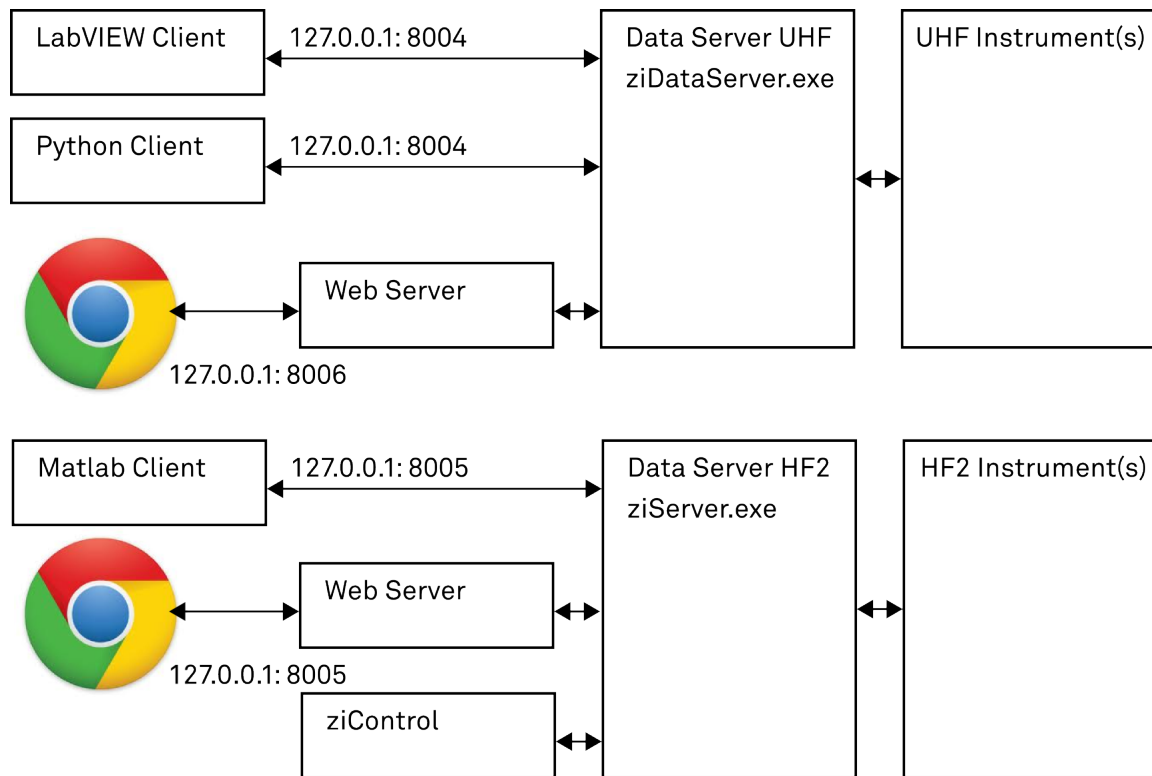


Figure 1.2. Server address and port handling for HF2 and UHF instruments for the case where the API client and data server are running on the same PC. In this case the server hostname is `localhost` and the default port value is 8004 for UHF Instruments and 8005 for HF2 Instruments.

Data Server Hostname (MF instruments)

In the case of MF instruments the data server is running on the instrument itself and as such an API client from a PC is always accessing the data server remotely. Thus, in this case the data server hostname is the value of the instrument's hostname. This will be the same hostname (but not port) that is used to run the LabOne User Interface in a web browser (when the Web Server is running on the MF instrument), see [Figure 1.3](#). Please see the Getting Started chapter of the MFLI User Manual for more details.

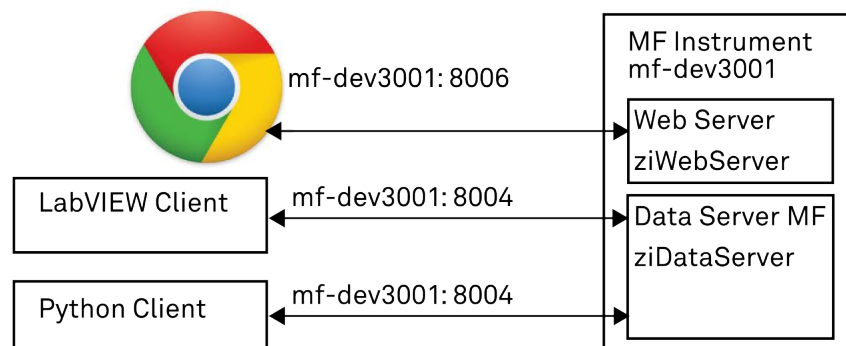


Figure 1.3. Server address and port handling on MF Instruments. Note, the data server is running on the instrument and the server hostname is the same as the instrument's hostname. The default data server port value is 8004 for MF Instruments. In this example, the MF has device serial dev3001.

API Level and Connectivity Examples

The last parameter to specify, the API level, specifies the version of the API to use for the session. In short, an API Level of 1 must be used for HF2 devices and an API Level of 5 is recommended for other instruments. Since the default API Level is 1, it is necessary to specify this parameter for UHF and MF instruments. A more detailed explanation of API Levels is provided in [Section 1.3.2](#).

For example, to initialize a session to the HF2's data server running on the `localhost` with the LabOne Python API, the following commands should be used:

```
>>> import zhinst.ziPython
>>> daq = zhinst.ziPython.ziDAQServer('localhost', 8005, 1)
```

and in order to connect to the data server running on the MF instrument with device serial 'dev3001' with the LabOne Matlab API:

```
>> ziDAQ('connect', 'mf-dev3001', 8004, 5)
```

1.3.2. LabOne API Levels

All of the LabOne APIs are based on an internal core API. Needless to say, we try as hard as possible to make any improvements in our core API backwards compatible for the convenience of our users. We take care that existing programs do not need to be changed upon a new software release. Occasionally, however, we do have to make a breaking change in our API by removing some old functionality. This old functionality is, however, phased out over several software releases. First, the functionality is marked as deprecated and the user is informed via a depreciation warning (this can be turned off). This indicator warns that this function may be unsupported in the future. If we have to break some functionality we use a so-called **API level**.

With support of new devices and features we need to break functionality on the ziAPI.h e.g. data returned by poll commands. In order to still support the old functionality we introduced API levels. If a program is only using old functionality the API level 1 (default) can be used. If a user needs new functionality, they need to use a higher API level. This will usually need some changes on the existing code.

Available API levels as of LabOne Software Release 15.01 are:

- API Level 1: HF2 support, basic UHF support.
- API Level 4: UHF support with timestamps and PWA, name clean-up.
- API Level 5: Introduction of scope offset for extended (non-hardware) scope inputs (UHF, MF Instruments).

Note that Levels 2 and 3 are used only internally and are not available to the general public.

Note

The HF2 Series only supports API Level 1.

Note

New UHFLI and MFLI API users are recommended to use API Level 5.

API Level 4 Features

The new features in API Level 4 are:

- Timestamps are available for any settings or data node.
- Greatly improved Scope data transfer rates (and new Scope data structure).
- Greatly improved UHF Boxcar and PWA support.

API Level 5 Features

API Level 5 was introduced in LabOne Release 15.01 to accommodate a necessary change in the Scope data structure:

- The Scope data structure was extended with the new field "channeloffset" which contains the offset value that must be added to the scaled wave value in order to obtain the physical value recorded by the scope. For previous hardware scope "inputselects" there is essentially no change, since their offset is always zero. However, for the extended values of "inputselects", such as PID Out value, (available with the DIG option) the offset is determined by the values of "limitlower" and "limitupper" configured by the user.

1.4. Configuring and Obtaining Data from an Instrument

1.4.1. Finding settings: The Node Hierarchy

In order to communicate with an Zurich Instruments device via text-based commands in an API, it is necessary to understand how the settings and measurement data of the instrument are accessed. All the settings and data of the instrument are organized in a file-system-like hierarchical structure. The features of the instrument, such as demodulators, are accessed as branches in this tree and their individual settings are leaves of these branches. It is also possible to browse branches inside the tree as if the user were navigating in a file-system. This hierarchy is used, no matter which interface you use when performing measurements.

An example demonstrating the hierarchy is the representation of the first demodulator on the device, given by the node:

```
/devX/demods/0
```

which, as we've already noted, is very similar to a **path** on a computer's file-system. Note that, the top level of the path is the device that you are connected to. The demodulators are then given as a top-level **node** under your device-node and the node of the first demodulator is indexed by 0. This path represents a branch in the node hierarchy which, in this case, if we explore further, has the following nodes:

```
/devX/demods/0/adcselect  
/devX/demods/0/order  
/devX/demods/0/timeconstant  
/devX/demods/0/rate  
/devX/demods/0/trigger  
/devX/demods/0/oscselect  
/devX/demods/0/harmonic  
/devX/demods/0/phaseshift  
/devX/demods/0/sinc  
/devX/demods/0/sample
```

These nodes are **leaves**, the most bottom-level nodes which represent a setting of an instrument or a field that can be read to retrieve measurement data. For example, `/devX/demods/0/adcselect` is the leaf that controls the setting corresponding to the choice of signal input for the first demodulator. To set the index of the signal input the user writes to this node. The leaf `/devX/demods/0/sample` is the leaf where the demodulator's output (timestamp, demodulated x-value, demodulated y-value) are written at the frequency specified by `/devX/demods/0/rate`. In order to obtain the demodulator output you read the values from this node by **polling** this node. Polling a node sends a request from the client to ziServer to obtain the data from the node at that particular point in time.

Note

The numbering on the front panel of the UHFLI, MFLI and HF2 Instruments and the block numbering in the graphical user interfaces generally start with 1 (1-based indexing). Note, that when accessing settings and data via a software interface, the numbering starting with 0 (0-based indexing).

Note

A useful method to learn about the nodes of your instrument is to look at the output of the history in the bottom of the graphical user interface. The status line always shows the last applied command and you can view the entire history by clicking the "Show Log" or "Show History" button. You will find paths like

```
/devx/sigins/0/ac = 1
```

after you switched on the AC mode for signal input 1, or

```
/devx/demods/1/rate = 7200.000000
```

after changing the readout rate of demodulator 2 to 7.2 kHz.

1.4.2. Obtaining Data from the Instrument

The subscribe and poll commands

The easiest way to obtain data from an instrument is via the `poll` command, available in all of the LabOne API interfaces. The `poll` command is a function for synchronous data recording from specified [nodes](#) of an instrument. Synchronous means that the interface is blocked during execution of the command, see [Section 2.1.4](#) for asynchronous alternatives. `poll` takes two obligatory input arguments **recording time** and **timeout**.

The `subscribe` and `unsubscribe` commands are used to select the [nodes](#) from which data should be recorded. After subscribing to the node, the Data Server's internal data buffer will start filling with data from the subscribed nodes. The `poll` command will return the data that was recorded for the specified recording time (obligatory input argument) and any data that was already in the buffer since the last `poll`. To get rid of the data from earlier measurements it's possible to clear the buffer before polling by using the `flush` command.

In order to avoid losing data (the Data Server has a finite amount of memory available for its data buffers), long recording times (> 20s, depending on sampling rates and available memory) should be avoided. However, since internal data buffering on the Data Server ensures that no data is lost between `poll` commands, it's possible to record for longer periods of time by using the `poll` command inside a loop. In order to check that no data has been lost during a poll, the demodulator sample's `time` flags can be checked, see [the section called "Demodulator Sample Data Structure"](#).

If no data was stored in the Data Server's data buffer after issuing a `poll`, the command will wait for the data until the timeout time. If the buffer is empty after timeout time passed, `poll` will throw an error.

Note

One of the LabOne [ziCore Modules](#) could be a more efficient choice for data retrieval than the comparably low-level `poll` command, see the [Section ziCore Modules](#) in [Part II](#).

Demodulator Sample Data Structure

An instrument's demodulator data is returned as a data structure (typically a `struct`) with the following fields (regardless of which API Level is used):

<code>timestamp</code>	The instrument's timestamp of the measured demodulator data <code>uint64</code> . Divide by the instrument's clockbase (<code>/dev123/clockbase</code>) to obtain the time in seconds.
<code>x</code>	The demodulator <code>x</code> value in Volts [<code>double</code>].
<code>y</code>	The demodulator <code>y</code> value in Volts [<code>double</code>].
<code>frequency</code>	The current frequency used by the demodulator in Hertz [<code>double</code>].
<code>phase</code>	The oscillator's phase in Radians (not the demodulator phase) [<code>double</code>].
<code>auxin0</code>	The auxiliary input channel 0 value in Volts [<code>double</code>].
<code>auxin1</code>	The auxiliary input channel 1 value in Volts [<code>double</code>].
<code>bits</code>	The value of the digital input/output (DIO) connector. [<code>integer</code>].
<code>time.dataloss</code>	Indicator of sample loss (including block loss) [<code>bool</code>].
<code>time.blockloss</code>	Indication of data block loss over the socket connection. This may be the result of a too long break between subsequent poll commands [<code>bool</code>].
<code>time.invalidtimestamp</code>	Indication of invalid time stamp data as a result of a sampling rate change during the measurement [<code>bool</code>].

Note

[Chapter 6](#) contains some details of other data structures.

1.5. Instrument-Specific Considerations

This section describes some instrument-specific considerations when programming with the LabOne APIs.

1.5.1. UHF-Specific Considerations

UHF Lock-in Amplifiers perform an automatic calibration 10 minutes after power-up of the Instrument. This internal calibration is necessary to achieve the specifications of the system. However, if necessary, it can be ran manually by setting the device node `/devN/system/calib/calibrate` to 1 and then disabled using the `/devN/system/calib/auto` node.

The calibration routine takes about 200 ms and during that time the transfer of measurement data will be stopped on the Data Server level. If a [ziAPI](#) (LabOne C API) or [LabVIEW](#) client is polling data during this time, the user will experience data loss; ziAPI has no functionality to deal with such a streaming interrupt. Clients polling data from [ziCore-based APIs](#) (i.e. Matlab or Python APIs) will be informed of data loss, which allows the user to ignore this data.

Please see the UHF User Manual for more information about device calibration.

Chapter 2. ziCore Programming Overview

The LabOne APIs provide interfaces to configure, acquire data from, and run integral functionality of your Zurich Instruments device. These high-level interfaces are, however, just thin application layers based on a shared core API, `ziCore`. This chapter aims to describe the common functionality that's available to any of the interfaces (Matlab, Python, C, LabVIEW) based on `ziCore`.

Refer to:

- [Section 2.1](#) for [An Introduction to ziCore-based APIs](#).
- [Section 2.2](#) for the [Sweeper Module](#).
- [Section 2.3](#) for the [zoomFFT Module](#).
- [Section 2.4](#) for the [Software Trigger \(Recorder\) Module](#).
- [Section 2.5](#) for the [Device Settings Module](#).
- [Section 2.6](#) for the [PLL Advisor Module](#).
- [Section 2.7](#) for some `ziCore` programming [Tips and Tricks](#).

2.1. An Introduction to ziCore-based APIs

All of the LabOne APIs are based on a central API called `ziCore`. This allows them to share a common structure which provides a uniform interface for programming Zurich Instruments devices. The aim of this section is to familiarize the user with the key `ziCore` programming concepts which can then be used in any of the LabOne APIs (Matlab, Python, LabVIEW and C).

2.1.1. Software Architecture

Each of the `ziCore`-based APIs are designed to have a minimal code footprint: They are simply small interface layers that use the functionality derived from `ziCore`, a central C++ API. The derived API interfaces (Matlab, Python, LabVIEW and C) provide a familiar interface to the user and allow the user to receive and manipulate data from their instrument using the API language's native data types and formats. See [Section 1.1](#) for an overview of the LabOne software architecture.

2.1.2. ziCore Modules

In addition to the usual API commands available for instrument configuration and data retrieval, e.g., `setInt`, `poll`, `ziCore`-based APIs also provide a number of so-called **Modules**: high-level interfaces that perform common tasks such as sweeping data or performing **FFTs**.

The Module's functionality is implemented in `ziCore` and each derived high-level API simply provides an interface to that module from the API's native environment. This design ensures that the user can expect the same behavior from each module irrespective of which API is being used; if the user is familiar with a module available in one high-level programming API, it is quick and easy to start using the module in a different API. In particular, the LabOne User Interface is also based on `ziCore` and as such, the user can expect the same behavior using a `ziCore`-based API that is experienced in the LabOne User Interface, see [Figure 2.1](#).

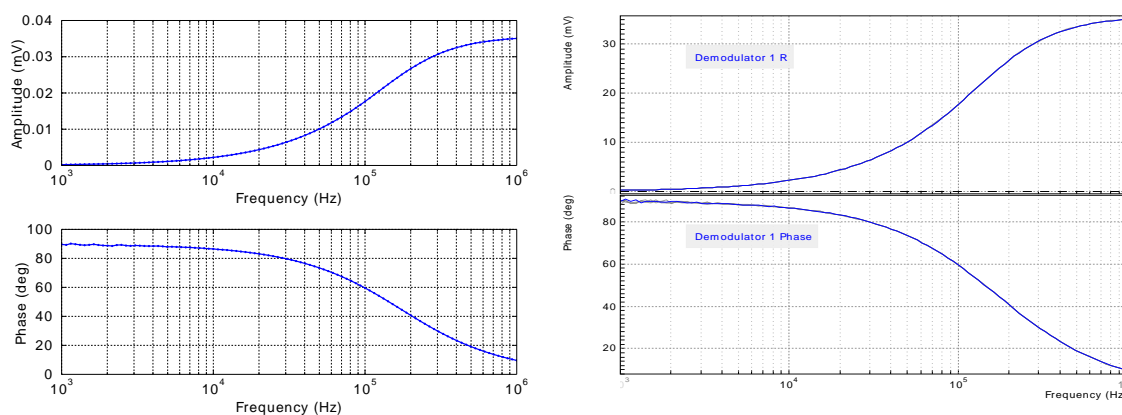


Figure 2.1. The same results and behavior can be obtained from Modules in any `ziCore`-based interface; **Sweeper Module** results from the LabOne Matlab API (left) and the LabOne User Interface (right) using the same Sweeper and instrument settings.

The modules currently available in `ziCore` are:

- The **Sweeper Module** for obtaining data whilst performing a sweep of one of the instrument's setting, e.g., measuring a frequency response.
- The **zoomFFT Module** for calculating the FFT of demodulator output.

- The [Software Trigger \(Recorder\) Module](#) for recording instrument data **asynchronously** based upon user-defined triggers.
- The [Device Settings Module](#) for saving and loading instrument settings to and from (XML) files.
- The [PLL Advisor Module](#) for modeling/simulating the PLL (phase-locked loop) incorporated in the instrument (available for UHF Lock-in Amplifiers only).

In addition to providing a unified-interface between APIs, modules also provide a uniform work-flow regardless of the functionality the module performs (e.g., sweeping, recording data), see [Section 2.1.3](#).

An important difference to low-level ziCore API commands is that Modules execute their commands **asynchronously**, see [Section 2.1.4](#).

Note

The LabOne User Interface Command Log can be set to store commands in either Matlab or Python formats which can then be used to start writing custom programs, see [Section 2.7](#).

Note

Much of the same functionality is provided in ziControl, but ziControl UI is not based on ziCore.

2.1.3. ziCore Module Work-Flow

Regardless of the Module's function, all ziCore Modules follow same work flow in all of the derived interfaces:

- create (instantiate) an instance of the module,
- **set** the module's parameters using `path`, `value` pairs,
- **subscribe** to instrument nodes from which to obtain data (note, this is a module subscribe, which is different from a normal API session subscribe command),
- **execute** the module (this starts the module's thread),
- wait until the module has **finished** executing; intermediate reading of data is possible,
- **read** the module's data,
- **clear** the module to remove it from memory.

The highlighted words above are commands for all the Modules. For interface-specific concepts when using Modules see the following Sections:

- [Using ziCore Modules in the LabOne Matlab API](#),
- [Using ziCore Modules in the LabOne Python API](#),
- [Using ziCore Modules in the LabOne LabVIEW API](#),
- [Using ziCore Modules in the LabOne C API](#).

2.1.4. Synchronous versus Asynchronous Commands

The low-level API commands such as `setInt` and `poll` are **synchronous** commands, that is the interface will be blocked until that command has finished executing; the user can not run any

commands in the meantime. Another feature of ziCore's Modules is that each instantiation of a Module creates a new [Thread](#) and, as such, the commands executed by a Module are performed **asynchronously**. Asynchronous means that the task is performed in the background and the interface's process is available to perform other tasks in the meantime, i.e., Module commands are non-blocking for the user.

2.2. Sweeper Module

The Sweeper Module allows the user to perform sweeps as in the Sweeper Tab of the LabOne User Interface. In general, the Sweeper can be used to obtain data when measuring a DUT's response to varying (or **sweeping**) one instrument setting while other instrument settings are kept constant.

2.2.1. Configuring the Sweeper

In the following we briefly describe how to configure the Sweeper Module. See [Table 2.1](#) for a full list of the Sweeper's parameters and [Table 2.2](#) for a description of the Sweeper's outputs.

Specifying the Instrument Setting to Sweep

The Sweeper's `sweep/gridnode` parameter, the so-called **sweep parameter**, specifies the instrument's setting to be swept, specified as a path to an instrument's [node](#). This is typically an oscillator frequency in a [Frequency Response Analyzer](#), e.g., `/dev123/oscs/0/freq`, but a wide range of instrument settings can be chosen, such as a signal output amplitude or a PID controller's setpoint.

Specifying the Range of Values for the Sweep Parameter

The Sweeper will change the sweep parameter's value `sweep/samplecount` times within the **range** of values specified by `sweep/start` and `sweep/stop`. The `sweep/xmapping` parameter specifies whether the spacing between two sequential values in the range is linear (`=0`) or logarithmic (`=1`).

Controlling the Scan mode: The Selection of Range Values

The `sweep/scan` parameter defines the **order** that the values in the specified range are written to the sweep parameter. In sequential scan mode (`=0`), the sweep parameter's values change incrementally from smaller to larger values, see [Figure 2.4](#). In order to scan the sweep parameter's in the opposite direction, i.e., from larger to smaller values, reverse scan mode (`=3`) can be used.

In binary scan mode (`=1`) the first sweep parameter's value is taken as the value in the middle of the range, then the range is split into two halves and the next two values for the sweeper parameter are the values in the middle of those halves. This process continues until all the values in the range were assigned to the sweeper parameter, see [Figure 2.6](#). Binary scan mode ensures that the sweep parameter uses values from the entire range near the beginning of a measurement, which allows the user to get feedback quickly about the measurement's entire range. Since the Sweeper Module is an [asynchronous](#) interface, it's possible to continuously read and plot data whilst the sweep measurement is ongoing and update points in a graph dynamically.

In bidirectional scan mode (`=2`) the sweeper parameter's values are first set from smaller to larger values as in sequential mode, but are then set in reverse order from larger to smaller values, see [Figure 2.5](#). This allows for effects in the sweep parameter to be observed that depend on the order of changes in the sweep parameter's values.

Controlling how the Sweeper sets the Demodulator's Time Constant

The `sweep/bandwidthcontrol` parameter specifies which demodulator filter bandwidth (equivalently time constant) the Sweeper should set for the current measurement point. The user can either specify the bandwidth manually (`=0`), in which case the value of the current demodulator filter's bandwidth is simply used for all measurement points; specify a fixed bandwidth (`=1`), specified by `sweep/bandwidth`, for all measurement points; or specify that

the Sweeper sets the demodulator's bandwidth automatically (=2). Note, to use either Fixed or Manual mode, `sweep/bandwidth` must be set to a value > 0 (even though in manual mode it is ignored).

Specifying the Sweeper's Settling Time

For each change in the sweep parameter that takes effect on the instrument the Sweeper waits before recording measurement data in order to allow the measured signal to settle. This behavior is configured by two parameters in the `sweep/settling/` branch: `sweep/settling/time` and `sweep/settling/inaccuracy`.

The `sweep/settling/time` parameter specifies the minimum time in seconds to wait before recording measurement data for that sweep point. This can be used to specify the settling time required by the user's experimental setup before measuring the response in their system.

The `sweep/settling/inaccuracy` parameter is used to derive the settling time to allow for the lock-in amplifier's demodulator filter response to settle following a change of value in the sweep parameter. More precisely, the `sweep/settling/inaccuracy` parameter specifies the amount of settling time as the time required to attain the specified remaining proportion [1e-13, 0.1] of an incoming step function. Based upon the value of `sweep/settling/inaccuracy` and the demodulator filter order, the number of demodulator filter time constants to wait is calculated and written to `sweep/settling/tc` (upon calling the module's `execute()` command) which can then be read back by the user. See [Table 2.1](#) for recommended values of `sweep/settling/inaccuracy`. The relationship between `sweep/settling/inaccuracy` and `sweep/settling/tc` is plotted in [Figure 2.2](#).

The actual amount of time the Sweeper Module will wait after setting a new sweep parameter value before recording measurement data is defined in [Equation 2.1](#). For a frequency sweep, the `sweep/settling/inaccuracy` parameter will tend to influence the settling time at lower frequencies, whereas `sweep/settling/time` will tend to influence the settling time at higher frequencies.

$$t_s = \max(\text{sweep_settling_tc} \times tc, \text{sweep_settling_time})$$

Equation 2.1. The settling time t_s used by the Sweeper for each measurement point; the amount of time between setting the sweep parameter and recording measurement data is determined by the `sweep/settling/tc` and `sweep/settling/time`.

Note, although it is recommended to use `sweep/settling/inaccuracy`, it is still possible to set the settling time via `sweep/settling/tc` instead of `sweep/settling/inaccuracy` (the parameter applied will be simply the last one that is set by the user).

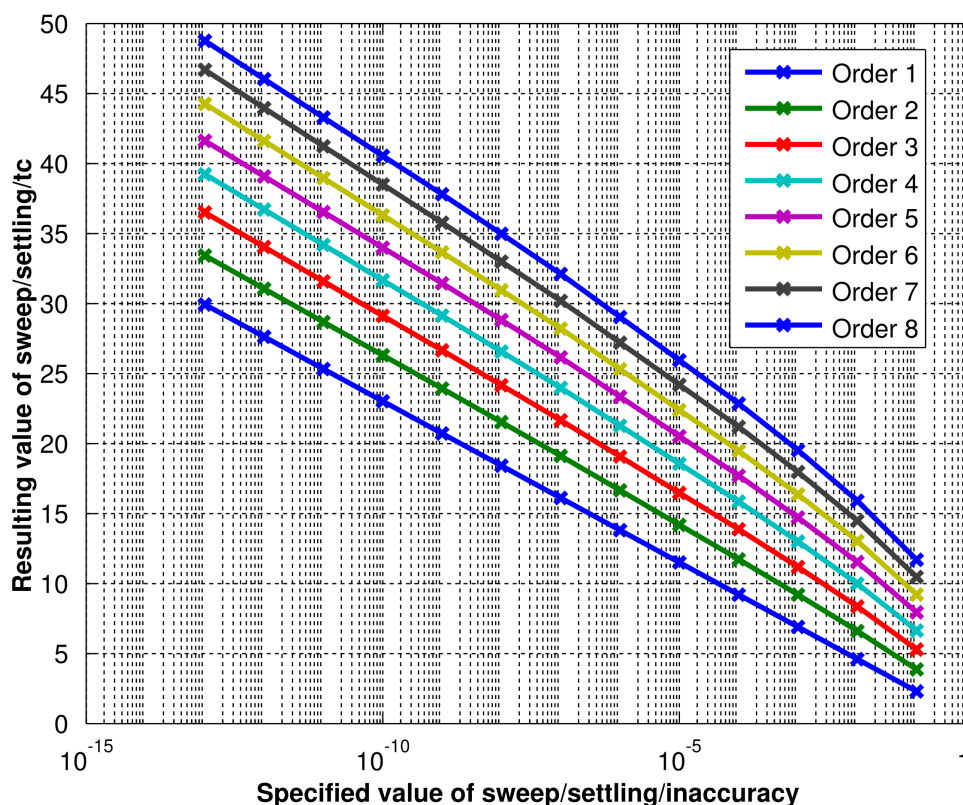


Figure 2.2. A plot showing the values of the Sweeper's sweep/settling/tc as calculated from sweep/settling/inaccuracy parameter and their dependency on demodulator order filter.

Specifying which Data to Measure

Which measurement data is actually returned by the Sweeper's read command is configured by subscribing to `node paths` using the Sweeper Module's subscribe command.

Specifying how the Measurement Data is Averaged

One Sweeper measurement point is obtained by averaging recorded data which is configured via the parameters in the `sweep/averaging/` branch.

The `sweep/averaging/tc` parameter specifies the minimum time window in factors of demodulator filter time constants during which samples will be recorded in order to average for one returned sweeper measurement point. The `sweep/averaging/sample` parameter specifies the minimum number of data samples that should be recorded and used for the average. The Sweeper takes both these settings into account for the measurement point's average according to Equation 2.2.

$$N = \max(\text{sweep_averaging_tc} \times \text{tc} \times \text{sampling_rate}, \text{sweep_averaging_sample})$$

Equation 2.2. The number of samples N used to average one sweeper measurement point is determined by the parameters `sweep/averaging/tc` and `sweep/averaging/sample`.

Note, the value of the demodulator filter's time constant may be controlled by the Sweeper depending on the value of `sweep/bandwidthcontrol` and `sweep/bandwidth`, see [above](#), [Controlling how the Sweeper sets the Demodulator's Time Constant](#). For a frequency sweep, the

`sweep/averaging/tc` parameter will tend to influence the number of samples recorded at lower frequencies, whereas `sweep/averaging/sample` will influence averaging behavior at higher frequencies.

An Explanation of Settling and Averaging Times in a Frequency Sweep

Figure 2.3 shows which demodulator samples are used in order to calculate an averaged measurement point in a frequency sweep. This explanation of the Sweeper's parameters is specific to the following commonly-used Sweeper settings:

- `sweep/gridnode` is set to an oscillator frequency, e.g., `/dev123/oscs/0/freq`.
- `sweep/bandwidthcontrol` is set to 2, corresponding to automatic bandwidth control, i.e., the Sweeper will set the demodulator's filter bandwidth settings optimally for each frequency used.
- `sweep/scan` is set to 0, corresponding to sequential scan mode for the range of frequency values swept, i.e, the frequency is increasing for each measurement point made.

Each one of the three red segments in the demodulator data correspond to the data used to calculate one single Sweeper measurement point. The light blue bars correspond to the time the sweeper should wait as indicated by `sweep/settling/tc` (this is calculated by the Sweeper Module from the specified the `sweep/settling/inaccuracy` parameter). The purple bars correspond to the time specified by the `sweep/settling/time` parameter. The sweeper will wait for the maximum of these two times according to Equation 2.1. When measuring at lower frequencies the Sweeper sets a smaller demodulator filter bandwidth (due to automatic `sweep/bandwidthcontrol`) corresponding to a larger demodulator filter time constant. Therefore, the `sweep/settling/tc` parameter dominates the settling time used by the Sweeper at low frequencies and at high frequencies the `sweep/settling/time` parameter takes effect. Note, that the light blue bars corresponding to the value of `sweep/settling/tc` get shorter for each measurement point (larger frequency used \rightarrow shorter time constant required), whereas the purple bars corresponding to `sweep/settling/time` stay a constant length for each measurement point. Similarly, the `sweep/averaging/tc` parameter (yellow bars) dominates the Sweeper's averaging behavior at low frequencies, whereas `sweep/averaging/samples` (green bars) specifies the behavior at higher frequencies, see also Equation 2.2.

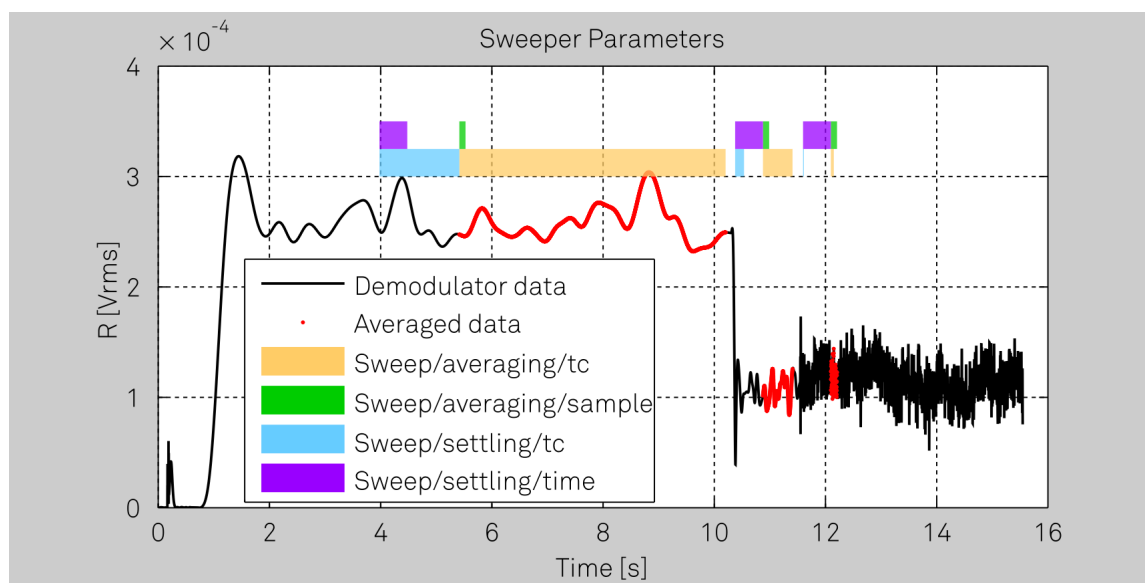


Figure 2.3. Plot demonstrating how the Sweeper records three measurement points from demodulator data when using automatic bandwidth control in a frequency sweep. Please see [An Explanation of Settling and Averaging Times in a Frequency Sweep](#), above, for a detailed explanation.

Average Power and Standard Deviation of the Measured Data

The Sweeper returns measurement data upon calling the Sweeper's `read()` function. This returns not only the averaged measured samples (e.g. `r`) but also their average power (`rwr`) and standard deviation (`rstddev`). In order to obtain reliable values from this statistical data, please ensure that the `sweep/averaging` branch parameters are configured correctly. It's recommended to use at least a value of 12 for `sweep/averaging/sample` to ensure enough values are used to calculate the standard deviation and 5 for `sweep/averaging/tc` in order to prevent aliasing effects from influencing the result.

Table 2.1. Sweeper Parameters

Setting/Path	Type	Unit	Description
<code>sweep/device</code>	byte array	-	The device ID to perform the sweep on, e.g., <code>dev123</code> (compulsory parameter).
<code>sweep/gridnode</code>	byte array	Node	The device parameter (specified by node) to be swept, e.g., <code>"oscs/0/freq"</code> .
<code>sweep/start</code>	double	Many	The start value of the sweep parameter.
<code>sweep/stop</code>	double	Many	The stop value of the sweep parameter.
<code>sweep/samplecount</code>	uint64	-	The number of measurement points to set the sweep on.
<code>sweep/endless</code>	bool	-	Enable Endless mode; run the sweeper continuously.
<code>sweep/remainingtime</code>	double	Seconds	Read only: Reports the remaining time of the current sweep. A valid number is only displayed once the sweeper has been started. An undefined sweep time is indicated as NAN.
<code>sweep/averaging/sample</code>	uint64	Samples	Sets the number of data samples per sweeper parameter point that is considered in the measurement. The maximum of this value and <code>sweep/averaging/tc</code> is taken as the effective calculation time. See Figure 2.3 .
<code>sweep/averaging/tc</code>	double	Seconds	Sets the effective measurement time per sweeper parameter point that is considered in the measurement. The maximum between of this value and <code>sweep/averaging/sample</code> is taken as the effective calculation time. See Figure 2.3 .
<code>sweep/bandwidthcontrol</code>	uint64	-	Specify how the sweeper should specify the bandwidth of each measurement point, Automatic is recommended, in particular for logarithmic sweeps and assures the whole spectrum is covered. 0=Manual (the sweeper module leaves the demodulator bandwidth settings entirely untouched); 1=Fixed (use the value from <code>sweep/bandwidth</code>); 2=Automatic. Note, to use either Fixed or Manual mode, <code>sweep/bandwidth</code> must be set to a value > 0 (even though in manual mode it is ignored).

Setting/Path	Type	Unit	Description
sweep/ bandwidthoverlap	bool	-	If enabled the bandwidth of a sweep point may overlap with the frequency of neighboring sweep points. The effective bandwidth is only limited by the maximal bandwidth setting and omega suppression. As a result, the bandwidth is independent of the number of sweep points. For frequency response analysis bandwidth overlap should be enabled to achieve maximal sweep speed.
sweep/bandwidth	double	Hz	Defines the measurement bandwidth when using Fixed bandwidth mode (sweep/bandwidthcontrol=1), and corresponds to the noise equivalent power bandwidth (NEP).
sweep/order	uint64	-	Defines the filter roll off to use in Fixed bandwidth mode (sweep/bandwidthcontrol=1). Valid values are between 1 (6 dB/octave) and 8 (48 dB/octave).
sweep/maxbandwidth	double	Hz	Specifies the maximum bandwidth used when in Auto bandwidth mode (sweep/bandwidthcontrol=2) (sweep/bandwidthcontrol=2). The default is 1.25 MHz.
sweep/ omegasuppression	double	dB	Damping of omega and 2omega components when in Auto bandwidth mode (sweep/bandwidthcontrol=2). Default is 40dB in favor of sweep speed. Use a higher value for strong offset values or 3omega measurement methods.
sweep/loopcount	uint64	-	The number of sweeps to perform.
sweep/phaseunwrap	bool	-	Enable unwrapping of slowly changing phase evolutions around the +/-180 degree boundary.
sweep/sincfilter	bool	-	Enables the sinc filter if the sweep frequency is below 50 Hz. This will improve the sweep speed at low frequencies as omega components do not need to be suppressed by the normal low pass filter.
sweep/scan	uint64	-	Selects the scanning type: 0=Sequential (incremental scanning from start to stop value, see Figure 2.4); 1=Binary (Non-sequential sweep continues increase of resolution over entire range, see Figure 2.6), 2=Bidirectional (Sequential sweep from Start to Stop value and back to Start again, Figure 2.5), 3=Reverse (reverse sequential scanning from stop to start value).
sweep/settling/time	double	Seconds	Minimum wait time in seconds between setting the new sweep parameter value and the start of the measurement. The maximum between this value and sweep/settling/

Setting/Path	Type	Unit	Description
			t_c is taken as effective settling time. See Figure 2.3 .
sweep/settling/inaccuracy	double	-	Demodulator filter settling inaccuracy defining the wait time between a sweep parameter change and recording of the next sweep point. The settling time is calculated as the time required to attain the specified remaining proportion [1e-13, 0.1] of an incoming step function. Typical inaccuracy values: 10m for highest sweep speed for large signals, 100u for precise amplitude measurements, 100n for precise noise measurements. Depending on the order of the demodulator filter the settling inaccuracy will define the number of filter time constants the sweeper has to wait. The maximum between this value and the settling time is taken as wait time until the next sweep point is recorded. The relationship between sweep/settling/inaccuracy and sweep/settling/tc is plotted in Figure 2.2 .
sweep/settling/tc	double	TC	Minimum wait time in factors of the time constant (TC) between setting the new sweep parameter value and the start of the measurement. This filter settling time is preferably configured via the sweep/settling/inaccuracy (see discussion in Section 2.2.1 and Figure 2.2). The maximum between this value and sweep/settling/time is taken as effective settling time. See Figure 2.3 .
sweep/xmapping	uint64	-	Selects the spacing of the grid used by sweep/gridnode (the sweep parameter): 0=linear and 1=logarithmic distribution of sweep parameter values.
sweep/historylength	uint64		Maximum number of entries stored in the measurement history.
sweep/clearhistory	bool	-	Remove all records from the history list.
sweep/directory	byte array	-	The directory to which sweeper measurements are saved to via <code>save()</code> .
sweep/savepath	byte array	-	This parameter is deprecated, see sweep/directory.
sweep/fileformat	byte array	-	The format of the file for saving sweeper measurements. 0=Matlab, 1=CSV.

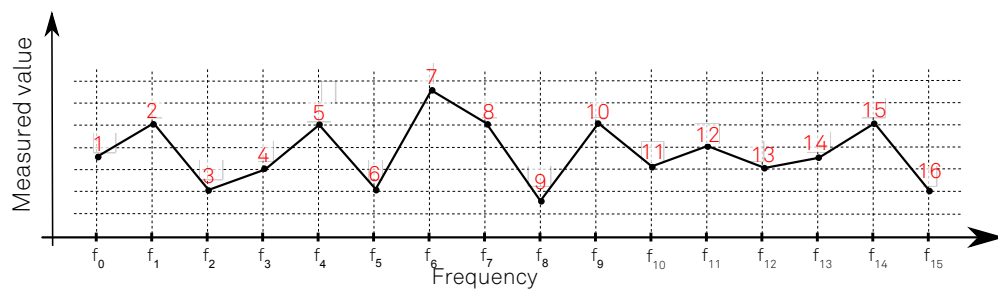


Figure 2.4. Sweeper scanning modes: Sequential ($\text{sweep}/\text{scan} = 0$).

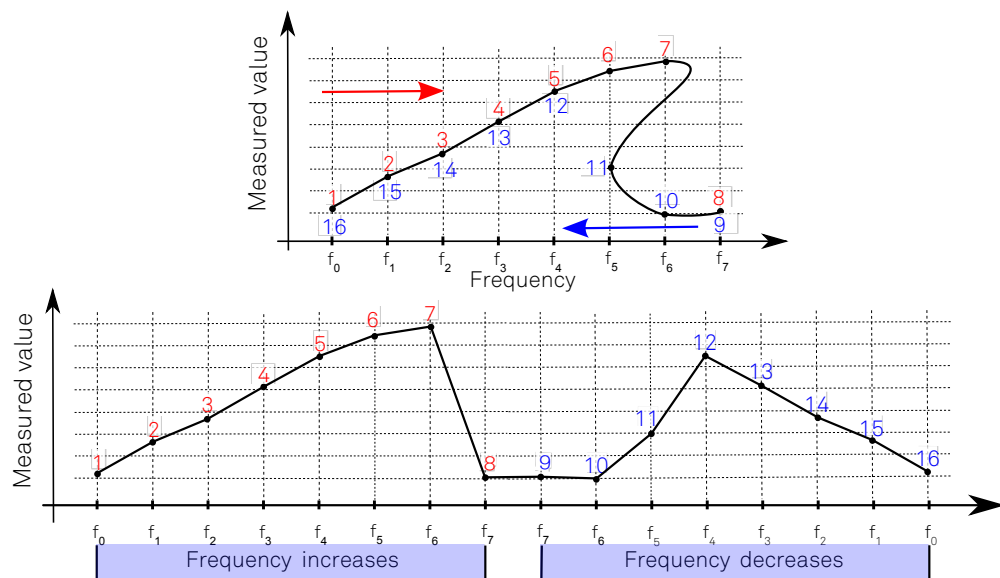


Figure 2.5. Sweeper scanning modes: Bidirectional ($\text{sweep}/\text{scan} = 2$).

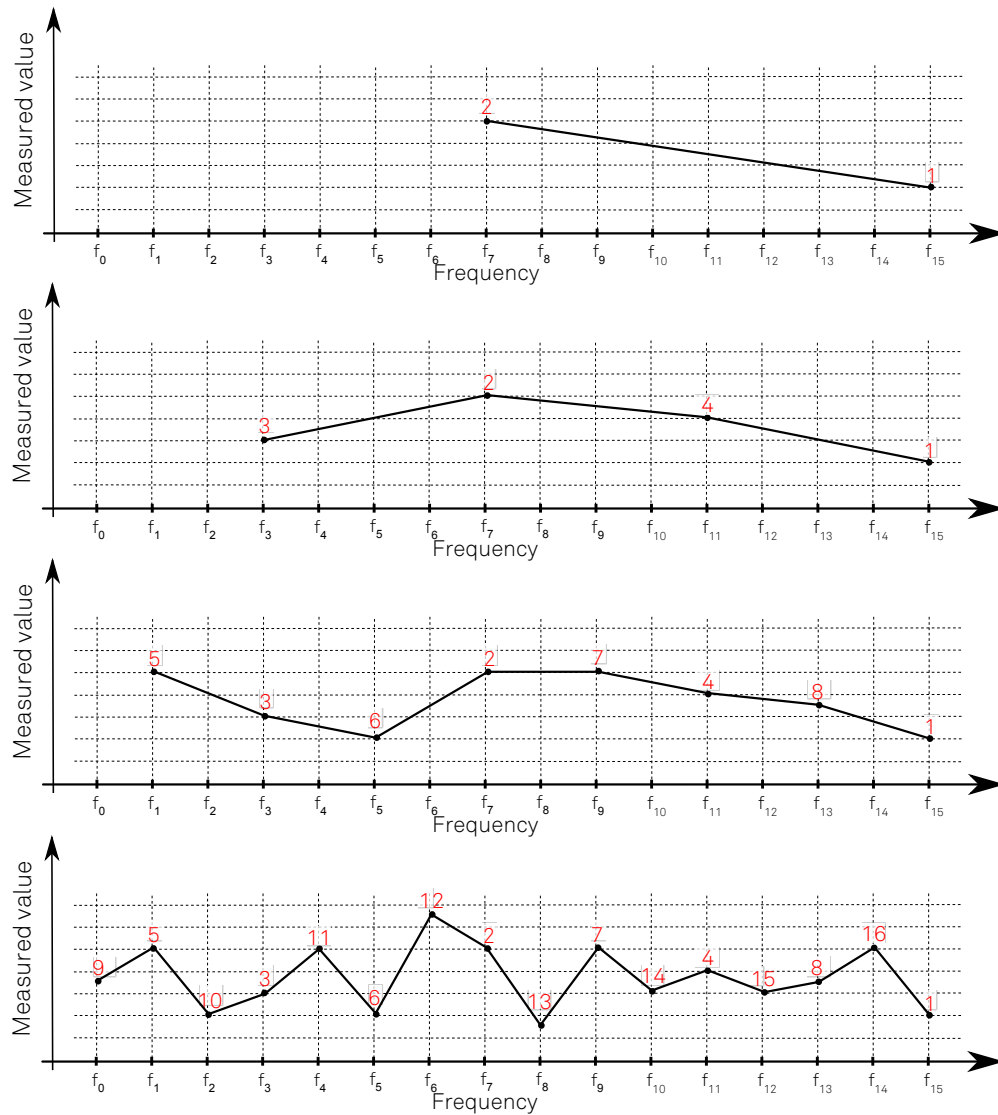


Figure 2.6. Sweeper scanning modes: Binary (sweep/scan = 1).

Table 2.2. Sweeper Output Values

Name	Type	Unit	Description
auxin0	double	Volts	Auxiliary Input 1 value.
auxin1	double	Volts	Auxiliary Input 2 value.
auxin0pwr	double	Volts ²	Average power of Auxiliary Input 1 value.
auxin1pwr	double	Volts ²	Average power of Auxiliary Input 2 value.
auxin0stddev	double	Volts	Standard deviation of Auxiliary Input 1 value.
auxin1stddev	double	Volts	Standard deviation of Auxiliary Input 2 value.
frequency	double	Hz	The oscillator frequency for each measurement point (for a frequency sweep this is the same as grid).
frequencypwr	double	Hz ²	Average power of the oscillator frequency.
frequencystddev	double	Hz	Standard deviation of the oscillator frequency.

Name	Type	Unit	Description
phase	double	Radians	Demodulator phase value.
phasestddev	double	Radians	Standard deviation of demodulator phase value (phase noise).
phaserpwr	double	Radians ²	Average power of demodulator phase value (phase noise).
r	double	VoltsRMS	Demodulator R value.
rstddev	double	VoltsRMS	Standard deviation of demodulator R value.
rpwr	double	Volts ²	Average power of demodulator x value.
x	double	Volts	Demodulator x value.
xstddev	double	Volts	Standard deviation of demodulator x value.
xpwr	double	Volts ²	Average power of demodulator x value.
y	double	Volts	Demodulator y value.
ystddev	double	Volts	Standard deviation of demodulator y value.
ypwr	double	Volts ²	Average power of demodulator y value.
bandwidth	double	Hz	Demodulator filter's bandwidth as calculated from sweep/tc (if performing a frequency sweep).
bandwidthmode	integer	-	The value of the sweep/bandwidthcontrol used for the sweep.
count	integer	-	The number of measurement points actually used by the sweeper when averaging the data. This depends on the values of the parameters in the sweep/averaging/branch.
grid	double	Many	Values of sweeping setting (frequency values at which demodulator samples where recorded).
flags	integer	-	Reserved for future use.
settling	double	Seconds	The waiting time for each measurement point.
samplecount	uint64	-	The number of swept measurement points (the value of sweep/samplecount).
sampleformat	integer	-	Reserved for future use.
sweepmode	integer	-	The value of the sweep/scan used for the sweep.
tc	double	Seconds	Demodulator's filter time constant as set for each measurement point.
tcmeas	double	Seconds	Reserved for future use.
timestamp	uint64	Ticks	A timestamp that gets updated each time a new measurement point has been recorded by the sweeper (divide by the device's clockbase to obtain seconds). It is not part of the sweeper's measurement data and only relevant for intermediate reads of sweeper data (before the current sweep has finished).

Name	Type	Unit	Description
settimestamp	uint64	Ticks	The timestamp at which we verify that the frequency for the current measurement point was set on the device (by reading back demodulator data).
nexttimestamp	uint64	Ticks	The timestamp at which we can obtain the data for that measurement point, i.e., <code>nexttimestamp - settimestamp</code> corresponds roughly to the demodulator filter settling time.

2.3. zoomFFT Module

The zoomFFT Module corresponds to the Spectrum Tab of the LabOne User Interface. It allows the user to perform Fast Fourier Transforms (FFT) on a specified demodulator's output.

See [Table 2.3](#) for the input parameters to configure the ZoomFFT Module and [Table 2.4](#) for a description of the ZoomFFT's outputs.

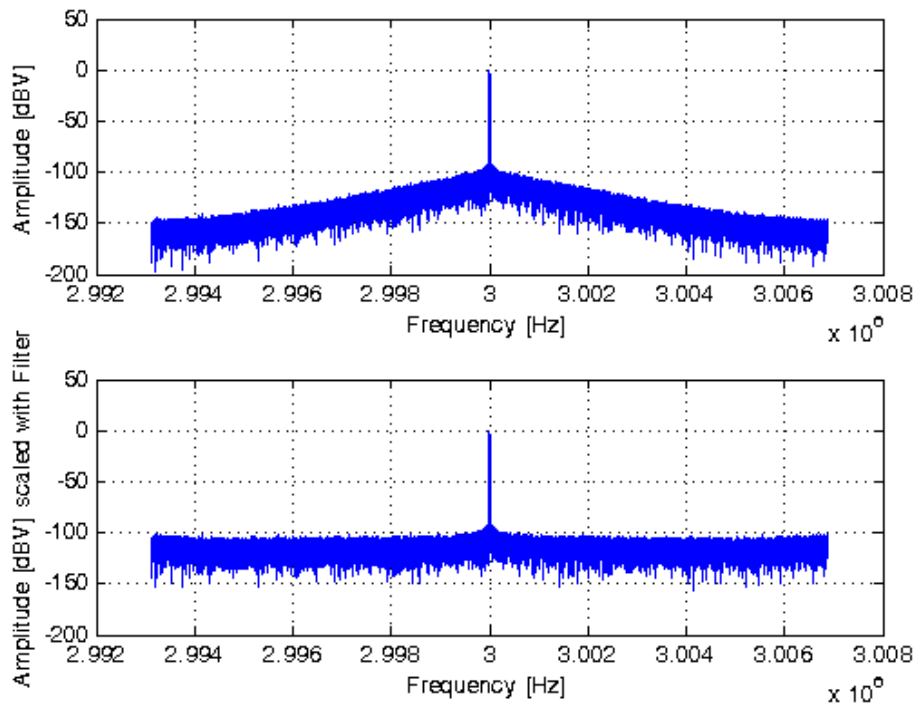


Figure 2.7. A plot of an FFT created by one of the LabOne Matlab API examples.

Table 2.3. ZoomFFT Input Parameters

Setting/Path	Type	Unit	Description
zoomFFT/device	byte array	-	The device ID to perform the FFT on, e.g., dev123 (compulsory parameter).
zoomFFT/absolute	bool	-	Shifts the frequencies so that the center frequency becomes the demodulation frequency rather than 0 Hz.
zoomFFT/bit	uint64	-	Number of lines of the FFT spectrum (powers of 2). Increasing the bits increases the frequency resolution of the spectrum.
zoomFFT/endless	bool	-	Enable Endless mode; run the zoomFFT continuously.
zoomFFT/loopcount	uint64	-	The number of FFTs to perform.
zoomFFT/mode	uint64	-	Select the source signal for the FFT. 0=FFT(x+iy), 1=FFT(R), 2=FFT(phase), 3=FFT(Freq)
zoomFFT/overlap	double	-	Overlap of the demodulator data used for the FFT. Use 0 for no overlap and 0.99 for maximal overlap.

Setting/Path	Type	Unit	Description
zoomFFT/settling/tc	double	TC	Minimum wait time in factors of the time constant (TC) before starting the measurement. The maximum between this value and zoomFFT/settling/time is taken as effective settling time.
zoomFFT/settling/time	double	Seconds	Minimum wait time in seconds before starting the measurement. The maximum between this value and zoomFFT/settling/tc is taken as effective settling time.
zoomFFT/window	uint64	-	The type of FFT window to use. 0=Rectangular, 1=Hann, 2=Hamming, 3=Blackman Harris.

Table 2.4. ZoomFFT Output Values

Name	Type	Unit	Description
x	double	Volts	The real part, x, of the complex FFT result.
y	double	Volts	The imaginary part, y, of the complex FFT result.
r	double	VoltsRMS	The absolute value, R, of the complex FFT result.
timestamp	uint64	Ticks	Demodulator timestamp of the measurement (divide by the device's clockbase to obtain seconds)
center	double	Hz	The center frequency (corresponds to the demodulation frequency).
rate	double	-	Sampling rate of the demodulator.
filter	double	-	The filter envelope; the filter compensation value for each gridnode.
bandwidth	double	Hz	The bandwidth of the demodulator
grid	double	Hz	The frequency grid.
nenbw	double	-	The normalized equivalent noise bandwidth.
resolution	double	Hz	FFT resolution: Spectral resolution defined by the reciprocal acquisition time (sample rate, number of samples recorded).
aliasingreject	double	dB	How much damping is present at the border of your spectrum.

2.4. Software Trigger (Recorder) Module

The Recorder Module corresponds to the Software Trigger Tab of the LabOne User Interface. It allows the user to record bursts of instrument data based upon pre-defined trigger criteria similar to that of a laboratory oscilloscope, see [Figure 2.8](#) for an example. The types of trigger available are listed in [Table 2.5](#).

Table 2.5. Overview of the trigger types available in the Software Trigger Module.

Trigger Type	Description	trigger/N/type
Manual	For simple recording.	0
Edge	Edge trigger with level hysteresis and noise rejection, see Figure 2.9 .	1
Digital	Digital trigger with bit masking.	2
Pulse	Pulse width trigger with level hysteresis and noise reduction, see Figure 2.10 and Figure 2.11 .	3
Tracking (edge or pulse)	Level tracking trigger to compensate signal drift, see Figure 2.12 .	4
Hardware Trigger	UHFLI and MFLI only. Trigger on one of the instrument's hardware trigger channels.	6

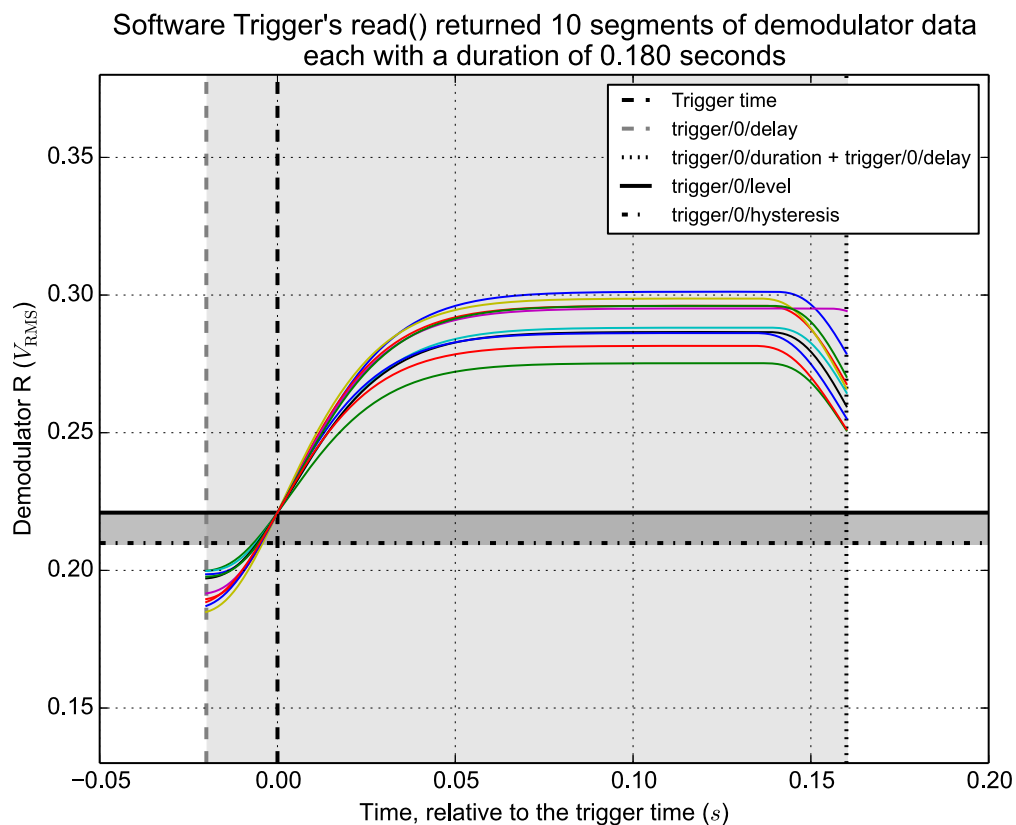


Figure 2.8. The plot produced by `example_swtrigger_edge.py`, an example distributed with the LabOne Python API. The plot shows 10 bursts of data from a single demodulator; each burst was recorded when the demodulator's R value exceeded a specified threshold using a positive edge trigger. See [Section 4.2.3](#) for help getting started with the Python examples.

See Table 2.6 for the input parameters to configure the Software Trigger's Module. Note that some parameters effect all triggers, e.g., `trigger/endless`, whereas some are configured on a per-trigger basis, e.g., `trigger/N/duration`, where N is the index of the trigger, starting at zero. The data output when using the Software Trigger's `read` command has the same format as returned by `ziCore's poll` command.

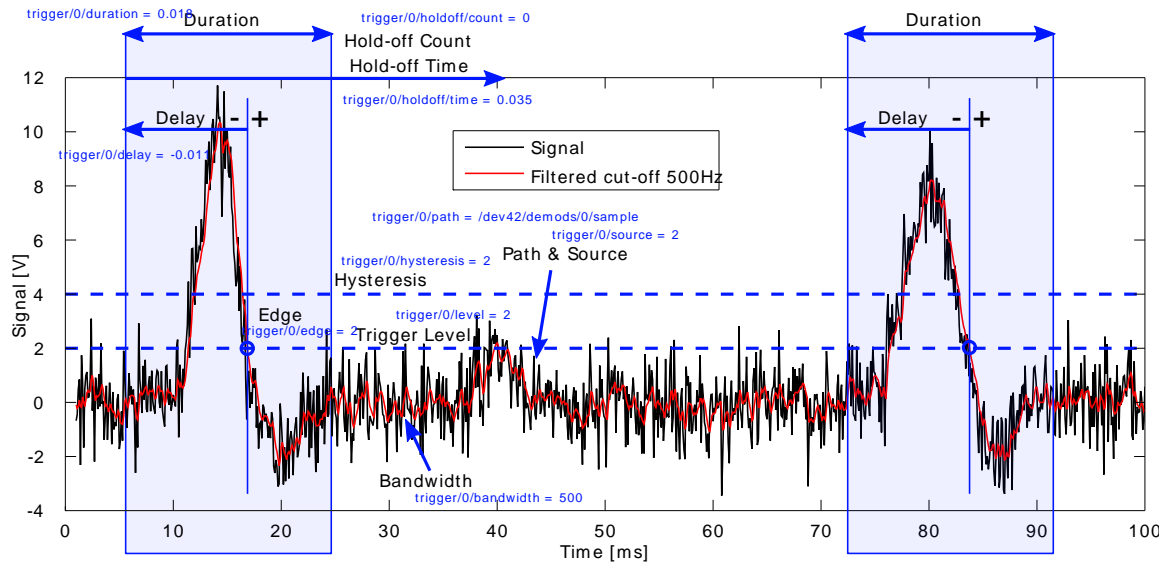


Figure 2.9. Explanation of the Software Trigger Module's parameters for an Edge Trigger.

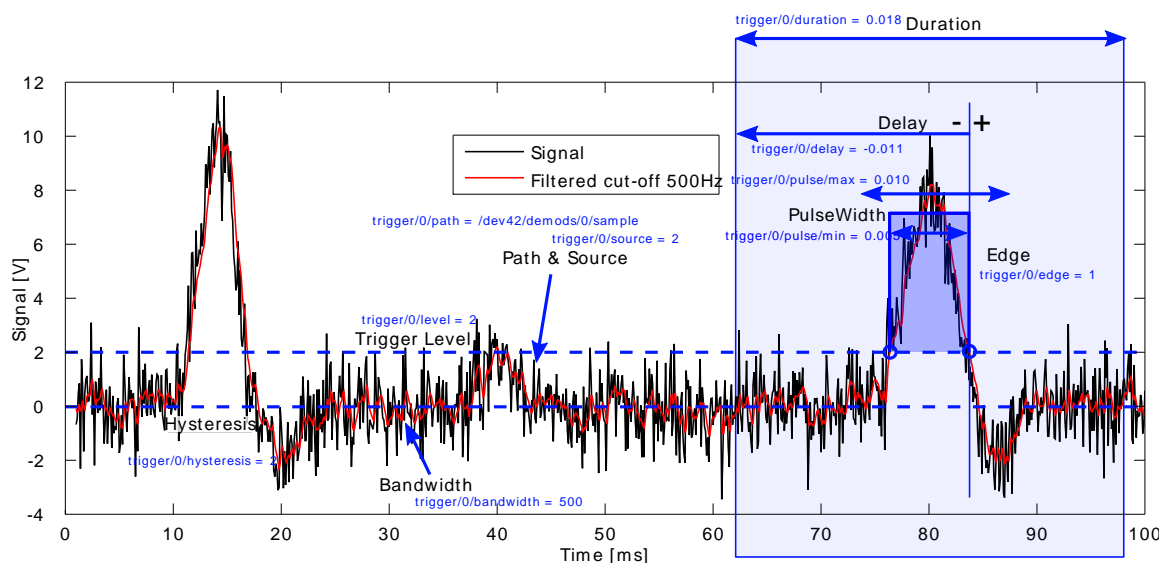


Figure 2.10. Explanation of the Software Trigger Module's parameters for a positive Pulse Trigger.

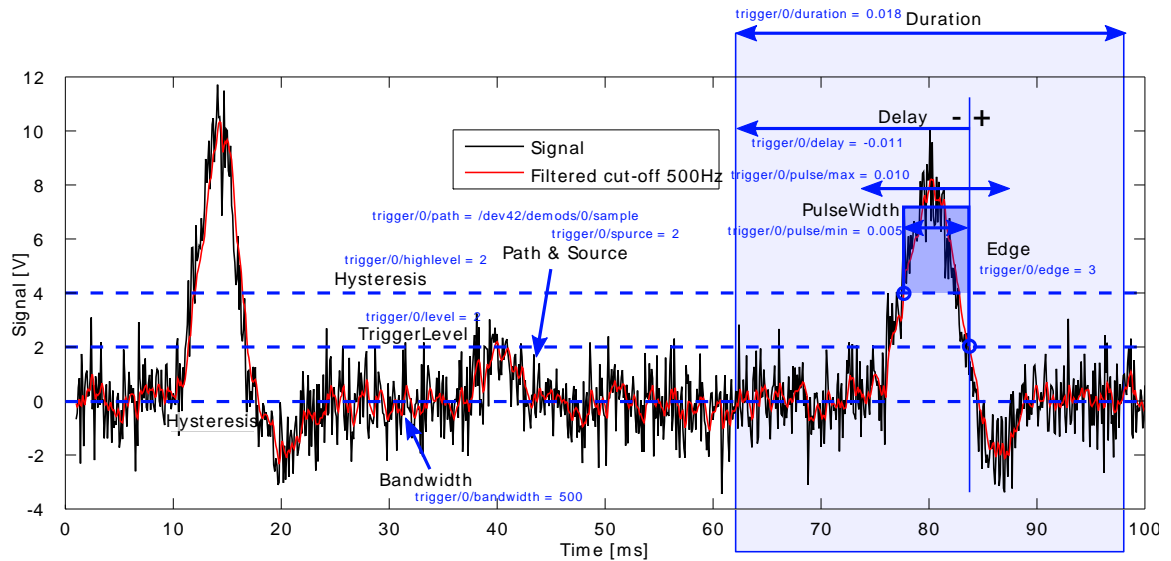


Figure 2.11. Explanation of the Software Trigger parameters for a positive or negative Pulse Trigger.

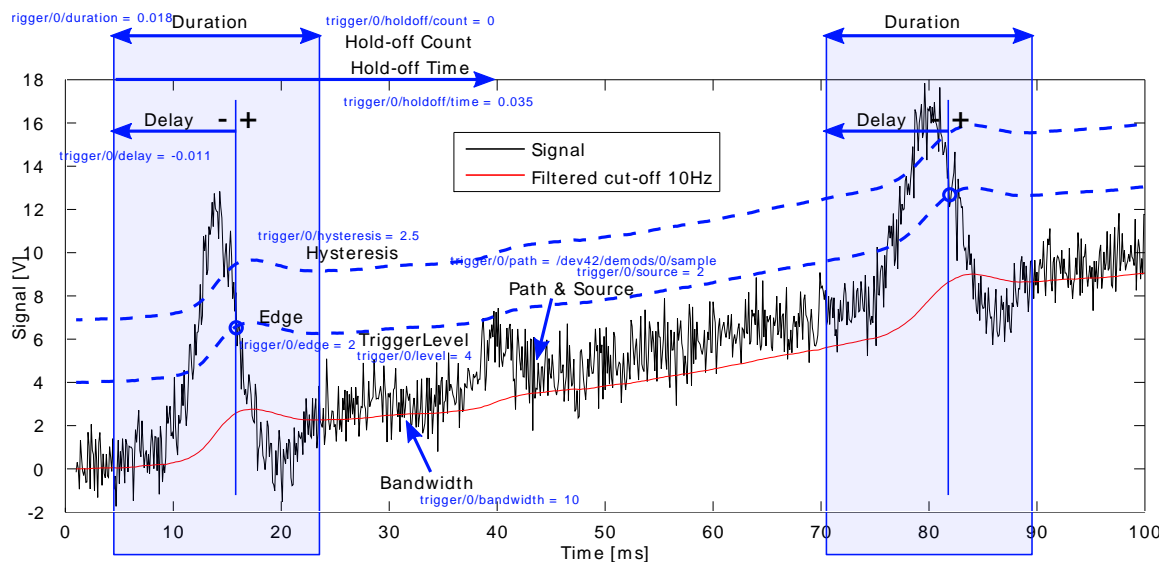


Figure 2.12. Explanation of the Software Trigger Module's parameters for a Tracking Trigger.

Table 2.6. Software Trigger Input Parameters.

Setting/Path	Type	Unit	Description
trigger/device	byte array	-	The device ID to execute the software trigger, e.g., dev123 (compulsory parameter).
trigger/buffersize	double	Seconds	Set the buffer size of the trigger object. The recommended buffer size is $2 \times \text{trigger}/N/\text{duration}$.
trigger/flags	uint64	-	Define the SW Trigger's behaviour if sampleloss is encountered: Fill holes (=0x01), align data that contains a timestamp (=0x02), throw EOFError if sampleloss is detected.

Setting/Path	Type	Unit	Description
trigger/endless	uint64	-	Enable endless triggering 1=enable; 0=disable.
trigger/forcetrigger	uint64	-	Force a trigger.
trigger/filename	byte array	-	This parameter is deprecated. If specified, i.e. not empty, it enables automatic saving of data in single triggering mode (<code>trigger/endless = 0</code>).
trigger/savepath	byte array	-	The directory where files are saved when saving data.
trigger/fileformat	byte array	-	The format of the file for saving data. 0=Matlab, 1=CSV.
trigger/historylength	uint64	-	Maximum number of entries stored in the measurement history.
trigger/clearhistory	uint64	-	Clear the measurement history
trigger/triggered	uint64	-	Has the software trigger triggered? 1=Yes, 0=No (read only).
trigger/N/bandwidth	double	Hz	Only for Tracking Triggers. The bandwidth used in the calculation of the exponential running average of the source signal.
trigger/N/bitmask	uint64	-	Only for Digital triggers. Specify the bitmask used with <code>trigger/N/bits</code> . The trigger value is bits AND bit mask (bitwise).
trigger/N/bits	uint64	-	Only for Digital triggers. Specify the bits used for the Digital trigger value. The trigger value is bits AND bit mask (bitwise)
trigger/N/count	uint64	-	The number of triggers to save.
trigger/N/delay	uint64	Seconds	The amount of time to record data before the trigger was activated, Delay: Time delay of trigger frame position (left side) relative to the trigger edge. For delays smaller than 0, trigger edge inside trigger frame (pre trigger). For delays greater than 0, trigger edge before trigger frame (post trigger), see Figure 2.9 .
trigger/N/duration	double	Seconds	The length of time to record data for, see Figure 2.9 .
trigger/N/edge	uint64	-	Define on which signal edge to trigger. Triggers when the trigger input signal crosses the trigger level from either low to high (edge=1), high to low (edge=2) or both (edge=3). Used for Trigger Type edge, pulse, tracking edge and tracking pulse. In the case of pulse trigger, the value specifies a positive (edge=1) or negative (edge=2) pulse relative to the trigger level (edge=3 specifies either positive or negative).
trigger/N/findlevel	uint64	-	Automatically find the value of <code>trigger/N/level</code> based on the current signal value.

Setting/Path	Type	Unit	Description
trigger/N/level	uint64	Many	Specify the main trigger level value.
trigger/N/holdoff/count	uint64	-	The holdoff count, the number of skipped triggers until the next trigger is recorded again.
trigger/N/holdoff/time	double	Seconds	The holdoff time, the amount of time until the next trigger is recorded again. A hold off time smaller than @trigger/0/duration@ will produce overlapping trigger frames.
trigger/N/hysteresis	double	Many	Specify the hysteresis value (the trigger is re-armed after the signal exceeds trigger/N/level and then falls below trigger/N/hysteresis, if using positive edge).
trigger/N/pulse/max	double	-	Only for Pulse triggers: The maximum pulse width to trigger on. See Figure 2.10 .
trigger/N/pulse/min	double	-	Only for Pulse triggers: The minimum pulse width to trigger on. See Figure 2.10 .
trigger/N/retrigger	uint64	-	1=enable, 0=disable. Enable to allow re-triggering within one trigger duration. If enabled continue recording data in one segment if another trigger comes within the previous trigger's duration. If disabled the triggers will be recorded as separate events.
trigger/N/triggernode	byte array	-	<p>Path and signal of the node that should be used for triggering, separated by a dot (.), e.g. /devN/demods/0/sample.x.</p> <p>SAMPLE.X Demodulator X value SAMPLE.Y Demodulator Y value SAMPLE.R Demodulator Magnitude SAMPLE.THETA Demodulator Phase SAMPLE.AUXIN0 Auxiliary Input 1 value SAMPLE.AUXIN1 Auxiliary Input 2 value SAMPLE.DIO Digital I/O value</p> <p>Over HW Trigger paths may also be specified (device-class dependent). Overrides values from trigger/0/path and trigger/0/source.</p>
trigger/N/type	uint64	-	The trigger type, see Table 2.5
trigger/0/grid/mode	int	-	Enable grid mode. In grid mode a matrix instead of a vector is returned by read(). Each trigger becomes a row in the matrix and each trigger's data is interpolated onto a new grid defined by the number of columns: 0: Disable, 1: Enable grid mode with nearest neighbour interpolation, 2: Enable grid mode with linear interpolation.
trigger/0/grid/operation	int	-	If running in endless mode, either replace or average the data in the grid's matrix.

Setting/Path	Type	Unit	Description
trigger/0/grid/cols	int	-	Specify the number of columns in the grid's matrix. The data from each row is interpolated onto a grid with the specified number of columns.
trigger/0/grid/rows	int	-	Specify the number of rows in the grid's matrix. Each row is the data recorded from one trigger interpolated onto the columns.
trigger/0/grid/direction	int	-	The direction to organize data in the grid's matrix: 0: Forward. 1: Reverse. 2: Bidirectional. Forward - the data in each row is ordered chronologically, e.g., the first data point in each row corresponds to the first timestamp in the trigger data. Reverse - the data in each row is ordered reverse chronologically, e.g., the first data point in each row corresponds to the last timestamp in the trigger data. Bidirectional - the ordering of the data alternates between Forward and Backward ordering from row-to-row, the first row is Forward ordered.
trigger/N/path	byte array	-	This parameter is deprecated, see the trigger/N/triggernode parameter.
trigger/N/source	uint64	-	This parameter is deprecated, see the trigger/N/triggernode parameter.
trigger/N/hwtrigsource	uint64	-	This parameter is deprecated, see the trigger/N/triggernode parameter.

Note

For the pulse trigger type, there is a subtle difference between the way the trigger level and the hysteresis are used for positive/negative pulse triggering (trigger/N/edge= 1 or 2) and both (trigger/N/edge= 3). The difference can be seen in [Figure 2.10](#) and [Figure 2.11](#).

2.5. Device Settings Module

The Device Settings Module provides functionality for saving and loading device settings to and from file. The file is saved in [XML](#) format.

In general, users are recommended to use the utility functions provided by the APIs instead of using the Device Settings module directly. The Matlab API provides `ziSaveSettings()` and `ziLoadSettings()` and the Python API provides `zhinst.utils.save_settings()` and `zhinst.utils.load_settings`. These are convenient wrappers to the Device Settings module for loading settings asynchronously, i.e., these functions block until loading or saving has completed, the desired behavior in most cases. Advanced users can use the Device Settings module directly if they need to implement loading or saving a synchronously (non-blocking).

See [Table 2.7](#) for the input parameters to configure the Device Settings Module.

Table 2.7. Device Settings Input Parameters

Setting/Path	Type	Description
<code>deviceSettings/device</code>	byte array	The device ID to save the settings for, e.g., <code>dev123</code> (compulsory parameter).
<code>deviceSettings/command</code>	byte array	The command to issue: "load" (load settings from file); "save" (read device settings and save to file) or "read" (just read the device settings) (compulsory parameter).
<code>deviceSettings/filename</code>	byte array	The name of the file to load or save to.
<code>deviceSettings/path</code>	byte array	The path containing the file to load from or save to.

Table 2.8. Device Settings Parameters for use only by the LabOne Web Server.

Setting/Path	Type	Description
<code>deviceSettings/throwonerror</code>	uint64	Throw an exception if there was error executing the command.
<code>deviceSettings/errortext</code>	byte array	The error text used in error messages.
<code>deviceSettings/finished</code>	uint64	The status of the command (read-only).

2.6. PLL Advisor Module

The PLL Advisor Module corresponds to the PLL Advisor section of the LabOne User Interface PLL tab. The PLL Advisor is a mathematical model of the PLL incorporated in the instrument and provides a convenient way to tune parameters to obtain an optimal feedback loop performance for the desired application.

Note

Note the PLL Advisor Module is only available for UHF Lock-in Amplifiers.

Table 2.9. PLL Advisor Parameters.

Setting/Path	Type	Unit	Description
<code>pllAdvisor/bode</code>	struct	-	Output parameter. Contains the resulting bode plot of the PLL simulation.
<code>pllAdvisor/calculate</code>	uint64	-	Issues a command for the PLL Advisor to calculate values. Set the value to 1 to start the calculation.
<code>pllAdvisor/center</code>	double	Hz	Center frequency of the PLL oscillator. The PLL frequency shift is relative to this center frequency.
<code>pllAdvisor/d</code>	double	Hz/deg s	The PID differential gain.
<code>pllAdvisor/demodbw</code>	double	Hz	The demodulator bandwidth to use for the PLL loop filter.
<code>pllAdvisor/i</code>	double	Hz/deg/s	The PID integral gain
<code>pllAdvisor/mode</code>	uint64	-	Sets the PLL operating mode. Currently only open-loop mode is supported.
<code>pllAdvisor/order</code>	uint64	-	Demodulator filter order to use for the PLL loop filter.
<code>pllAdvisor/p</code>	double	Hz/deg	The PID proportional gain.
<code>pllAdvisor/pllbw</code>	double	Hz	The demodulator bandwidth to use for the PLL loop filter.
<code>pllAdvisor/pm</code>	double	deg	Output parameter. Simulated phase margin of the PLL with the current settings. The phase margin should be greater than 45 deg and preferably greater than 65 deg for stable conditions.
<code>pllAdvisor/pmfreq</code>	double	-	Output parameter. Simulated phase margin frequency.
<code>pllAdvisor/q</code>	double	-	Quality factor. Currently not used.
<code>pllAdvisor/rate</code>	double	Hz	PLL Advisor sampling rate of the PLL control loop.
<code>pllAdvisor/stable</code>	bool	-	Output parameter. When 1, the PLL Advisor found a stable solution with the given settings. When 0, revise your settings and rerun the PLL Advisor.

Setting/Path	Type	Unit	Description
pllAdvisor/targetbw	double	Hz	Requested PLL bandwidth. Higher frequencies may need manual tuning.
pllAdvisor/targetfail	bool	-	Output parameter. 1 indicates the simulated PLL BW is smaller than the Target BW.

2.7. Tips and Tricks

Use the LabOne User Interface's Command Log to start programming

If you use the LabOne User Interface to perform a measurement, you can obtain the commands sent to your instrument in the "Command Log" by clicking the "Show Log" button in the status bar at the bottom of the User Interface. Be sure to set the "Log Format" of the Command Log in the "User Interface" section of the Config Tab first: The log is available in Matlab and Python formats and can be used as a starting point for your own custom program.

Use the included examples to get started programming

Both the [LabOne Matlab API](#) and the [LabOne Python API](#) come with examples to help you get started programming. In particular, both APIs have at least one example for each of the [ziCore modules](#).

Load LabOne User Interface settings files from the APIs.

The [XML](#) files used for device settings can be loaded and saved from the LabOne User Interface or from any of the [ziCore](#)-based APIs. This means that an instrument can be conveniently configured via the LabOne User Interface and then its settings saved to file. This settings file can then be loaded via an API in order to configure an instrument for a script. See the [Section 2.5](#).

Use the API's logging capabilities.

All of the LabOne APIs can write a log which can contain useful debugging or status information. See the relevant section in the API's chapter for more details:

- [Enabling Logging in the LabOne Matlab API](#),
- [Enabling Logging in the LabOne Python API](#),
- [Error Handling and Logging in the LabOne C API](#).

Part II. LabOne APIs

This part of the Programming Manual documents language-specific installation and usage for each of the LabOne APIs. For details of common functionality and features that are shared by all the LabOne APIs please refer to [Part I](#).

Refer to:

- [Chapter 3](#) for the [LabOne Matlab API](#) (`ziDAQ`).
- [Chapter 4](#) for the [LabOne Python API](#) (`ziPython`).
- [Chapter 5](#) for the [LabOne LabVIEW API](#).
- [Chapter 6](#) for the [LabOne C API](#) (`ziAPI`).

Chapter 3. Matlab Programming

The Mathworks' numerical computing environment [Matlab®](#) has powerful tools for data analysis and visualization that can be used to create graphical user interfaces or automatically generate reports of experimental results in various formats. LabOne's Matlab API, also known as `ziDAQ`, "Zurich Instruments Data Acquisition", enables the user to stream data from their instrument directly into Matlab allowing them to take full advantage of this powerful environment.

This chapter aims to help you get started using Zurich Instruments LabOne's Matlab API, `ziDAQ`, to control your instrument, please refer to:

- [Section 3.1](#) for help [Installing the LabOne Matlab API](#).
- [Section 3.2](#) for help [Getting Started with the LabOne Matlab API](#) and [Running the Examples](#).
- [Section 3.3](#) for some [LabOne Matlab API Tips and Tricks](#).
- [Section 3.4](#) for help [Troubleshooting the LabOne Matlab API](#).
- [Section 3.5](#) for [LabOne Matlab API \(ziDAQ\) Command Reference](#).

Note

This section and the provided examples are no substitute for a Matlab tutorial. See either Mathworks' online [Documentation Center](#) or one of the many online resources, for example, the [Matlab Programming Wikibook](#) for help to get started programming with Matlab.

3.1. Installing the LabOne Matlab API

3.1.1. Requirements

To use LabOne's Matlab API, `ziDAQ`, a Matlab installation and license on either Windows or Linux is required. On Windows, Matlab R2009b (or newer) is required, both 32-bit and 64-bit platforms are supported. On Linux, Matlab R2014b (or newer) and a 64-bit platform is required. No additional Matlab Toolboxes are required to use `ziDAQ`.

The LabOne Matlab API `ziDAQ` is included in a standard LabOne installation and is also available as a separate package (see below, [Separate Matlab Package](#)). No installation as such is required, only a few configuration steps must be performed to use `ziDAQ` in Matlab. Both the main LabOne installer and the separate LabOne Matlab API package are available from Zurich Instruments' [download page](#).

Separate Matlab Package

The separate Matlab API package should be used if you would like to:

1. Use the Matlab API on Mac OS X (the main LabOne installer is not available for Mac OS X).
2. Use the Matlab API to work with an instrument remotely (i.e., on a separate PC from where the Data Server is running) and you do not require a full LabOne installation. This is the case, for example, with MF Instruments.
3. Use the Matlab API on a PC where you do not have administrator rights.

3.1.2. Windows or Linux

No additional installation steps are required to use `ziDAQ` on either Windows or Linux; it's only necessary to add the folder containing LabOne's Matlab API library to Matlab's search path. This is done as following:

1. Start Matlab and either set the "Current Folder" (current working directory) to the Matlab API folder in your LabOne installation or the extracted zip archive of the separate Matlab API package (see above, [Separate Matlab Package](#)) as appropriate.

If using a LabOne installation on Windows this is typically:

```
C:\Program Files\Zurich Instruments\LabOne\API\MATLAB2012\
```

and on Linux this is the location where you unpacked the LabOne `.tar.gz` file:

```
[PATH]/LabOne64/API/MATLAB2012/
```

2. In the Matlab Command Window, run the Matlab script `ziAddPath` located in the `MATLAB2012` directory:

```
>> ziAddPath;
```

On Windows (similar for Linux) you should see the following output in Matlab's Command Window:

```
Added ziDAQ's Driver, Utilities and Examples directories to Matlab's path  
for this session.
```

To make this configuration persistent across Matlab sessions either:

1. Run the 'pathtool' command in the Matlab Command Window and add the following paths WITH SUBFOLDERS to the Matlab search path:

```
C:\Program Files\Zurich Instruments\LabOne\API\MATLAB2012\
```

or

2. Add the following line to your Matlab startup.m file:

```
run('C:\Program Files\Zurich Instruments\LabOne\API\MATLAB2012\ziAddPath');
```

This is sufficient configuration if you would only like to use ziDAQ in the current Matlab session.

3. To make this configuration persistent between Matlab sessions do either one of the next two steps (as also indicated by the output of ziAddPath):
 - a. Run the pathtool and click "Add with Subfolders". Browse to the "MATLAB2012" directory that was located above in Step 1 and click "OK".
 - b. Edit your startup.m to contain the line indicated in the output from Step 2 above. For more help on Matlab's startup.m file, type the following in Matlab's Command Window:

```
>> docsearch('startup.m')
```

4. Verify your Matlab configuration as described in [Section 3.1.3](#).

3.1.3. Verifying Successful Matlab Configuration

In order to verify that Matlab is correctly configured to use ziDAQ please perform the following steps:

1. Ensure that the correct Data Server is running for your HF2 or UHF Instrument (the Data Server on MF Instruments starts when the device is powered on). The quickest way to check is to start the User Interface for your device, see [Section 1.1](#) for more details.
2. Proceed either of the following two ways:
 - a. The easiest way to verify correct configuration is run one of the Matlab API's examples. In the Matlab command Window run, for example, example_poll with your device ID as the input argument:

```
>> example_poll('dev123'); % Replace with your device ID.
```

If this fails, please try issuing the connect command, as described in the next method.

- b. If a device is not currently available, correct Matlab API configuration can be checked by initializing a API session to the Data Server without device communication.

An API session with the Data Server is created using ziDAQ's connect (the port specifies which Data Server to connect to on the localhost) cf. [Section 1.3.1](#)). In the Matlab command window type one of the following:

```
■ >> ziDAQ('connect', 'localhost', 8005) % 8005 for HF2 Series
```

```
■ >> ziDAQ('connect', 'localhost', 8004, 5) % 8004 for UHFLI
```

```
■ >> ziDAQ('connect', mf-hostname, 8004, 5) % 8004 for MFLI (see below)
```

Note, using 'localhost' above assumes that the Data Server is running on the same computer from which you are using Matlab. See [Section 1.3.1](#) for information about port choice and connecting to the Data Server. For MFLI instruments the hostname/IP address

of the MFLI instrument must be provided (the value of `mf-hostname`), see [Section 1.3.1](#) and the Getting Started chapter of the MFLI User Manual for more information.

3. If no error is reported then Matlab is correctly configured to use `ziDAQ` - congratulations! Otherwise, please try the steps listed in [Troubleshooting the LabOne Matlab API](#).

3.2. Getting Started with the LabOne Matlab API

This section introduces the user to the LabOne Matlab API.

3.2.1. Contents of the LabOne Matlab API

Alongside the driver for interfacing with your Zurich Instruments device, the LabOne Matlab API includes many files for documentation, utility functions and examples. See the `Contents.m` file located in a LabOne Matlab API directory (see Step 1 in [Section 3.1.2](#) for its typical location) for a description of the API's sub-folders and files. Run the command:

```
>> doc('Contents')
```

in the Matlab Command Window in the LabOne Matlab API directory to access the following contents interactively in Matlab.

```
% ziDAQ : The LabOne Matlab API for interfacing with Zurich Instruments Devices
%
% FILES
%   ziAddPath - add the LabOne Matlab API drivers, utilities and examples to
%               Matlab's Search Path for the current session
%   README.txt - a README briefly describing how to get started with ziDAQ
%
% DIRECTORIES
%   Driver/    - contains Matlab driver for interfacing with Zurich Instruments
%               devices
%   Utils/     - contains some utility functions for common tasks
%   Examples/  - contains examples for performing measurements on Zurich
%               Instruments devices
%
% DRIVER
%   Driver/ziDAQ.m          - ziDAQ command reference documentation.
%   Driver/ziDAQ.mex*      - ziDAQ API driver
%
% UTILS
%   ziAutoConnect          - Create a connection to a Zurich Instruments
%                           server (Deprecated: See ziCreateAPISession).
%   ziAutoDetect           - Return the ID of a connected device (if only one
%                           device is connected)
%   ziBW2TC                - Convert demodulator 3dB bandwidth to timeconstant
%   ziCheckPathInData      - Check whether a node is present in data and non-empty
%   ziCreateAPISession     - Create an API session for the specified device with
%                           the correct Data Server.
%   ziDevices              - Return a cell array of connected Zurich Instruments
%                           devices
%   ziGetDefaultSettingsPath - Get the default settings file path from the
%                           ziDeviceSettings ziCore module
%   ziGetDefaultSigoutMixerChannel - return the default output mixer channel
%   ziLoadSettings         - Load instrument settings from file
%   ziSaveSettings         - Save instrument settings to file
%   ziSiginAutorange       - Activate the device's autorange functionality
%   ziTC2BW                - Convert demodulator timeconstants to 3 dB Bandwidth
%
% EXAMPLES/COMMON - Examples that will run on any Zurich Instruments Device
%   example_connect        - A simple example to demonstrate how to
%                           connect to a Zurich Instruments device
%   example_connect_config - Connect to and configure a Zurich
%                           Instruments device
%   example_pid_advisor    - Setup and optimize a PID for internal
%                           PLL mode
%   example_poll           - Record demodulator data using
%                           ziDAQServer's synchronous poll function
%   example_record_async   - Record data asynchronously using ziDAQ's
%                           record module
```

```
% example_save_device_settings_simple - Save and load device settings
%                                     synchronously using ziDAQ's utility
%                                     functions
% example_save_device_settings_expert - Save and load device settings
%                                     asynchronously with ziDAQ's
%                                     devicesettings module
% example_scope - Record scope data using ziDAQServer's
%                                     synchronous poll function
% example_spectrum - Perform an FFT using ziDAQ's zoomFFT
%                                     module (Spectrum Tab of the LabOne UI)
% example_sweeper - Perform a frequency sweep using ziDAQ's
%                                     sweep module
% example_sweeper_rstddev_fixedbw - Perform a frequency sweep plotting the
%                                     stddev in demodulator output R using
%                                     ziDAQ's sweep module
% example_sweeper_two_demods - Perform a frequency sweep saving data
%                                     from 2 demodulators using ziDAQ's sweep
%                                     module
% example_swtrigger_edge - Record demodulator data upon a rising
%                                     edge trigger via ziDAQ's SW Trigger
%                                     module
% example_swtrigger_digital - Record data using a digital trigger via
%                                     ziDAQ's SW Trigger module
% example_swtrigger_grid - Record demodulator data, interpolated
%                                     on a grid from multiple triggers
%                                     using the SW Trigger's Grid Mode.
%
% EXAMPLES/UHF - Examples specific to the UHF Series
% uhf_example_boxcar - Record boxcar data using ziDAQServer's
%                                     synchronous poll function
% uhf_example_scope_offset - Record scope/digitizer data using
%                                     ziDAQServer's synchronous poll function
%
% EXAMPLES/HF2 - Examples specific to the HF2 Series
% hf2_example_autorange - determine and set an appropriate range
%                                     for a sign channel
% hf2_example_poll_hardware_trigger - Poll demodulator data in combination
%                                     with a HW trigger
% hf2_example_scope - Record scope data using ziDAQServer's
%                                     synchronous poll function
% hf2_example_zsync_poll - Synchronous demodulator sample timestamps
%                                     from multiple HF2s via the Zsync feature
```

Note

On Windows the MEX-file is called either `ziDAQ.mexw64` or `ziDAQ.mexw32` for 64-bit and 32-bit platforms respectively and on Linux it is called `ziDAQ.mexa64` or `ziDAQ.mexa32`. If more than one MEX-file is present, Matlab automatically selects the correct MEX-file for the current platform.

3.2.2. Using the Built-in Documentation

To access `ziDAQ`'s documentation within Matlab, type either of the following in the Matlab Command Window:

```
>> help ziDAQ
```

```
>> doc ziDAQ
```

This documentation is located in the file `MATLAB2012/Driver/ziDAQ.m`. See [Section 3.5, LabOne Matlab API \(ziDAQ\) Command Reference](#) for a printer friendly version.

3.2.3. Running the Examples

Prerequisites for running the Matlab examples:

1. Matlab is configured for ziDAQ as described above in [Section 3.1](#).
2. The Data Server program is running and the instrument is discoverable, this is the case if the instrument can be seen in the User Interface.
3. Signal Output 1 of the instrument is connected to Signal Input 1 via a BNC cable; many of the Matlab examples measure on this hardware channel.

See [Section 3.2.1](#) for a list of available examples bundled with the LabOne Matlab API. All the examples follow the same structure and take one input argument: the device ID of the instrument they are to be ran with. For example:

```
>> example_sweeper('dev123');
```

The example should produce some output in the Matlab Command Window, such as:

```
ziDAQ version Jul 7 2015 accessing server localhost 8005.  
Will run the example on `dev123`, an `HF2LI` with options `MFK|PLL|MOD|RTK|PID`.  
Sweep progress 9%  
Sweep progress 19%  
Sweep progress 30%  
Sweep progress 42%  
Sweep progress 52%  
Sweep progress 58%  
Sweep progress 68%  
Sweep progress 79%  
Sweep progress 91%  
Sweep progress 100%  
ziDAQ: AtExit called
```

Most examples will also plot some data in a Matlab figure, see [Figure 3.1](#) for an example. If you encounter an error message please ensure that the [above prerequisites](#) are fulfilled and see [Section 3.4](#) for help troubleshooting the error.

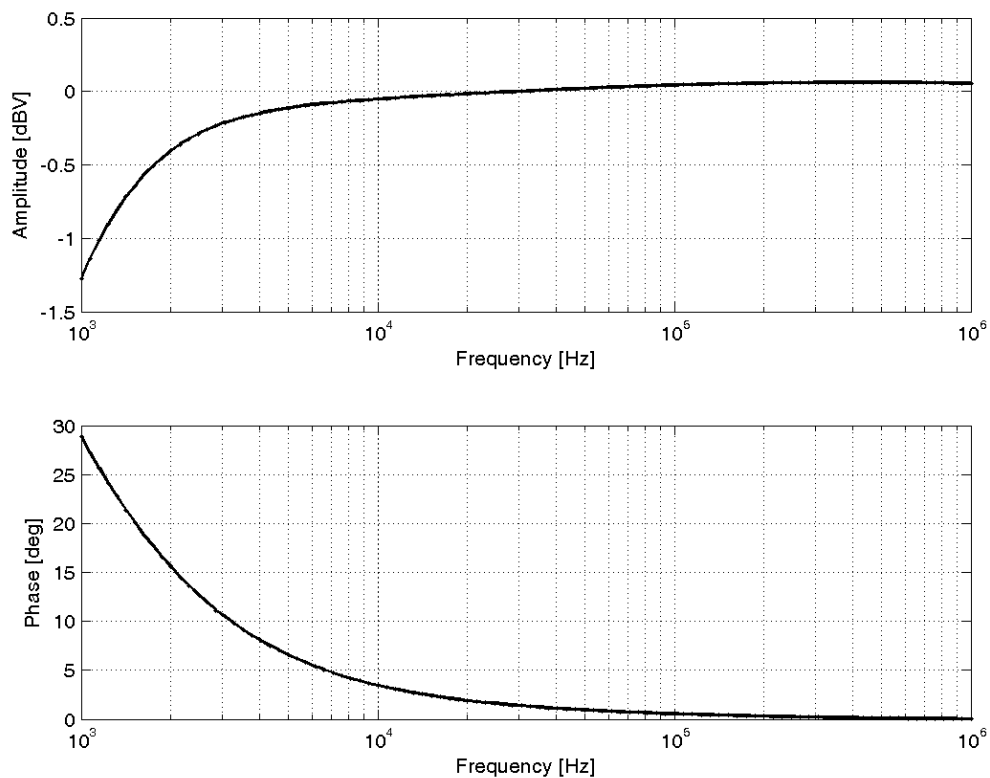


Figure 3.1. The plot produced by the LabOne Matlab API example `example_sweeper.m`; the plots show the instruments demodulator output when performing a frequency sweep over a simple feedback cable.

Note

The examples serve as a starting point for your own measurement needs. However, before editing the m-files, be sure to copy them to your own user space (they could be overwritten upon updating your LabOne installation) and give them a unique name to avoid name conflicts in Matlab.

3.2.4. Using ziCore Modules in the LabOne Matlab API

In the LabOne Matlab API [ziCore Modules](#) are configured and controlled via Matlab "handles". For example, in order to use the [Sweeper Module](#) a handle is created via:

```
>> h = ziDAQ('sweep');
```

and the Module's parameters are configured using the `set` command and specifying the Module's handle with a `path, value` pair, for example:

```
>> ziDAQ('set', h, 'sweep/start', 1.2e5);
```

The parameters can be read-back using the `get` command, which supports wildcards, for example:

```
>> sweep_params = ziDAQ('get', h, 'sweep/*');
```

The variable `sweep_params` now contains a struct of all the Sweeper's parameters. The other main Module commands are used similarly, e.g., `ziDAQ('execute', h)` to start the sweeper. See [Section 2.1.2](#) for more help with Modules and a description of their parameters.

3.2.5. Enabling Logging in the LabOne Matlab API

Logging from the API is not enabled by default upon initializing a server session with `ziDAQ`, it must be enabled (after using `connect`) with the `setDebugLevel` command. For example,

```
>> ziDAQ('setDebugLevel', 0);
```

sets the API's logging level to 0, which provides the most verbose logging output. The other log levels are defined as following:

```
trace:0, info:1, debug:2, warning:3, error:4, fatal:5, status:6.
```

It is also possible for the user to write their own messages directly to `ziDAQ`'s log using the `writeDebugLog` command. For example to write a log message of `info` severity level:

```
>> ziDAQ('writeDebugLog', 1, 'Hello log!');
```

On Windows the logs can be found by navigating to the Zurich Instruments Logs" folder entry in the Windows Start Menu: Programs → Zurich Instruments → LabOne Servers → Logs. This will open an Explorer window displaying folders containing log files from various LabOne components, in particular, the `ziDAQLog` folder contains logs from the LabOne Matlab API. On Linux, the logs can be found at `"/tmp/ziDAQLog_USERNAME"`, where `"USERNAME"` is the same as the output of the `"whoami"` command.

3.3. LabOne Matlab API Tips and Tricks

In this section some tips and tricks for working with the LabOne Matlab API are provided.

The structure of **ziDAQ** commands.

All LabOne Matlab API commands are based on a call to the Matlab function `ziDAQ()`. The first argument to `ziDAQ()` specifies the API command to be executed and is an obligatory argument. For example, a session is instantiated between the API and the Data Server with the Matlab command `ziDAQ('connect')`. Depending on the type of command specified, optional arguments may be required. For example, to obtain an integer node value, the node path must be specified as a second argument to the `'getInt'` command:

```
s = ziDAQ('getInt','/dev123/sigouts/0/on');
```

where the output argument contains the current value of the specified node.

To set an integer node value, both the node path and the value to be set must be specified as the second and third arguments:

```
ziDAQ('setInt','/dev123/sigouts/0/on', 1);
```

See the [LabOne Matlab API \(ziDAQ\) Command Reference](#) for a list of all available commands.

Data Structures returned by **ziDAQ**.

The output arguments that `ziDAQ` returns are designed to use the native data structures that Matlab users are familiar with and that reflect the data's location in the instruments node hierarchy. For example, when the `poll` command returns data from the instruments fourth demodulator (located in the node hierarchy as `/dev123/demods/3/sample`), the output argument contains a nested `struct` in which the data can be accessed by

```
data = ziDAQ('poll', poll_length, poll_timeout);  
x = data.dev123.demods(4).sample.x;  
y = data.dev123.demods(4).sample.y;
```

The instrument's node tree uses zero-based indexing; Matlab uses one-based indexing.

See the tip [Data Structures returned by ziDAQ](#): The **fourth** demodulator sample located at `/dev123/demods/3/sample`, is indexed in the data structure returned by `poll` as `data.dev123.demods(4).sample`.

Explicitly convert **uint64** data types to **double**.

Matlab's native data type is double-precision floating point and doesn't support performing calculations with other data types such as 64-bit unsigned integers, for example:

```
>> a = uint64(2); b = uint64(1); a - b  
? Undefined function or method 'minus' for input arguments of type 'uint64'.
```

Due to this limitation, be sure to convert demodulator timestamps to `double` before performing calculations. For example, in the following, both `clockbase` and `timestamp` (both 64-bit unsigned

integers) need to be converted to double before converting the timestamps from the instrument's native "ticks" to seconds via the instrument's clockbase:

```
data = ziDAQ('poll', 1.0, 500);           % poll data
sample = data.(device).demods(0).sample; % get the sample from the zeroth demod
% convert timestamps from ticks to seconds via the device's clockbase
% (the ADC's sampling rate), specify reference start time via t0.
clockbase = double(ziDAQ('getInt',['/' device '/clockbase']));
t = (double(sample.timestamp) - double(sample.timestamp(1)))/clockbase;
```

Use the utility function **ziCheckPathInData**.

Checking that a sub-structure in the nested data structure returned by `poll` actually exists can be cumbersome and can require multiple nested `if` statements; this can be avoided by using the utility function `ziCheckPathInData`. For example, the code:

```
data = ziDAQ('poll', poll_length, poll_timeout );
if isfield(data,device)
    if isfield(data.(device),'demods')
        if length(data.(device).demods) >= channel
            if ~isempty(data.(device).demods(channel).sample)
                % do something with the demodulator sample...
```

can be replaced by:

```
data = ziDAQ('poll', poll_length, poll_timeout );
if ziCheckPathInData( data, ['/' device '/demods/' demod_c '/sample']);
    % do something with the demodulator sample...
```

3.4. Troubleshooting the LabOne Matlab API

This section intends to solve possible error messages than can occur when using `ziDAQ` in Matlab.

Error message: "Undefined function or method 'ziDAQ' for input arguments of type '*'"

Matlab can not find the LabOne Matlab API library. Check whether the `MATLAB2012/Driver` subfolder of your LabOne installation is in the Matlab Search Path by using the command:

```
>> path
```

and repeating the steps to configure Matlab's search path in [Section 3.1.2](#).

Error message: "Undefined function or method 'example_sweeper'"

Matlab can not find the example. Check whether the `MATLAB2012/Examples/Common` subfolder (respectively `MATLAB2012/Examples/UHF` or `MATLAB2012/Examples/HF2`) of your LabOne installation are in the Matlab Search Path by using the command:

```
>> path
```

and repeating the steps to configure Matlab's search path in [Section 3.1.2](#).

Error message: "Error using: ziDAQ ZIAPIException with status code: 32870. Connection invalid."

The Matlab API can not connect to the Data Server. Please check that the correct port was used; that the correct server is running for your device and that the device is connected to the server, see [Section 1.3.1](#).

Error Message: "Error using: ziAutoConnect at 63 ziAutoConnect(): failed to find a running server or failed to find a connected a device..."

The utility function `ziAutoConnect()` located in `MATLAB2012/Utils/` tries to determine which Data Server is running and whether any devices are connected to that Data Server. It is only supported by UHFLI and HF2 Series instruments, MFLI instruments are not supported. Some suggestions to verify the problem:

- Please verify in the User Interface, whether a device is connected to the Data Server running on your computer.
- If the Data Server is running on a different computer, connect manually to the Data Server via `ziDAQ`'s `connect` function:

```
>> ziDAQ('connect', hostname, port);
```

where `hostname` should be replaced by the IP of the computer the Data Server is running on and `port` is specified as in [Section 1.3.1](#).

Error Message: "Error using: ziDAQ ZIAPException on path /dev123/sigins/0/imp50 with status code: 16387. Value or Node not found"

The API is connected to the Data Server, but the command failed to find the specified node. Please:

- Check whether your instrument is connected to the Data Server in the User Interface; if it is not connected the instruments device node tree, e.g., /dev123/, will not be constructed by the Data Server.
- Check whether the node path is spelt correctly.
- Explore the node tree to verify the node actually exists with the `listNodes` command:

```
>> ziDAQ('listNodes', '/dev123/sigins/0', 3)
```

Error Message: "using: ziDAQ Server not connected. Use 'ziDAQ('connect', ...) first."

A `ziDAQ` command was issued before initializing a connection to the Data Server. First use the `connect` command:

```
>> ziDAQ('connect', hostname, port);
```

where `hostname` should be replaced by the IP address of the computer the Data Server is running on and `port` is specified as in [Section 1.3.1](#). If the Data Server is running on the same computer, use `'localhost'` as the `hostname`.

Error Message: "Attempt to execute SCRIPT ziDAQ as a function: ziDAQ.m"

There could be a problem with your LabOne Matlab API installation. The call to `ziDAQ()` is trying to call the help file `ziDAQ.m` as a function instead of calling the `ziDAQ()` function defined in the MEX-file. In this case you need to ensure that the `ziDAQ` MEX-file is in your search path as described in [Section 3.1](#) and navigate away from the `Driver` directory. Secondly, ensure that the LabOne Matlab MEX-file is in the `Driver` folder as described in [Section 3.2.1](#).

3.5. LabOne Matlab API (ziDAQ) Command Reference

```
%
% Copyright 2009-2016, Zurich Instruments Ltd, Switzerland
% This software is a preliminary version. Function calls and
% parameters may change without notice.
%
% This version of ziDAQ is linked against:
% * Matlab 7.9.0.529, R2009b, Windows,
% * Matlab 8.4.0.145, R2014b, Linux64.
% You can check which version of Matlab you are using Matlab's `ver` command.
% A list of compatible Matlab and ziDAQ versions is available here:
% www.zhinst.com/labone/compatibility
%
% ziDAQ is an interface for communication with Zurich Instruments Data Servers.
%
% Usage: ziDAQ(command, [option1], [option2])
%       command = 'clear', 'connect', 'connectDevice',
%               'disconnectDevice', 'finished', 'flush', 'get',
%               'getAsEvent', 'getAuxInSample', 'getBytes',
%               'getDIO', 'getDouble', 'getInt',
%               'getSample', 'listNodes', 'logOn', 'logOff',
%               'poll', 'pollEvent', 'programRT', 'progress', 'read',
%               'record', 'setByte', 'setDouble', 'syncSetDouble',
%               'setInt', 'syncSetInt', 'subscribe',
%               'sweep', 'trigger', 'unsubscribe', 'update',
%               'zoomFFT', 'deviceSettings'
%
% Preconditions: ZI Server must be running (check task manager)
%
%       ziDAQ('connect', [host = '127.0.0.1'], [port = 8005], [apiLevel = 1]);
%       [host] = Server host string (default is localhost)
%       [port] = Port number (double)
%               Use port 8005 to connect to the HF2 Data Server
%               Use port 8004 to connect to the MF or UHF Data Server
%       [apiLevel] = Compatibility mode of the API interface (int64)
%               Use API level 1 to use code written for HF2.
%               Higher API levels are currently only supported
%               for MF and UHF devices. To get full functionality for
%               MF and UHF devices use API level 5.
%       To disconnect use 'clear ziDAQ'
%
%       result = ziDAQ('getConnectionAPILevel');
%               Returns ziAPI level used for the active connection.
%
%       ziDAQ('connectDevice', device, interface);
%       device (string) = Device serial to connect (e.g. 'DEV2000')
%       interface (string) = Interface, e.g., 'USB', '1GbE', '10GbE'.
%       Connect with the data server to a specified device over the
%       specified interface. The device must be visible to the server.
%       If the device is already connected the call will be ignored.
%       The function will block until the device is connected and
%       the device is ready to use. This method is useful for UHF
%       devices offering several communication interfaces.
%
%       ziDAQ('disconnectDevice', device);
%       device (string) = Device serial of device to disconnect.
%       This function will return immediately. The disconnection of
%       the device may not yet finished.
%
%       result = ziDAQ('listNodes', path, flags);
%       path (string) = Node path or partial path, e.g.,
%                       '/dev100/demods/'.
```

```

%           flags (int64) = Define which nodes should be returned, set the
%           following bits to obtain the described behavior:
%           int64(0) -> ZI_LIST_NODES_NONE 0x00
%           The default flag, returning a simple
%           listing of the given node
%           int64(1) -> ZI_LIST_NODES_RECURSIVE 0x01
%           Returns the nodes recursively
%           int64(2) -> ZI_LIST_NODES_ABSOLUTE 0x02
%           Returns absolute paths
%           int64(4) -> ZI_LIST_NODES_LEAFONLY 0x04
%           Returns only nodes that are leafs,
%           which means the they are at the
%           outermost level of the tree.
%           int64(8) -> ZI_LIST_NODES_SETTINGSONLY 0x08
%           Returns only nodes which are marked
%           as setting
%           Flags may also be combined, e.g., set flags to bitor(1, 2)
%           to return paths recursively and printed as absolute paths.
%
% result = ziDAQ('getSample', path);
% path (string) = Node path
% Returns a single demodulator sample (including
% DIO and AuxIn). For more efficient data recording
% use the subscribe and poll functions.
%
% result = ziDAQ('getAuxInSample', path);
% path (string) = Node path
% Returns a single auxin sample. Note, the auxin data
% is averaged in contrast to the auxin data embedded
% in the demodulator sample.
%
% result = ziDAQ('getDIO', path);
% path (string) = Node path.
% Returns a single DIO sample.
%
% result = ziDAQ('getDouble', path);
% path (string) = Node path
%
% result = ziDAQ('getInt', path);
% path (string) = Node path
%
% result = ziDAQ('getByte', path);
% path (string) = Node path
%
% ziDAQ('setDouble', path, value);
% path (string) = Node path
% value (double) = Setting value
%
% ziDAQ('syncSetDouble', path, value);
% Deprecated, see the 'sync' command.
% path (string) = Node path
% value (double) = Setting value
%
% ziDAQ('setInt', path, value);
% path (string) = Node path
% value (int64) = Setting value
%
% ziDAQ('syncSetInt', path, value);
% Deprecated, see the 'sync' command.
% path (string) = Node path
% value (int64) = Setting value
%
% ziDAQ('setByte', path, value);
% path (string) = Node path
% value (double) = Setting value
%
% ziDAQ('vectorWrite', path, value);

```

```

%           path (string) = Vector node path
%           value (vector of (u)int8, (u)int16, (u)int32, (u)int64,
%               float, double; or string) = Setting value
%
%   ziDAQ('subscribe', path);
%       path (string) = Node path
%       Subscribe to the specified path to receive streaming data
%       or setting data if changed. Use either 'poll' command to
%       obtain the subscribed data.
%
%   ziDAQ('unsubscribe', path);
%       path (string) = Node path
%       Unsubscribe from the node paths specified via 'subscribe'.
%       Use a wildcard ('*') to unsubscribe from all data.
%
%   ziDAQ('getAsEvent', path);
%       path (string) = Node path
%       Triggers a single event on the path to return the current
%       value. The result can be fetched with the 'poll' or 'pollEvent'
%       command.
%
%   ziDAQ('update');
%       Detect HF2 devices connected to the USB. On Windows this
%       update is performed automatically.
%
%   ziDAQ('get', path, [settingsOnly]);
%       path (string) = Node path
%       Gets a structure of the node data from the specified
%       branch. High-speed streaming nodes (e.g. /devN/demods/0/sample)
%       are not returned. Wildcards (*) may be used, in which case
%       read-only nodes are ignored.
%       [settingsOnly] (uint32) = Specify which type of nodes to include
%       in the result. Allowed:
%           ZI_LIST_NODES_SETTINGSONLY = 8 (default)
%           ZI_LIST_NODES_NONE = 0 (all nodes)
%
%   ziDAQ('flush');
%       Deprecated, see the 'sync' command.
%       Flush all data in the socket connection and API buffers.
%       Call this function before a subscribe with subsequent poll
%       to get rid of old streaming data that might still be in
%       the buffers.
%
%   ziDAQ('echoDevice', device);
%       Deprecated, see the 'sync' command.
%       device (string) = device serial, e.g. 'dev100'.
%       Sends an echo command to a device and blocks until
%       answer is received. This is useful to flush all
%       buffers between API and device to enforce that
%       further code is only executed after the device executed
%       a previous command.
%
%   ziDAQ('sync');
%       Synchronize all data paths. Ensures that get and poll
%       commands return data which was recorded after the
%       setting changes in front of the sync command. This
%       sync command replaces the functionality of all 'syncSet*',
%       'flush', and 'echoDevice' commands.
%
%   ziDAQ('programRT', device, filename);
%       device (string) = device serial, e.g. 'dev100'.
%       filename (string) = filename of RT program.
%       HF2 devices only; writes down a real-time program. Requires
%       the Real time Option must be available for the specified
%       HF2 device.
%
%   result = ziDAQ('secondsTimeStamp', [timestamps]);

```

```

%           timestamps (uint64) = vector of uint64 device ticks
%           Deprecated. In order to convert timestamps to seconds divide the
%           timestamps by the value instrument's clockbase device node,
%           e.g., /dev99/clockbase.
%           [Converts a timestamp vector of uint64 ticks
%           into a double vector of timestamps in seconds (HF2 Series).]
%
% Synchronous Interface
%
%           ziDAQ('poll', duration, timeout, [flags]);
%           duration (double) = Recording time in [s]
%           timeout (int64) = Poll timeout in [ms]
%           [flags] (uint32) = Flags that specify data polling properties
%           Bit[0] FILL : Fill data loss holes
%           Bit[1] ALIGN : Align data of several demodulators
%           Bit[2] THROW : Throw if data loss is detected
%           Records data for the specified time. This function call
%           is blocking. Use ziDAQRecorder's asynchronous interface for
%           long recording durations.
%
%           result = ziDAQ('pollEvent', timeout);
%           timeout (int64) = Poll timeout in [ms]
%           Execute a single poll command. This is a low-level
%           function. The poll function is better suited for most
%           cases.
%
% Asynchronous Interface
%
% Trigger Parameters
%   trigger/buffersize      double Set the buffersize [s] of the trigger
%                               object. The recommended buffer size is
%                               2*trigger/0/duration.
%   trigger/flags           int    Record flags.
%                               FILL = 0x0001 : Fill holes.
%                               ALIGN = 0x0002 : Align data that contains a
%                               timestamp.
%                               THROW = 0x0004 : Throw if sample loss
%                               is detected.
%   trigger/device          string The device serial to use the software trigger
%                               with, e.g. dev123 (compulsory parameter).
%   trigger/endless         bool   Enable endless triggering 1=enable; 0=disable.
%   trigger/forcetrigger    bool   Force a trigger.
%   trigger/0/triggernode   string Path and signal of the node that should be
%                               used for triggering, separated by a dot (.),
%                               e.g. /devN/demods/0/sample.x
%                               Overrides values from trigger/0/path and
%                               trigger/0/source.
%   trigger/0/path          string The path to the demod sample to trigger on,
%                               e.g. demods/3/sample, see also trigger/0/source
%                               DEPRECATED - use trigger/0/triggernode instead
%   trigger/0/source        int    Signal that is used to trigger on.
%                               0 = x
%                               1 = y
%                               2 = r
%                               3 = angle
%                               4 = frequency
%                               5 = phase
%                               6 = auxiliary input 0 / parameter 0
%                               7 = auxiliary input 1 / parameter 1
%                               DEPRECATED - use trigger/0/triggernode instead
%   trigger/0/count         int    Number of trigger edges to record.
%   trigger/0/type          int    Trigger type used. Some parameters are
%                               only valid for special trigger types.
%                               0 = trigger off
%                               1 = analog edge trigger on source
%                               2 = digital trigger mode on DIO source
%                               3 = analog pulse trigger on source

```

%			4 = analog tracking trigger on source
%			5 = hardware trigger on trigger line source
%			6 = tracking edge trigger on source
%			7 = event count trigger on counter source
%	trigger/0/edge	int	Trigger edge
%			1 = rising edge
%			2 = falling edge
%			3 = both
%	trigger/0/findlevel	bool	Automatically find the value of trigger/0/level based on the current signal value.
%	trigger/0/bits	int	Digital trigger condition.
%	trigger/0/bitmask	int	Bit masking for bits used for triggering. Used for digital trigger.
%	trigger/0/delay	double	Trigger frame position [s] (left side) relative to trigger edge.
%			delay = 0 -> trigger edge at left border.
%			delay < 0 -> trigger edge inside trigger frame (pretrigger).
%			delay > 0 -> trigger edge before trigger frame (posttrigger).
%	trigger/0/duration	double	Recording frame length [s]
%	trigger/0/level	double	Trigger level voltage [V].
%	trigger/0/hysteresis	double	Trigger hysteresis [V].
%	trigger/0/retrigger	int	Record more than one trigger in a trigger frame. If a trigger event is currently being recorded and another trigger event is detected within the duration of the current trigger event, extend the size of the trigger frame to include the duration of the new trigger event.
%	trigger/triggered	bool	Has the software trigger triggered? 1=Yes, 0=No (read only).
%	trigger/0/bandwidth	double	Filter bandwidth [Hz] for pulse and tracking triggers.
%	trigger/0/holdoff/count	int	Number of skipped triggers until the next trigger is recorded again.
%	trigger/0/holdoff/time	double	Hold off time [s] before the next trigger is recorded again. A hold off time smaller than the duration will produce overlapped trigger frames.
%	trigger/0/hwtrigsource	int	Only available for devices that support hardware triggering. Specify the channel to trigger on.
%			DEPRECATED - use trigger/0/triggernode instead
%	trigger/0/pulse/min	double	Minimal pulse width [s] for the pulse trigger.
%	trigger/0/pulse/max	double	Maximal pulse width [s] for the pulse trigger.
%	trigger/0/grid/mode	int	Enable grid mode. In grid mode a matrix instead of a vector is returned. Each trigger becomes a row in the matrix and each trigger's data is interpolated onto a new grid defined by the number of columns:
%			0: Disable
%			1: Enable with nearest neighbour interpolation
%			2: Enable with linear interpolation.
%	trigger/0/grid/operation	int	If running in endless mode, either replace or average the data in the grid's matrix.
%	trigger/0/grid/cols	int	Specify the number of columns in the grid's matrix. The data from each row is interpolated onto a grid with the specified number of columns.
%	trigger/0/grid/rows	int	Specify the number of rows in the grid's matrix. Each row is the data recorded from one trigger interpolated onto the columns.
%	trigger/0/grid/direction	int	The direction to organize data in the grid's

```

matrix:
0: Forward.
    The data in each row is ordered chronologically, e.g., the first data point in each row corresponds to the first timestamp in the trigger data.
1: Reverse.
    The data in each row is ordered reverse chronologically, e.g., the first data point in each row corresponds to the last timestamp in the trigger data.
2: Bidirectional.
    The ordering of the data alternates between Forward and Backward ordering from row-to-row. The first row is Forward ordered.

trigger/filename      string This parameter is deprecated. If specified, i.e., not empty, it enables automatic saving of data in single trigger mode (trigger/endsless = 0).
trigger/directory     string The directory where files are saved.
trigger/fileformat    int    The format of the file for saving data.
                                0 = Matlab,
                                1 = CSV,
                                2 = ZView (Impedance data only).
trigger/historylength bool    Maximum number of entries stored in the measurement history.
trigger/clearhistory  bool    Remove all records from the history list.

handle = ziDAQ('record' duration, timeout);
    duration (double) = The module's internal buffersize to use when recording data [s]. The recommended size is 2*trigger/0/duration parameter. Note that this can be modified via the trigger/buffersize parameter.
                        DEPRECATED, set 'buffersize' param instead.
    timeout (int64) = Poll timeout [ms]. - DEPRECATED, ignored
Create an instance of the ziDAQRecorder class (note that the module's thread is not yet started) and return a Matlab handle with which to access it.
Before the thread can actually be started (via 'execute'):
- the desired data to record must be specified via the module's 'subscribe' command,
- the device serial (e.g., dev100) that will be used must be set.
The real measurement is started upon calling the 'execute' function. After that the trigger will start recording data and verifying for incoming triggers.

result = ziDAQ('listNodes', handle, path, flags);
    handle = Matlab handle (reference) specifying an instance of the ziDAQRecorder class.
    path (string) = Module parameter path
    flags (int64) = Define which module parameters paths should be returned, set the following bits to obtain the described behaviour:
    flags = int64(0) -> ZI_LIST_NODES_NONE 0x00
                The default flag, returning a simple listing of the given path
    int64(1) -> ZI_LIST_NODES_RECURSIVE 0x01
                Returns the paths recursively
    int64(2) -> ZI_LIST_NODES_ABSOLUTE 0x02
                Returns absolute paths
    int64(4) -> ZI_LIST_NODES_LEAFONLY 0x04
                Returns only paths that are leafs, which means the they are at the outermost level of the tree.
    int64(8) -> ZI_LIST_NODES_SETTINGSONLY 0x08

```

```
%           Returns only paths which are marked
%           as setting
%           Flags may also be combined, e.g., set flags to bitor(1, 2)
%           to return paths recursively and printed as absolute paths.
%
%           ziDAQ('subscribe', handle, path);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           path (string) = Node path to record data from.
%           Subscribe to device nodes. Call multiple times to
%           subscribe to multiple node paths. After subscription the
%           recording process can be started with the 'execute'
%           command. During the recording process paths can not be
%           subscribed or unsubscribed.
%
%           ziDAQ('unsubscribe', handle, path);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           path (string) = Node path to record data from. Use wildcard
%           ('*') to select all.
%           Unsubscribe from one or several nodes. During the
%           recording process paths can not be subscribed or
%           unsubscribed.
%
%           ziDAQ('get', handle, path);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           path (string) = Path string of the module parameter. Must
%           start with 'sweep/'.
%           Get module parameters. Wildcards are supported, e.g. 'sweep/*'.
%
%           ziDAQ('set', handle, path, value);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           path (string) = Path string of the module parameter. Must
%           start with 'sweep/'.
%           value = The value to set the module parameter to, see the list
%           of module parameters for the correct type.
%           Set the specified module parameter value.
%
%           ziDAQ('execute', handle);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           Start the recorder. After that command any trigger will
%           start the measurement. Subscription or unsubscription
%           is not possible until the recording is finished.
%
%           ziDAQ('trigger', handle);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           Force a trigger to manually record one duration of the
%           subscribed data.
%
%           result = ziDAQ('finished', handle);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           Returns 1 if the recording is finished, otherwise 0.
%
%           result = ziDAQ('read', handle);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           Read out the recorded data; transfer the recorded data to
%           Matlab.
%
%           ziDAQ('finish', handle);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
```



```

%           Stop recording data. The recording may be restarted by
%           calling 'execute' again.
%
% result = ziDAQ('progress', handle);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           Report the progress of the measurement with a number
%           between 0 and 1.
%
%           ziDAQ('clear', handle);
%           handle = Matlab handle (reference) specifying an instance of
%           the ziDAQRecorder class.
%           Stop the module's thread.
%
% Sweep Module
%
% Sweep Parameters
%   sweep/device          string  Device that should be used for
%                               the parameter sweep, e.g. 'dev99'.
%   sweep/start           double  Sweep start frequency [Hz]
%   sweep/stop            double  Sweep stop frequency [Hz]
%   sweep/gridnode        string  Path of the node that should be
%                               used for sweeping. For frequency
%                               sweep applications this will be e.g.
%                               'oscs/0/freq'. The device name of
%                               the path can be omitted and is given
%                               by sweep/device.
%   sweep/loopcount       int     Number of sweep loops (default 1)
%   sweep/endless         int     Endless sweeping (default 0)
%                               0 = Use loopcount value
%                               1 = Endless sweeping enabled, ignore
%                               loopcount
%   sweep/samplecount     int     Number of samples per sweep
%   sweep/settling/time   double  Settling time before measurement is
%                               performed, in [s]
%   sweep/settling/tc     double  Settling precision
%                               5 ~ low precision
%                               15 ~ medium precision
%                               50 ~ high precision
%   sweep/settling/inaccuracy int  Demodulator filter settling inaccuracy
%                               that defines the wait time between a
%                               sweep parameter change and recording of
%                               the next sweep point. The settling time
%                               is calculated as the time required to
%                               attain the specified remaining proportion
%                               [1e-13, 0.1] of an incoming step
%                               function. Typical inaccuracy
%                               values:
%                               - 10m for highest sweep speed for large
%                               signals,
%                               - 100u for precise amplitude measurements,
%                               - 100n for precise noise measurements.
%                               Depending on the order the settling
%                               accuracy will define the number of filter
%                               time constants the sweeper has to
%                               wait. The maximum between this value and
%                               the settling time is taken as wait time
%                               until the next sweep point is recorded.
%   sweep/xmapping        int     Sweep mode
%                               0 = linear
%                               1 = logarithmic
%   sweep/scan            int     Scan type
%                               0 = sequential
%                               1 = binary
%                               2 = bidirectional
%                               3 = reverse
%   sweep/bandwidth       double  Fixed bandwidth [Hz]

```

```

%
% sweep/bandwidthcontrol int      0 = Automatic calculation (obsolete)
%                               Sets the bandwidth control mode (default 2)
%                               0 = Manual (user sets bandwidth and order)
%                               1 = Fixed (uses fixed bandwidth value)
%                               2 = Auto (calculates best bandwidth value)
%                               Equivalent to the obsolete bandwidth = 0
%                               setting
% sweep/bandwidthoverlap bool     Sets the bandwidth overlap mode (default 0). If
%                               enabled the bandwidth of a sweep point may
%                               overlap with the frequency of neighboring sweep
%                               points. The effective bandwidth is only limited
%                               by the maximal bandwidth setting and omega
%                               suppression. As a result, the bandwidth is
%                               independent of the number of sweep points. For
%                               frequency response analysis bandwidth overlap
%                               should be enabled to achieve maximal sweep
%                               speed (default: 0).
%                               0 = Disable
%                               1 = Enable
% sweep/order int                Defines the filter roll off to use in Fixed
%                               bandwidth selection.
%                               Valid values are between 1 (6 dB/octave)
%                               and 8 (48 dB/octave). An order of 0
%                               triggers a read-out of the order from the
%                               selected demodulator.
% sweep/maxbandwidth double      Maximal bandwidth used in auto bandwidth
%                               mode in [Hz]. The default is 1.25MHz.
% sweep/omegasuppression double   Damping in [dB] of omega and 2omega components.
%                               Default is 40dB in favor of sweep speed.
%                               Use higher value for strong offset values or
%                               3omega measurement methods.
% sweep/averaging/tc double       Min averaging time [tc]
%                               0 = no averaging (see also time!)
%                               5 ~ low precision
%                               15 ~ medium precision
%                               50 ~ high precision
% sweep/averaging/sample int      Min samples to average
%                               1 = no averaging (if averaging/tc = 0)
% sweep/phaseunwrap bool          Enable unwrapping of slowly changing phase
%                               evolutions around the +/-180 degree boundary.
% sweep/sincfilter bool           Enables the sinc filter if the sweep frequency
%                               is below 50 Hz. This will improve the sweep
%                               speed at low frequencies as omega components
%                               do not need to be suppressed by the normal
%                               low pass filter.
% sweep/filename string           This parameter is deprecated. If specified,
%                               i.e. not empty, it enables automatic saving of
%                               data in single sweep mode (sweep/endsless = 0).
% sweep/directory string          The directory where files are located when
%                               saving sweeper measurements.
% sweep/fileformat int            The format of the file for saving sweeper
%                               measurements:
%                               0 = Matlab,
%                               1 = CSV,
%                               2 = ZView (Impedance data only).
% sweep/historylength bool        Maximum number of entries stored in the
%                               measurement history.
% sweep/clearhistory bool         Remove all records from the history list.
%
% Note:
% Settling time = max(settling.tc * tc, settling.time)
% Averaging time = max(averaging.tc * tc, averaging.sample / sample-rate)
%
% handle = ziDAQ('sweep', timeout);
%           timeout = Poll timeout in [ms] - DEPRECATED, ignored
%           Creates a sweep class. The thread is not yet started.
%           Before the thread start subscribe and set command have

```

```
%
%      to be called. To start the real measurement use the
%      execute function.
%
%      result = ziDAQ('listNodes', handle, path, flags);
%      path (string) = Module parameter path
%      flags (int64) = Define which module parameters paths should be
%                      returned, set the following bits to obtain the
%                      described behaviour:
%                      int64(0) -> ZI_LIST_NODES_NONE 0x00
%                          The default flag, returning a simple
%                          listing of the given path
%                      int64(1) -> ZI_LIST_NODES_RECURSIVE 0x01
%                          Returns the paths recursively
%                      int64(2) -> ZI_LIST_NODES_ABSOLUTE 0x02
%                          Returns absolute paths
%                      int64(4) -> ZI_LIST_NODES_LEAFSONLY 0x04
%                          Returns only paths that are leafs,
%                          which means the they are at the
%                          outermost level of the tree.
%                      int64(8) -> ZI_LIST_NODES_SETTINGSONLY 0x08
%                          Returns only paths which are marked
%                          as setting
%      Flags may also be combined, e.g., set flags to bitor(1, 2)
%      to return paths recursively and printed as absolute paths.
%
%      ziDAQ('subscribe', handle, path);
%      Subscribe to one or several nodes. After subscription
%      the recording process can be started with the 'execute'
%      command. During the recording process paths can not be
%      subscribed or unsubscribed.
%      handle = Reference to the ziDAQSweeper class.
%      path = Path string of the node. Use wild card to
%      select all. Alternatively also a list of path
%      strings can be specified.
%
%      ziDAQ('unsubscribe', handle, path);
%      Unsubscribe from one or several nodes. During the
%      recording process paths can not be subscribed or
%      unsubscribed.
%      handle = Reference to the ziDAQSweeper class.
%      path = Path string of the node. Use wild card to
%      select all. Alternatively also a list of path
%      strings can be specified.
%
%      ziDAQ('execute', handle);
%      Start the sweep. Subscription or unsubscription
%      is no more possible until the sweep is finished.
%
%      result = ziDAQ('finished', handle);
%      handle = Handle of the sweep session.
%      Returns 1 if the sweep is finished, otherwise 0.
%
%      result = ziDAQ('read', handle);
%      handle = Handle of the sweep session.
%      Transfer the sweep data to Matlab.
%
%      result = ziDAQ('progress', handle);
%      Report the progress of the measurement with a number
%      between 0 and 1.
%
%      ziDAQ('finish', handle);
%      Stop the sweep. The sweep may be restarted by
%      calling 'execute' again.
%
%      ziDAQ('clear', handle);
%      handle = Handle of the sweep session.
%      Stop the current sweep.
```

```

%
%      ziDAQ('save', handle);
%      Save the measured data to a file.
%      handle = Handle of the sweep session.
%      [filename] = File in which to store the data.
%
% Zoom FFT Module
%
% Zoom FFT Parameters
%      zoomFFT/device      string  Device that should be used for
%                                the zoom FFT, e.g. 'dev99'.
%      zoomFFT/bit         int     Number of FFT points 2^bit
%      zoomFFT/mode        int     Zoom FFT mode
%                                0 = Perform FFT on X+iY
%                                1 = Perform FFT on R
%                                2 = Perform FFT on Phase
%      zoomFFT/loopcount   int     Number of zoom FFT loops (default 1)
%      zoomFFT/endless     int     Perform endless zoom FFT (default 0)
%                                0 = Use loopcount value
%                                1 = Endless zoom FFT enabled, ignore
%                                    loopcount
%      zoomFFT/overlap     double  FFT overlap 0 = none, [0..1]
%      zoomFFT/settling/time double Settling time before measurement is performed
%      zoomFFT/settling/tc double Settling time in time constant units before
%                                the FFT recording is started.
%                                5 ~ low precision
%                                15 ~ medium precision
%                                50 ~ high precision
%      zoomFFT/window      int     FFT window (default 1 = Hann)
%                                0 = Rectangular
%                                1 = Hann
%                                2 = Hamming
%                                3 = Blackman Harris 4 term
%      zoomFFT/absolute    bool    Shifts the frequencies so that the center
%                                frequency becomes the demodulation frequency
%                                rather than 0 Hz.
%
%      handle = ziDAQ('zoomFFT', timeout);
%      timeout = Poll timeout in [ms] - DEPRECATED, ignored
%      Creates a zoom FFT class. The thread is not yet started.
%      Before the thread start subscribe and set command have
%      to be called. To start the real measurement use the
%      execute function.
%
%      result = ziDAQ('listNodes', handle, path, flags);
%      path (string) = Module parameter path
%      flags (int64) = Define which module parameters paths should be
%                      returned, set the following bits to obtain the
%                      described behaviour:
%                      int64(0) -> ZI_LIST_NODES_NONE 0x00
%                          The default flag, returning a simple
%                          listing of the given path
%                      int64(1) -> ZI_LIST_NODES_RECURSIVE 0x01
%                          Returns the paths recursively
%                      int64(2) -> ZI_LIST_NODES_ABSOLUTE 0x02
%                          Returns absolute paths
%                      int64(4) -> ZI_LIST_NODES_LEAFSONLY 0x04
%                          Returns only paths that are leafs,
%                          which means the they are at the
%                          outermost level of the tree.
%                      int64(8) -> ZI_LIST_NODES_SETTINGSONLY 0x08
%                          Returns only paths which are marked
%                          as setting
%      Flags may also be combined, e.g., set flags to bitor(1, 2)
%      to return paths recursively and printed as absolute paths.
%
%      ziDAQ('subscribe', handle, path);

```

```

%          Subscribe to one or several nodes. After subscription
%          the recording process can be started with the 'execute'
%          command. During the recording process paths can not be
%          subscribed or unsubscribed.
%          handle = Reference to the ziDAQZoomFFT class.
%          path = Path string of the node. Use wild card to
%          select all. Alternatively also a list of path
%          strings can be specified.
%
%          ziDAQ('unsubscribe', handle, path);
%          Unsubscribe from one or several nodes. During the
%          recording process paths can not be subscribed or
%          unsubscribed.
%          handle = Reference to the ziDAQZoomFFT class.
%          path = Path string of the node. Use wild card to
%          select all. Alternatively also a list of path
%          strings can be specified.
%
%          ziDAQ('execute', handle);
%          Start the zoom FFT. Subscription or unsubscription
%          is no more possible until the zoomFFT is finished.
%
%          result = ziDAQ('finished', handle);
%          handle = Handle of the zoom FFT session.
%          Returns 1 if the zoom FFT is finished, otherwise 0.
%
%          result = ziDAQ('read', handle);
%          handle = Handle of the zoom FFT session.
%          Transfer the zoomFFT data to Matlab.
%
%          result = ziDAQ('progress', handle);
%          Report the progress of the measurement with a number
%          between 0 and 1.
%
%          ziDAQ('finish', handle);
%          Stop the zoomFFT. The zoom FFT may be restarted by
%          calling 'execute' again.
%
%          ziDAQ('clear', handle);
%          handle = Handle of the zoom FFT session.
%          Stop the current zoom FFT.
%
% Device Settings Module
%
% Device Settings Parameters
%   deviceSettings/device      string  Device whose settings are to be
%                                     saved/loaded, e.g. 'dev99'.
%   deviceSettings/path        string  Path where the settings files are to
%                                     be located. If not set, the default
%                                     settings location of the LabOne
%                                     software is used.
%   deviceSettings/filename    string  The file to which the settings are to
%                                     be saved/loaded.
%   deviceSettings/command      string  The save/load command to execute.
%                                     'save' = Read device settings and save
%                                     to file.
%                                     'load' = Load settings from file and
%                                     write to device.
%                                     'read' = Read device settings only
%                                     (no save).
%
%   handle = ziDAQ('deviceSettings', timeout);
%   timeout = Poll timeout in [ms] - DEPRECATED, ignored
%   Creates a device settings class for saving/loading device
%   settings to/from a file. Before the thread start, set the path,
%   filename and command parameters. To run the command, use the
%   execute function.

```

```

%
% result = ziDAQ('listNodes', handle, path, flags);
%           path (string) = Module parameter path
%           flags (int64) = Define which module parameters paths should be
%                           returned, set the following bits to obtain the
%                           described behaviour:
%                           int64(0) -> ZI_LIST_NODES_NONE 0x00
%                               The default flag, returning a simple
%                               listing of the given path
%                           int64(1) -> ZI_LIST_NODES_RECURSIVE 0x01
%                               Returns the paths recursively
%                           int64(2) -> ZI_LIST_NODES_ABSOLUTE 0x02
%                               Returns absolute paths
%                           int64(4) -> ZI_LIST_NODES_LEAFSONLY 0x04
%                               Returns only paths that are leafs,
%                               which means the they are at the
%                               outermost level of the tree.
%                           int64(8) -> ZI_LIST_NODES_SETTINGSONLY 0x08
%                               Returns only paths which are marked
%                               as setting
%                           Flags may also be combined, e.g., set flags to bitor(1, 2)
%                           to return paths recursively and printed as absolute paths.
%
%           ziDAQ('subscribe', handle, path);
%           Not relevant for the device settings module.
%
%           ziDAQ('unsubscribe', handle, path);
%           Not relevant for the device settings module.
%
%           ziDAQ('execute', handle);
%           Execute the command.
%
% result = ziDAQ('finished', handle);
%           handle = Handle of the device settings session.
%           Returns 1 if the command is finished, otherwise 0.
%
% result = ziDAQ('read', handle);
%           handle = Handle of the device settings session.
%           Transfer the device settings to Matlab.
%           Not relevant since device settings are saved to a file.
%
% result = ziDAQ('progress', handle);
%           Report the progress of the command with a number
%           between 0 and 1.
%
%           ziDAQ('finish', handle);
%           Stop the device settings module. The module may be restarted by
%           calling 'execute' again.
%
%           ziDAQ('clear', handle);
%           handle = Handle of the device settings session.
%           End the current device settings thread.
%
% PLL Advisor Module
%
% PLL Advisor Parameters
%   pllAdvisor/bode      struct  Output parameter. Contains the resulting bode
%                               plot of the PLL simulation.
%   pllAdvisor/calculate int     Command to calculate values. Set to 1 to start
%                               the calculation.
%   pllAdvisor/center    double  Center frequency of the PLL oscillator. The PLL
%                               frequency shift is relative to this center
%                               frequency.
%   pllAdvisor/d          int     Differential gain.
%   pllAdvisor/demodbw    int     Demodulator bandwidth used for the PLL loop
%                               filter.
%   pllAdvisor/i          double  Integral gain.

```

```

%     pllAdvisor/mode         double  Select PLL Advisor mode. Currently only one mode
%                               (open loop) is supported.
%     pllAdvisor/order        double  Demodulator order used for the PLL loop filter.
%     pllAdvisor/p             int     Proportional gain.
%     pllAdvisor/pllbw         int     Demodulator bandwidth used for the PLL loop
%                                     filter.
%     pllAdvisor/pm            int     Output parameter. Simulated phase margin of the
%                                     PLL with the current settings. The phase margin
%                                     should be greater than 45 deg and preferably
%                                     greater than 65 deg for stable conditions.
%     pllAdvisor/pmfreq        int     Output parameter. Simulated phase margin
%                                     frequency.
%     pllAdvisor/q             int     Quality factor. Currently not used.
%     pllAdvisor/rate          int     PLL Advisor sampling rate of the PLL control
%                                     loop.
%     pllAdvisor/stable        int     Output parameter. When 1, the PLL Advisor found
%                                     a stable solution with the given settings. When
%                                     0, revise your settings and rerun the PLL
%                                     Advisor.
%     pllAdvisor/targetbw      int     Requested PLL bandwidth. Higher frequencies may
%                                     need manual tuning.
%     pllAdvisor/targetfail    int     Output parameter. 1 indicates the simulated PLL
%                                     BW is smaller than the Target BW.
%
%     handle = ziDAQ('pllAdvisor', timeout);
%                               timeout = Poll timeout in [ms] - DEPRECATED, ignored
%                               Creates a PLL Advisor class for simulating the PLL in the
%                               device. Before the thread start, set the command parameters,
%                               call execute() and then set the "calculate" parameter to start
%                               the simulation.
%
%     result = ziDAQ('listNodes', handle, path, flags);
%                               path (string) = Module parameter path
%                               flags (int64) = Define which module parameters paths should be
%                                               returned, set the following bits to obtain the
%                                               described behaviour:
%                               int64(0) -> ZI_LIST_NODES_NONE 0x00
%                                       The default flag, returning a simple
%                                       listing of the given path
%                               int64(1) -> ZI_LIST_NODES_RECURSIVE 0x01
%                                       Returns the paths recursively
%                               int64(2) -> ZI_LIST_NODES_ABSOLUTE 0x02
%                                       Returns absolute paths
%                               int64(4) -> ZI_LIST_NODES_LEAFSONLY 0x04
%                                       Returns only paths that are leafs,
%                                       which means the they are at the
%                                       outermost level of the tree.
%                               int64(8) -> ZI_LIST_NODES_SETTINGSONLY 0x08
%                                       Returns only paths which are marked
%                                       as setting
%                               Flags may also be combined, e.g., set flags to bitor(1, 2)
%                               to return paths recursively and printed as absolute paths.
%
%     ziDAQ('subscribe', handle, path);
%                               Subscribe to one or several nodes.
%
%     ziDAQ('unsubscribe', handle, path);
%                               Unsubscribe from one or several nodes..
%
%     ziDAQ('execute', handle);
%                               Start the PLL Advisor.
%
%     result = ziDAQ('finished', handle);
%                               handle = Handle of the PLL Advisor session.
%                               Returns 1 if the command is finished, otherwise 0.
%
%     result = ziDAQ('read', handle);

```

```

%             handle = Handle of the PLL Advisor session.
%             Read pllAdvisor data. If the simulation is still ongoing only a
%             subset of the data is returned.
%
% result = ziDAQ('progress', handle);
%             Report the progress of the command with a number
%             between 0 and 1.
%
%             ziDAQ('finish', handle);
%             Stop the PLL Advisor module.
%
%             ziDAQ('clear', handle);
%             handle = Handle of the PLL Advisor session.
%             End the current PLL Advisor thread.
%
% PID Advisor Module
%
% PID Advisor Parameters
%     pidAdvisor/advancedmode    int    Disable automatic calculation of the
%                                     start and stop value.
%     pidAdvisor/auto            int    Automatic response calculation triggered
%                                     by parameter change.
%     pidAdvisor/bode            struct  Output parameter. Contains the resulting
%                                     bode plot of the PID simulation.
%     pidAdvisor/bw              double  Output parameter. Calculated system
%                                     bandwidth.
%     pidAdvisor/calculate       int    In/Out parameter. Command to calculate
%                                     values. Set to 1 to start the
%                                     calculation.
%     pidAdvisor/display/freqstart double  Start frequency for Bode plot.
%                                     For disabled advanced mode the start
%                                     value is automatically derived from the
%                                     system properties.
%     pidAdvisor/display/freqstop double  Stop frequency for Bode plot.
%     pidAdvisor/display/timestart double  Start time for step response.
%     pidAdvisor/display/timestop double  Stop time for step response.
%     pidAdvisor/dut/bw          double  Bandwidth of the DUT (device under test).
%     pidAdvisor/dut/damping     double  Damping of the second order
%                                     low pass filter.
%     pidAdvisor/dut/delay       double  IO Delay of the feedback system
%                                     describing the earliest response for
%                                     a step change.
%     pidAdvisor/dut/fcenter     double  Resonant frequency of the of the modelled
%                                     resonator.
%     pidAdvisor/dut/gain        double  Gain of the DUT transfer function.
%     pidAdvisor/dut/q           double  quality factor of the modelled resonator.
%     pidAdvisor/dut/source      int    Type of model used for the external
%                                     device to be controlled by the PID.
%                                     source = 1: Low-pass first order
%                                     source = 2: Low-pass second order
%                                     source = 3: Resonator frequency
%                                     source = 4: Internal PLL
%                                     source = 5: VCO
%                                     source = 6: Resonator amplitude
%     pidAdvisor/impulse         struct  Output parameter. Impulse response
%                                     (not yet supported).
%     pidAdvisor/index           int    PID index for parameter detection.
%     pidAdvisor/pid/autobw      int    Adjusts the demodulator bandwidth to fit
%                                     best to the specified target bandwidth
%                                     of the full system.
%     pidAdvisor/pid/d           double  In/Out parameter. Differential gain.
%     pidAdvisor/pid/dlimittimeconstant double  In/Out parameter. Differential filter
%                                     timeconstant.
%     pidAdvisor/pid/i           double  In/Out parameter. Integral gain.
%     pidAdvisor/pid/mode        double  Select PID Advisor mode. Mode value is
%                                     bit coded, bit 0: P, bit 1: I, bit 2: D,

```



```

% bit 3: D filter limit.
% pidAdvisor/pid/p double In/Out parameter. Proportional gain.
% pidAdvisor/pid/rate double In/Out parameter. PID Advisor sampling
% rate of the PID control loop.
% pidAdvisor/pid/targetbw double PID system target bandwidth.
% pidAdvisor/pm double Output parameter. Simulated phase margin
% of the PID with the current settings.
% The phase margin should be greater than
% 45 deg and preferably greater than 65 deg
% for stable conditions.
% pidAdvisor/pmfreq double Output parameter. Simulated phase margin
% frequency.
% pidAdvisor/stable int Output parameter. When 1, the PID Advisor
% found a stable solution with the given
% settings. When 0, revise your settings
% and rerun the PID Advisor.
% pidAdvisor/step struct Output parameter. Contains the resulting
% step response plot of the PID simulation.
% pidAdvisor/targetbw double Requested PID bandwidth. Higher
% frequencies may need manual tuning.
% pidAdvisor/targetfail int Output parameter. 1 indicates the
% simulated PID BW is smaller than the
% Target BW.
% pidAdvisor/tf/closedloop int Switch the response calculation mode
% between closed or open loop.
% pidAdvisor/tf/input int Start point for the plant response
% simulation for open or closed loops.
% pidAdvisor/tf/output int End point for the plant response
% simulation for open or closed loops.
% pidAdvisor/tune int Optimize the PID parameters so that
% the noise of the closed-loop
% system gets minimized.
%
% handle = ziDAQ('pidAdvisor', timeout);
% timeout = Poll timeout in [ms] - DEPRECATED, ignored
% Creates a PID Advisor class for simulating the PID in the
% device. Before the thread start, set the command parameters,
% call execute() and then set the "calculate" parameter to start
% the simulation.
%
% result = ziDAQ('listNodes', handle, path, flags);
% handle = Handle of the PID Advisor session.
% path (string) = Module parameter path
% flags (int64) = Define which module parameters paths should be
% returned, set the following bits to obtain the
% described behaviour:
% int64(0) -> ZI_LIST_NODES_NONE 0x00
% The default flag, returning a simple
% listing of the given path
% int64(1) -> ZI_LIST_NODES_RECURSIVE 0x01
% Returns the paths recursively
% int64(2) -> ZI_LIST_NODES_ABSOLUTE 0x02
% Returns absolute paths
% int64(4) -> ZI_LIST_NODES_LEAFSONLY 0x04
% Returns only paths that are leafs,
% which means the they are at the
% outermost level of the tree.
% int64(8) -> ZI_LIST_NODES_SETTINGSONLY 0x08
% Returns only paths which are marked
% as settings
%
% Flags may also be combined, e.g., set flags to bitor(1, 2)
% to return paths recursively and printed as absolute paths.
%
% ziDAQ('subscribe', handle, path);
% handle = Handle of the PID Advisor session.
% Subscribe to one or several nodes.
%

```

```

%      ziDAQ('unsubscribe', handle, path);
%          handle = Handle of the PID Advisor session.
%          Unsubscribe from one or several nodes..
%
%      ziDAQ('get', handle, path);
%          handle = Handle of the PID Advisor session.
%          path = Path string of the node.
%
%      ziDAQ('execute', handle);
%          handle = Handle of the PID Advisor session.
%          Starts the pidAdvisor if not yet running.
%
%      ziDAQ('trigger', handle);
%          Not applicable to this module.
%
%      result = ziDAQ('finished', handle);
%          handle = Handle of the PID Advisor session.
%          Returns 1 if the command is finished, otherwise 0.
%
%      result = ziDAQ('read', handle);
%          handle = Handle of the PID Advisor session.
%          Read pidAdvisor data. If the simulation is still ongoing only a
%          subset of the data is returned.
%
%      result = ziDAQ('progress', handle);
%          handle = Handle of the PID Advisor session.
%          Report the progress of the command with a number
%          between 0 and 1.
%
%      ziDAQ('finish', handle);
%          handle = Handle of the PID Advisor session.
%          Stop the PID Advisor module.
%
%      ziDAQ('clear', handle);
%          handle = Handle of the PID Advisor session.
%          End the current PID Advisor thread.
%
%      ziDAQ('save', handle);
%          Save the measured data to a file.
%          handle = Handle of the PID Advisor session.
%          [filename] = File name string (without extension)..
%
% Debugging Functions
%
%      ziDAQ('setDebugLevel', debuglevel);
%          debuglevel (int) = Debug level (trace:0, info:1, debug:2,
%          warning:3, error:4, fatal:5, status:6).
%          Enables debug log and sets the debug level.
%
%      ziDAQ('writeDebugLog', severity, message);
%          severity (int) = Severity (trace:0, info:1, debug:2, warning:3,
%          error:4, fatal:5, status:6).
%          message (str) = Message to output to the log.
%          Outputs message to the debug log (if enabled).
%
%      ziDAQ('logOn', flags, filename, [style]);
%          flags = LOG_NONE:          0x00000000
%                  LOG_SET_DOUBLE:   0x00000001
%                  LOG_SET_INT:       0x00000002
%                  LOG_SET_BYTE:      0x00000004
%                  LOG_SYNC_SET_DOUBLE: 0x00000010
%                  LOG_SYNC_SET_INT:  0x00000020
%                  LOG_SYNC_SET_BYTE: 0x00000040
%                  LOG_GET_DOUBLE:    0x00000100
%                  LOG_GET_INT:       0x00000200
%                  LOG_GET_BYTE:      0x00000400
%                  LOG_GET_DEMOD:     0x00001000

```

```
%
%          LOG_GET_DIO:          0x00002000
%          LOG_GET_AUXIN:       0x00004000
%          LOG_LISTNODES:       0x00010000
%          LOG_SUBSCRIBE:       0x00020000
%          LOG_UNSUBSCRIBE:     0x00040000
%          LOG_GET_AS_EVENT:    0x00080000
%          LOG_UPDATE:         0x00100000
%          LOG_POLL_EVENT:     0x00200000
%          LOG_POLL:           0x00400000
%          LOG_ALL :           0xffffffff
%
filename = Log file name
[style] = LOG_STYLE_TELNET: 0 (default)
%          LOG_STYLE_MATLAB: 1
%          LOG_STYLE_PYTHON: 2
%
Log all API commands sent to the Data Server. This is useful
%
% for debugging.
%
ziDAQ('logOff');
%
Turn of message logging.
%
```

Chapter 4. Python Programming

Python is open source software, freely available for download from [Python's official website](#). Python is a high-level programming language with an extensive standard library renowned for its "batteries included" approach. Combined with the [NumPy](#) package for scientific computing, Python is a powerful computational tool for scientists that does not require expensive software licenses. The Zurich Instruments LabOne Python API, also known as `ziPython` enables the user to configure and stream data from their instrument directly into Python.

This chapter aims to help you get started using the Zurich Instruments LabOne Python API, `ziPython`, to control your instrument, please refer to:

- [Section 4.1](#) for help Installing the LabOne Python API.
- [Section 4.2](#) for help Getting Started with the LabOne Python API and Running the Examples.
- [Section 4.3](#) for LabOne Python API Tips and Tricks.
- [Section 4.4](#) for the LabOne Python API (`ziPython`) Command Reference.

Note

This chapter and the provided examples are not intended to be a Python tutorial. For help getting started with Python itself, see either the [Python Tutorial](#) or one of the many online resources, for example, the [learnpython.org](#). The [Interactive Python Course](#) is an interesting resource for those already familiar with Python basics.

4.1. Installing the LabOne Python API

4.1.1. Requirements

In order to install and use the LabOne Python API you require:

1. Either a Python 2.7 or a Python 3.5 installation on either Windows or Linux.
2. The [NumPy](#) python package installed for your Python installation.
3. The correct version of [ziPython](#) for your Python version and platform, available from the Zurich Instruments [download page](#).

Note

Linux users must also ensure they download the version of [ziPython](#) that is Unicode compatible with their Linux distribution's Python installation, see [Section 4.1.4](#) for help determining which version is required.

Note

Important: If you your system already has an existing [ziPython](#) installation older than version 14.08, please be sure to either manually uninstall [ziPython](#) or manually remove the existing [zhinst](#) installation folder. This is due to improvements in the [zhinst](#) package structure in 14.08 (examples for different device classes are now organized in separate module/sub-directories) and the Python installer simply overwrites the existing installation, leading to a duplication of some files. For help locating `[PYTHONROOT]\lib\site-packages\zhinst\` on your system, please see the section called "[Locating the zhinst Installation Folder and Examples](#)".

4.1.2. Recommended Python Packages

The following Python packages can additionally be useful for programming with the LabOne Python API:

1. [Matplotlib](#) - recommended to plot the output from many of [ziPython](#)'s examples.
2. [SciPy](#) - recommended to load data saved from the LabOne UI in binary Matlab format (.mat).

Note

Unofficial pre-compiled 32-bit and 64-bit Windows binaries of NumPy, SciPy and matplotlib are available from Christoph Gohlke's [pythonlibs page](#).

4.1.3. Windows Installation

To install ziPython on Windows execute the .msi installer available from the Zurich Instruments [download page](#). It will guide you through the installation process as displayed in the following screenshots.

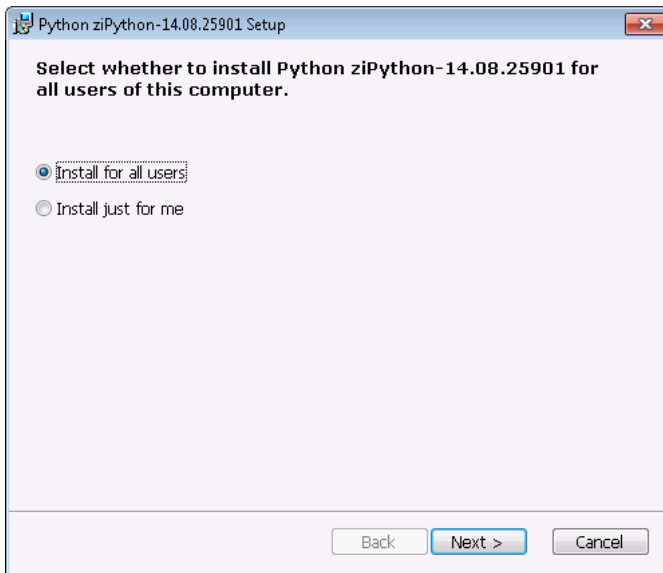


Figure 4.1. Windows ziPython installation: Step 1.

If multiple Python Installations are available on your system, the installer will ask which Python version the ziPython package should be installed. The ziPython package will be installed in selected versions in the folder `[PYTHONROOT]\lib\site-packages\zhinst\`.

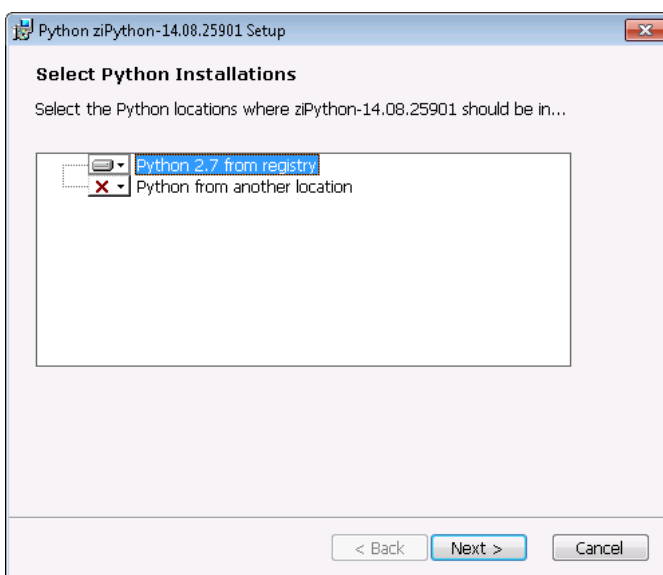


Figure 4.2. Windows ziPython installation: Step 2.

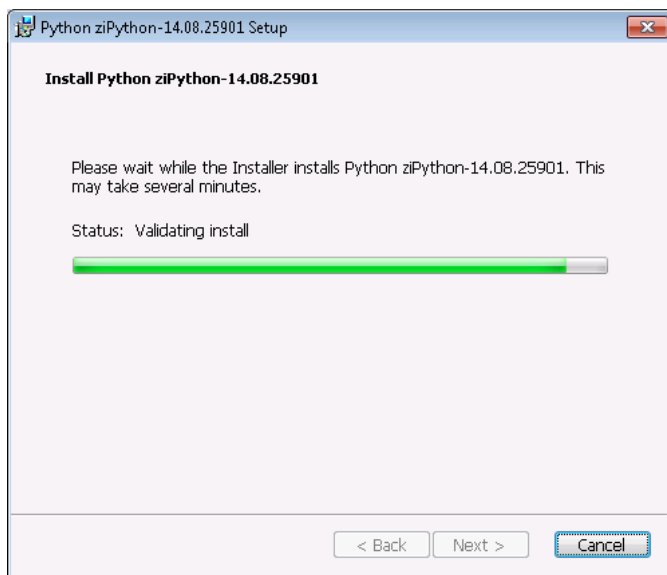


Figure 4.3. Windows ziPython installation: Step 3.

4.1.4. Linux Installation

In addition to the [requirements above](#) and selecting the correct version of ziPython for your Python distribution (2.7 or 3.5) and platform (32-bit or 64-bit), on Linux the correct Unicode version must also be installed. This is because some Python distributions on Linux are compiled to use UCS-2 character encoding, whereas some use UCS-4.

Determining the correct Unicode version of ziPython for your Python distribution

In order to determine which version of Unicode your Python distribution uses, please type the following commands in the interactive shell of your target Python distribution:

```
>>> import sys
>>> print sys.maxunicode
```

If the last command prints:

- 65535, use the UCS-2 version of ziPython,
- 1114111, use the UCS-4 version of ziPython.

Note

The installation needs root access rights. If you do not have these permissions, ask your system administrator for help.

To install ziPython on a Debian-derived distribution such as Ubuntu perform the following steps:

1. If required, install Python, NumPy and matplotlib (with elevated access rights):

```
$ sudo apt-get install python python-numpy python-matplotlib
```

2. Unpack the ziPython software bundle:

```
$ tar xzf ziPython-[version]-[build]-[linux32|linux64].tar.gz
```

3. Change directory into the unpacked folder and run the setup script `setup.py` as following:

```
$ cd ziPython-[version]-[build]-(linux32|linux64)
$ python setup.py build
$ sudo python setup.py install --install-layout=deb # Elevated access rights
```

It's possible to skip the `build` step (`install` will automatically perform this step), but splitting the steps avoids creating a directory in your user space which is owned by `root`.

4.2. Getting Started with the LabOne Python API

This section introduces the user to the LabOne Python API.

4.2.1. Contents of the LabOne Python API

Alongside the driver for interfacing with your Zurich Instruments device, the LabOne Python API includes utility functions and examples. See:

- [Section 4.4.1](#) to see which examples are available in `ziPython`.
- [Section 4.4.2](#) to see which utility functions are available in `ziPython`.

4.2.2. Using the Built-in Documentation

`ziPython`'s built-in documentation can be accessed using the `help` command in a python interactive shell:

- On module level:

```
>>> import zhinst.ziPython as ziPython
>>> help(ziPython)
```

- On class level, for example, for the Sweeper Module:

```
>>> import zhinst.ziPython as ziPython
>>> help(ziPython.ziDAQSweeper)
```

- On function level, for example, for the `ziDAQServer.poll` method:

```
>>> import zhinst.ziPython as ziPython
>>> help(ziPython.ziDAQServer.poll)
```

See [Section 4.4, LabOne Python API \(ziPython\) Command Reference](#) for a printer friendly version of the built-in documentation.

4.2.3. Running the Examples

Prerequisites for running the Python examples:

1. The `zhinst` package is installed as described above in [Section 4.1](#).
2. The Data Server program is running and the instrument is discoverable, this is the case if the instrument can be seen in the User Interface.
3. Signal Output 1 of the instrument is connected to Signal Input 1 via a BNC cable; many of the Python examples measure on this hardware channel.

It's also recommended to install the [Matplotlib](#) Python package in order to plot the data obtained in many of the examples, see [Section 4.1.2](#).

The API examples are available in the module `zhinst.examples`, which is organized into sub-modules according to the target Instrument class:

- `zhinst.examples.common`: examples compatible with any class of instrument,
- `zhinst.examples.uhf`: examples only compatible with the UHF Lock-in Amplifier,
- `zhinst.examples.hf2`: examples only compatible with HF2 Series Instruments.

All the examples follow the same structure and take one input argument: The device ID of the instrument to run the example with. The recommended way to run a `ziPython` example is to

import the example's module in an interactive shell and call the `run_example()` function. For example, to run the `zoomFFT` Module example:

```
>>> import zhinst.examples
>>> # Use do_plot=False if matplotlib is unavailable
>>> zhinst.examples.common.example_spectrum.run_example('dev123', do_plot=True);
```

The example should produce some output in the Python shell, such as:

```
Will perform 1 zoomFFTs
Individual zoomFFT 100.00 complete.
sample contains 1 zoomFFTs
Number of lines in first zoomFFT: 65535
```

Most examples will also plot the retrieved data using matplotlib, see [Figure 4.4](#) for an example. If you encounter an error message please ensure that the [above prerequisites](#) are fulfilled.

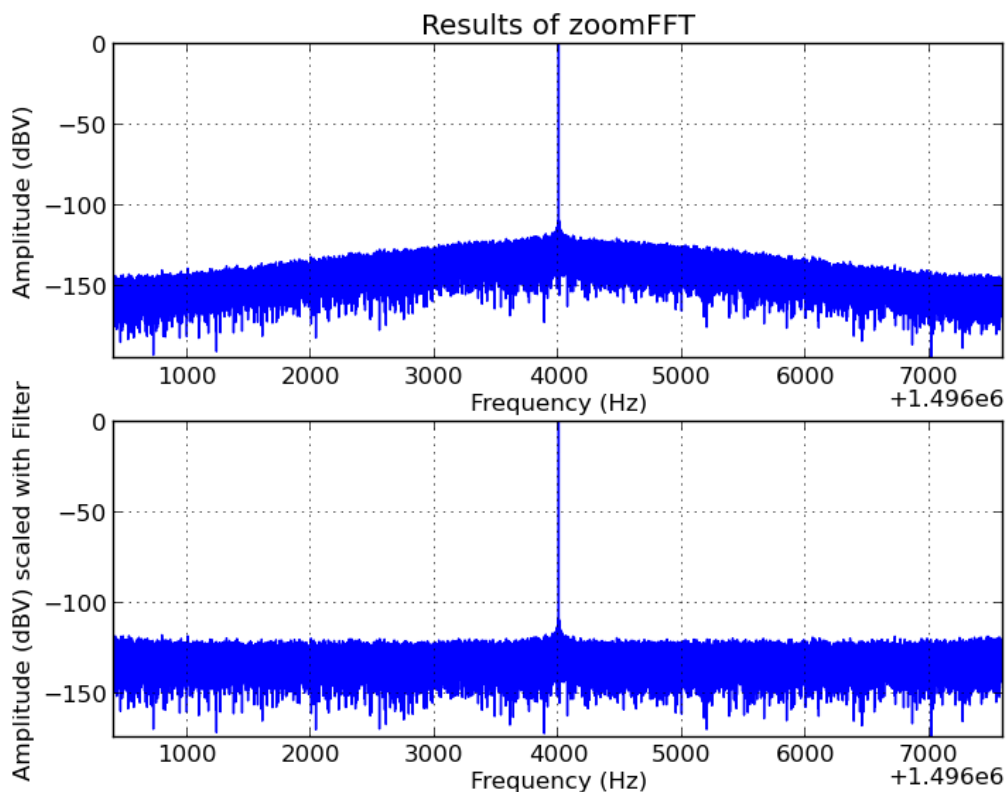


Figure 4.4. The plot produced by the LabOne Python API example `example_spectrum.py`; the plots show the results of an FFT performed with `ziCore`'s `zoomFFT` module on demodulator output obtained over a simple feedback cable.

Exploring which Examples are available

Python's help system can be used to see which examples are available for a particular device class; when help is called on the module the available examples are listed under the "Package Contents" section. For example, for the `zhinst.examples.common` package:

```
>>> help('zhinst.examples.common')

Help on package zhinst.examples.common in zhinst.examples:

NAME
```

```
zhinst.examples.common - Zurich Instruments LabOne Python API Examples (for any
instrument class).

PACKAGE CONTENTS
  example_connect
  example_connect_config
  example_pid_advisor
  example_poll
  example_record_edge_trigger
  example_save_device_settings_expert
  example_save_device_settings_simple
  example_scope
  example_spectrum
  example_sweeper
  example_swtrigger_edge
  example_swtrigger_grid
  example_swtrigger_trackingedge
  example_zoomfft

DATA
  __all__ = ['example_connect', 'example_connect_config', 'example_pid_a...

FILE
  /home/ci/.pyenv/versions/3.5.1/lib/python3.5/site-packages/zhinst/examples/
  common/__init__.py
```

Locating the **zhinst** Installation Folder and Examples

The examples distributed with the **zhinst** package can serve as a starting point to program your own measurement needs. The example python files, however, are generally not installed in user space. In order to ensure that you have sufficient permission to edit the examples and that your modifications are not overwritten by a later upgrade of the **zhinst** package, please copy them to your own user space before editing them.

The examples are contained in a subfolder of the **zhinst** package installation folder

```
[PYTHONROOT]\lib\site-packages\zhinst\
```

If you are unsure about the location of your **PYTHONROOT**, the `__path__` attribute of the **zhinst** module can be used in order to determine its location, for example,

```
>>> import zhinst
>>> print zhinst.__path__
```

will output something similar to:

```
C:\Python27\lib\site-packages\zhinst
```

4.2.4. Using **ziCore** Modules in the LabOne Python API

In the LabOne Python API **ziCore** [Modules](#) are configured and controlled by instantiating an object of the Module's class. For example, in order to use the [Sweeper Module](#) a sweeper object is created as following:

```
>>> daq = ziPython.ziDAQServer('localhost', 8004, 5) # Create a connection to the
  Data Server
  (running on the same PC as the API client)
>>> sweeper = daq.sweep();
```

Note, that since creating a Module object without an API connection to the Data Server does not make sense, the Sweeper object is instantiated via the `sweep` method of the `ziDAQServer` class, not directly from the `ziDAQSweeper` class.

The Module's parameters are configured using the Module's `set` method and specifying a path, value pair, for example:

```
>>> sweeper.set('sweep/start', 1.2e5);
```

The parameters can be read-back using the `get` method, which supports wildcards, for example:

```
>>> sweep_params = sweeper.get('sweep/*');
```

The variable `sweep_params` now contains a dictionary of all the Sweeper's parameters. The other main Module commands are similarly used, e.g., `sweeper.execute()`, to start the sweeper. See [Section 2.1.2](#) for more help with Modules and a description of their parameters.

4.2.5. Enabling Logging in the LabOne Python API

Logging from the API is not enabled by default upon initializing a server session with `ziPython`, it must be enabled (after using `connect`) with the `setDebugLevel` command. For example,

```
>>> daq.setDebugLevel(0)
```

sets the API's logging level to 0, which provides the most verbose logging output. The other log levels are defined as following:

```
trace:0, info:1, debug:2, warning:3, error:4, fatal:5, status:6.
```

It is also possible for the user to write their own messages directly to `ziPython`'s log using the `writeDebugLog` command. For example to write a log message of `info` severity level:

```
>>> daq.writeDebugLog(1, 'Hello log!')
```

On Windows the logs can be found by navigating to the Zurich Instruments Logs" folder entry in the Windows Start Menu: Programs → Zurich Instruments → LabOne Servers → Logs. This will open an Explorer window displaying folders containing log files from various LabOne components, in particular, the `ziPythonLog` folder contains logs from the LabOne Python API. On Linux, the logs can be found at `"/tmp/ziPythonLog_USERNAME"`, where `"USERNAME"` is the same as the output of the `"whoami"` command.

4.3. LabOne Python API Tips and Tricks

In this section some tips and tricks for working with the LabOne Python API are provided.

Data Structures returned by **ziPython**.

The output arguments that **ziPython** returns are designed to use the native data structures that Python users are familiar with and that reflect the data's location in the instruments node hierarchy. For example, when the `poll` command returns data from the instruments fourth demodulator (located in the node hierarchy as `/dev123/demods/3/sample`), the output argument contains a tree of nested dictionaries in which the data can be accessed by

```
data = daq.poll( poll_length, poll_timeout);
x = data['dev123']['demods']['4']['sample']['x'];
y = data['dev123']['demods']['4']['sample']['y'];
```

Tell **poll** to return a flat dictionary

By default, the data returned by `poll` is contained in a tree of nested dictionaries that closely mimics the tree structure of the instrument node hierarchy. By setting the optional fifth argument of `poll` to `True`, the data will be a flat dictionary. This can help avoid many nested `if` statements in order to check that the expected data was returned by `poll`. For example:

```
daq.subscribe('/dev123/demods/0/sample')
flat_dictionary_key = False
data = daq.poll(0.1, 200, 1, flat_dictionary_key)
if 'dev123' in data:
    if 'demods' in data['device']:
        if '0' in data['device']['demods']:
            # access the demodulator data:
            x = data['dev123']['demods']['0']['sample']['x']
            y = data['dev123']['demods']['0']['sample']['y']
```

Could be rewritten more concisely as:

```
daq.subscribe('/dev123/demods/0/sample')
flat_dictionary_key = True
data = daq.poll(0.1, 200, 1, flat_dictionary_key)
if '/dev123/demods/0/sample' in data:
    # access the demodulator data:
    x = data['/dev123/demods/0/sample']['x']
    y = data['/dev123/demods/0/sample']['y']
```

Use the Utility Routines to load Data saved from the LabOne UI and ziControl in Python.

The utilities package `zhinst.utils` contains several routines to help loading `.csv` or `.mat` files saved from either the LabOne User Interface or `ziControl` into Python. These functions are generally minimal wrappers around `NumPy` (`genfromtxt()`) or `SciPy` (`loadmat()`) routines. However, the function `load_labone_demod_csv()` is optimized to load demodulator data saved in `.csv` format by the LabOne UI (since it specifies the `.csv` columns' dtypes explicitly) and the function `load_zicontrol_zibin()` can directly load data saved in binary format from `ziControl`. See [Section 4.4.2](#) for reference documentation on these commands.

4.4. LabOne Python API (ziPython) Command Reference

The following reference documentation for ziPython is available in from within a python session using python's help (see [Section 4.2.2](#)) command; It is included here for convenience.

The documentation is grouped by module and class as following:

- [Help for the zhinst Python Package](#)
- [Help for zhinst's Utility Functions](#)
- [Help for ziPython's ziDAQServer class](#)
- [Help for ziPython's ziDeviceSettings class](#)
- [Help for ziPython's ziDAQSweeper class](#)
- [Help for ziPython's ziDAQZoomFFT class](#)
- [Help for ziPython's ziDAQRecorder class](#)
- [Help for ziPython's ziPllAdvisor class](#)
- [Help for ziPython's ziPidAdvisor class](#)

4.4.1. Help for the **zhinst** Python Package

```
>>> help('zhinst')

Help on package zhinst:

NAME
    zhinst - Zurich Instruments LabOne Python API

DESCRIPTION
    Contains the API driver, utility functions and examples for Zurich Instruments
    devices.

PACKAGE CONTENTS
    examples (package)
    utils
    ziPython

DATA
    __all__ = ['ziPython', 'utils']

FILE
    /home/ci/.pyenv/versions/3.5.1/lib/python3.5/site-packages/zhinst/__init__.py
```

4.4.2. Help for zhinst's Utility Functions

```
>>> help('zhinst.utils')

Help on module zhinst.utils in zhinst:

NAME
    zhinst.utils - Zurich Instruments LabOne Python API Utility Functions.

DESCRIPTION
    This module provides basic utility functions for:

    - Creating an API session by connecting to an appropriate Data Server.
```

- Detecting devices.
- Loading and saving device settings.
- Loading data saved by either the Zurich Instruments LabOne User Interface or ziControl into Python as numpy structured arrays.

FUNCTIONS

`autoConnect(default_port=None, api_level=None)`

Try to connect to a Zurich Instruments Data Server with an attached available UHF or HF2 device.

Important: `autoConnect()` does not support MFLI devices.

Args:

`default_port` (int, optional): The default port to use when connecting to the Data Server (specify 8005 for the HF2 Data Server and 8004 for the UHF Data Server).

`api_level` (int, optional): The API level to use, either 1, 4 or 5. HF2 only supports Level 1, Level 5 is recommended for UHF and MFLI devices.

Returns:

`ziDAQServer`: An instance of the `ziPython.ziDAQServer` class that is used for communication to the Data Server.

Raises:

`RuntimeError`: If no running Data Server is found or no device is found that is attached to a Data Server.x

If `default_port` is not specified (`=None`) then first try to connect to a HF2, if no server devices are found then try to connect to an UHF. This behaviour is useful for the API examples. If we cannot connect to a server and/or detect a connected device raise a `RuntimeError`.

If `default_port` is 8004 try to connect to a UHF; if it is 8005 try to connect to an HF2. If no server and device is detected on this port raise a `RuntimeError`.

`autoDetect(daq, exclude=None)`

Return a string containing the first device ID (not in the exclude list) that is attached to the Data Server connected via `daq`, an instance of the `ziPython.ziDAQServer` class.

Args:

`daq` (`ziDAQServer`): An instance of the `ziPython.ziDAQServer` class (representing an API session connected to a Data Server).

`exclude` (list of str, optional): A list of strings specifying devices to exclude. `autoDetect()` will not return the name of a device in this list.

Returns:

A string specifying the first device ID not in exclude.

Raises:

`RuntimeError`: If no device was found.

`RuntimeError`: If `daq` is not an instance of `ziPython.ziDAQServer`.

Example:

```

zhinst.utils
daq = zhinst.utils.autoConnect()
device = zhinst.utils.autoDetect(daq)

bw2tc(bandwidth, order)
Convert the demodulator 3 dB bandwidth to its equivalent timeconstant for the
specified demodulator order.

Inputs:

    bandwidth (double): The demodulator 3dB bandwidth to convert.

    order (int): The demodulator order (1 to 8) for which to convert the
    bandwidth.

Output:

    timeconstant (double): The equivalent demodulator timeconstant.

bwtc_scaling_factor(order)
Return the appropriate scaling factor for bandwidth to timeconstant
conversion for the provided demodulator order.

check_for_sampleloss(timestamps)
Check whether timestamps are equidistantly spaced, if not, it is an
indication that sampleloss has occurred whilst recording the demodulator
data.

This function assumes that the timestamps originate from continuously saved
demodulator data, during which the demodulator sampling rate was not
changed.

Arguments:

    timestamp (numpy array): a 1-dimensional array containing
    demodulator timestamps

Returns:

    idx (numpy array): a 1-dimensional array indicating the indices in
    timestamp where sampleloss has occurred. An empty array is returned if no
    sampleloss was present.

create_api_session(device_serial, maximum_supported_apilevel,
required_devtype='.*', required_options=None, required_err_msg='')
Create an API session for the specified device.

Args:

    device_serial (str): A string specifying the device serial number. For
    example, 'uhf-dev2123' or 'dev2123'.

    maximum_supported_apilevel (int): The maximum API Level that is supported
    by the code where the returned API session will be used. The maximum API
    Level you may use is defined by the device class. HF2 only supports API
    Level 1 and other devices support API Level 5. You should try to use the
    maximum level possible to enable extended API features.

    required_devtype (str): The required device type, e.g., 'HF2LI' or
    'MFLI'. This is given by the value of the device node
    '/devX/features/devtype' or the 'devicetype' discovery property. Raise an
    exception if the specified device_serial's devtype does not match the
    'required_devtype'.

    required_options (list of str|None): The required device option set. E.g.,
    ['MF', 'PID']. This is given by the value of the device node
    '/devX/features/options' or the 'options' discovery property. Raise an

```


exception if the specified `device_serial`'s option set does contain the ``required_options``.

`required_error_msg (str)` : An additional error message to print if either the device specified by the ``device_serial`` is not the ``required_devtype`` or does not have the ``required_options``.

Returns:

`daq (ziDAQServer)`: An instance of the `ziPython.ziDAQServer` class (representing an API session connected to a Data Server).

`device (str)`: The device's ID, this is the string that specifies the device's node branch in the data server's node tree.

`props (dict)`: The device's discovery properties as returned by the `ziDiscovery.get()` method.

`default_output_mixer_channel(discovery_props, output_channel=0)`
Return an instrument's default output mixer channel based on the specified ``devicetype`` and ``options`` discovery properties and the hardware output channel.

This utility function is used by the `ziPython` examples and returns a node available under the `/devX/sigouts/0/{amplitudes,enables}/` branches.

Args:

`discovery_props (dict)`: A device's discovery properties as returned by `ziDiscovery`'s `get()` method.

`output_channel (int, optional)`: The zero-based index of the hardware output channel for which to return an output mixer channel.

Returns:

`output_mixer_channel (int)`: The zero-based index of an available signal output mixer channel.

Raises:

Exception: If an invalid signal input index was provided.

`devices(daq)`

Return a list of strings containing the device IDs that are attached to the Data Server connected via `daq`, an instance of the `ziPython.ziDAQServer` class. Returns an empty list if no devices are found.

Args:

`daq (ziDAQServer)`: An instance of the `ziPython.ziDAQServer` class (representing an API session connected to a Data Server).

Returns:

A list of strings of connected device IDs. The list is empty if no devices are detected.

Raises:

`RuntimeError`: If `daq` is not an instance of `ziPython.ziDAQServer`.

Example:

```
import zhinst.utils
```

```
    daq = zhinst.utils.autoConnect() # autoConnect not supported for MFLI
devices
    device = zhinst.utils.autoDetect(daq)

get_default_settings_path(daq)
    Return the default path used for settings by the ziDeviceSettings module.
```

Arguments:

daq (instance of ziDAQServer): A ziPython API session.

Returns:

settings_path (str): The default ziDeviceSettings path.

```
load_labone_csv(fname)
    Load a CSV file containing generic data as saved by the LabOne User
    Interface into a numpy structured array.
```

Arguments:

filename (str): The filename of the CSV file to load.

Returns:

sample (numpy ndarray): A numpy structured array of shape (num_points,) whose field names correspond to the column names in the first line of the CSV file. num_points is the number of lines in the CSV file - 1.

Example:

```
import zhinst.utils
# Load the CSV file of PID error data (node: /dev2004/pids/0/error)
data = zhinst.utils.load_labone_csv('dev2004_pids_0_error_00000.csv')
import matplotlib.pyplot as plt
# Plot the error
plt.plot(data['timestamp'], data['value'])

load_labone_demod_csv(fname, column_names=('chunk', 'timestamp', 'x', 'y',
'freq', 'phase', 'dio', 'trigger', 'auxin0', 'auxin1'))
    Load a CSV file containing demodulator samples as saved by the LabOne User
    Interface into a numpy structured array.
```

Arguments:

fname (file or str): The file or filename of the CSV file to load.

column_names (list or tuple of str, optional): A list (or tuple) of column names to load from the CSV file. Default is to load all columns.

Returns:

sample (numpy ndarray): A numpy structured array of shape (num_points,) whose field names correspond to the column names in the first line of the CSV file. num_points is the number of lines in the CSV file - 1.

Example:

```
import zhinst.utils
sample =
zhinst.utils.load_labone_demod_csv('dev2004_demods_0_sample_00000.csv',
('timestamp', 'x', 'y'))
import matplotlib.pyplot as plt
import numpy as np
plt.plot(sample['timestamp'], np.abs(sample['x'] + 1j*sample['y']))

load_labone_mat(filename)
```

A wrapper function for loading a MAT file as saved by the LabOne User Interface with `scipy.io's loadmat()` function. This function is included mainly to document how to work with the data structure return by `scipy.io.loadmat()`.

Arguments:

`filename (str)`: the name of the MAT file to load.

Returns:

`data (dict)`: a nested dictionary containing the instrument data as specified in the LabOne User Interface. The nested structure of ``data`` corresponds to the path of the data's node in the instrument's node hierarchy.

Further comments:

The MAT file saved by the LabOne User Interface (UI) is a Matlab V5.0 data file. The LabOne UI saves the specified data using native Matlab data structures in the same format as are returned by commands in the LabOne Matlab API. More specifically, these data structures are nested Matlab structs, the nested structure of which correspond to the location of the data in the instrument's node hierarchy.

Matlab structs are returned by `scipy.io.loadmat()` as dictionaries, the name of the struct becomes a key in the dictionary. However, as for all objects in MATLAB, structs are in fact arrays of structs, where a single struct is an array of shape (1, 1). This means that each (nested) dictionary that is returned (corresponding to a node in node hierarchy) is loaded by `scipy.io.loadmat` as a 1-by-1 array and must be indexed as such. See the ``Example`` section below.

For more information please refer to the following link:

<http://docs.scipy.org/doc/scipy/reference/tutorial/io.html#matlab-structs>

Example:

```
device = 'dev88'
# See ``Further explanation`` above for a comment on the indexing:
timestamp = data[device][0,0]['demods'][0,0]['sample'][0,0]['timestamp'][0]
x = data[device][0,0]['demods'][0,0]['sample'][0,0]['x'][0]
y = data[device][0,0]['demods'][0,0]['sample'][0,0]['y'][0]
import matplotlib.pyplot as plt
import numpy as np
plt.plot(timestamp, np.abs(x + 1j*y))

# If multiple demodulator's are saved, data from the second demodulator,
# e.g., is accessed as following:
x = data[device][0,0]['demods'][0,1]['sample'][0,0]['x'][0]
```

`load_settings(daq, device, filename)`
Load a LabOne settings file to the specified device. This function is synchronous; it will block until loading the settings has finished.

Arguments:

`daq (instance of ziDAQServer)`: A ziPython API session.

`device (str)`: The device ID specifying where to load the settings, e.g., 'dev123'.

`filename (str)`: The filename of the xml settings file to load. The filename can include a relative or full path.

Raises:

RuntimeError: If loading the settings times out.

Examples:

```
import zhinst.utils as utils
daq = utils.autoConnect()
dev = utils.autoDetect(daq)

# Then, e.g., load settings from a file in the current directory:
utils.load_settings(daq, dev, 'my_settings.xml')
# Then, e.g., load settings from the default LabOne settings path:
filename = 'default_ui.xml'
path = utils.get_default_settings_path(daq)
utils.load_settings(daq, dev, path + os.sep + filename)

load_zicontrol_csv(filename, column_names=('t', 'x', 'y', 'freq', 'dio',
'auxin0', 'auxin1'))
Load a CSV file containing demodulator samples as saved by the ziControl
User Interface into a numpy structured array.
```

Arguments:

filename (str): The file or filename of the CSV file to load.

column_names (list or tuple of str, optional): A list (or tuple) of column names (demodulator sample field names) to load from the CSV file. Default is to load all columns.

Returns:

sample (numpy ndarray): A numpy structured array of shape (num_points,) whose field names correspond to the field names of a ziControl demodulator sample. num_points is the number of lines in the CSV file - 1.

Example:

```
import zhinst.utils
sample = zhinst.utils.load_labone_csv('Freq1.csv', ('t', 'x', 'y'))
import matplotlib.pyplot as plt
import numpy as np
plt.plot(sample['t'], np.abs(sample['x'] + 1j*sample['y']))

load_zicontrol_zibin(filename, column_names=('t', 'x', 'y', 'freq', 'dio',
'auxin0', 'auxin1'))
Load a ziBin file containing demodulator samples as saved by the ziControl
User Interface into a numpy structured array. This is for data saved by
ziControl in binary format.
```

Arguments:

filename (str): The filename of the .ziBin file to load.

column_names (list or tuple of str, optional): A list (or tuple) of column names to load from the CSV file. Default is to load all columns.

Returns:

sample (numpy ndarray): A numpy structured array of shape (num_points,) whose field names correspond to the field names of a ziControl demodulator sample. num_points is the number of sample points saved in the file.

Further comments:

Specifying a fewer names in ``column_names`` will not result in a speed-up as all data is loaded from the binary file by default.

Example:

```
import zhinst.utils
sample = zhinst.utils.load_zicontrol_zibin('Freq1.ziBin')
import matplotlib.pyplot as plt
import numpy as np
plt.plot(sample['t'], np.abs(sample['x'] + 1j*sample['y']))
```

`save_settings(daq, device, filename)`
Save settings from the specified device to a LabOne settings file. This function is synchronous; it will block until saving the settings has finished.

Arguments:

`daq` (instance of `ziDAQServer`): A ziPython API session.

`device` (str): The device ID specifying where to load the settings, e.g., 'dev123'.

`filename` (str): The filename of the LabOne xml settings file. The filename can include a relative or full path.

Raises:

`RuntimeError`: If saving the settings times out.

Examples:

```
import zhinst.utils as utils
daq = utils.autoConnect()
dev = utils.autoDetect(daq)

# Then, e.g., save settings to a file in the current directory:
utils.save_settings(daq, dev, 'my_settings.xml')

# Then, e.g., save settings to the default LabOne settings path:
filename = 'my_settings_example.xml'
path = utils.get_default_settings_path(daq)
utils.save_settings(daq, dev, path + os.sep + filename)
```

`signin_aurange(daq, device, in_channel)`
Perform an automatic adjustment of the signal input range based on the measured input signal. This utility function starts the functionality implemented in the device's firmware and waits until it has completed. The range is set by the firmware based on the measured input signal's amplitude measured over approximately 100 ms.

Requirements:

A devtype that supports autorange functionality on the firmware level, e.g., UHFLI, MFLI, MFIA.

Arguments:

`daq` (instance of `ziDAQServer`): A ziPython API session.

`device` (str): The device ID on which to perform the signal input autorange.

`in_channel` (int): The index of the signal input channel to autorange.

Raises:

`AssertionError`: If the functionality is not supported by the device or an invalid `in_channel` was specified.

`RuntimeError`: If autorange functionality does not complete within the timeout.

Example:

```
import zhinst.utils
device_serial = 'dev2006'
(daq, _, _) = zhinst.utils.create_api_session(device_serial, 5)
input_channel = 0
zhinst.utils.signin_autorange(daq, device_serial, input_channel)

tc2bw(timeconstant, order)
Convert the demodulator timeconstant to its equivalent 3 dB bandwidth for the
specified demodulator order.
```

Inputs:

timeconstant (double): The equivalent demodulator timeconstant.

order (int): The demodulator order (1 to 8) for which to convert the bandwidth.

Output:

bandwidth (double): The demodulator 3dB bandwidth to convert.

DATA

```
LABONE_DEMOD_DTYPE = [('chunk', 'u8'), ('timestamp', 'u8'), ('x', 'f8'...
LABONE_DEMOD_FORMATS = ('u8', 'u8', 'f8', 'f8', 'f8', 'f8', 'u4', 'u4'...
LABONE_DEMOD_NAMES = ('chunk', 'timestamp', 'x', 'y', 'freq', 'phase',...
ZICONTROL_DTYPE = [('t', 'f8'), ('x', 'f8'), ('y', 'f8'), ('freq', 'f8'...
ZICONTROL_FORMATS = ('f8', 'f8', 'f8', 'f8', 'u4', 'f8', 'f8')
ZICONTROL_NAMES = ('t', 'x', 'y', 'freq', 'dio', 'auxin0', 'auxin1')
logger = <logging.Logger object>
print_function = _Feature((2, 6, 0, 'alpha', 2), (3, 0, 0, 'alpha', 0)...
```

FILE

```
/home/ci/.pyenv/versions/3.5.1/lib/python3.5/site-packages/zhinst/utils.py
```

4.4.3. Help for ziPython's ziDAQServer class

```
>>> help('zhinst.ziPython.ziDAQServer')
```

Help on class ziDAQServer in zhinst.ziPython:

```
zhinst.ziPython.ziDAQServer = class ziDAQServer(Boost.Python.instance)
|   Class to connect with a Zurich Instruments data server.
|
|   Method resolution order:
|       ziDAQServer
|       Boost.Python.instance
|       builtins.object
|
|   Methods defined here:
|
|   __init__(...)
|       __init__( (object)arg1) -> None
|
|       __init__( (object)arg1, (str)arg2, (int)arg3) -> None :
|           Connect to the server by using host address and port number.
|           arg1: Reference to the ziDAQServer class.
|           arg2: Host string e.g. '127.0.0.1' for localhost.
|           arg3: Port number e.g. 8004 for the ziDataServer.
|
|       __init__( (object)arg1, (str)arg2, (int)arg3, (int)arg4) -> None :
|           Connect to the server by using host address and port number.
```

```

    |         arg1: Reference to the ziDAQServer class.
    |         arg2: Host string e.g. '127.0.0.1' for localhost.
    |         arg3: Port number e.g. 8004 for the ziDataServer.
    |         arg4: API level number.
    |
    | __reduce__ = <unnamed Boost.Python function>(...)
    |
    | awgModule(...)
    |     awgModule( (ziDAQServer)arg1) -> ziAwgModule :
    |         Create a awgModule class. This will start a thread for running an
    |         asynchronous awgModule.
    |         arg1: Reference to the ziDAQServer class.
    |
    | connect(...)
    |     connect( (ziDAQServer)arg1) -> None
    |
    | connectDevice(...)
    |     connectDevice( (ziDAQServer)arg1, (str)arg2, (str)arg3, (str)arg4) -> None :
    |         Connect with the data server to a specified device over the specified
    |         interface. The device must be visible to the server. If the device is
    |         already connected the call will be ignored. The function will block
    |         until the device is connected and the device is ready to use. This
    |         method is useful for UHF devices offering several communication
    |         interfaces.
    |         arg1: Reference to the ziDAQServer class.
    |         arg2: Device serial.
    |         arg3: Device interface.
    |         arg4: Optional interface parameters string.
    |
    |         connectDevice( (ziDAQServer)arg1, (str)arg2, (str)arg3) -> None
    |
    | deviceSettings(...)
    |     deviceSettings( (ziDAQServer)arg1) -> ziDeviceSettings :
    |         Create a deviceSettings class. This will start a thread for running an
    |         asynchronous deviceSettings.
    |         arg1: Reference to the ziDAQServer class.
    |         arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED,
ignored
    |
    |         deviceSettings( (ziDAQServer)arg1, (int)arg2) -> ziDeviceSettings
    |
    | disconnect(...)
    |     disconnect( (ziDAQServer)arg1) -> None
    |
    | disconnectDevice(...)
    |     disconnectDevice( (ziDAQServer)arg1, (str)arg2) -> None :
    |         Disconnect a device on the data server. This function will return
    |         immediately. The disconnection of the device may not yet finished.
    |         arg1: Reference to the ziDAQServer class.
    |         arg2: Device serial string of device to disconnect.
    |
    | echoDevice(...)
    |     echoDevice( (ziDAQServer)arg1, (str)arg2) -> None :
    |         Sends an echo command to a device and blocks until
    |         answer is received. This is useful to flush all
    |         buffers between API and device to enforce that
    |         further code is only executed after the device executed
    |         a previous command.
    |         arg1: Reference to the ziDAQServer class.
    |         arg2: Device string e.g. 'dev100'.
    |
    | flush(...)
    |     flush( (ziDAQServer)arg1) -> None :
    |         Flush all data in the socket connection and API buffers.
    |         Call this function before a subscribe with subsequent poll
    |         to get rid of old streaming data that might still be in
    |         the buffers.

```

```

        arg1: Reference to the ziDAQServer class.

get(...)
    get( (ziDAQServer)arg1, (str)arg2, (bool)arg3, (int)arg4) -> object :
    Return a dict with all nodes from the specified sub-tree.
    High-speed streaming nodes (e.g. /devN/demods/0/sample)
    are not returned. Wildcards (*) may be used, in which case
    read-only nodes are ignored.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string of the node. Use wild card to
            select all.
        arg3[optional]: Specify which type of data structure to return.
            Return data either as a flat dict (True) or as a nested
            dict tree (False). Default = False.
        arg4[optional]: Specify which type of nodes to include in the
            result. Allowed:
            ZI_LIST_NODES_SETTINGSONLY = 8 (default)
            ZI_LIST_NODES_NONE = 0 (all nodes)

    get( (ziDAQServer)arg1, (str)arg2 [, (bool)arg3]) -> object

getAsEvent(...)
    getAsEvent( (ziDAQServer)arg1, (str)arg2) -> None :
    Trigger an event on the specified node. The node data is returned by a
    subsequent poll command.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string of the node.

getAuxInSample(...)
    getAuxInSample( (ziDAQServer)arg1, (str)arg2) -> object :
    Returns a single auxin sample. The auxin data is averaged in contrast to
    the auxin data embedded in the demodulator sample.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string

getByte(...)
    getByte( (ziDAQServer)arg1, (str)arg2) -> object :
    Get a byte array (string) value from the specified node.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string of the node.

getConnectionAPILevel(...)
    getConnectionAPILevel( (ziDAQServer)arg1) -> int :
    Returns ziAPI level used for the active connection.

getDIO(...)
    getDIO( (ziDAQServer)arg1, (str)arg2) -> object :
    Returns a single DIO sample.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string

getDouble(...)
    getDouble( (ziDAQServer)arg1, (str)arg2) -> float :
    Get a double value from the specified node.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string of the node.

getInt(...)
    getInt( (ziDAQServer)arg1, (str)arg2) -> int :
    Get a integer value from the specified node.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string of the node.

getList(...)
    getList( (ziDAQServer)arg1, (str)arg2) -> object :
    Return a list with all nodes from the specified sub-tree.
        arg1: Reference to the ziDAQServer class.

```



```

        arg2: Path string of the node. Use wild card to
            select all.

getSample(...)
    getSample( (ziDAQServer)arg1, (str)arg2) -> object :
        Returns a single demodulator sample (including DIO and AuxIn). For more
        efficient data recording use subscribe and poll methods.
        arg1: Reference to the ziDAQServer class.
        arg2: Path string

impedanceModule(...)
    impedanceModule( (ziDAQServer)arg1) -> ziImpedanceModule :
        Create a impedanceModule class. This will start a thread for running an
        asynchronous impedanceModule.
        arg1: Reference to the ziDAQServer class.

listNodes(...)
    listNodes( (ziDAQServer)arg1, (str)arg2, (int)arg3) -> list :
        This function returns a list of node names found at the specified path.
        arg1: Reference to the ziDAQRecorder class.
        arg2: Path for which the nodes should be listed. The path may
            contain wildcards so that the returned nodes do not
            necessarily have to have the same parents.
        arg3: Enum that specifies how the selected nodes are listed.
            ziPython.ziListEnum.none -> 0x00
                The default flag, returning a simple
                listing of the given node
            ziPython.ziListEnum.recursive -> 0x01
                Returns the nodes recursively
            ziPython.ziListEnum.absolute -> 0x02
                Returns absolute paths
            ziPython.ziListEnum.leafonly -> 0x04
                Returns only nodes that are leafs,
                which means the they are at the
                outermost level of the tree.
            ziPython.ziListEnum.settingsonly -> 0x08
                Returns only nodes which are marked
                as setting
            Or any combination of flags can be used.

logOff(...)
    logOff( (ziDAQServer)arg1) -> None :
        Disables logging of commands sent to a server.
        arg1: Reference to the ziDAQServer class.

logOn(...)
    logOn( (ziDAQServer)arg1, (int)arg2, (str)arg3, (int)arg4) -> None :
        Enables logging of commands sent to a server.
        arg1: Reference to the ziDAQServer class.
        arg2: Flags (LOG_NONE: 0x00000000
            LOG_SET_DOUBLE: 0x00000001
            LOG_SET_INT: 0x00000002
            LOG_SET_BYTE: 0x00000004
            LOG_SYNC_SET_DOUBLE: 0x00000010
            LOG_SYNC_SET_INT: 0x00000020
            LOG_SYNC_SET_BYTE: 0x00000040
            LOG_GET_DOUBLE: 0x00000100
            LOG_GET_INT: 0x00000200
            LOG_GET_BYTE: 0x00000400
            LOG_GET_DEMOD: 0x00001000
            LOG_GET_DIO: 0x00002000
            LOG_GET_AUXIN: 0x00004000
            LOG_LISTNODES: 0x00010000
            LOG_SUBSCRIBE: 0x00020000
            LOG_UNSUBSCRIBE: 0x00040000
            LOG_GET_AS_EVENT: 0x00080000
            LOG_UPDATE: 0x00100000

```

```

|                                     LOG_POLL_EVENT:      0x00200000
|                                     LOG_POLL:             0x00400000
|                                     LOG_ALL :             0xffffffff)
|
|     arg3: Log file name.
|     arg4: Log style (LOG_STYLE_TELNET: 0 (default),
|            LOG_STYLE_MATLAB: 1, LOG_STYLE_PYTHON: 2).
|
|     logOn( (ziDAQServer)arg1, (int)arg2, (str)arg3) -> None
|
| pidAdvisor(...)
|     pidAdvisor( (ziDAQServer)arg1) -> ziPidAdvisor :
|         Create a pidAdvisor class. This will start a thread for running an
|         asynchronous pidAdvisor.
|         arg1: Reference to the ziDAQServer class.
|         arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED,
ignored
|
|     pidAdvisor( (ziDAQServer)arg1, (int)arg2) -> ziPidAdvisor
|
| pllAdvisor(...)
|     pllAdvisor( (ziDAQServer)arg1) -> ziPllAdvisor :
|         Create a pllAdvisor class. This will start a thread for running an
|         asynchronous pllAdvisor.
|         arg1: Reference to the ziDAQServer class.
|         arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED,
ignored
|
|     pllAdvisor( (ziDAQServer)arg1, (int)arg2) -> ziPllAdvisor
|
| poll(...)
|     poll( (ziDAQServer)arg1, (float)arg2, (int)arg3, (int)arg4, (bool)arg5) ->
object :
|         This function returns subscribed data previously in the API's buffers or
|         obtained during the specified time. It returns a dict tree containing
|         the recorded data. This function blocks until the recording time is
|         elapsed.     arg1: Reference to the ziDAQServer class.
|         arg2: Recording time in [s]. The function will block during that.
|             time.
|         arg3: Poll timeout in [ms]. Recommended value is 500ms.
|         arg4[optional]: Poll flags.
|             FILL = 0x0001 : Fill holes.
|             ALIGN = 0x0002 : Align data that contains a
|                 timestamp.
|             THROW = 0x0004 : Throw EOFError exception if sample
|                 loss is detected.
|         arg5[optional]: Specify which type of data structure to return.
|             Return data either as a flat dict (True) or as a nested
|             dict tree (False). Default = False.
|
|     poll( (ziDAQServer)arg1, (float)arg2, (int)arg3 [, (int)arg4]) -> object
|
| pollEvent(...)
|     pollEvent( (ziDAQServer)arg1, (int)arg2) -> object :
|         Execute a single poll command. Note: only one data packet will be
|         fetched. To get all data waiting in the buffers this command should be
|         executed continuously until nothing is returned anymore. This is a low
|         level command. Use the poll command or asynchronous recording instead.
|         arg1: Reference to the ziDAQServer class.
|         arg2: Poll timeout in [ms]. Recommended value is 500ms.
|
| programRT(...)
|     programRT( (ziDAQServer)arg1, (str)arg2, (str)arg3) -> None :
|         Program RT.
|         arg1: Device identifier e.g. 'dev99'.
|         arg2: File name of the RT program.
|
| record(...)

```

```
| record( (ziDAQServer)arg1) -> ziDAQRecorder :  
|     Create a recording class. This will start a thread for asynchronous  
|     recording.  
|     arg1: Reference to the ziDAQServer class.  
|     arg2: Maximum recording time for single triggers in [s]. -  
DEPRECATED, set 'bufferize' param instead  
|     arg3: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED,  
ignored  
|     arg4[optional]: Record flags. - DEPRECATED, set 'flags' param instead  
|         FILL = 0x0001 : Fill holes.  
|         ALIGN = 0x0002 : Align data that contains a  
|             timestamp.  
|         THROW = 0x0004 : Throw EOFError exception if  
|             sample loss is detected.  
|  
| record( (ziDAQServer)arg1, (float)arg2, (int)arg3 [, (int)arg4]) ->  
ziDAQRecorder  
|  
| revision(...)  
|     revision( (ziDAQServer)arg1) -> int :  
|         Get the revision number of the Python interface of Zurich Instruments.  
|         arg1: Reference to the ziDAQServer class.  
|  
| saveEngine(...)  
|     saveEngine( (ziDAQServer)arg1) -> ziSaveEngine :  
|         Create a saveEngine class. This will start a thread for running an  
|         asynchronous saveEngine.  
|         arg1: Reference to the ziDAQServer class.  
|  
|     saveEngine( (ziDAQServer)arg1, (int)arg2) -> ziSaveEngine  
|  
| set(...)  
|     set( (ziDAQServer)arg1, (object)arg2) -> None :  
|         arg1: Reference to the ziDAQServer class.  
|         arg2: A list of path/value pairs.  
|  
| setByte(...)  
|     setByte( (ziDAQServer)arg1, (str)arg2, (object)arg3) -> None :  
|         arg1: Reference to the ziDAQServer class.  
|         arg2: Path string of the node.  
|  
| setDebugLevel(...)  
|     setDebugLevel( (ziDAQServer)arg1, (int)arg2) -> None :  
|         Enables debug log and sets the debug level.  
|         arg1: Reference to the ziDAQServer class.  
|         arg2: Debug level (trace:0, info:1, debug:2, warning:3, error:4,  
|             fatal:5, status:6).  
|  
| setDouble(...)  
|     setDouble( (ziDAQServer)arg1, (str)arg2, (float)arg3) -> None :  
|         arg1: Reference to the ziDAQServer class.  
|         arg2: Path string of the node.  
|  
| setInt(...)  
|     setInt( (ziDAQServer)arg1, (str)arg2, (int)arg3) -> None :  
|         arg1: Reference to the ziDAQServer class.  
|         arg2: Path string of the node.  
|  
| subscribe(...)  
|     subscribe( (ziDAQServer)arg1, (object)arg2) -> None :  
|         Subscribe to one or several nodes. Fetch data with the poll  
|         command. In order to avoid fetching old data that is still in the  
|         buffer execute a flush command before subscribing to data streams.  
|         arg1: Reference to the ziDAQServer class.  
|         arg2: Path string of the node. Use wild card to  
|             select all. Alternatively also a list of path
```

```

|                                     strings can be specified.
|
| sweep(...)
|     sweep( (ziDAQServer)arg1) -> ziDAQSweeper :
|         Create a sweeper class. This will start a thread for asynchronous
|         sweeping.
|         arg1: Reference to the ziDAQServer class.
|         arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED,
ignored
|
|     sweep( (ziDAQServer)arg1, (int)arg2) -> ziDAQSweeper
|
| sync(...)
|     sync( (ziDAQServer)arg1) -> None :
|         Synchronize all data path. Ensures that get and poll
|         commands return data which was recorded after the
|         setting changes in front of the sync command. This
|         sync command replaces the functionality of all syncSet,
|         flush, and echoDevice commands.
|         arg1: Reference to the ziDAQServer class.
|
| syncSetDouble(...)
|     syncSetDouble( (ziDAQServer)arg1, (str)arg2, (float)arg3) -> float :
|         arg1: Reference to the ziDAQServer class.
|         arg2: Path string of the node.
|
| syncSetInt(...)
|     syncSetInt( (ziDAQServer)arg1, (str)arg2, (int)arg3) -> int :
|         arg1: Reference to the ziDAQServer class.
|         arg2: Path string of the node.
|
| unsubscribe(...)
|     unsubscribe( (ziDAQServer)arg1, (object)arg2) -> None :
|         Unsubscribe data streams. Use this command after recording to avoid
|         buffer overflows that may increase the latency of other command.
|         arg1: Reference to the ziDAQServer class.
|         arg2: Path string of the node. Use wild card to
|             select all. Alternatively also a list of path
|             strings can be specified.
|
| update(...)
|     update( (ziDAQServer)arg1) -> None :
|         Check if additional devices are attached. This function is not needed
|         for servers running under windows as devices will be detected
|         automatically.
|         arg1: Reference to the ziDAQServer class.
|
| vectorWrite(...)
|     vectorWrite( (ziDAQServer)arg1, (str)arg2, (object)arg3) -> None :
|         arg1: Reference to the ziDAQServer class.
|         arg2: Path string of the node.
|         arg3: Vector ((u)int8, (u)int16, (u)int32, (u)int64, float, double)
or string to write.
|
| version(...)
|     version( (ziDAQServer)arg1) -> str :
|         Get version string of the Python interface of Zurich Instruments.
|         arg1: Reference to the ziDAQServer class.
|
| writeDebugLog(...)
|     writeDebugLog( (ziDAQServer)arg1, (int)arg2, (str)arg3) -> None :
|         Outputs message to the debug log (if enabled).
|         arg1: Reference to the ziDAQServer class.
|         arg2: Severity (trace:0, info:1, debug:2, warning:3, error:4,
|             fatal:5, status:6).
|         arg3: Message to output to the log.

```

```
| zoomFFT(...)
|     zoomFFT( (ziDAQServer)arg1) -> ziDAQZoomFFT :
|         Create a zoomFFT class. This will start a thread for running an
|         asynchronous zoomFFT.
|         arg1: Reference to the ziDAQServer class.
|         arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED,
ignored
|
|         zoomFFT( (ziDAQServer)arg1, (int)arg2) -> ziDAQZoomFFT
|
| -----
| Data and other attributes defined here:
|
| __instance_size__ = 56
|
| -----
| Methods inherited from Boost.Python.instance:
|
| __new__(*args, **kwargs) from Boost.Python.class
|     Create and return a new object. See help(type) for accurate signature.
|
| -----
| Data descriptors inherited from Boost.Python.instance:
|
| __dict__
|
| __weakref__
```

4.4.4. Help for ziPython's ziDeviceSettings class

An instance of ziDeviceSettings is initialized using the deviceSettings method from ziDAQServer:

```
>>> help('zhinst.ziPython.ziDAQServer.deviceSettings')
```

Help on built-in function deviceSettings in zhinst.ziPython.ziDAQServer:

```
zhinst.ziPython.ziDAQServer.deviceSettings = deviceSettings(...)
deviceSettings( (ziDAQServer)arg1) -> ziDeviceSettings :
    Create a deviceSettings class. This will start a thread for running an
    asynchronous deviceSettings.
    arg1: Reference to the ziDAQServer class.
    arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED, ignored

deviceSettings( (ziDAQServer)arg1, (int)arg2) -> ziDeviceSettings
```

Reference help for the ziDeviceSettings class.

```
>>> help('zhinst.ziPython.ziDeviceSettings')
```

Help on class ziDeviceSettings in zhinst.ziPython:

```
zhinst.ziPython.ziDeviceSettings = class ziDeviceSettings(Boost.Python.instance)
| Method resolution order:
|     ziDeviceSettings
|     Boost.Python.instance
|     builtins.object
|
| Methods defined here:
|
| __init__(...)
|     Raises an exception
|     This class cannot be instantiated from Python
```

```

| __reduce__ = <unnamed Boost.Python function>(...)
|
| clear(...)
|     clear( (ziDeviceSettings)arg1) -> None :
|         End the deviceSettings thread.
|
| execute(...)
|     execute( (ziDeviceSettings)arg1) -> None :
|         Execute the save/load command.
|
| finish(...)
|     finish( (ziDeviceSettings)arg1) -> None :
|         Stop the load/save command. The command may be restarted by calling
|         'execute' again.
|
| finished(...)
|     finished( (ziDeviceSettings)arg1) -> bool :
|         Check if the command execution has finished. Returns True if finished.
|
| get(...)
|     get( (ziDeviceSettings)arg1, (str)arg2, (bool)arg3) -> object :
|         Return a dict with all nodes from the specified sub-tree.
|         arg1: Reference to the ziDeviceSettings class.
|         arg2: Path string of the node. Use wild card to
|             select all.
|         arg3[optional]: Specify which type of data structure to return.
|             Return data either as a flat dict (True) or as a nested
|             dict tree (False). Default = False.
|
|     get( (ziDeviceSettings)arg1, (str)arg2) -> object
|
| listNodes(...)
|     listNodes( (ziDeviceSettings)arg1, (str)arg2, (int)arg3) -> list :
|         This function returns a list of node names found at the specified path.
|         arg1: Reference to the ziDeviceSettings class.
|         arg2: Path for which the nodes should be listed. The path may
|             contain wildcards so that the returned nodes do not
|             necessarily have to have the same parents.
|         arg3: Enum that specifies how the selected nodes are listed.
|             ziPython.ziListEnum.none -> 0x00
|                 The default flag, returning a simple
|                 listing if the given node
|             ziPython.ziListEnum.recursive -> 0x01
|                 Returns the nodes recursively
|             ziPython.ziListEnum.absolute -> 0x02
|                 Returns absolute paths
|             ziPython.ziListEnum.leafsonly -> 0x04
|                 Returns only nodes that are leafs,
|                 which means the they are at the
|                 outermost level of the tree.
|             ziPython.ziListEnum.settingsonly -> 0x08
|                 Returns only nodes which are marked
|                 as setting
|         Or combinations of flags can be used.
|
| progress(...)
|     progress( (ziDeviceSettings)arg1) -> object :
|         Reports the progress of the command with a number between
|         0 and 1.
|
| read(...)
|     read( (ziDeviceSettings)arg1, (bool)arg2) -> object :
|         Read device settings. Only relevant for the save command.
|         arg1[optional]: Specify which type of data structure to return.
|             Return data either as a flat dict (True) or as a nested
|             dict tree (False). Default = False.

```

```

    read( (ziDeviceSettings)arg1) -> object

save(...)
    save( (ziDeviceSettings)arg1, (str)arg2) -> None :
        Not relevant for the deviceSettings module.

set(...)
    set( (ziDeviceSettings)arg1, (str)arg2, (float)arg3) -> None :
        Device Settings Parameters
        Path name      Type      Description
        deviceSettings/device    string  Device that should be used for
                                         loading/saving device settings,
                                         e.g. 'dev99'.
        deviceSettings/path      string  Directory where settings files should be
                                         located. If not set, the default settings
                                         location of the LabOne software is used.
        deviceSettings/filename string  Name of settings file to use
        deviceSettings/command   string  The command to execute
                                         'save' = Read device settings and save to
                                         file.
                                         'load' = Load settings from file and
                                         write to device.
                                         'read' = Read device settings only
                                         (no save).

    set( (ziDeviceSettings)arg1, (str)arg2, (int)arg3) -> None

    set( (ziDeviceSettings)arg1, (str)arg2, (str)arg3) -> None

    set( (ziDeviceSettings)arg1, (object)arg2) -> None :
        arg1: Reference to the ziDeviceSettings class.
        arg2: A list of path/value pairs.

subscribe(...)
    subscribe( (ziDeviceSettings)arg1, (str)arg2) -> None :
        Not relevant for the deviceSettings module.

trigger(...)
    trigger( (ziDeviceSettings)arg1) -> None :
        Not applicable to this module.

unsubscribe(...)
    unsubscribe( (ziDeviceSettings)arg1, (str)arg2) -> None :
        Not relevant for the deviceSettings module.

-----
Methods inherited from Boost.Python.instance:

__new__(*args, **kwargs) from Boost.Python.class
    Create and return a new object.  See help(type) for accurate signature.

-----
Data descriptors inherited from Boost.Python.instance:

__dict__
__weakref__

```

4.4.5. Help for ziPython's ziDAQSweeper class

An instance of ziDAQSweeper is initialized using the sweep method from ziDAQServer:

```
>>> help('zhinst.ziPython.ziDAQServer.sweep')
```

Help on built-in function sweep in zhinst.ziPython.ziDAQServer:

```
zhinst.ziPython.ziDAQServer.sweep = sweep(...)
sweep( (ziDAQServer)arg1) -> ziDAQSweeper :
    Create a sweeper class. This will start a thread for asynchronous
    sweeping.
    arg1: Reference to the ziDAQServer class.
    arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED, ignored

sweep( (ziDAQServer)arg1, (int)arg2) -> ziDAQSweeper
```

Reference help for the ziDAQSweeper class.

```
>>> help('zhinst.ziPython.ziDAQSweeper')
```

Help on class ziDAQSweeper in zhinst.ziPython:

```
zhinst.ziPython.ziDAQSweeper = class ziDAQSweeper(Boost.Python.instance)
|   Method resolution order:
|       ziDAQSweeper
|       Boost.Python.instance
|       builtins.object
|
|   Methods defined here:
|
|       __init__(...)
|           Raises an exception
|           This class cannot be instantiated from Python
|
|       __reduce__ = <unnamed Boost.Python function>(...)
|
|       clear(...)
|           clear( (ziDAQSweeper)arg1) -> None :
|               End the sweeper thread.
|
|       execute(...)
|           execute( (ziDAQSweeper)arg1) -> None :
|               Start the sweeper. Subscription or unsubscription is no more
|               possible until the sweep is finished.
|
|       finish(...)
|           finish( (ziDAQSweeper)arg1) -> None :
|               Stop sweeping. The sweeping may be restarted by calling
|               'execute' again.
|
|       finished(...)
|           finished( (ziDAQSweeper)arg1) -> bool :
|               Check if the sweep has finished. Returns True if finished.
|
|       get(...)
|           get( (ziDAQSweeper)arg1, (str)arg2, (bool)arg3) -> object :
|               Return a dict with all nodes from the specified sub-tree.
|               arg1: Reference to the ziDAQSweeper class.
|               arg2: Path string of the node. Use wild card to
|                   select all.
|               arg3[optional]: Specify which type of data structure to return.
|                   Return data either as a flat dict (True) or as a nested
|                   dict tree (False). Default = False.
|
|           get( (ziDAQSweeper)arg1, (str)arg2) -> object
|
|       listNodes(...)
|           listNodes( (ziDAQSweeper)arg1, (str)arg2, (int)arg3) -> list :
|               This function returns a list of node names found at the specified path.
|               arg1: Reference to the ziDAQRecorder class.
|               arg2: Path for which the nodes should be listed. The path may
|                   contain wildcards so that the returned nodes do not
```



```

necessarily have to have the same parents.
arg3: Enum that specifies how the selected nodes are listed.
ziPython.ziListEnum.none -> 0x00
    The default flag, returning a simple
    listing if the given node
ziPython.ziListEnum.recursive -> 0x01
    Returns the nodes recursively
ziPython.ziListEnum.absolute -> 0x02
    Returns absolute paths
ziPython.ziListEnum.leafonly -> 0x04
    Returns only nodes that are leafs,
    which means the they are at the
    outermost level of the tree.
ziPython.ziListEnum.settingsonly -> 0x08
    Returns only nodes which are marked
    as setting
Or combinations of flags can be used.

progress(...)
    progress( (ziDAQSweeper)arg1) -> object :
        Reports the progress of the measurement with a number between
        0 and 1.

read(...)
    read( (ziDAQSweeper)arg1, (bool)arg2) -> object :
        Read sweep data. If the sweeping is still ongoing only a subset
        of sweep data is returned. If huge data sets
        are recorded call this method to keep memory usage reasonable.
        arg1[optional]: Specify which type of data structure to return.
        Return data either as a flat dict (True) or as a nested
        dict tree (False). Default = False.

    read( (ziDAQSweeper)arg1) -> object

save(...)
    save( (ziDAQSweeper)arg1, (str)arg2) -> None :
        Save sweeper data to file.
        arg1: Reference to the ziDAQSweeper class.
        arg2: File name string (without extension).

set(...)
    set( (ziDAQSweeper)arg1, (str)arg2, (float)arg3) -> None :
        Sweep Parameters
        Path name      Type      Description
        sweep/device    string    Device that should be used for
                                the parameter sweep, e.g. 'dev99'.
        sweep/start     double   Sweep start frequency [Hz]
        sweep/stop      double   Sweep stop frequency [Hz]
        sweep/gridnode   string    Path of the node that should be
                                used for sweeping. For frequency
                                sweep applications this will be e.g.
                                'oscs/0/freq'. The device name of
                                the path can be omitted and is given
                                by sweep/device.
        sweep/loopcount  int      Number of sweep loops (default 1)
        sweep/endless    int      Sweep endless (default 0)
                                0 = endless off, use loopcount,
                                1 = endless on, ignore loopcount.
        sweep/samplecount  int      Number of samples per sweep.
        sweep/settling/time double    Settling time before measurement is
                                performed.
        sweep/settling/tc  double    Shows the approximate settling precision
                                in time constant units as specified by
                                setting/inaccuracy (calculated upon
                                execute()). Setting this parameter
                                directly is now deprecated and may not
                                be supported in future versions.

```

		5 ~ low precision 15 ~ medium precision 50 ~ high precision
sweep/settling/inaccuracy	int	Demodulator filter settling inaccuracy that defines the wait time between a sweep parameter change and recording of the next sweep point. The settling time is calculated as the time required to attain the specified remaining proportion [1e-13, 0.1] of an incoming step function. Typical inaccuracy values: 10m ~ for highest sweep speed for large signals 100u ~ for precise amplitude measurements 100n ~ for precise noise measurements. Depending on the order, the settling inaccuracy will define the number of filter time constants the sweeper has to wait. The maximum between this value and the settling time is taken as wait time until the next sweep point is recorded.
sweep/xmapping	int	Sweep mode: 0 = linear, 1 = logarithmic.
sweep/scan	int	Scan type: 0 = sequential, 1 = binary, 2 = bidirectional, 3 = reverse.
sweep/bandwidth	double	Fixed bandwidth [Hz], 0 = Automatic calculation.
sweep/bandwidthoverlap	int	Sets the bandwidth overlap mode, (default 0): 0 = disabled 1 = enabled If enabled the bandwidth of a sweep point may overlap with the frequency of neighboring sweep points. The effective bandwidth is only limited by the maximal bandwidth setting and omega suppression. As a result, the bandwidth is independent of the number of sweep points. For frequency response analysis bandwidth overlap should be enabled to achieve maximal sweep speed.
sweep/bandwidthcontrol	int	Sets the bandwidth control mode, (default 2): 0 = Manual (user sets bandwidth and order), 1 = Fixed (uses fixed bandwidth value), 2 = Auto (calculates best bandwidth value) Equivalent to the obsolete bandwidth = 0 setting.
sweep/order	int	Defines the filter roll off to use in Fixed bandwidth selection. Valid values are between 0 (6 dB/octave) and 8 (48 dB/octave). An order of 0 triggers a read-out of the order from the selected demodulator.
sweep/maxbandwidth	double	Maximal bandwidth used in auto bandwidth mode in [Hz]. The default is 1.25MHz.
sweep/omegasuppression	double	Damping in [dB] of omega and 2-omega components. Default is 40dB in favor of sweep speed.

		Use higher value for strong offset values or 3-omega measurement methods.
sweep/averaging/tc	double	Min averaging time [tc] 0 = no averaging (see also time!) 5 ~ low precision 15 ~ medium precision 50 ~ high precision
sweep/averaging/sample	int	Min samples to average 1 = no averaging (if averaging/tc = 0)
sweep/phaseunwrap	bool	Enable unwrapping of slowly changing phase evolutions around the +/-180 degree boundary.
sweep/sincfilter	bool	Enables the sinc filter if the sweep frequency is below 50 Hz. This will improve the sweep speed at low frequencies as omega components do not need to be suppressed by the normal low pass filter.
sweep/filename	string	This parameter is deprecated. If specified, i.e. not empty, it enables automatic saving of data in single sweep mode (sweep/endless = 0).
sweep/directory	string	The directory where files are located when saving sweeper measurements.
sweep/fileformat	string	The format of the file for saving sweeper measurements: 0 = Matlab, 1 = CSV, 2 = ZView (Impedance data only).
sweep/historylength	bool	Maximum number of entries stored in the measurement history.
sweep/clearhistory	bool	Remove all records from the history list.
set((ziDAQSweeper)arg1, (str)arg2, (int)arg3) -> None		
set((ziDAQSweeper)arg1, (str)arg2, (str)arg3) -> None		
set((ziDAQSweeper)arg1, (object)arg2) -> None :		
arg1: Reference to the ziDAQSweeper class.		
arg2: A list of path/value pairs.		
subscribe(...)		
subscribe((ziDAQSweeper)arg1, (str)arg2) -> None :		
Subscribe to one or several nodes. After subscription the sweep process can be started with the 'execute' command. During the sweep process paths can not be subscribed or unsubscribed.		
arg1: Reference to the ziDAQSweeper class.		
arg2: Path string of the node. Use wild card to select all. Alternatively also a list of path strings can be specified.		
trigger(...)		
trigger((ziDAQSweeper)arg1) -> None :		
Execute a manual trigger.		
unsubscribe(...)		
unsubscribe((ziDAQSweeper)arg1, (str)arg2) -> None :		
Unsubscribe from one or several nodes. During the sweep process paths can not be subscribed or unsubscribed.		
arg1: Reference to the ziDAQSweeper class.		
arg2: Path string of the node. Use wild card to select all. Alternatively also a list of path strings can be specified.		

Methods inherited from Boost.Python.instance:		

```
| __new__(*args, **kwargs) from Boost.Python.class
|     Create and return a new object.  See help(type) for accurate signature.
|
| -----
| Data descriptors inherited from Boost.Python.instance:
|
| __dict__
|
| __weakref__
```

4.4.6. Help for ziPython's ziDAQZoomFFT class

An instance of ziDAQZoomFFT is initialized using the zoomFFT method from ziDAQServer:

```
>>> help('zhinst.ziPython.ziDAQServer.zoomFFT')
```

Help on built-in function zoomFFT in zhinst.ziPython.ziDAQServer:

```
zhinst.ziPython.ziDAQServer.zoomFFT = zoomFFT(...)
zoomFFT( (ziDAQServer)arg1) -> ziDAQZoomFFT :
    Create a zoomFFT class. This will start a thread for running an
    asynchronous zoomFFT.
        arg1: Reference to the ziDAQServer class.
        arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED, ignored

zoomFFT( (ziDAQServer)arg1, (int)arg2) -> ziDAQZoomFFT
```

Reference help for the ziDAQZoomFFT class.

```
>>> help('zhinst.ziPython.ziDAQZoomFFT')
```

Help on class ziDAQZoomFFT in zhinst.ziPython:

```
zhinst.ziPython.ziDAQZoomFFT = class ziDAQZoomFFT(Boost.Python.instance)
|   Method resolution order:
|       ziDAQZoomFFT
|       Boost.Python.instance
|       builtins.object
|
|   Methods defined here:
|
|   __init__(...)
|       Raises an exception
|       This class cannot be instantiated from Python
|
|   __reduce__ = <unnamed Boost.Python function>(...)
|
|   clear(...)
|       clear( (ziDAQZoomFFT)arg1) -> None :
|           End the zoom FFT thread.
|
|   execute(...)
|       execute( (ziDAQZoomFFT)arg1) -> None :
|           Start the zoom FFT. Subscription or unsubscription is no more
|           possible until the zoom FFT is finished.
|
|   finish(...)
|       finish( (ziDAQZoomFFT)arg1) -> None :
|           Stop the zoom FFT. The zoom FFT may be restarted by calling
|           'execute' again.
|
|   finished(...)
|       finished( (ziDAQZoomFFT)arg1) -> bool :
|           Check if the zoom FFT has finished. Returns True if finished.
```

```

get(...)
    get( (ziDAQZoomFFT)arg1, (str)arg2, (bool)arg3) -> object :
        Return a dict with all nodes from the specified sub-tree.
        arg1: Reference to the ziDAQZoomFFT class.
        arg2: Path string of the node. Use wild card to
            select all.
        arg3[optional]: Specify which type of data structure to return.
            Return data either as a flat dict (True) or as a nested
            dict tree (False). Default = False.

    get( (ziDAQZoomFFT)arg1, (str)arg2) -> object

listNodes(...)
    listNodes( (ziDAQZoomFFT)arg1, (str)arg2, (int)arg3) -> list :
        This function returns a list of node names found at the specified path.
        arg1: Reference to the ziDAQZoomFFT class.
        arg2: Path for which the nodes should be listed. The path may
            contain wildcards so that the returned nodes do not
            necessarily have to have the same parents.
        arg3: Enum that specifies how the selected nodes are listed.
            ziPython.ziListEnum.none -> 0x00
                The default flag, returning a simple
                listing if the given node
            ziPython.ziListEnum.recursive -> 0x01
                Returns the nodes recursively
            ziPython.ziListEnum.absolute -> 0x02
                Returns absolute paths
            ziPython.ziListEnum.leafsonly -> 0x04
                Returns only nodes that are leafs,
                which means the they are at the
                outermost level of the tree.
            ziPython.ziListEnum.settingsonly -> 0x08
                Returns only nodes which are marked
                as setting
        Or combinations of flags can be used.

progress(...)
    progress( (ziDAQZoomFFT)arg1) -> object :
        Reports the progress of the measurement with a number between
        0 and 1.

read(...)
    read( (ziDAQZoomFFT)arg1, (bool)arg2) -> object :
        Read zoom FFT data. If the zoom FFT is still ongoing only a subset
        of zoom FFT data is returned.
        arg1[optional]: Specify which type of data structure to return.
            Return data either as a flat dict (True) or as a nested
            dict tree (False). Default = False.

    read( (ziDAQZoomFFT)arg1) -> object

save(...)
    save( (ziDAQZoomFFT)arg1, (str)arg2) -> None :
        Save zoom FFT data to file.
        arg1: Reference to the ziDAQZoomFFT class.
        arg2: File name string (without extension).

set(...)
    set( (ziDAQZoomFFT)arg1, (str)arg2, (float)arg3) -> None :
        Zoom FFT Parameters
        Path name      Type      Description
        zoomFFT/device  string    Device that should be used for
            the zoom FFT, e.g. 'dev99'.
        zoomFFT/bit     int       Number of FFT points 2^bit
        zoomFFT/mode     int       Zoom FFT mode
            0 = Perform FFT on X+iY
            1 = Perform FFT on R

```

		2 = Perform FFT on Phase
zoomFFT/loopcount	int	Number of zoomFFT loops (default 1)
zoomFFT/endless	int	Perform endless zoomFFT (default 0)
		0 = endless off, use loopcount
		1 = endless on, ignore loopcount
zoomFFT/overlap	double	FFT overlap 0 = none, [0..1]
zoomFFT/settling/time	double	Settling time before measurement is performed.
zoomFFT/settling/tc	double	Settling time in time constant units before the FFT recording is started.
		5 ~ low precision
		15 ~ medium precision
		50 ~ high precision
zoomFFT/window	int	FFT window (default 1 = Hann)
		0 = Rectangular
		1 = Hann
		2 = Hamming
		3 = Blackman Harris 4 term
zoomFFT/absolute	bool	Shifts the frequencies so that the center frequency becomes the demodulation frequency rather than 0 Hz.

set((ziDAQZoomFFT)arg1, (str)arg2, (int)arg3) -> None

set((ziDAQZoomFFT)arg1, (str)arg2, (str)arg3) -> None

set((ziDAQZoomFFT)arg1, (object)arg2) -> None :

 arg1: Reference to the ziDAQZoomFFT class.

 arg2: A list of path/value pairs.

subscribe(...)

 subscribe((ziDAQZoomFFT)arg1, (str)arg2) -> None :

 Subscribe to one or several nodes. After subscription the zoom FFT process can be started with the 'execute' command. During the zoom FFT process paths can not be subscribed or unsubscribed.

 arg1: Reference to the ziDAQZoomFFT class.

 arg2: Path string of the node. Use wild card to select all. Alternatively also a list of path strings can be specified.

trigger(...)

 trigger((ziDAQZoomFFT)arg1) -> None :

 Execute a manual trigger.

unsubscribe(...)

 unsubscribe((ziDAQZoomFFT)arg1, (str)arg2) -> None :

 Unsubscribe from one or several nodes. During the zoom FFT process paths can not be subscribed or unsubscribed.

 arg1: Reference to the ziDAQZoomFFT class.

 arg2: Path string of the node. Use wild card to select all. Alternatively also a list of path strings can be specified.

Methods inherited from Boost.Python.instance:

__new__(*args, **kwargs) from Boost.Python.class

 Create and return a new object. See help(type) for accurate signature.

Data descriptors inherited from Boost.Python.instance:

__dict__

__weakref__

4.4.7. Help for ziPython's ziDAQRecorder class

An instance of ziDAQRecorder is initialized using the record method from ziDAQServer:

```
>>> help('zhinst.ziPython.ziDAQServer.record')
```

Help on built-in function record in zhinst.ziPython.ziDAQServer:

```
zhinst.ziPython.ziDAQServer.record = record(...)
  record( (ziDAQServer)arg1) -> ziDAQRecorder :
    Create a recording class. This will start a thread for asynchronous
    recording.
    arg1: Reference to the ziDAQServer class.
    arg2: Maximum recording time for single triggers in [s]. - DEPRECATED,
    set 'bufferize' param instead
    arg3: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED, ignored
    arg4[optional]: Record flags. - DEPRECATED, set 'flags' param instead
        FILL = 0x0001 : Fill holes.
        ALIGN = 0x0002 : Align data that contains a
                        timestamp.
        THROW = 0x0004 : Throw EOFError exception if
                        sample loss is detected.

    record( (ziDAQServer)arg1, (float)arg2, (int)arg3 [, (int)arg4]) -> ziDAQRecorder
```

Reference help for the ziDAQRecorder class.

```
>>> help('zhinst.ziPython.ziDAQRecorder')
```

Help on class ziDAQRecorder in zhinst.ziPython:

```
zhinst.ziPython.ziDAQRecorder = class ziDAQRecorder(Boost.Python.instance)
|   Method resolution order:
|       ziDAQRecorder
|       Boost.Python.instance
|       builtins.object
|
|   Methods defined here:
|
|   __init__(...)
|       Raises an exception
|       This class cannot be instantiated from Python
|
|   __reduce__ = <unnamed Boost.Python function>(...)
|
|   clear(...)
|       clear( (ziDAQRecorder)arg1) -> None :
|           End the recording thread.
|
|   execute(...)
|       execute( (ziDAQRecorder)arg1) -> None :
|           Start the recorder. After that command any trigger will start
|           the measurement. Subscription or unsubscription is no more
|           possible until the recording is finished.
|
|   finish(...)
|       finish( (ziDAQRecorder)arg1) -> None :
|           Stop recording. The recording may be restarted by calling
|           'execute' again.
|
|   finished(...)
|       finished( (ziDAQRecorder)arg1) -> bool :
|           Check if the recording has finished. Returns True if finished.
```

```

| get(...)
|     get( (ziDAQRecorder)arg1, (str)arg2, (bool)arg3) -> object :
|         Return a dict with all nodes from the specified sub-tree.
|         arg1: Reference to the ziDAQRecorder class.
|         arg2: Path string of the node. Use wild card to
|             select all.
|         arg3[optional]: Specify which type of data structure to return.
|             Return data either as a flat dict (True) or as a nested
|             dict tree (False). Default = False.
|
|     get( (ziDAQRecorder)arg1, (str)arg2) -> object
|
| listNodes(...)
|     listNodes( (ziDAQRecorder)arg1, (str)arg2, (int)arg3) -> list :
|         This function returns a list of node names found at the specified path.
|         arg1: Reference to the ziDAQRecorder class.
|         arg2: Path for which the nodes should be listed. The path may
|             contain wildcards so that the returned nodes do not
|             necessarily have to have the same parents.
|         arg3: Enum that specifies how the selected nodes are listed.
|             ziPython.ziListEnum.none -> 0x00
|                 The default flag, returning a simple
|                 listing if the given node
|             ziPython.ziListEnum.recursive -> 0x01
|                 Returns the nodes recursively
|             ziPython.ziListEnum.absolute -> 0x02
|                 Returns absolute paths
|             ziPython.ziListEnum.leafsonly -> 0x04
|                 Returns only nodes that are leafs,
|                 which means the they are at the
|                 outermost level of the tree.
|             ziPython.ziListEnum.settingsonly -> 0x08
|                 Returns only nodes which are marked
|                 as setting
|         Or combinations of flags can be used.
|
| progress(...)
|     progress( (ziDAQRecorder)arg1) -> object :
|         Reports the progress of the measurement with a number between
|         0 and 1.
|
| read(...)
|     read( (ziDAQRecorder)arg1, (bool)arg2) -> object :
|         Read recorded data. If the recording is still ongoing only a subset
|         of recorded data is returned. If many triggers or huge data sets
|         are recorded call this method to keep memory usage reasonable.
|         arg1[optional]: Specify which type of data structure to return.
|             Return data either as a flat dict (True) or as a nested
|             dict tree (False). Default = False.
|
|     read( (ziDAQRecorder)arg1) -> object
|
| save(...)
|     save( (ziDAQRecorder)arg1, (str)arg2) -> None :
|         Save trigger data to file.
|         arg1: Reference to the ziDAQRecorder class.
|         arg2: File name string (without extension).
|
| set(...)
|     set( (ziDAQRecorder)arg1, (str)arg2, (float)arg3) -> None :
|         Trigger Parameters
|         Path name      Type      Description
|         trigger/buffersize  double  Overwrite the buffersize [s] of the
|                                     trigger object (set when it was
|                                     instantiated). Recommended buffer size
|                                     is 2*trigger/0/duration.
|         trigger/flags      int      Record flags.

```



```

|                                     FILL = 0x0001 : Fill holes.
|                                     ALIGN = 0x0002 : Align data that contains
a                                     timestamp.
|                                     THROW = 0x0004 : Throw EOFError exception
|
|                                     sample loss is detected.
|                                     trigger/device string The device ID to execute the software
|                                     trigger, e.g. dev123
|                                     DEPRECATED - use absolute path in
trigger/0/triggernode
|                                     trigger/endless bool Enable endless triggering:
|                                     1 = enable,
|                                     0 = disable.
|                                     trigger/forcetrigger bool Force a trigger.
|                                     trigger/0/triggernode string Path and signal of the node that should be
|                                     used for triggering, separated by a dot
(.),
|                                     e.g. /devN/demods/0/sample.x
|                                     Overrides values from trigger/0/path and
|                                     trigger/0/source.
|                                     trigger/0/path string The path to the demod sample to trigger
|                                     on, e.g. demods/3/sample, see also
|                                     trigger/0/source.
|                                     DEPRECATED - use trigger/0/triggernode
instead.
|                                     trigger/0/source int Signal that is used to trigger on.
|                                     0 = x [X_SOURCE]
|                                     1 = y [Y_SOURCE]
|                                     2 = r [R_SOURCE]
|                                     3 = angle [ANGLE_SOURCE]
|                                     4 = frequency [FREQUENCY_SOURCE]
|                                     5 = phase [PHASE_SOURCE]
|                                     [AUXINO_SOURCE / PARAM0_SOURCE]
|                                     6 = auxiliary input 0 / parameter 0
|                                     [AUXIN1_SOURCE / PARAM1_SOURCE]
|                                     7 = auxiliary input 1 / parameter 1
|                                     DEPRECATED - use trigger/0/triggernode
instead.
|                                     trigger/0/count int Number of trigger edges to record.
|                                     trigger/0/type int Trigger type used. Some parameters are
|                                     only valid for special trigger nodes
|                                     and/or types.
|                                     0 = trigger off
|                                     1 = analog edge trigger on source
|                                     2 = digital trigger mode on DIO
|                                     3 = analog pulse trigger on source
|                                     4 = analog tracking trigger on source
|                                     5 = hardware trigger on trigger line
source
|                                     6 = tracking edge trigger on source
|                                     7 = event count trigger on counter source
|                                     trigger/0/edge int Trigger edge
|                                     1 = rising edge
|                                     2 = falling edge
|                                     3 = both
|                                     trigger/0/findlevel bool Automatically find the value of
|                                     trigger/0/level based on
|                                     the current signal value.
|                                     trigger/0/bits int Digital trigger condition.
|                                     trigger/0/bitmask int Bit masking for bits used for
|                                     triggering. Used for digital trigger.
|                                     trigger/0/delay double Trigger frame position [s] (left side)
|                                     relative to trigger edge.
|                                     delay = 0 -> trigger edge at left
|                                     border.
|                                     delay < 0 -> trigger edge inside

```

```

|                                     trigger frame (pretrigger).
|                                     delay > 0 -> trigger edge before
|                                     trigger frame (posttrigger).
|
| trigger/0/duration      double Recording frame length [s].
| trigger/0/level        double Trigger level voltage [V].
| trigger/0/hysteresis   double Trigger hysteresis [V].
| trigger/0/retrigger    int    Record more than one trigger in a
|                             trigger frame.
|
| trigger/triggered      bool   Has the software trigger triggered?
|                             1=Yes, 0=No (read only).
|
| trigger/0/bandwidth     double Filter bandwidth [Hz] for pulse and
|                             tracking triggers.
|
| trigger/0/holdoff/count int    Number of skipped triggers until the
|                             next trigger is recorded again.
|
| trigger/0/holdoff/time  double Hold off time [s] before the next
|                             trigger is recorded again. A hold off
|                             time smaller than the duration will
|                             produce overlapped trigger frames.
|
| trigger/0/hwtrigsources int    Only available for devices that support
|                             hardware triggering. Specify the channel
|                             to trigger on.
|                             DEPRECATED - use trigger/0/triggernode
|
| instead.
|
| trigger/0/pulse/min     double Minimal pulse width [s] for the pulse
|                             trigger.
|
| trigger/0/pulse/max     double Maximal pulse width [s] for the pulse
|                             trigger.
|
| trigger/0/grid/mode     int    Enable grid mode. In grid mode a matrix
|                             instead of a vector is returned. Each
|                             trigger becomes a row in the matrix and
|
| each
|
|                             trigger's data is interpolated onto a new
|                             grid defined by the number of columns:
|                             0: Disable
|                             1: Enable with nearest neighbour
|
| interpolation
|
|                             2: Enable with linear interpolation.
|
| trigger/0/grid/operation int    If running in endless mode, either replace
|
| or
|
| trigger/0/grid/cols     int    average the data in the grid's matrix.
|                             Specify the number of columns in the
|
| grid's
|                             matrix. The data from each row is
|
| interpolated
|
|                             onto a grid with the specified number of
|                             columns.
|
| trigger/0/grid/rows     int    Specify the number of rows in the grid's
|                             matrix. Each row is the data recorded from
|
| one
|
|                             trigger interpolated onto the columns.
|
| trigger/0/grid/direction int    The direction to organize data in the
|
| grid's
|                             matrix:
|                             0: Forward.
|                             The data in each row is ordered chrono-
|                             logically, e.g., the first data point
|
| in
|
|                             each row corresponds to the first
|                             timestamp in the trigger data.
|
|                             1: Reverse.
|                             The data in each row is ordered reverse
|                             chronologically, e.g., the first data
|                             point in each row corresponds to the
|
| last
|
|                             timestamp in the trigger data.

```

```

|                                     2: Bidirectional.
|                                     The ordering of the data alternates
between
|                                     Forward and Backward ordering from row-
to-
|                                     row. The first row is Forward ordered.
|
| trigger/filename                    string This parameter is deprecated. If
|                                     specified, i.e. not empty, it enables
|                                     automatic saving of data in single
|                                     trigger mode (trigger/endless = 0).
|
| trigger/directory                  string The directory where files are saved.
| trigger/fileformat                  string The format of the file for saving data:
|                                     0 = Matlab,
|                                     1 = CSV,
|                                     2 = ZView (Impedance data only).
|
| trigger/historylength              bool   Maximum number of entries stored in the
|                                     measurement history.
|
| trigger/clearhistory               bool   Remove all records from the history list.
|
| set( (ziDAQRecorder)arg1, (str)arg2, (int)arg3) -> None
|
| set( (ziDAQRecorder)arg1, (str)arg2, (str)arg3) -> None
|
| set( (ziDAQRecorder)arg1, (object)arg2) -> None :
|     arg1: Reference to the ziDAQRecorder class.
|     arg2: A list of path/value pairs.
|
| subscribe(...)
|     subscribe( (ziDAQRecorder)arg1, (str)arg2) -> None :
|         Subscribe to one or several nodes. After subscription the recording
|         process can be started with the 'execute' command. During the
|         recording process paths can not be subscribed or unsubscribed.
|         arg1: Reference to the ziDAQRecorder class.
|         arg2: Path string of the node. Use wild card to
|             select all. Alternatively also a list of path
|             strings can be specified.
|
| trigger(...)
|     trigger( (ziDAQRecorder)arg1) -> None :
|         Execute a manual trigger.
|
| unsubscribe(...)
|     unsubscribe( (ziDAQRecorder)arg1, (str)arg2) -> None :
|         Unsubscribe from one or several nodes. During the
|         recording process paths can not be subscribed or unsubscribed.
|         arg1: Reference to the ziDAQRecorder class.
|         arg2: Path string of the node. Use wild card to
|             select all. Alternatively also a list of path
|             strings can be specified.
|
| -----
| Methods inherited from Boost.Python.instance:
|
| __new__(*args, **kwargs) from Boost.Python.class
|     Create and return a new object. See help(type) for accurate signature.
|
| -----
| Data descriptors inherited from Boost.Python.instance:
|
| __dict__
|
| __weakref__

```

4.4.8. Help for ziPython's ziP11Advisor class

An instance of `ziPllAdvisor` is initialized using the `pllAdvisor` method from `ziDAQServer`:

```
>>> help('zhinst.ziPython.ziDAQServer.pllAdvisor')
```

Help on built-in function `pllAdvisor` in `zhinst.ziPython.ziDAQServer`:

```
zhinst.ziPython.ziDAQServer.pllAdvisor = pllAdvisor(...)
  pllAdvisor( (ziDAQServer)arg1) -> ziPllAdvisor :
    Create a pllAdvisor class. This will start a thread for running an
    asynchronous pllAdvisor.
    arg1: Reference to the ziDAQServer class.
    arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED, ignored

  pllAdvisor( (ziDAQServer)arg1, (int)arg2) -> ziPllAdvisor
```

Reference help for the `ziPllAdvisor` class.

```
>>> help('zhinst.ziPython.ziPllAdvisor')
```

Help on class `ziPllAdvisor` in `zhinst.ziPython`:

```
zhinst.ziPython.ziPllAdvisor = class ziPllAdvisor(Boost.Python.instance)
| Method resolution order:
|   ziPllAdvisor
|   Boost.Python.instance
|   builtins.object
|
| Methods defined here:
|
|   __init__(...)
|       Raises an exception
|       This class cannot be instantiated from Python
|
|   __reduce__ = <unnamed Boost.Python function>(...)
|
|   clear(...)
|       clear( (ziPllAdvisor)arg1) -> None :
|           End the pllAdvisor thread.
|
|   execute(...)
|       execute( (ziPllAdvisor)arg1) -> None :
|           Starts the pllAdvisor if not yet running.
|
|   finish(...)
|       finish( (ziPllAdvisor)arg1) -> None :
|           Stop the pllAdvisor.
|
|   finished(...)
|       finished( (ziPllAdvisor)arg1) -> bool :
|           Check if the command execution has finished. Returns True if finished.
|
|   get(...)
|       get( (ziPllAdvisor)arg1, (str)arg2, (bool)arg3) -> object :
|           Return a dict with all nodes from the specified sub-tree.
|           arg1: Reference to the ziPllAdvisor class.
|           arg2: Path string of the node. Use wild card to
|               select all.
|           arg3[optional]: Specify which type of data structure to return.
|               Return data either as a flat dict (True) or as a nested
|               dict tree (False). Default = False.
|
|       get( (ziPllAdvisor)arg1, (str)arg2) -> object
|
|   listNodes(...)
|       listNodes( (ziPllAdvisor)arg1, (str)arg2, (int)arg3) -> list :
|           This function returns a list of node names found at the specified path.
```

```

|         arg1: Reference to the ziPllAdvisor class.
|         arg2: Path for which the nodes should be listed. The path may
|               contain wildcards so that the returned nodes do not
|               necessarily have to have the same parents.
|         arg3: Enum that specifies how the selected nodes are listed.
|               ziPython.ziListEnum.none -> 0x00
|                   The default flag, returning a simple
|                   listing if the given node
|               ziPython.ziListEnum.recursive -> 0x01
|                   Returns the nodes recursively
|               ziPython.ziListEnum.absolute -> 0x02
|                   Returns absolute paths
|               ziPython.ziListEnum.leafsonly -> 0x04
|                   Returns only nodes that are leafs,
|                   which means they are at the
|                   outermost level of the tree.
|               ziPython.ziListEnum.settingsonly -> 0x08
|                   Returns only nodes which are marked
|                   as setting
|               Or combinations of flags can be used.
|
|     progress(...)
|         progress( (ziPllAdvisor)arg1) -> object :
|             Reports the progress of the command with a number between
|             0 and 1.
|
|     read(...)
|         read( (ziPllAdvisor)arg1, (bool)arg2) -> object :
|             Read pllAdvisor data. If the simulation is still ongoing only a subset
|             of the data is returned.
|             arg1[optional]: Specify which type of data structure to return.
|             Return data either as a flat dict (True) or as a nested
|             dict tree (False). Default = False.
|
|         read( (ziPllAdvisor)arg1) -> object
|
|     save(...)
|         save( (ziPllAdvisor)arg1, (str)arg2) -> None :
|             Save PLL advisor data to file.
|             arg1: Reference to the ziPllAdvisor class.
|             arg2: File name string (without extension).
|
|     set(...)
|         set( (ziPllAdvisor)arg1, (str)arg2, (float)arg3) -> None :
|             PLL Advisor Parameters
|
|             Path name      Type      Description
|             pllAdvisor/bode      struct  Output parameter. Contains the resulting
|                                         bode plot of the PLL simulation.
|
|             pllAdvisor/calculate int      Command to calculate values. Set to 1 to
|                                         start the calculation.
|
|             pllAdvisor/center  double   Center frequency of the PLL oscillator.
|                                         The PLL frequency shift is relative to
|                                         this center frequency.
|
|             pllAdvisor/d        double   Differential gain.
|             pllAdvisor/demodbw  double   Demodulator bandwidth used for the PLL
|                                         loop filter
|
|             pllAdvisor/i        double   Integral gain.
|             pllAdvisor/mode     double   Select PLL Advisor mode. Currently only
|                                         one mode (open loop) is supported.
|
|             pllAdvisor/order    double   Demodulator order used for the PLL loop
|                                         filter.
|
|             pllAdvisor/p        double   Proportional gain.
|             pllAdvisor/pllbw    double   Demodulator bandwidth used for the PLL
|                                         loop filter.
|
|             pllAdvisor/pm       double   Output parameter. Simulated phase margin
|                                         of the PLL with the current settings. The
|                                         phase margin should be greater than 45 deg

```

```

|                                     and preferably greater than 65 deg for
|                                     stable conditions.
|         pllAdvisor/pmfreq           double Output parameter. Simulated phase margin
|                                     frequency.
|         pllAdvisor/q               double Quality factor. Currently not used.
|         pllAdvisor/rate            double PLL Advisor sampling rate of the PLL
|                                     control loop.
|         pllAdvisor/stable          int    Output parameter. When 1, the PLL Advisor
|                                     found a stable solution with the given
|                                     settings. When 0, revise your settings and
|                                     rerun the PLL Advisor.
|         pllAdvisor/targetbw       int    Requested PLL bandwidth. Higher
|                                     frequencies may need manual tuning.
|         pllAdvisor/targetfail     int    Output parameter. 1 indicates the
|                                     simulated PLL BW is smaller than the
|                                     Target BW.
|
|         set( (ziPllAdvisor)arg1, (str)arg2, (int)arg3) -> None
|
|         set( (ziPllAdvisor)arg1, (str)arg2, (str)arg3) -> None
|
|         set( (ziPllAdvisor)arg1, (object)arg2) -> None :
|             arg1: Reference to the ziPllAdvisor class.
|             arg2: A list of path/value pairs.
|
|         subscribe(...)
|             subscribe( (ziPllAdvisor)arg1, (str)arg2) -> None :
|                 Subscribe to one or several nodes.
|
|         trigger(...)
|             trigger( (ziPllAdvisor)arg1) -> None :
|                 Not applicable to this module.
|
|         unsubscribe(...)
|             unsubscribe( (ziPllAdvisor)arg1, (str)arg2) -> None :
|                 Unsubscribe from one or several nodes.
|
| -----
| Methods inherited from Boost.Python.instance:
|
| __new__(*args, **kwargs) from Boost.Python.class
|     Create and return a new object. See help(type) for accurate signature.
|
| -----
| Data descriptors inherited from Boost.Python.instance:
|
| __dict__
|
| __weakref__

```

4.4.9. Help for ziPython's ziPidAdvisor class

An instance of ziPidAdvisor is initialized using the pidAdvisor method from ziDAQServer:

```
>>> help('zhinst.ziPython.ziDAQServer.pidAdvisor')
```

Help on built-in function pidAdvisor in zhinst.ziPython.ziDAQServer:

```

zhinst.ziPython.ziDAQServer.pidAdvisor = pidAdvisor(...)
pidAdvisor( (ziDAQServer)arg1) -> ziPidAdvisor :
    Create a pidAdvisor class. This will start a thread for running an
    asynchronous pidAdvisor.
    arg1: Reference to the ziDAQServer class.
    arg2: Timeout in [ms]. Recommended value is 500ms. - DEPRECATED, ignored

```

```
pidAdvisor( (ziDAQServer)arg1, (int)arg2) -> ziPidAdvisor
```

Reference help for the ziPidAdvisor class.

```
>>> help('zhinst.ziPython.ziPidAdvisor')
```

Help on class ziPidAdvisor in zhinst.ziPython:

```
zhinst.ziPython.ziPidAdvisor = class ziPidAdvisor(Boost.Python.instance)
|   Method resolution order:
|       ziPidAdvisor
|       Boost.Python.instance
|       builtins.object
|
|   Methods defined here:
|
|   __init__(...)
|       Raises an exception
|       This class cannot be instantiated from Python
|
|   __reduce__ = <unnamed Boost.Python function>(...)
|
|   clear(...)
|       clear( (ziPidAdvisor)arg1) -> None :
|           End the pidAdvisor thread.
|
|   execute(...)
|       execute( (ziPidAdvisor)arg1) -> None :
|           Starts the pidAdvisor if not yet running.
|
|   finish(...)
|       finish( (ziPidAdvisor)arg1) -> None :
|           Stop the pidAdvisor.
|
|   finished(...)
|       finished( (ziPidAdvisor)arg1) -> bool :
|           Check if the command execution has finished. Returns True if finished.
|
|   get(...)
|       get( (ziPidAdvisor)arg1, (str)arg2, (bool)arg3) -> object :
|           Return a dict with all nodes from the specified sub-tree.
|           arg1: Reference to the ziPidAdvisor class.
|           arg2: Path string of the node. Use wild card to
|               select all.
|           arg3[optional]: Specify which type of data structure to return.
|               Return data either as a flat dict (True) or as a nested
|               dict tree (False). Default = False.
|
|       get( (ziPidAdvisor)arg1, (str)arg2) -> object
|
|   listNodes(...)
|       listNodes( (ziPidAdvisor)arg1, (str)arg2, (int)arg3) -> list :
|           This function returns a list of node names found at the specified path.
|           arg1: Reference to the ziPidAdvisor class.
|           arg2: Path for which the nodes should be listed. The path may
|               contain wildcards so that the returned nodes do not
|               necessarily have to have the same parents.
|           arg3: Enum that specifies how the selected nodes are listed.
|               ziPython.ziListEnum.none -> 0x00
|                   The default flag, returning a simple
|                   listing of the given node
|               ziPython.ziListEnum.recursive -> 0x01
|                   Returns the nodes recursively
|               ziPython.ziListEnum.absolute -> 0x02
|                   Returns absolute paths
|               ziPython.ziListEnum.leafsonly -> 0x04
```

```

|                                     Returns only leaf nodes,
|                                     which means the they are at the
|                                     outermost level of the tree.
|                                     ziPython.ziListEnum.settingsonly -> 0x08
|                                     Returns only nodes which are marked
|                                     as settings
|                                     Or combinations of flags can be used.
|
| progress(...)
|     progress( (ziPidAdvisor)arg1) -> object :
|         Reports the progress of the command with a number between
|         0 and 1.
|
| read(...)
|     read( (ziPidAdvisor)arg1, (bool)arg2) -> object :
|         Read pidAdvisor data. If the simulation is still ongoing, only a subset
|         of the data is returned.
|         arg1[optional]: Specify which type of data structure to return.
|         Return data either as a flat dict (True) or as a nested
|         dict tree (False). Default = False.
|
|     read( (ziPidAdvisor)arg1) -> object
|
| save(...)
|     save( (ziPidAdvisor)arg1, (str)arg2) -> None :
|         Save PID advisor data to file.
|         arg1: Reference to the ziPidAdvisor class.
|         arg2: File name string (without extension).
|
| set(...)
|     set( (ziPidAdvisor)arg1, (str)arg2, (float)arg3) -> None :
|         PID Advisor Parameters
|
|         Path name      Type      Description
|         pidAdvisor/advancedmode    int      Disable automatic calculation of
|                                           the start and stop value.
|
|         pidAdvisor/auto             int      Automatic response calculation
|                                           triggered by parameter change.
|
|         pidAdvisor/bode             struct   Output parameter. Contains the
|                                           resulting bode plot of the PID
|                                           simulation.
|
|         pidAdvisor/bw              double   Output parameter. Calculated system
|                                           bandwidth.
|
|         pidAdvisor/calculate        int      In/Out parameter. Command to
|                                           calculate values. Set to 1 to start
|                                           the calculation.
|
|         pidAdvisor/display/freqstart double   Start frequency for Bode plot.
|                                           For disabled advanced mode the
|                                           start value is automatically
|                                           derived from the system properties.
|
|         pidAdvisor/display/freqstop double   Stop frequency for Bode plot.
|         pidAdvisor/display/timestart double   Start time for step response.
|         pidAdvisor/display/timestop double   Stop time for step response.
|         pidAdvisor/dut/bw           double   Bandwidth of the DUT (device under
|                                           test).
|
|         pidAdvisor/dut/damping      double   Damping of the second order low
|                                           pass filter.
|
|         pidAdvisor/dut/delay        double   IO Delay of the feedback system
|                                           describing the earliest response
|                                           for a step change.
|
|         pidAdvisor/dut/fcenter      double   Resonant frequency of the of the
|                                           modelled resonator.
|
|         pidAdvisor/dut/gain         double   Gain of the DUT transfer function.
|         pidAdvisor/dut/q            double   quality factor of the modelled
|                                           resonator.
|
|         pidAdvisor/dut/source       int      Type of model used for the external
|                                           device to be controlled by the PID.
|                                           source = 1: Low-pass first order

```


		source = 2: Low-pass second order
		source = 3: Resonator frequency
		source = 4: Internal PLL
		source = 5: VCO
		source = 6: Resonator amplitude
pidAdvisor/impulse	struct	Output parameter. Impulse response (not yet supported).
pidAdvisor/index	int	PID index for parameter detection.
pidAdvisor/pid/autobw	int	Adjusts the demodulator bandwidth to fit best to the specified target bandwidth of the full system.
pidAdvisor/pid/d	double	In/Out parameter. Differential gain.
pidAdvisor/pid/dlimittimeconstant	double	In/Out parameter. Differential filter timeconstant.
pidAdvisor/pid/i	double	In/Out parameter. Integral gain.
pidAdvisor/pid/mode	double	Select PID Advisor mode. Mode value is bit coded, bit 0: P, bit 1: I, bit 2: D, bit 3: D filter limit.
pidAdvisor/pid/p	double	In/Out parameter. Proportional gain.
pidAdvisor/pid/rate	double	In/Out parameter. PID Advisor sampling rate of the PID control loop.
pidAdvisor/pid/targetbw	double	PID system target bandwidth.
pidAdvisor/pm	double	Output parameter. Simulated phase margin of the PID with the current settings. The phase margin should be greater than 45 deg and preferably greater than 65 deg for stable conditions.
pidAdvisor/pmfreq	double	Output parameter. Simulated phase margin frequency.
pidAdvisor/stable	int	Output parameter. When 1, the PID Advisor found a stable solution with the given settings. When 0, revise your settings and rerun the PID Advisor.
pidAdvisor/step	struct	Output parameter. Contains the resulting step response plot of the PID simulation.
pidAdvisor/targetbw	double	Requested PID bandwidth. Higher frequencies may need manual tuning.
pidAdvisor/targetfail	int	Output parameter. 1 indicates the simulated PID BW is smaller than the Target BW.
pidAdvisor/tf/closedloop	int	Switch the response calculation mode between closed or open loop.
pidAdvisor/tf/input	int	Start point for the plant response simulation for open or closed loops.
pidAdvisor/tf/output	int	End point for the plant response simulation for open or closed loops.
pidAdvisor/tune	int	Optimize the PID parameters so that the noise of the closed-loop system gets minimized.
set((ziPidAdvisor)arg1, (str)arg2, (int)arg3) -> None		
set((ziPidAdvisor)arg1, (str)arg2, (str)arg3) -> None		
set((ziPidAdvisor)arg1, (object)arg2) -> None :		
arg1: Reference to the ziPidAdvisor class.		
arg2: A list of path/value pairs.		

```
| subscribe(...)
|     subscribe( (ziPidAdvisor)arg1, (str)arg2) -> None :
|         Subscribe to one or several nodes.
|
| trigger(...)
|     trigger( (ziPidAdvisor)arg1) -> None :
|         Not applicable to this module.
|
| unsubscribe(...)
|     unsubscribe( (ziPidAdvisor)arg1, (str)arg2) -> None :
|         Unsubscribe from one or several nodes.
|
| -----
| Methods inherited from Boost.Python.instance:
|
| __new__(*args, **kwargs) from Boost.Python.class
|     Create and return a new object.  See help(type) for accurate signature.
|
| -----
| Data descriptors inherited from Boost.Python.instance:
|
| __dict__
|
| __weakref__
```

Chapter 5. LabVIEW Programming

Interfacing with your Zurich Instruments device via National Instruments' [LabVIEW®](#) is an efficient choice in terms of development time and run-time performance. LabVIEW is a graphical programming language designed to interface with laboratory equipment via so-called VIs ("virtual instruments"), whose key strength is the ease of displaying dynamic signals obtained from your instrument.

This chapter aims to help you get started using the Zurich Instruments LabOne LabVIEW API to control your instrument, please refer to:

- [Section 5.1](#) for help [Installing the LabOne LabVIEW API](#).
- [Section 5.2](#) for help [Getting Started with the LabOne LabVIEW API](#) and [running the examples](#).
- [Section 5.3](#) for some [LabVIEW Programming Tips and Tricks](#).

Note

This section and the provided examples are no substitute for a general LabVIEW tutorial. See, for example, the National Instruments [website](#) for help to get started programming with LabVIEW.

5.1. Installing the LabOne LabVIEW API

5.1.1. Requirements

A LabVIEW 2009 (or higher) installation is required on either Windows or Linux in order to use the LabOne LabVIEW API.

The LabOne LabVIEW API is included in a standard LabOne installation and is also available as a separate package (see below, [Separate LabVIEW Package](#)). In order to make the LabOne LabVIEW API available for use within LabVIEW, a directory needs to be copied to a specific directory of your LabVIEW installation. Both the main LabOne installer and the separate LabOne LabVIEW API package are available from Zurich Instruments' [download page](#).

Separate LabVIEW Package

The separate LabVIEW API package should be used if you would like to either:

1. Use the LabVIEW API on Mac OS X (the main LabOne installer is not available for Mac OS X).
2. Use the LabVIEW API to work with an instrument remotely (i.e., on a separate PC from where the Data Server is running) and you do not require a full LabOne installation. This is the case, for example, with MF Instruments.

5.1.2. Windows Installation

1. Locate the `instr.lib` directory in your LabVIEW installation and delete any previous Zurich Instruments API directories. The `instr.lib` directory is typically located at:

```
C:\Program Files\National Instruments\LabVIEW 201x\instr.lib\
```

Previous Zurich Instruments installations will be directories located in the `instr.lib` directory that are named either:

- Zurich Instruments HF2, or
- Zurich Instruments LabOne.

These folders may simply be deleted (administrator rights required).

2. On Windows, either navigate to the `API\LabVIEW` subdirectory of your LabOne installation or, in the case of the separate installer (see [Separate LabVIEW Package](#)), the directory of the unzipped LabOne LabVIEW package, and copy the subdirectory

```
Zurich Instruments LabOne
```

to the `instr.lib` directory in your LabVIEW installation as located in Step 1. Note, you will need administrator rights to copy to this directory.

In the case of copying from a LabOne installation, this folder is typically located at:

```
C:\Program Files\Zurich Instruments\LabOne\API\LabVIEW\
```

3. Restart LabVIEW and verify your installation as described in [Section 5.1.4](#).

5.1.3. Linux Installation

1. Locate the `instr.lib` directory in your LabVIEW installation and remove any previous Zurich Instruments API installations. The `instr.lib` directory is typically located at:

```
/usr/local/natinst/LabVIEW-201x/instr.lib/
```

Previous Zurich Instruments installations will be folders located in the `instr.lib` directory that are named either:

- Zurich Instruments HF2, or
- Zurich Instruments LabOne.

These folders may simply be deleted (administrator rights required).

2. Navigate to the path where you unpacked LabOne or the separate LabVIEW package and copy the subdirectory

```
Zurich Instruments LabOne/
```

to the `instr.lib` directory in your LabVIEW installation as located in Step 1. Note, you will need administrator rights to copy to this directory.

Note, when copying from the main LabOne tarball, the `Zurich Instruments LabOne/` directory is located in

```
[PATH]/LabOneLinux64/API/LabVIEW/
```

3. Restart LabVIEW and verify your installation as described in [Section 5.1.4](#).

5.1.4. Verifying your Installation

If the LabOne LabVIEW API palette can be accessed from within LabVIEW, the LabOne LabVIEW API is correctly installed. See [Section 5.2.1](#) for help finding the palette.

5.2. Getting Started with the LabOne LabVIEW API

5.2.1. Locating the LabOne LabVIEW VI Palette

In order to locate the LabOne LabVIEW VIs start LabVIEW and create a new VI. In the VI's "Block Diagram" (CTRL-e) you can to access the LabOne LabVIEW API palette with a mouse right-click and browsing the tree under "Instrument I/O" → "Instr. Drivers", see [Figure 5.1](#).

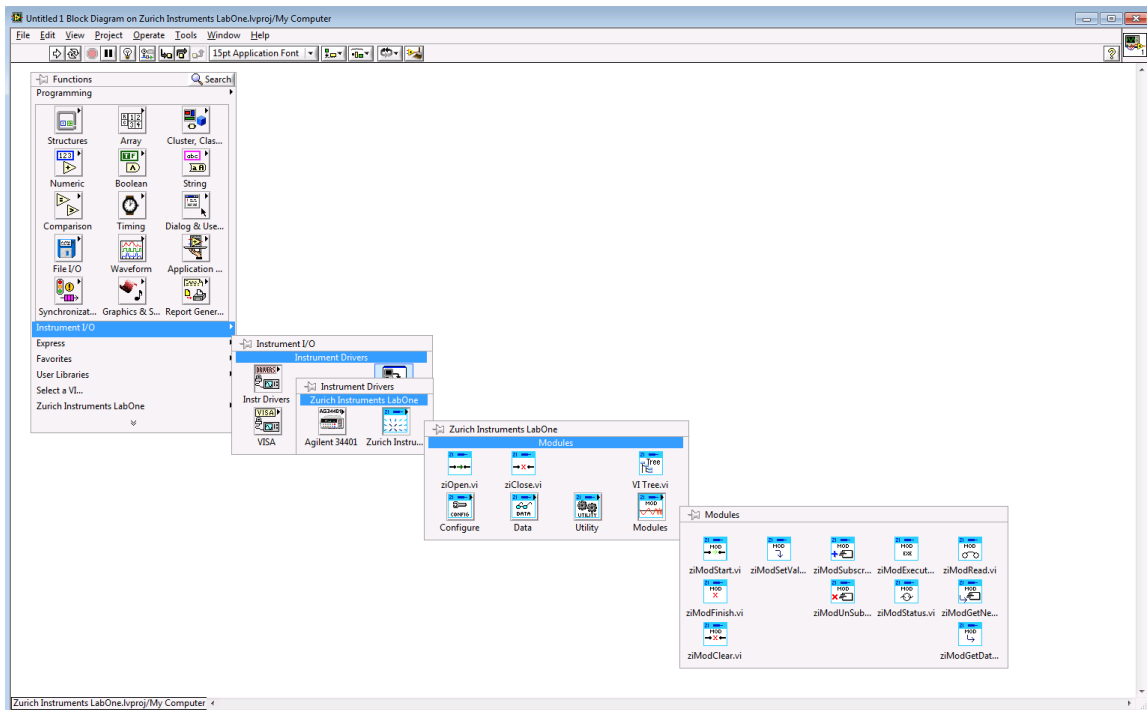


Figure 5.1. Locating the LabOne LabVIEW Palette

5.2.2. LabOne LabVIEW Programming Concepts

As described in [Section 1.1](#) a LabVIEW program communicates to a Zurich Instrument device via a software program running on the PC called the data server. In general, the outline of the instruction flow for a LabVIEW virtual instrument is as following:

1. Initialization: Open a connection from the API to the data server program.
2. Configuration: Perform the instrument's settings. For example, using the virtual instrument `ziSetValueDouble.vi`.
3. Data: Read data from the instrument.
4. Utility: Perform data analysis on the read data, potentially repeating Step 2 and/or Step 3.
5. Close: Terminate the API's connection to the data server program.

The `VI Tree.vi` included the LabOne LabVIEW API demonstrates this flow and lists common VIs used for working with a Zurich Instruments device, see [Figure 5.2](#). The `VI Tree.vi` can be found either via the LabOne VI palette, see [Section 5.2.1](#), or by opening the file in the `Public` folder of your LabOne LabVIEW installation, typically located at:

```
C:\Program Files\National Instruments\LabVIEW 2012\instr.lib\Zurich
Instruments LabOne\Public\VI Tree.vi
```

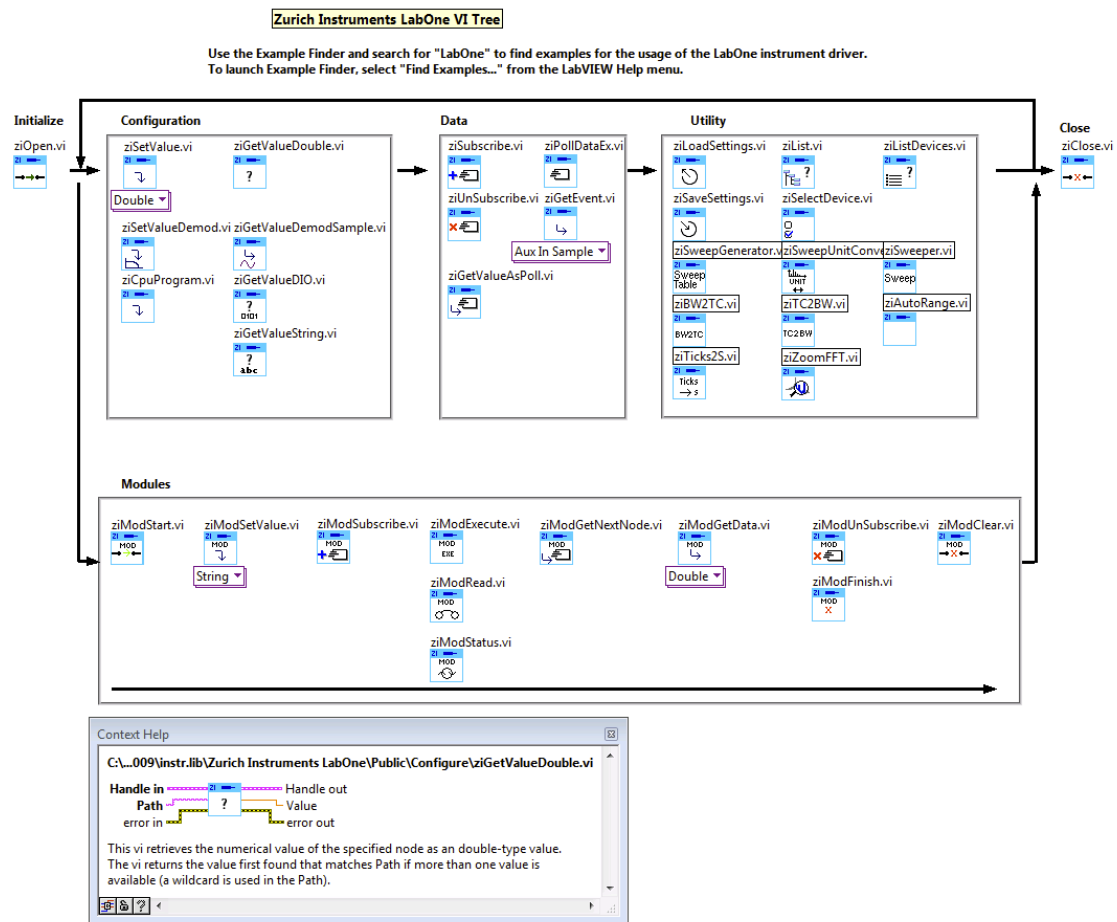


Figure 5.2. An overview of the LabOne LabVIEW VIs is given in VI Tree.vi. Press CTRL-h after selecting one of the VIs to obtain help.

5.2.3. Using ziCore Modules in the LabOne LabVIEW API

LabOne [ziCore Modules](#) (e.g. Sweeper) enable high-level measurement tools to use with your Zurich instrument device in LabVIEW. The outline of the instruction flow for a LabVIEW Module is as following:

1. Initialization: Create a `ziModHandle` from a `ziHandle` `ziModStart.vi`.
2. Configuration: Perform the module's settings. For example, using the virtual instrument `ziModSetValue.vi`.
3. Subscribe: Define the recorded data node `ziModSubscribe.vi`.
4. Execute: Start the operation of the module `ziModExecute.vi`.
5. Data: Read data from the module. For example, using the `ziModGetNextNode.vi` and `ziModGetData.vi`.
6. Utility: Perform data analysis on the read data, potentially repeating Step 2, Step 3 and/or Step 4.
7. Clear: Terminate the API's connection to the module `ziModClear.vi`.

5.2.4. Finding help for the LabOne VIs from within LabVIEW

As is customary for LabVIEW, built-in help for LabOne's VIs can be obtained by selecting the VI with the mouse in a block diagram and pressing CTRL-h to view the VI's context help. See [Figure 5.2](#) for an example.

5.2.5. Finding the LabOne LabVIEW API Examples

Many examples come bundled with the LabOne LabVIEW API which demonstrate the most important concepts of working with Zurich Instrument devices. The easiest way to browse the list of available examples is via the NI Example Finder: In LabVIEW select "Find Examples..." from the "Help" menu-bar and search for "LabOne", see [Figure 5.3](#).

The examples are located in the directory `instr.lib/Zurich Instruments LabOne/Examples` found in LabVIEW installation directory. In order to modify an example for your needs, please copy it to your local workspace.

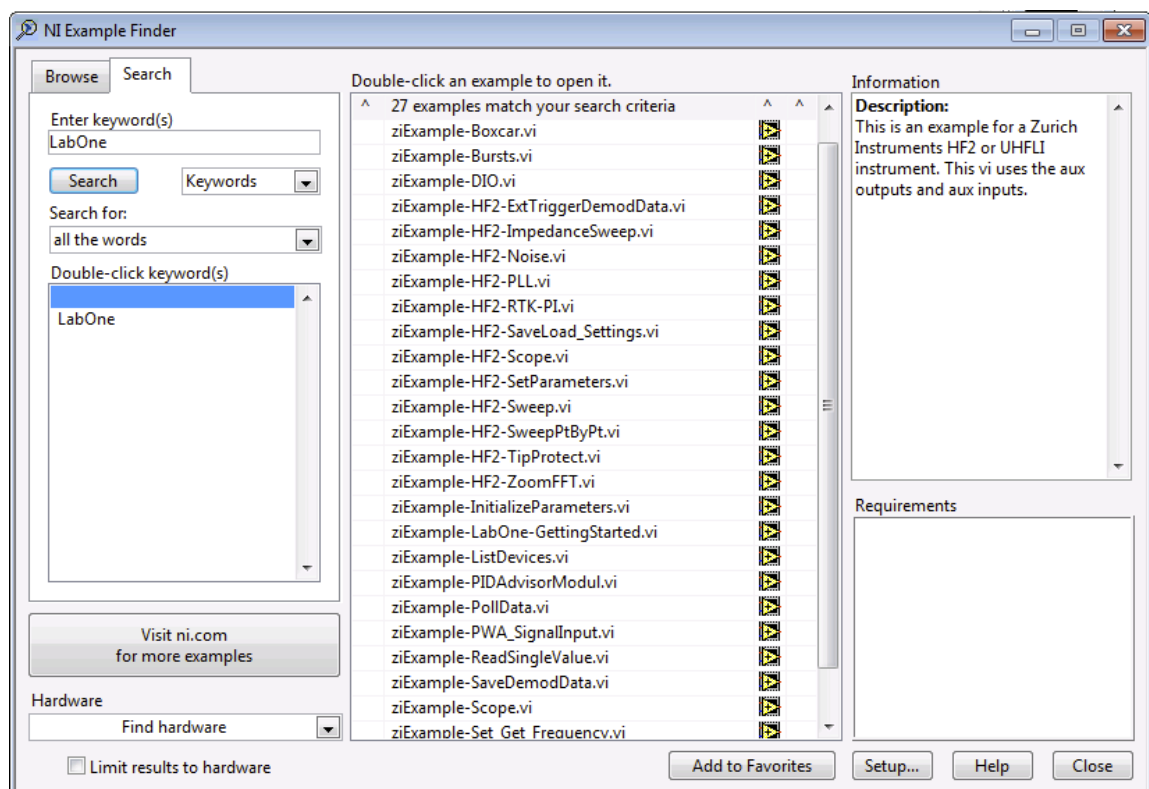


Figure 5.3. Search for "LabOne" in NI's Example Finder to find examples to run with your instrument.

5.2.6. Running the LabOne Example VIs

This section describes how to run a LabOne LabVIEW example on your instrument.

Note

Please ensure that the example you would like to run is supported by your instrument class and its options set. For example, examples for HF2 Instruments can be found in the Example Finder (see [Section 5.2.5](#)) by searching for "HF2", examples for the UHFLI by searching for "UHFLI" and examples for the MFLI by searching for "MFLI".

Device Connection

After opening one of the LabOne LabVIEW examples, please ensure that the example is configured to run on the desired instrument type. `ziOpen.vi` establishes a connection to a Data Server. The address is of the format `{<host>}{:<port>}:{<Device ID>}`. Usually it is sufficient to provide the Device ID only highlighted in [Figure 5.4](#). The Device ID corresponds to the serial number (S/N) found on the instrument rear panel. The host and port are then determined by network discovery. Should the discovery not work, prepend `<host>:<port>::` to the Device ID. Examples are "myhf2.company.com:8004::dev466" or "myhf2.company.com:8004". In the latter case the first found instrument on the data server listening on "myhf2.company.com:8004" will be selected.

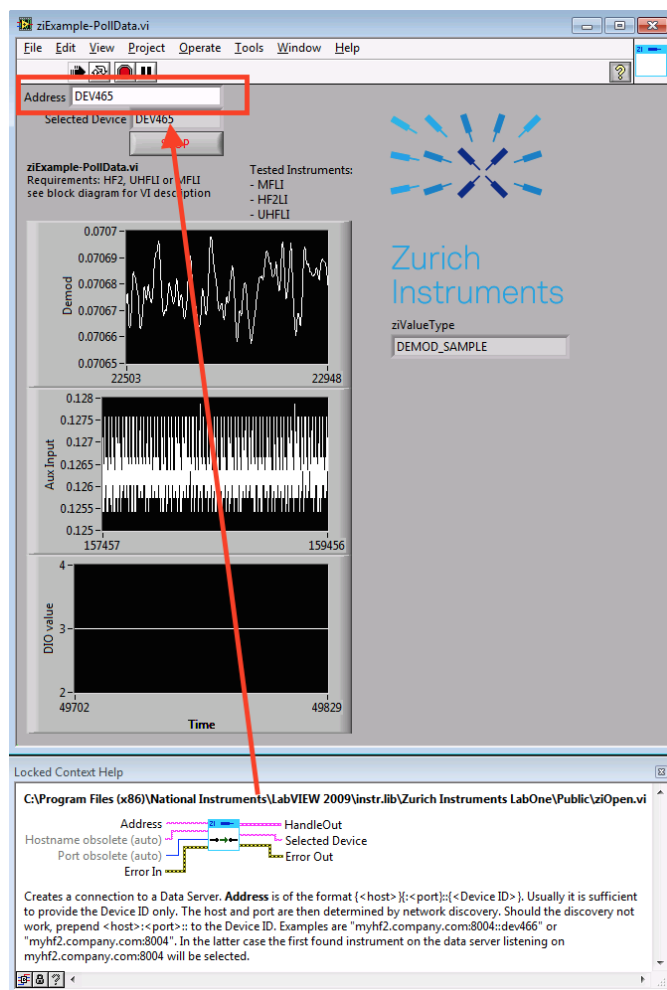


Figure 5.4. LabOne LabVIEW Example Poll Data: Device selection.

Running the VI and Block Diagram

The example can be ran as any LabVIEW program; by clicking the "Run" icon in the icon bar. Be sure to check the example's code and explanation by pressing CTRL-e to view the example's block diagram, see [Figure 5.5](#).

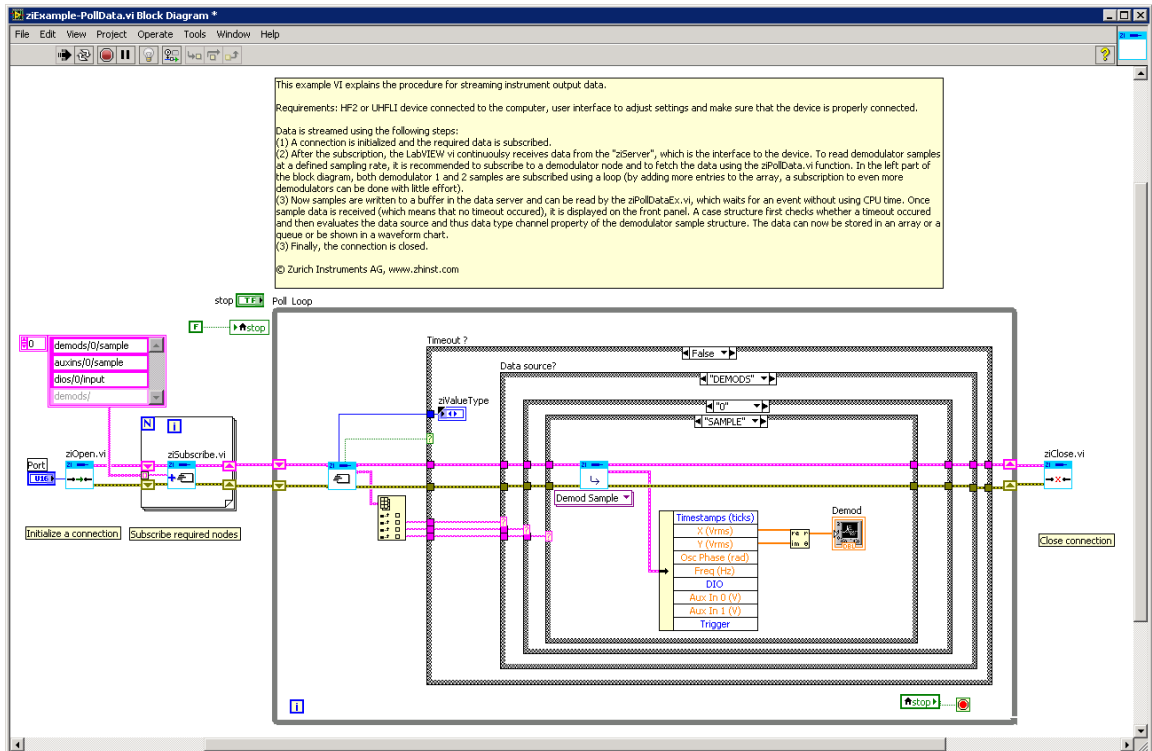


Figure 5.5. LabOne LabVIEW Example Poll Data: Block Diagram.

5.3. LabVIEW Programming Tips and Tricks

Use the User Interface's command log or Server's text interface while programming with LabVIEW

As with all other interfaces, LabVIEW uses the "path" and "nodes" concept to address settings on an instrument, see [Section 1.1](#). In order to learn about or verify the nodes available it can be very helpful to view the command log in the User Interface (see the bar in the bottom of the screen) to see which node has been configured during a previous setting change. The text interface (HF2 Series) provides a convenient way to explore the node hierarchy.

Always close ziHandles and ziModHandles or LabVIEW runs out of memory

If you use the "Abort Execution" button of LabVIEW, your LabVIEW program will not close any existing connections to the ziServer. Any open connection inside of LabVIEW will persist and continue to consume about 12 MB of RAM so that with time you will run out of memory. Completely exit LabVIEW in order to release the memory again.

Use shift registers

The structure of efficient LabVIEW code is distinguished by signals being "piped through" by use of shift registers in loops and by the absence of object replication. Using shift registers in LabVIEW avoids copying of data and, more important, running the garbage collector frequently.

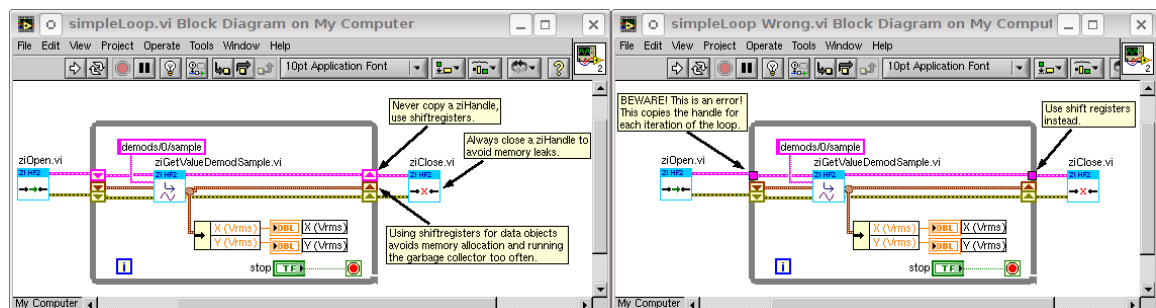


Figure 5.6. Examples of simple LabVIEW programs for the Zurich Instruments HF2 Series. Left: A well implemented loop, Right: An example for-loop gone wrong.

Chapter 6. C Programming

The LabOne C API, also known as ziAPI, provides a simple and robust way to communicate with the Data Server. It enables you to get or set parameters and receive streaming data.

6.1. Getting Started

After installing the LabOne software package and relevant drivers for your instrument you are ready start programming with ziAPI. All you need is a C compiler, linker and editor.

The structure of a program using ziAPI can be split into three parts: initialization/connection, data manipulation and disconnection/cleanup. The basic object that is always used is the ziConnection data structure. First, ziConnection has to be initialized by calling [ziAPIInit](#). After initialization ziConnection is ready to connect to a ziServer by calling [ziAPIConnect](#). Then ziConnection is ready to be used for getting and setting parameters and streaming data. When ziConnection is not needed anymore the established connection to the ziServer has to be hung up using [ziAPIDisconnect](#) before cleaning it up by calling [ziAPIDestroy](#).

6.1.1. Example

Below you find a simple program, which sets the demodulator rate of all demods for all devices.

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

int main() {
    ZIResult_enum retVal;
    char* errBuffer;

    ZIConnection conn;

    // Initialize ZIConnection.
    if ((retVal = ziAPIInit(&conn)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(retVal, &errBuffer, NULL);
        fprintf(stderr, "Can't init Connection: %s\n", errBuffer);
        return 1;
    }

    // Connect to the Data Server: Use port 8005 for the HF2 Data Server, use
    // 8004 for the UHF and MF Data Servers. HF2 only support ZI_API_VERSION_1,
    // see the LabOne Programming Manual for an explanation of API Levels.
    char serverAddress[] = "localhost";
    if ((retVal = ziAPIConnectEx(conn, serverAddress, 8004, ZI_API_VERSION_5, NULL)) !=
        ZI_INFO_SUCCESS) {
        ziAPIGetError(retVal, &errBuffer, NULL);
        fprintf(stderr, "Error, can't connect to the Data Server: `%s`.\n", errBuffer);
    } else {
        // Set all demodulator rates of all devices to 150 Hz
        if ((retVal = ziAPISetValueD(conn,
                                     "/dev1046/demods/*/rate",
                                     150)) != ZI_INFO_SUCCESS) {
            ziAPIGetError(retVal, &errBuffer, NULL);
            fprintf(stderr, "Can't set parameter: %s\n", errBuffer);
        }

        // Disconnect from the Data Server. Since ZIAPIDisconnect always returns
        // ZI_INFO_SUCCESS no error handling is required.
        ziAPIDisconnect(conn);
    }

    // Destroy the ZIConnection. Since ZIAPIDestroy always returns
    // ZI_INFO_SUCCESS, no error handling is required.
    ziAPIDestroy(conn);

    return 0;
}
```

}

6.2. Module Documentation

6.2.1. Connecting to Data Server

This section describes how to initialize the `ZIConnection` and establish a connection to Data Server as well as how to disconnect after all data handling is done and cleanup the `ZIConnection`.

Typedefs

- `typedef ZIConnection`
The `ZIConnection` is a connection reference; it holds information and helper variables about a connection to the Data Server. There is nothing in this reference which the user may use, so it is hidden and instead a dummy pointer is used. See [ziAPIInit](#) for how to create a `ZIConnection`.

Enumerations

- `enum ZIAPIVersion_enum { ZI_API_VERSION_0, ZI_API_VERSION_1, ZI_API_VERSION_4, ZI_API_VERSION_5 }`

Functions

- `ZIResult_enum ziAPIInit (ZIConnection* conn)`
Initializes a `ZIConnection` structure.
- `ZIResult_enum ziAPIDestroy (ZIConnection conn)`
Destroys a `ZIConnection` structure.
- `ZIResult_enum ziAPIConnect (ZIConnection conn, const char* hostname, uint16_t port)`
Connects the `ZIConnection` to Data Server.
- `ZIResult_enum ziAPIDisconnect (ZIConnection conn)`
Disconnects an established connection.
- `ZIResult_enum ziAPIListImplementations (char* implementations, uint32_t bufferSize)`
Returns the list of supported implementations.
- `ZIResult_enum ziAPIConnectEx (ZIConnection conn, const char* hostname, uint16_t port, ZIAPIVersion_enum apiLevel, const char* implementation)`
Connects to Data Server and enables extended `ziAPI`.
- `ZIResult_enum ziAPIGetConnectionAPILevel (ZIConnection conn, ZIAPIVersion_enum* apiLevel)`
Returns `ziAPI` level used for the connection `conn`.
- `ZIResult_enum ziAPIGetRevision (unsigned int* revision)`
Retrieves the revision of `ziAPI`.

Detailed Description

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

int main() {
    ZIResult_enum retVal;
    ZIConnection conn;
    char* errBuffer;

    // Initialize ZIConnection.
    if ((retVal = ziAPIInit(&conn)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(retVal, &errBuffer, NULL);
        fprintf(stderr, "Can't init Connection: %s\n", errBuffer);
        return 1;
    }

    // Connect to the Data Server: Use port 8005 for the HF2 Data Server, use
    // 8004 for the UHF and MF Data Servers. HF2 only support ZI_API_VERSION_1,
    // see the LabOne Programming Manual for an explanation of API Levels.
    char serverAddress[] = "localhost";
    if ((retVal = ziAPIConnectEx(conn, serverAddress, 8004, ZI_API_VERSION_5, NULL)) !=
        ZI_INFO_SUCCESS) {
        ziAPIGetError(retVal, &errBuffer, NULL);
        fprintf(stderr, "Error, can't connect to the Data Server: `%s`.\n", errBuffer);
    } else {
        /*
         * Do something using ZIConnection here.
         */

        // Since ZIAPIDisconnect always returns ZI_INFO_SUCCESS
        // no error handling is required.
        ziAPIDisconnect(conn);
    }

    // Since ZIAPIDestroy always returns ZI_INFO_SUCCESS
    // no error handling is required.
    ziAPIDestroy(conn);

    return 0;
}
```


Enumeration Type Documentation

enum ZIAPIVersion_enum

Enumerator:

- ZI_API_VERSION_0
- ZI_API_VERSION_1
- ZI_API_VERSION_4
- ZI_API_VERSION_5

Function Documentation

ziAPIInit

ZIResult_enum ziAPIInit (**ZIConnection*** conn)

Initializes a **ZIConnection** structure.

This function initializes the structure so that it is ready to connect to Data Server. It allocates memory and sets up the infrastructure needed.

Parameters:

[out] conn

Pointer to **ZIConnection** that is to be initialized

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_MALLOC on memory allocation failure
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIDestroy](#), [ziAPIConnect](#), [ziAPIDisconnect](#)

See [Connection](#) for an example

ziAPIDestroy

ZIResult_enum ziAPIDestroy (**ZIConnection** conn)

Destroys a **ZIConnection** structure.

This function frees all memory that has been allocated by **ziAPIInit**. If it is called with an uninitialized **ZIConnection** struct it may result in segmentation faults as well when it is called with a struct for which **ziAPIDestroy** already has been called.

Parameters:

[in] conn

Pointer to **ZIConnection** struct that has to be destroyed

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

ziAPIInit, **ziAPIConnect**, **ziAPIDisconnect**

See **Connection** for an example

ziAPIConnect

ZIResult_enum ziAPIConnect (**ZIConnection** conn, const char* hostname, uint16_t port)

Connects the ZIConnection to Data Server.

Connects to Data Server using a **ZIConnection** and prepares for data exchange. For most cases it is enough to just give a reference to the connection and give NULL for hostname and 0 for the port, so it connects to localhost on the default port.

Parameters:

[in] conn

Pointer to **ZIConnection** with which the connection should be established

[in] hostname

Name of the Host to which it should be connected, if NULL "localhost" will be used as default

[in] port

The Number of the port to connect to. If 0, default port of the local Data Server will be used (8005)

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_HOSTNAME if the given host name could not be found
- ZI_ERROR_SOCKET_CONNECT if no connection could be established
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_SOCKET_INIT if initialization of the socket failed
- ZI_ERROR_CONNECTION when the Data Server didn't return the correct answer
- ZI_ERROR_TIMEOUT when initial communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIDisconnect](#), [ziAPIInit](#), [ziAPIDestroy](#)

See [Connection](#) for an example

ziAPIDisconnect

ZIResult_enum `ziAPIDisconnect (ZIConnection conn)`

Disconnects an established connection.

Disconnects from Data Server. If the connection has not been established and the function is called it returns without doing anything.

Parameters:

[in] conn
Pointer to ZIConnection to be disconnected

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIConnect](#), [ziAPIInit](#), [ziAPIDestroy](#)

See [Connection](#) for an example

ziAPIListImplementations

ZIResult_enum **ziAPIListImplementations** (char* implementations, uint32_t bufferSize)

Returns the list of supported implementations.

Returned names are defined by implementations in the linked library and may change depending on software version.

Parameters:

[out] implementations

Pointer to a buffer receiving a newline-delimited list of the names of all the supported ziAPI implementations. The string is zero-terminated.

[in] bufferSize

The size of the buffer assigned to the implementations parameter

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_LENGTH if the length of the char-buffer given by MaxLen is too small for all elements
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIConnectEx](#)

ziAPIConnectEx

ZIResult_enum ziAPIConnectEx (**ZIConnection** conn, const char* hostname, uint16_t port, ZIAPIVersion_enum apiLevel, const char* implementation)

Connects to Data Server and enables extended ziAPI.

With apiLevel=ZI_API_VERSION_1 and implementation=NULL, this call is equivalent to plain [ziAPIConnect](#). With other version and implementation values enables corresponding ziAPI extension and connection using different implementation.

Parameters:

[in] conn

Pointer to the ZIConnection with which the connection should be established

[in] hostname

Name of the host to which it should be connected, if NULL "localhost" will be used as default

[in] port

The number of the port to connect to. If 0 the port of the local Data Server will be used

[in] apiLevel

Specifies the ziAPI compatibility level to use for this connection (1 or 4).

[in] implementation

Specifies implementation to use for a connection, must be one of the returned by ziAPIListImplementations or NULL to select default implementation

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_HOSTNAME if the given host name could not be found
- ZI_ERROR_SOCKET_CONNECT if no connection could be established
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_SOCKET_INIT if initialization of the socket failed
- ZI_ERROR_CONNECTION when the Data Server didn't return the correct answer or requested implementation is not found or doesn't support requested ziAPI level
- ZI_ERROR_TIMEOUT when initial communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIListImplementations](#), [ziAPIConnect](#), [ziAPIDisconnect](#), [ziAPIInit](#), [ziAPIDestroy](#), [ziAPIGetConnectionVersion](#)

See [Connection](#) for an example

ziAPIGetConnectionAPILevel

ZIResult_enum ziAPIGetConnectionAPILevel (**ZIConnection** conn, ZIAPIVersion_enum* apiLevel)

Returns ziAPI level used for the connection conn.

Parameters:

[in] conn

Pointer to ZIConnection

[out] apiLevel

Pointer to preallocated ZIAPIVersion_enum, receiving the ziAPI level

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION if level can not be determined due to conn is not connected

See Also:

[ziAPIConnectEx](#), [ziAPIGetVersion](#)

ziAPIGetRevision

ZIResult_enum ziAPIGetRevision (unsigned int* revision)

Retrieves the revision of ziAPI.

Sets an unsigned int with the revision (build number) of the ziAPI you are using.

Parameters:

[in] revision

Pointer to an unsigned int to fill up with the revision.

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

6.2.2. Tree

All parameters and streams are organized in a tree. You can list the whole tree, parts of it or single items using [ziAPIListNodes](#) or you may update the tree with nodes of newly connected devices by using [ziAPIUpdateDevices](#).

Enumerations

- `enum ZIListNodes_enum { ZI_LIST_NODES_NONE, ZI_LIST_NODES_RECURSIVE, ZI_LIST_NODES_ABSOLUTE, ZI_LIST_NODES_LEAFONLY, ZI_LIST_NODES_SETTINGSONLY, ZI_LIST_NONE, ZI_LIST_RECURSIVE, ZI_LIST_ABSOLUTE, ZI_LIST_LEAFONLY, ZI_LIST_SETTINGSONLY }`

Defines the values of the flags used in [ziAPIListNodes](#).

Functions

- `ZIResult_enum ziAPIListNodes (ZIConnection conn, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)`
Returns all child nodes found at the specified path.
- `ZIResult_enum ziAPIUpdateDevices (ZIConnection conn)`
Search for the newly connected devices and update the tree.
- `ZIResult_enum ziAPIConnectDevice (ZIConnection conn, const char* deviceSerial, const char* deviceInterface, const char* interfaceParams)`
Connect a device to the server.
- `ZIResult_enum ziAPIDisconnectDevice (ZIConnection conn, const char* deviceSerial)`
Disconnect a device from the server.

Detailed Description

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

void PrintChildren(ZIConnection Conn,
                  char* Path) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    char NodesBuffer[8192];

    if ((RetVal = ziAPIListNodes(Conn,
                                Path,
                                NodesBuffer,
                                8192,
                                ZI_LIST_NODES_NONE)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Can't List Nodes: %s\n", ErrBuffer);
    } else {
        char* Ptr = NodesBuffer;
        char* LastPtr = Ptr;
```

```
// print out each node on a separate line with dash as prefix
for (; *Ptr != 0; Ptr++) {
    if (*Ptr == '\n') {
        *Ptr = 0;
        printf("- %s\n", LastPtr);
        LastPtr = Ptr + 1;
    }
}

// print out the last node
if (Ptr != LastPtr) {
    printf("- %s\n", LastPtr);
}
}
```

Enumeration Type Documentation

enum ZIListNodes_enum

Defines the values of the flags used in [ziAPIListNodes](#).

Enumerator:

- **ZI_LIST_NODES_NONE**
Default, return a simple listing of the given node immediate descendants.
- **ZI_LIST_NODES_RECURSIVE**
List the nodes recursively.
- **ZI_LIST_NODES_ABSOLUTE**
Return absolute paths.
- **ZI_LIST_NODES_LEAFONLY**
Return only leaf nodes, which means the nodes at the outermost level of the tree.
- **ZI_LIST_NODES_SETTINGSONLY**
Return only nodes which are marked as setting.
- **ZI_LIST_NONE**
Default, return a simple listing of the given node immediate descendants.
- **ZI_LIST_RECURSIVE**
List the nodes recursively.
- **ZI_LIST_ABSOLUTE**
Return absolute paths.
- **ZI_LIST_LEAFONLY**
Return only leaf nodes, which means the nodes at the outermost level of the tree.
- **ZI_LIST_SETTINGSONLY**
Return only nodes which are marked as setting.

Function Documentation

ziAPIListNodes

ZIResult_enum ziAPIListNodes (**ZIConnection** conn, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)

Returns all child nodes found at the specified path.

This function returns a list of node names found at the specified path. The path may contain wildcards so that the returned nodes do not necessarily have to have the same parents. The list is returned in a null-terminated char-buffer, each element delimited by a newline. If the maximum length of the buffer (bufferSize) is not sufficient for all elements, nothing will be returned and the return value will be **ZI_LENGTH**.

Parameters:

[in] conn

Pointer to the ZIConnection for which the node names should be retrieved.

[in] path

Path for which all children will be returned. The path may contain wildcard characters.

[out] nodes

Upon call filled with newline-delimited list of the names of all the children found. The string is zero-terminated.

[in] bufferSize

The length of the buffer used for the nodes output parameter.

[in] flags

A combination of flags (applied bitwise) as defined in **ZIListNodes_enum**.

Returns:

- **ZI_INFO_SUCCESS** on success
- **ZI_ERROR_CONNECTION** when the connection is invalid (not connected) or when a communication error occurred
- **ZI_ERROR_LENGTH** if the path's length exceeds **MAX_PATH_LEN** or the length of the char-buffer for the nodes given by bufferSize is too small for all elements
- **ZI_WARNING_OVERFLOW** when a FIFO overflow occurred
- **ZI_ERROR_COMMAND** on an incorrect answer of the server
- **ZI_ERROR_SERVER_INTERNAL** if an internal error occurred in Data Server
- **ZI_WARNING_NOTFOUND** if the given path could not be resolved
- **ZI_ERROR_TIMEOUT** when communication timed out
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See [Tree Listing](#) for an example

See Also:

ziAPIUpdate

ziAPIUpdateDevices

ZIResult_enum ziAPIUpdateDevices (**ZIConnection** conn)

Search for the newly connected devices and update the tree.

This function forces the Data Server to search for newly connected devices and to connect to run them

Parameters:

[in] conn
Pointer to ZIConnection

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIListNodes](#)

ziAPIConnectDevice

ZIResult_enum ziAPIConnectDevice (**ZIConnection** conn, const char* deviceSerial, const char* deviceInterface, const char* interfaceParams)

Connect a device to the server.

This function connects a device with deviceSerial via the specified deviceInterface for use with the server.

Parameters:

[in] conn

Pointer to the ZIConnection with which the connection should be established

[in] deviceSerial

The serial of the device to connect to, e.g., dev2100

[in] deviceInterface

The interface to use for the connection, e.g., USB|1GbE

[in] interfaceParams

Parameters for interface configuration

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIDisconnectDevice](#), [ziAPIConnect](#), [ziAPIDisconnect](#), [ziAPIInit](#)

ziAPIDisconnectDevice

ZIResult_enum `ziAPIDisconnectDevice (ZIConnection conn, const char* deviceSerial)`

Disconnect a device from the server.

This function disconnects a device specified by deviceSerial from the server.

Parameters:

[in] conn

Pointer to the [ZIConnection](#) with which the connection should be established

[in] deviceSerial

The serial of the device to connect to, e.g., dev2100

Returns:

- [ZI_INFO_SUCCESS](#) on success
- [ZI_ERROR_TIMEOUT](#) when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIConnectDevice](#), [ziAPIConnect](#), [ziAPIDisconnect](#), [ziAPIInit](#)

6.2.3. Set and Get Parameters

This section describes several functions for getting and setting parameters of different datatypes.

Functions

- `ZIResult_enum` `ziAPIGetValueD (ZIConnection conn, const char* path, ZIDoubleData* value)`
gets the double-type value of the specified node
- `ZIResult_enum` `ziAPIGetValueI (ZIConnection conn, const char* path, ZIntegerData* value)`
gets the integer-type value of the specified node
- `ZIResult_enum` `ziAPIGetDemodSample (ZIConnection conn, const char* path, ZIDemodSample* value)`
Gets the demodulator sample value of the specified node.
- `ZIResult_enum` `ziAPIGetDIOSample (ZIConnection conn, const char* path, ZIDIOSample* value)`
Gets the Digital I/O sample of the specified node.
- `ZIResult_enum` `ziAPIGetAuxInSample (ZIConnection conn, const char* path, ZIAuxInSample* value)`
gets the AuxIn sample of the specified node
- `ZIResult_enum` `ziAPIGetValueB (ZIConnection conn, const char* path, unsigned char* buffer, unsigned int* length, unsigned int bufferSize)`
gets the Bytearray value of the specified node
- `ZIResult_enum` `ziAPISetValueD (ZIConnection conn, const char* path, ZIDoubleData value)`
asynchronously sets a double-type value to one or more nodes specified in the path
- `ZIResult_enum` `ziAPISetValueI (ZIConnection conn, const char* path, ZIntegerData value)`
asynchronously sets an integer-type value to one or more nodes specified in a path
- `ZIResult_enum` `ziAPISetValueB (ZIConnection conn, const char* path, unsigned char* buffer, unsigned int length)`
asynchronously sets the binary-type value of one ore more nodes specified in the path
- `ZIResult_enum` `ziAPISyncSetValueD (ZIConnection conn, const char* path, ZIDoubleData* value)`
synchronously sets a double-type value to one or more nodes specified in the path
- `ZIResult_enum` `ziAPISyncSetValueI (ZIConnection conn, const char* path, ZIntegerData* value)`
synchronously sets an integer-type value to one or more nodes specified in a path

- `ZIResult_enum` `ziAPISyncSetValueB (ZIConnection conn, const char* path, uint8_t* buffer, uint32_t* length, uint32_t bufferSize)`
Synchronously sets the binary-type value of one ore more nodes specified in the path.
- `ZIResult_enum` `ziAPISync (ZIConnection conn)`
Synchronizes the session by dropping all pending data.
- `ZIResult_enum` `ziAPIEchoDevice (ZIConnection conn, const char* deviceSerial)`
Sends an echo command to a device and blocks until answer is received.
- `__inline ZIResult_enum` `ziAPIGetValueS (ZIConnection conn, char* path, DemodSample* value)`
- `__inline ZIResult_enum` `ziAPIGetValueDIO (ZIConnection conn, char* path, DIOSample* value)`
- `__inline ZIResult_enum` `ziAPIGetValueAuxIn (ZIConnection conn, char* path, AuxInSample* value)`

Function Documentation

ziAPIGetValueD

ZIResult_enum ziAPIGetValueD (**ZIConnection** conn, const char* path, ZIDoubleData* value)

gets the double-type value of the specified node

This function retrieves the numerical value of the specified node as an double-type value. The value first found is returned if more than one value is available (a wildcard is used in the path).

Parameters:

[in] conn

Pointer to ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to a double in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the path's length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIDoubleData ValueD;

    if ((RetVal = ziAPISetValueI(Conn,
                                "/dev1046/demods/*/rate",
                                100)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }
}
```

```
if ((RetVal = ziAPIGetValueD(Conn,
                            "/dev1046/demods/0/rate",
                            &ValueD)) != ZI\_INFO\_SUCCESS) {
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
} else {
    printf("Value = %f\n", ValueD);
}
}
```

See Also:

[ziAPISetValueD](#), [ziAPIGetValueAsPollData](#)

ziAPIGetValueI

ZIResult_enum ziAPIGetValueI (**ZIConnection** conn, const char* path, ZIIntegerData* value)

gets the integer-type value of the specified node

This function retrieves the numerical value of the specified node as an integer-type value. The value first found is returned if more than one value is available (a wildcard is used in the path).

Parameters:

[in] conn

Pointer to ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to an 64bit integer in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the path's length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIIntegerData ValueI;

    if ((RetVal = ziAPISetValueD(Conn,
                                "/dev1046/demods/*/rate",
                                5.53)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }

    if ((RetVal = ziAPIGetValueI(Conn,
                                "/dev1046/demods/0/rate",
```

```
                                &ValueI)) != ZI_INFO_SUCCESS) {
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
} else {
    printf("Value = %f\n", (float)ValueI);
}
}
```

See Also:

[ziAPISetValueI](#), [ziAPIGetValueAsPollData](#)

ziAPIGetDemodSample

ZIResult_enum ziAPIGetDemodSample (**ZIConnection** conn, const char* path, **ZIDemodSample*** value)

Gets the demodulator sample value of the specified node.

This function retrieves the value of the specified node as an **DemodSample** struct. The value first found is returned if more than one value is available (a wildcard is used in the path). This function is only applicable to paths matching DEMODS/[0-9]+/SAMPLE.

Parameters:

[in] conn

Pointer to ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to a **ZIDemodSample** struct in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the path's length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void GetSample(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    ZIDemodSample DemodSample;

    if ((RetVal = ziAPIGetDemodSample(Conn,
                                     "/dev1046/demods/0/sample",
                                     &DemodSample)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
```



```
        printf("TS = %f, X=%f, Y=%f\n",  
              (float)DemodSample.timeStamp,  
              DemodSample.x,  
              DemodSample.y);  
    }  
}
```

See Also:

[ziAPIGetValueAsPollData](#)

ziAPIGetDIOSample

ZIResult_enum ziAPIGetDIOSample (**ZIConnection** conn, const char* path, **ZIDIOSample*** value)

Gets the Digital I/O sample of the specified node.

This function retrieves the newest available DIO sample from the specified node. The value first found is returned if more than one value is available (a wildcard is used in the path). This function is only applicable to nodes ending in "/DIOS/[0-9]+/INPUT".

Parameters:

[in] conn

Pointer to the ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to a **ZIDIOSample** struct in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void GetSample(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    ZIDIOSample DIOSample;

    if ((RetVal = ziAPIGetDIOSample(Conn,
                                    "/dev1046/dios/0/output",
                                    &DIOSample)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
```

```
        printf("TS = %f, bits=%08x\n",  
              (float)DIOSample.timeStamp,  
              DIOSample.bits);  
    }  
}
```

See Also:

[ziAPIGetValueAsPollData](#)

ziAPIGetAuxInSample

ZIResult_enum ziAPIGetAuxInSample (**ZIConnection** conn, const char* path, **ZIAuxInSample*** value)

gets the AuxIn sample of the specified node

This function retrieves the newest available AuxIn sample from the specified node. The value first found is returned if more than one value is available (a wildcard is used in the path). This function is only applicable to nodes ending in "/AUXINS/[0-9]+/SAMPLE".

Parameters:

[in] conn

Pointer to the ziConnection with which the Value should be retrieved

[in] path

Path to the Node holding the value

[out] value

Pointer to an **ZIAuxInSample** struct in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void GetSample(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    ZIAuxInSample AuxInSample;

    if ((RetVal = ziAPIGetAuxInSample(Conn,
                                     "/dev1046/auxins/0/sample",
                                     &AuxInSample)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s\n", ErrBuffer);
    }
}
```

```
    } else {  
        printf("TS = %f, ch0=%f, ch1=%f\n",  
              (float)AuxInSample.timeStamp,  
              AuxInSample.ch0,  
              AuxInSample.ch1);  
    }  
}
```

See Also:

[ziAPIGetValueAsPollData](#)

ziAPIGetValueB

ZIResult_enum ziAPIGetValueB (**ZIConnection** conn, const char* path, unsigned char* buffer, unsigned int* length, unsigned int bufferSize)

gets the Bytearray value of the specified node

This function retrieves the newest available DIO sample from the specified node. The value first found is returned if more than one value is available (a wildcard is used in the path).

Parameters:

[in] conn

Pointer to the ziConnection with which the value should be retrieved

[in] path

Path to the Node holding the value

[out] buffer

Pointer to a buffer to store the retrieved data in

[out] length

Pointer to an unsigned int to store the length of data in. if an error occurred or the length of the passed buffer doesn't reach a zero will be returned

[in] bufferSize

The length of the passed buffer

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void PrintVersion(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
```

```
const char* Path = "ZI/ABOUT/VERSION";
unsigned char Buffer[0xff];
unsigned int Length;

if ((RetVal = ziAPIGetValueB(Conn,
                             Path,
                             Buffer,
                             &Length,
                             sizeof(Buffer) - 1)) != ZI_INFO_SUCCESS) {
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't get value: %s.\n", ErrBuffer);
} else {
    Buffer[Length] = 0;
    printf("%s=\"%s\"\n", Path, Buffer);
}
}
```

See Also:

[ziAPISetValueB](#), [ziAPIGetValueAsPollData](#)

ziAPISetValueD

ZIResult_enum ziAPISetValueD (**ZIConnection** conn, const char* path, ZIDoubleData value)

asynchronously sets a double-type value to one or more nodes specified in the path

This function sets the values of the nodes specified in path to Value. More than one value can be set if a wildcard is used. The function sets the value asynchronously which means that after the function returns you have no security to which value it is finally set nor at what point in time it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set.

[in] path

Path to the Node(s) for which the value(s) will be set to Value.

[in] value

The double-type value that will be written to the node(s).

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred.
- ZI_ERROR_READONLY on attempt to set a read-only node.
- ZI_ERROR_COMMAND on an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIIntegerData ValueI;

    if ((RetVal = ziAPISetValueD(Conn,
                                "/dev1046/demods/*/rate",
                                5.53)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }
}
```



```
    }

    if ((RetVal = ziAPIGetValueI(Conn,
                                "/dev1046/demods/0/rate",
                                &ValueI)) != ZI\_INFO\_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
        printf("Value = %f\n", (float)ValueI);
    }
}
```

See Also:

[ziAPIGetValueD](#). [ziAPISyncSetValueD](#)

ziAPISetValueI

ZIResult_enum ziAPISetValueI (**ZIConnection** conn, const char* path, ZIIntegerData value)

asynchronously sets an integer-type value to one or more nodes specified in a path

This function sets the values of the nodes specified in path to Value. More than one value can be set if a wildcard is used. The function sets the value asynchronously which means that after the function returns you have no security to which value it is finally set nor at what point in time it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set

[in] value

The int-type value that will be written to the node(s)

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred.
- ZI_ERROR_READONLY on attempt to set a read-only node.
- ZI_ERROR_COMMAND on an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIDoubleData ValueD;

    if ((RetVal = ziAPISetValueI(Conn,
                                "/dev1046/demods/*/rate",
                                100)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }
}
```

```
    }

    if ((RetVal = ziAPIGetValueD(Conn,
                                "/dev1046/demods/0/rate",
                                &ValueD)) != ZI\_INFO\_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
        printf("Value = %f\n", ValueD);
    }
}
```

See Also:

[ziAPIGetValueI](#). [ziAPISyncSetValueI](#)

ziAPISetValueB

ZIResult_enum ziAPISetValueB (**ZIConnection** conn, const char* path, unsigned char* buffer, unsigned int length)

asynchronously sets the binary-type value of one ore more nodes specified in the path

This function sets the values at the nodes specified in a path. More than one value can be set if a wildcard is used. The function sets the value asynchronously which means that after the function returns you have no security to which value it is finally set nor at what point in time it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set

[in] buffer

Pointer to the byte array with the data

[in] length

Length of the data in the buffer

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred.
- ZI_ERROR_READONLY on attempt to set a read-only node.
- ZI_ERROR_COMMAND on an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values.
- ZI_ERROR_TIMEOUT when communication timed out.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void ProgramCPU(ZIConnection Conn,
                unsigned char* Buffer,
                int Len) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    if ((RetVal = ziAPISetValueB(Conn,
```

```
        "/dev1046/cpus/0/program",
        Buffer,
        Len)) != ZI_INFO_SUCCESS) {
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
}
}
```

See Also:

[ziAPIGetValueB](#). [ziAPISyncSetValueB](#)

ziAPISyncSetValueD

ZIResult_enum ziAPISyncSetValueD (**ZIConnection** conn, const char* path, ZIDoubleData* value)

synchronously sets a double-type value to one or more nodes specified in the path

This function sets the values of the nodes specified in path to Value. More than one value can be set if a wildcard is used. The function sets the value synchronously. After returning you know that it is set and to which value it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set to value

[in] value

Pointer to a double-type containing the value to be written. When the function returns value holds the effectively written value.

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_READONLY on attempt to set a read-only node
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIGetValueD](#), [ziAPISetValueD](#)

ziAPISyncSetValueI

ZIResult_enum ziAPISyncSetValueI (**ZIConnection** conn, const char* path, **ZIIntegerData*** value)

synchronously sets an integer-type value to one or more nodes specified in a path

This function sets the values of the nodes specified in path to value. More than one value can be set if a wildcard is used. The function sets the value synchronously. After returning you know that it is set and to which value it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the node(s) for which the value(s) will be set

[in] value

Pointer to a int-type containing then value to be written. when the function returns value holds the effectively written value.

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_READONLY on attempt to set a read-only node
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIGetValueI](#), [ziAPISetValueI](#)

ziAPISyncSetValueB

ZIResult_enum ziAPISyncSetValueB (**ZIConnection** conn, const char* path, uint8_t* buffer, uint32_t* length, uint32_t bufferSize)

Synchronously sets the binary-type value of one or more nodes specified in the path.

This function sets the values at the nodes specified in a path. More than one value can be set if a wildcard is used. This function sets the value synchronously. After returning you know that it is set and to which value it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set

[in] buffer

Pointer to the byte array with the data

[in] length

Length of the data in the buffer

[in] bufferSize

Length of the data in the buffer

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_READONLY on attempt to set a read-only node
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIGetValueB](#), [ziAPISetValueB](#)

ziAPISync

ZIResult_enum ziAPISync (**ZIConnection** conn)

Synchronizes the session by dropping all pending data.

This function drops any data that is pending for transfer. Any data (including poll data) retrieved afterwards is guaranteed to be produced not earlier than the call to ziAPISync. This ensures in particular that any settings made prior to the call to ziAPISync have been propagated to the device, and the data retrieved afterwards is produced with the new settings already set to the hardware. Note, however, that this does not include any required settling time.

Parameters:

[in] conn

Pointer to the ZIConnection that is to be synchronized

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIEchoDevice

ZIResult_enum ziAPIEchoDevice (**ZIConnection** conn, const char* deviceSerial)

Sends an echo command to a device and blocks until answer is received.

This is useful to flush all buffers between API and device to enforce that further code is only executed after the device executed a previous command. Per device echo is only implemented for HF2. For other device types it is a synonym to ziAPISync, and deviceSerial parameter is ignored.

Parameters:

[in] conn

Pointer to the ZIConnection that is to be synchronized

[in] deviceSerial

The serial of the device to get the echo from, e.g., dev2100

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIGetValueS

```
__inline ZIResult_enum ziAPIGetValueS ( ZIConnection conn, char* path,  
DemodSample* value )
```

ziAPIGetValueDIO

```
__inline ZIResult_enum ziAPIGetValueDIO ( ZIConnection conn, char* path,  
DIOSample* value )
```

ziAPIGetValueAuxIn

```
__inline ZIResult_enum ziAPIGetValueAuxIn ( ZIConnection conn, char* path,  
AuxInSample* value )
```

6.2.4. Data Streaming

This section describes how to perform data streaming. It allows for recording at high data rates without sample loss.

Data Structures

- struct [ZIEvent](#)
This struct holds event data forwarded by the Data Server.
- struct [ziEvent](#)
This struct holds event data forwarded by the Data Server.
Deprecated: See [ZIEvent](#).

Functions

- [ZIEvent*](#) [ziAPIAllocateEventEx](#) ()
Allocates [ZIEvent](#) structure and returns the pointer to it.
Attention!!! It is the client code responsibility to deallocate the structure by calling [ziAPIDeallocateEventEx](#)!
- void [ziAPIDeallocateEventEx](#) ([ZIEvent*](#) ev)
Deallocates [ZIEvent](#) structure created with [ziAPIAllocateEventEx\(\)](#).
- [ZIResult_enum](#) [ziAPISubscribe](#) ([ZIConnection](#) conn, const char* path)
subscribes the nodes given by path for [ziAPIPollDataEx](#)
- [ZIResult_enum](#) [ziAPIUnSubscribe](#) ([ZIConnection](#) conn, const char* path)
unsubscribes to the nodes given by path
- [ZIResult_enum](#) [ziAPIPollDataEx](#) ([ZIConnection](#) conn, [ZIEvent*](#) ev, uint32_t timeOutMilliseconds)
checks if an event is available to read
- [ZIResult_enum](#) [ziAPIGetValueAsPollData](#) ([ZIConnection](#) conn, const char* path)
triggers a value request, which will be given back on the poll event queue
- [__inline ZIResult_enum](#) [ziAPIPollData](#) ([ZIConnection](#) conn, [ziEvent*](#) ev, int timeOut)
Checks if an event is available to read. Deprecated: See [ziAPIPollDataEx\(\)](#).

Detailed Description

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>
#include <stdlib.h>

#include "ziAPI.h"

void EventLoop(ZIConnection Conn) {
    ZIResult_enum RetVal;
```

```
char* ErrBuffer;

ZIEvent* Event;
unsigned int Cnt = 0;

/*
   Allocate ZIEvent in heap memory instead of getting it from stack will
   secure against stack overflows especially in windows.
*/
if ((Event = ziAPIAllocateEventEx()) == NULL) {
    fprintf(stderr, "Can't allocate memory\n");
    return;
}

// Subscribe to a node, e.g., a demodulator sample.
if ((RetVal = ziAPISubscribe(Conn, "/dev1024/demod/0/sample")) != ZI_INFO_SUCCESS)
{
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't subscribe: %s\n", ErrBuffer);

    ziAPIDeallocateEventEx(Event);

    return;
}

// loop 1000 times
while (Cnt < 1000) {
    // get all demod rates from all devices every 10th cycle
    if (++Cnt % 10 == 0) {
        if ((RetVal =
            ziAPIGetValueAsPollData(
                Conn, "/dev1046/demods/*/rate")) != ZI_INFO_SUCCESS) {
            ziAPIGetError(RetVal, &ErrBuffer, NULL);
            fprintf(stderr, "Error, can't get value as poll data: %s.\n",
                ErrBuffer);

            break;
        }
    }
}

// Poll data until no more data is available.
while (1) {
    if ((RetVal = ziAPIPollDataEx(
        Conn, Event, 0)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't poll data: %s.\n", ErrBuffer);

        break;
    } else {
        // The field Count of the Event struct is zero when no data has been
        // polled
        if (Event->valueType != ZI_VALUE_TYPE_NONE && Event->count > 0) {
            /*
               process the received event here
            */
        } else {
            // no more data is available so go on
            break;
        }
    }
}

if (ziAPIUnsubscribe(Conn, "") != ZI_INFO_SUCCESS) {
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't unsubscribe: %s.\n", ErrBuffer);
}
```

```
    ziAPIDeallocateEventEx(Event);  
}
```


Data Structure Documentation

struct ZIEvent

This struct holds event data forwarded by the Data Server.

```
#include "ziAPI.h"

typedef struct ZIEvent {
    uint32_t valueType;
    uint32_t count;
    uint8_t path[256];
    void* untyped;
    ZIDoubleData* doubleData;
    ZIDoubleDataTS* doubleDataTS;
    ZIIntegerData* integerData;
    ZIIntegerDataTS* integerDataTS;
    ZIByteArray* byteArray;
    ZIByteArrayTS* byteArrayTS;
    ZICntSample* cntSample;
    ZITreeChangeData* treeChangeData;
    TreeChange* treeChangeDataOld;
    ZIDemodSample* demodSample;
    ZIAuxInSample* auxInSample;
    ZIDIOSample* dioSample;
    ZIScopeWave* scopeWave;
    ZIScopeWaveEx* scopeWaveEx;
    ScopeWave* scopeWaveOld;
    ZIPWAWave* pwaWave;
    ZISweeperWave* sweeperWave;
    ZISpectrumWave* spectrumWave;
    ZIAdvisorWave* advisorWave;
    ZIAsyncReply* asyncReply;
    ZIVectorData* vectorData;
    ZIImpedanceSample* impedanceSample;
    uint64_t alignment;
    union ZIEvent::@6 value;
    uint8_t data[0x400000];
} ZIEvent;
```

Data Fields

- `uint32_t valueType`
Specifies the type of the data held by the [ZIEvent](#), see [ZIValueType_enum](#).
- `uint32_t count`
Number of values available in this event.
- `uint8_t path`
The path to the node from which the event originates.
- `void* untyped`
For convenience. The void field doesn't have a corresponding data type.
- `ZIDoubleData* doubleData`
when `valueType == ZI_VALUE_TYPE_DOUBLE_DATA`
- `ZIDoubleDataTS* doubleDataTS`
when `valueType == ZI_VALUE_TYPE_DOUBLE_DATA_TS`

- `ZIntegerData*` `integerData`
when `valueType == ZI_VALUE_TYPE_INTEGER_DATA`
- `ZIntegerDataTS*` `integerDataTS`
when `valueType == ZI_VALUE_TYPE_INTEGER_DATA_TS`
- `ZByteArray*` `byteArray`
when `valueType == ZI_VALUE_TYPE_BYTE_ARRAY`
- `ZByteArrayTS*` `byteArrayTS`
when `valueType == ZI_VALUE_TYPE_BYTE_ARRAY_TS`
- `ZCntSample*` `cntSample`
when `valueType == ZI_VALUE_TYPE_CNT_SAMPLE`
- `ZTreeChangeData*` `treeChangeData`
when `valueType == ZI_VALUE_TYPE_TREE_CHANGE_DATA`
- `TreeChange*` `treeChangeDataOld`
when `valueType == ZI_VALUE_TYPE_TREE_CHANGE_DATA_OLD`
- `ZIDemodSample*` `demodSample`
when `valueType == ZI_VALUE_TYPE_DEMOD_SAMPLE`
- `ZAuxInSample*` `auxInSample`
when `valueType == ZI_VALUE_TYPE_AUXIN_SAMPLE`
- `ZIDIOSample*` `dioSample`
when `valueType == ZI_VALUE_TYPE_DIO_SAMPLE`
- `ZIScopeWave*` `scopeWave`
when `valueType == ZI_VALUE_TYPE_SCOPE_WAVE`
- `ZIScopeWaveEx*` `scopeWaveEx`
when `valueType == ZI_VALUE_TYPE_SCOPE_WAVE_EX`
- `ScopeWave*` `scopeWaveOld`
when `valueType == ZI_VALUE_TYPE_SCOPE_WAVE_OLD`
- `ZIPWAWave*` `pwaWave`
when `valueType == ZI_VALUE_TYPE_PWA_WAVE`
- `ZISweeperWave*` `sweeperWave`
when `valueType == ZI_VALUE_TYPE_SWEEPER_WAVE`
- `ZISpectrumWave*` `spectrumWave`
when `valueType == ZI_VALUE_TYPE_SPECTRUM_WAVE`
- `ZIAdvisorWave*` `advisorWave`
when `valueType == ZI_VALUE_TYPE_ADVISOR_WAVE`
- `ZIAsyncReply*` `asyncReply`
when `valueType == ZI_VALUE_TYPE_ASYNC_REPLY`
- `ZIVectorData*` `vectorData`

- when valueType == ZI_VALUE_TYPE_VECTOR_DATA
- [ZIImpedanceSample](#)* impedanceSample
when valueType == ZI_VALUE_TYPE_IMPEDANCE_SAMPLE
- uint64_t alignment
ensure union size is 8 bytes
- union ZIEvent::@6 value
Convenience pointer to allow for access to the first entry in Data using the correct type according to [ZIEvent.valueType](#) field.
- uint8_t data
The raw value data.

Detailed Description

[ZIEvent](#) is used to give out events like value changes or errors to the user. Event handling functionality is provided by [ziAPISubscribe](#) and [ziAPIUnSubscribe](#) as well as [ziAPIPollDataEx](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

void ProcessEvent(ZIEvent* Event) {
    unsigned int j;

    switch (Event->valueType) {
    case ZI_VALUE_TYPE_DOUBLE_DATA:

        printf("%u elements of double data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", Event->value.doubleData[j]);

        break;

    case ZI_VALUE_TYPE_INTEGER_DATA:

        printf("%u elements of integer data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", (float)Event->value.integerData[j]);

        break;

    case ZI_VALUE_TYPE_DEMOD_SAMPLE:

        printf("%u elements of sample data %s\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("TS=%f, X=%f, Y=%f.\n",
                (float)Event->value.demodSample[j].timeStamp,
                Event->value.demodSample[j].x,
                Event->value.demodSample[j].y);
```

```
        break;

    case ZI_VALUE_TYPE_TREE_CHANGE_DATA:

        printf("%u elements of tree-changed data, %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++) {
            switch (Event->value.treeChangeDataOld[j].Action) {
                case ZI_TREE_ACTION_REMOVE:
                    printf("Tree removed: %s\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_ADD:
                    printf("treeChangeDataOld added: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_CHANGE:
                    printf("treeChangeDataOld changed: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;
            }
        }

        break;

    default:

        printf("Unexpected event value type: %d.\n", Event->valueType);
        break;
    }
}
```

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

struct `ziEvent`

This struct holds event data forwarded by the Data Server. Deprecated: See [ZIEvent](#).

```
#include "ziAPI.h"

typedef struct ziEvent {
    uint32_t Type;
    uint32_t Count;
    unsigned char Path[256];
    union ziEvent::Val Val;
    unsigned char Data[0x400000];
} ziEvent;
```

Data Structures

- union `ziEvent::Val`

Data Fields

- `uint32_t` Type
- `uint32_t` Count
- unsigned char Path
- union `ziEvent::Val` Val
- unsigned char Data

Detailed Description

[ziEvent](#) is used to give out events like value changes or errors to the user. Event handling functionality is provided by `ziAPISubscribe` and `ziAPIUnSubscribe` as well as `ziAPIPollDataEx`.

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

void ProcessEvent(ZIEvent* Event) {
    unsigned int j;

    switch (Event->valueType) {
    case ZI_VALUE_TYPE_DOUBLE_DATA:

        printf("%u elements of double data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", Event->value.doubleData[j]);
    }
```

```
        break;

    case ZI_VALUE_TYPE_INTEGER_DATA:

        printf("%u elements of integer data: %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", (float)Event->value.integerData[j]);

        break;

    case ZI_VALUE_TYPE_DEMOD_SAMPLE:

        printf("%u elements of sample data %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++)
            printf("TS=%f, X=%f, Y=%f.\n",
                   (float)Event->value.demodSample[j].timeStamp,
                   Event->value.demodSample[j].x,
                   Event->value.demodSample[j].y);

        break;

    case ZI_VALUE_TYPE_TREE_CHANGE_DATA:

        printf("%u elements of tree-changed data, %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++) {
            switch (Event->value.treeChangeDataOld[j].Action) {
                case ZI_TREE_ACTION_REMOVE:
                    printf("Tree removed: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_ADD:
                    printf("treeChangeDataOld added: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_CHANGE:
                    printf("treeChangeDataOld changed: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;
            }
        }

        break;

    default:

        printf("Unexpected event value type: %d.\n", Event->valueType);
        break;
    }
}
```

Data Structure Documentation

union ziEvent::Val

```
typedef union ziEvent::Val {  
    void* Void;  
    DemodSample* SampleDemod;  
    AuxInSample* SampleAuxIn;  
    DIOSample* SampleDIO;  
    ziDoubleType* Double;  
    ziIntegerType* Integer;  
    TreeChange* Tree;  
    ByteArrayData* ByteArray;  
    ScopeWave* Wave;  
    uint64_t alignment;  
} ziEvent::Val;
```

Data Fields

- void* Void
- DemodSample* SampleDemod
- AuxInSample* SampleAuxIn
- DIOSample* SampleDIO
- ziDoubleType* Double
- ziIntegerType* Integer
- TreeChange* Tree
- ByteArrayData* ByteArray
- ScopeWave* Wave
- uint64_t alignment

Function Documentation

ziAPIAllocateEventEx

ZIEvent* ziAPIAllocateEventEx ()

Allocates [ZIEvent](#) structure and returns the pointer to it. Attention!!! It is the client code responsibility to deallocate the structure by calling [ziAPIDeallocateEventEx](#)!

This function allocates a [ZIEvent](#) structure and returns the pointer to it. Free the memory using [ziAPIDeallocateEventEx](#).

See Also:

[ziAPIDeallocateEventEx](#)

ziAPIDeallocateEventEx

void ziAPIDeallocateEventEx ([ZIEvent*](#) ev)

Deallocates [ZIEvent](#) structure created with [ziAPIAllocateEventEx\(\)](#).

Parameters:

[in] ev
Pointer to [ZIEvent](#) structure to be deallocated..

See Also:

[ziAPIAllocateEventEx](#)

This function is the compliment to [ziAPIAllocateEventEx\(\)](#)

ziAPISubscribe

ZIResult_enum ziAPISubscribe (**ZIConnection** conn, const char* path)

subscribes the nodes given by path for [ziAPIPollDataEx](#)

This function subscribes to nodes so that whenever the value of the node changes the new value can be polled using [ziAPIPollDataEx](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one leaf can be subscribed to with one function call.

Parameters:

[in] conn

Pointer to the ziConnection for which to subscribe for

[in] path

Path to the nodes to subscribe

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPIUnSubscribe](#), [ziAPIPollDataEx](#), [ziAPIGetValueAsPollData](#)

ziAPIUnSubscribe

ZIResult_enum ziAPIUnSubscribe (**ZIConnection** conn, const char* path)

unsubscribes to the nodes given by path

This function is the complement to [ziAPISubscribe](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one node can be unsubscribed with one function call.

Parameters:

[in] conn

Pointer to the ziConnection for which to unsubscribe for

[in] path

Path to the Nodes to unsubscribe

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIPollDataEx](#), [ziAPIGetValueAsPollData](#)

ziAPIPollDataEx

ZIResult_enum ziAPIPollDataEx (**ZIConnection** conn, **ZIEvent*** ev, uint32_t
timeOutMilliseconds)

checks if an event is available to read

This function returns immediately if an event is pending. Otherwise it waits for an event for up to timeOutMilliseconds. All value changes that occur in nodes that have been subscribed to or in children of nodes that have been subscribed to are sent from the Data Server to the ziAPI session. For a description of how the data are available in the struct, refer to the documentation of struct [ziEvent](#). When no event was available within timeOutMilliseconds, the ziEvent::Type field will be ZI_DATA_NONE and the ziEvent::Count field will be zero. Otherwise these fields hold the values corresponding to the event that occurred.

Parameters:

[in] conn

Pointer to the [ZIConnection](#) for which events should be received

[out] ev

Pointer to a [ZIEvent](#) struct in which the received event will be written

[in] timeOutMilliseconds

Time to wait for an event in milliseconds. If -1 it will wait forever, if 0 the function returns immediately.

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIGetValueAsPollData](#), [ziEvent](#)

ziAPIGetValueAsPollData

ZIResult_enum ziAPIGetValueAsPollData ([ZIConnection](#) conn, const char* path)

triggers a value request, which will be given back on the poll event queue

Use this function to receive the value of one or more nodes as one or more events using [ziAPIPollDataEx](#), even when the node is not subscribed or no value change has occurred.

Parameters:

[in] conn

Pointer to the [ZIConnection](#) with which the value should be retrieved

[in] path

Path to the Node holding the value

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

ziAPIPollData

`__inline ZIResult_enum ziAPIPollData (ZIConnection conn, ziEvent* ev, int timeOut)`

Checks if an event is available to read. Deprecated: See [ziAPIPollDataEx\(\)](#).

Parameters:

[in] conn

Pointer to the [ZIConnection](#) for which events should be received

[out] ev

Pointer to a [ziEvent](#) struct in which the received event will be written

[in] timeOut

Time to wait for an event in milliseconds. If -1 it will wait forever, if 0 the function returns immediately.

Returns:

- ZI_SUCCESS On success.
- ZI_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_OVERFLOW When a FIFO overflow occurred.

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIGetValueAsPollData](#), [ziEvent](#)

6.2.5. API for fast asynchronous operation

Functions in this group are non-blocking, and on return only report errors that can be identified directly on a client side (e.g. not connected). Any further results (including errors like node not found) of the command processing is returned as a special event in poll data. Tags are used to match the asynchronous replies with the sent commands.

Functions

- [ZIResult_enum](#) `ziAPIAsyncSetDoubleData (ZIConnection conn, const char* path, ZIDoubleData value)`
- [ZIResult_enum](#) `ziAPIAsyncSetIntegerData (ZIConnection conn, const char* path, ZIIntegerData value)`
- [ZIResult_enum](#) `ziAPIAsyncSetByteArray (ZIConnection conn, const char* path, uint8_t* buffer, uint32_t length)`
- [ZIResult_enum](#) `ziAPIAsyncSubscribe (ZIConnection conn, const char* path, ZIAsyncTag tag)`
- [ZIResult_enum](#) `ziAPIAsyncUnSubscribe (ZIConnection conn, const char* path, ZIAsyncTag tag)`
- [ZIResult_enum](#) `ziAPIAsyncGetValueAsPollData (ZIConnection conn, const char* path, ZIAsyncTag tag)`

Function Documentation

ziAPIAsyncSetDoubleData

[ZIResult_enum](#) ziAPIAsyncSetDoubleData ([ZIConnection](#) conn, const char* path, ZIDoubleData value)

ziAPIAsyncSetIntegerData

[ZIResult_enum](#) ziAPIAsyncSetIntegerData ([ZIConnection](#) conn, const char* path, ZIIntegerData value)

ziAPIAsyncSetByteArray

ZIResult_enum ziAPIAsyncSetByteArray (**ZIConnection** conn, const char* path, uint8_t* buffer, uint32_t length)

ziAPIAsyncSubscribe

[ZIResult_enum](#) ziAPIAsyncSubscribe ([ZIConnection](#) conn, const char* path, ZIAsyncTag tag)

ziAPIAsyncUnSubscribe

[ZIResult_enum](#) ziAPIAsyncUnSubscribe ([ZIConnection](#) conn, const char* path, ZIAsyncTag tag)

ziAPIAsyncGetValueAsPollData

[ZIResult_enum](#) ziAPIAsyncGetValueAsPollData ([ZIConnection](#) conn, const char* path, ZIAsyncTag tag)

6.2.6. Error Handling and Logging in the LabOne C API

This section describes how to get more information when an error occurs.

Functions

- `ZIResult_enum` `ziAPIGetError (ZIResult_enum result, char** buffer, int* base)`
Returns a description and the severity for a `ZIResult_enum`.
- `ZIResult_enum` `ziAPIGetLastError (ZIConnection conn, char* buffer, uint32_t bufferSize)`
Returns the message from the last error that occurred.
- `void` `ziAPISetDebugLevel (int32_t debugLevel)`
Enable ziAPI's log and set the severity level of entries to be included in the log.
- `void` `ziAPIWriteDebugLog (int32_t debugLevel, const char* message)`
Write a message to ziAPI's log with the specified severity.

Detailed Description

In general, two types of errors can occur when using ziAPI. The two types are distinguished by the origin of the error: Whether it occurred within ziAPI itself or whether it occurred internally in the Zurich Instruments Core library.

All ziAPI functions (apart from a very few exceptions) return an exit code `ZIResult_enum`, which will be non-zero if the function call was not entirely successful. If the error originated in ziAPI itself, the exit code describes precisely the type of error that occurred (in other words, the exit code is not `ZI_ERROR_GENERAL`). In this case the error message corresponding to the exit code can be obtained with the function `ziAPIGetError`.

However, if the error has occurred internally, the exit code will be `ZI_ERROR_GENERAL`. In this case, the exit code does not describe the type of error precisely, instead a detailed error message is available to the user which can be obtained with the function `ziAPIGetLastError`. The function `ziAPIGetLastError` may be used with any function that takes a `ZIConnection` as an input argument (with the exception of `ziAPIInit`, `ziAPIDestroy`, `ziAPIConnect`, `ziAPIConnectEx`) and is the recommended function to use, if applicable, otherwise `ziAPIGetError` should be used.

The function `ziAPIGetLastError` was introduced in LabOne 15.11 due to the availability of `ziCoreModules` in ziAPI - its not desirable in general to map every possible error to an exit code in ziAPI; what is more relevant is the associated error message.

In addition to these two functions, ziAPI's log can be very helpful whilst debugging ziAPI-based programs. The log is not enabled by default; it's enabled by specifying a logging level with `ziAPISetDebugLevel`.

Function Documentation

ziAPIGetError

ZIResult_enum ziAPIGetError (**ZIResult_enum** result, char** buffer, int* base)

Returns a description and the severity for a **ZIResult_enum**.

This function returns a static char pointer to a description string for the given **ZIResult_enum** error code. It also provides a parameter returning the severity (info, warning, error). If the given error code does not exist a description for an unknown error and the base for an error will be returned. If a description or the base is not needed NULL may be passed. In general, it's recommended to use **ziAPIGetLastError** instead to get detailed error messages.

Parameters:

[in] result

A **ZIResult_enum** for which the description or base will be returned

[out] buffer

A pointer to a char array to return the description. May be NULL if no description is needed.

[out] base

The severity for the provided Status parameter:

- ZI_INFO_BASE For infos.
- ZI_WARNING_BASE For warnings.
- ZI_ERROR_BASE For errors.

Returns:

- ZI_INFO_SUCCESS Upon success.

ziAPIGetLastError

ZIResult_enum ziAPIGetLastError (**ZIConnection** conn, char* buffer, uint32_t bufferSize)

Returns the message from the last error that occurred.

This function can be used to obtain the error message from the last error that occurred associated with the provided **ZIConnection**. If the last ziAPI call is successful, then the last error message returned by ziAPIGetError is empty. Only ziAPI function calls that take **ZIConnection** as an input argument influence the message returned by ziAPIGetLastError, if they do not take **ZIConnection** as an input argument the last error message will neither be reset to be empty or set to an error message (in the case of the error). There are some exceptions to this rule, ziAPIGetLastError can also not be used with **ziAPIInit**, **ziAPIConnect**, **ziAPIConnectEx** and **ziAPIDestroy**. Note, a call to ziAPIGetLastError will also reset the last error message to empty if its call was successful. Since the buffer is left unchanged in the case of an error occurring in the call to ziAPIGetLastError it is safest to initialize the buffer with a known value, for example, "ziAPIGetLastError was not successful".

Parameters:

[in] conn

The **ZIConnection** from which to get the error message.

[out] buffer

A pointer to a char array to return the message.

[in] bufferSize

The length of the provided buffer.

Returns:

- **ZI_INFO_SUCCESS** Upon success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred. In this case the provided buffer is left unchanged.
- **ZI_ERROR_LENGTH** If the message's length exceeds the provided bufferSize, the message is truncated and written to buffer.

ziAPISetDebugLevel

void ziAPISetDebugLevel (int32_t debugLevel)

Enable ziAPI's log and set the severity level of entries to be included in the log.

Calling this function enables ziAPI's log at the specified severity level. On Windows the logs can be found by navigating to the Zurich Instruments "Logs" folder entry in the Windows Start Menu: Programs -> Zurich Instruments -> LabOne Servers -> Logs. This will open an Explorer window displaying folders containing log files from various LabOne components, in particular, the `ziAPILog` folder contains logs from ziAPI. On Linux, the logs can be found at `/tmp/ziAPILog_USERNAME`, where `USERNAME` is the same as the output of the `whoami` command.

Parameters:

[in] `debugLevel`

An integer specifying the log's severity level:

- `trace: 0,`
- `info: 1,`
- `debug: 2,`
- `warning: 3,`
- `error: 4,`
- `fatal: 5,`
- `status: 6.`

See Also:

[ziAPIWriteDebugLog](#)

ziAPIWriteDebugLog

void ziAPIWriteDebugLog (int32_t debugLevel, const char* message)

Write a message to ziAPI's log with the specified severity.

This function may be used to write a message to ziAPI's log from client code to assist with debugging. Note, this function is only available if the implementation used in [ziAPIConnectEx](#) is "ziAPI_Core" (the default implementation). Also logging must be first enabled using [ziAPISetDebugLevel](#).

Parameters:

[in] debugLevel

An integer specifying the severity of the message to write in the log:

- trace: 0,
- info: 1,
- debug: 2,
- warning: 3,
- error: 4,
- fatal: 5,
- status: 6.

[in] message

A character array comprising of the message to be written.

See Also:

[ziAPISetDebugLevel](#)

6.2.7. Using ziCore Modules in the LabOne C API

This sections describes ziAPI's interface for working with ziCore Modules. Modules provide a high-level interface for performing common measurement tasks such as sweeping data (Sweeper Module) or recording bursts of when certain trigger criteria have been fulfilled (Software Trigger Module). For an introduction to working with Modules please see the "ziCore Modules" section in the LabOne Programming Manual: .

Data Structures

- struct [ZISWTriggerHeader](#)
Structure to hold information about data returned from the SW Trigger Module.
- struct [ZIModuleHeaderSweeper](#)
Structure to hold information about data returned from the Sweep Module.
- struct [ZIModuleHeader](#)
Module-specific event header.
- struct [ZIModuleEvent](#)
This struct holds data of a single chunk from module lookup.

Typedefs

- typedef [ZIModuleEventPtr](#)
The pointer to a Module's data chunk to read out, updated via [ziAPIModGetChunk](#).

Enumerations

- enum [ZIModuleHeaderType_enum](#)
{ [ZI_MODULE_HEADER_TYPE_NONE](#),
[ZI_MODULE_HEADER_TYPE_SWTRIGGER](#),
[ZI_MODULE_HEADER_TYPE_SWEEPER](#) }
Enumerates all module header types.

Functions

- [ZIResult_enum](#) [ziAPIModCreate](#) ([ZIConnection](#) conn, [ZIModuleHandle](#)* handle, const char* moduleId)
Create a [ZIModuleHandle](#) that can be used for asynchronous measurement tasks.
- [ZIResult_enum](#) [ziAPIModSetDoubleData](#) ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, const char* path, [ZIDoubleData](#) value)
Sets a module parameter to the specified double type.
- [ZIResult_enum](#) [ziAPIModSetIntegerData](#) ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, const char* path, [ZIIntegerData](#) value)
Sets a module parameter to the specified integer type.

- `ZIResult_enum` `ziAPIModSetByteArray (ZIConnection conn, ZIModuleHandle handle, const char* path, uint8_t* buffer, uint32_t length)`
Sets a module parameter to the specified byte array.
- `ZIResult_enum` `ziAPIModListNodes (ZIConnection conn, ZIModuleHandle handle, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)`
Returns all child parameter node paths found under the specified parent module parameter path.
- `ZIResult_enum` `ziAPIModSubscribe (ZIConnection conn, ZIModuleHandle handle, const char* path)`
Subscribes to the nodes specified by path, these nodes will be recorded during module execution.
- `ZIResult_enum` `ziAPIModUnSubscribe (ZIConnection conn, ZIModuleHandle handle, const char* path)`
Unsubscribes to the nodes specified by path.
- `ZIResult_enum` `ziAPIModExecute (ZIConnection conn, ZIModuleHandle handle)`
Starts the module's thread and its associated measurement task.
- `ZIResult_enum` `ziAPIModTrigger (ZIConnection conn, ZIModuleHandle handle)`
Manually issue a trigger forcing data recording (SW Trigger Module only).
- `ZIResult_enum` `ziAPIModProgress (ZIConnection conn, ZIModuleHandle handle, ZIDoubleData* progress)`
Queries the current state of progress of the module's measurement task.
- `ZIResult_enum` `ziAPIModFinished (ZIConnection conn, ZIModuleHandle handle, ZIIntegerData* finished)`
Queries whether the module has finished its measurement task.
- `ZIResult_enum` `ziAPIModFinish (ZIConnection conn, ZIModuleHandle handle)`
Stops the module performing its measurement task.
- `ZIResult_enum` `ziAPIModSave (ZIConnection conn, ZIModuleHandle handle, const char* fileName)`
Saves the currently accumulated data to file.
- `ZIResult_enum` `ziAPIModRead (ZIConnection conn, ZIModuleHandle handle, const char* path)`
Make the currently accumulated data available for use in the C program.
- `ZIResult_enum` `ziAPIModNextNode (ZIConnection conn, ZIModuleHandle handle, char* path, uint32_t bufferSize, ZIValueType_enum* valueType, uint64_t* chunks)`

Make the data for the next node available for reading with [ziAPIModGetChunk](#).

- [ZIResult_enum](#) [ziAPIModGetChunk](#) ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, [uint64_t](#) chunkIndex, [ZIModuleEventPtr*](#) ev)
Get the specified data chunk from the current node.
- [ZIResult_enum](#) [ziAPIModEventDeallocate](#) ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, [ZIModuleEventPtr](#) ev)
Deallocate the [ZIModuleEventPtr](#) being used by the module.
- [ZIResult_enum](#) [ziAPIModClear](#) ([ZIConnection](#) conn, [ZIModuleHandle](#) handle)
Terminates the module's thread and destroys the module.

Data Structure Documentation

struct ZISWTriggerHeader

Structure to hold information about data returned from the SW Trigger Module.

```
#include "ziAPI.h"

typedef struct ZISWTriggerHeader {
    ZITimeStamp triggerStart;
    uint64_t triggerNumber;
    uint32_t cols;
    uint32_t rows;
    uint32_t repetitions;
    uint32_t flags;
    uint8_t reserved0[32];
} ZISWTriggerHeader;
```

Data Fields

- ZITimeStamp triggerStart
Trigger timestamp.
- uint64_t triggerNumber
Trigger counter since execution start.
- uint32_t cols
Number of columns.
- uint32_t rows
Number of rows.
- uint32_t repetitions
Number of repetitions in grid mode.
- uint32_t flags
Flags Bit (0): Finished (all repetitions recorded)
- uint8_t reserved0
Reserved space for future use.

struct ZIModuleHeaderSweeper

Structure to hold information about data returned from the Sweep Module.

```
#include "ziAPI.h"

typedef struct ZIModuleHeaderSweeper {
    char traceName[256];
} ZIModuleHeaderSweeper;
```

Data Fields

- char traceName

struct ZIModuleHeader

Module-specific event header.

```
#include "ziAPI.h"

typedef struct ZIModuleHeader {
    ZIModuleHeaderType_enum
        type;
    void* untyped;
    ZISWTriggerHeader* swTrigger;
    ZISweeperHeader* sweeper;
    union ZIModuleHeader::@7 ptr;
} ZIModuleHeader;
```

Data Fields

- `ZIModuleHeaderType_enum` type
- `void*` `untyped`
- `ZISWTriggerHeader*` `swTrigger`
- `ZISweeperHeader*` `sweeper`
- `union ZIModuleHeader::@7` `ptr`

struct ZIModuleEvent

This struct holds data of a single chunk from module lookup.

```
#include "ziAPI.h"

typedef struct ZIModuleEvent {
    uint64_t allocatedSize;
    ZIModuleHeader header;

    ZIEvent value[0];
} ZIModuleEvent;
```

Data Fields

- `uint64_t allocatedSize`
For internal use - never modify!
- `ZIModuleHeader header`
Module-specific event header.
- `ZIEvent value`
Defines location of stored `ZIEvent`.

Enumeration Type Documentation

enum ZIModuleHeaderType_enum

Enumerates all module header types.

Enumerator:

- ZI_MODULE_HEADER_TYPE_NONE
- ZI_MODULE_HEADER_TYPE_SWTRIGGER
- ZI_MODULE_HEADER_TYPE_SWEEPER

Function Documentation

ziAPIModCreate

ZIResult_enum ziAPIModCreate (**ZIConnection** conn, **ZIModuleHandle*** handle, const char* moduleId)

Create a **ZIModuleHandle** that can be used for asynchronous measurement tasks.

This function initializes a ziCore module and provides a pointer (handle) with which to access and work with it. Note that this function does not start the module's thread. Before the thread can be started (with **ziAPIModExecute**):

- the device serial (e.g., "dev100") to be used with module must be specified via **ziAPIModSetByteArray**.
- the desired data (node paths) to record during the measurement must be specified via **ziAPIModSubscribe**. The module's thread is stopped with **ziAPIModClear**.

Parameters:

[in] conn

The **ZIConnection** which should be used to initialize the module.

[out] handle

Pointer to the initialized **ZIModuleHandle**, which from then on can be used to reference the module.

[in] moduleId

The name specifying the type the module to create (only the following ziCore Modules are currently supported in ziAPI):

- "sweep" to initialize an instance of the Sweeper Module.
- "record" to initialize an instance of the Software Trigger (Recorder) Module.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** when the connection is invalid (not connected) or when a communication error occurred.
- **ZI_WARNING_NOTFOUND** if the provided moduleId was invalid.
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

ziAPIModExecute, **ziAPIModClear**

ziAPIModSetDoubleData

ZIResult_enum ziAPIModSetDoubleData (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, **ZIDoubleData** value)

Sets a module parameter to the specified double type.

This function is used to configure (set) module parameters which have double types.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to set data on.

[in] path

Path to the module parameter path.

[in] value

The double data to write to the path.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSetIntegerData](#), [ziAPIModSetByteArray](#)

ziAPIModSetIntegerData

ZIResult_enum ziAPIModSetIntegerData (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, **ZIntegerData** value)

Sets a module parameter to the specified integer type.

This function is used to configure (set) module parameters which have integer types.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to set data on.

[in] path

Path to the module parameter path.

[in] value

The integer data to write to the path.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_GENERAL** If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSetDoubleData](#), [ziAPIModSetByteArray](#)

ziAPIModSetByteArray

ZIResult_enum ziAPIModSetByteArray (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, uint8_t* buffer, uint32_t length)

Sets a module parameter to the specified byte array.

This function is used to configure (set) module parameters which have byte array types.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to set data on.

[in] path

Path to the module parameter path.

[in] buffer

Pointer to the byte array with the data.

[in] length

Length of the data in the buffer.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSetDoubleData](#), [ziAPIModSetIntegerData](#)

ziAPIModListNodes

ZIResult_enum ziAPIModListNodes (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)

Returns all child parameter node paths found under the specified parent module parameter path.

This function returns a list of parameter names found at the specified path. The path may contain wildcards. The list is returned in a null-terminated char-buffer, each element delimited by a newline. If the maximum length of the buffer (bufferSize) is not sufficient for all elements, nothing will be returned and the return value will be **ZI_ERROR_LENGTH**. Note, the provided path must match the module being addressed, i.e., path must exactly start with "sweep/" for the Sweeper Module.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** from which the parameter names should be retrieved.

[in] path

Path for which all children will be returned. The path may contain wildcard characters.

[out] nodes

Upon call filled with newline-delimited list of the names of all the children found. The string is zero-terminated.

[in] bufferSize

The length of the buffer specified as the nodes output parameter.

[in] flags

A combination of flags (applied bitwise) as defined in **ZIListNodes_enum**.

Returns:

- **ZI_INFO_SUCCESS** On success
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_LENGTH** If the path's length exceeds **MAX_PATH_LEN** or the length of the char-buffer for the nodes given by bufferSize is too small for all elements.
- **ZI_WARNING_OVERFLOW** When a FIFO overflow occurred.
- **ZI_ERROR_COMMAND** On an incorrect answer of the server.
- **ZI_ERROR_SERVER_INTERNAL** If an internal error occurred in Data Server.
- **ZI_WARNING_NOTFOUND** If the given path could not be resolved.
- **ZI_ERROR_TIMEOUT** When communication timed out.

- `ZI_ERROR_GENERAL` If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIModSubscribe

ZIResult_enum ziAPIModSubscribe (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path)

Subscribes to the nodes specified by path, these nodes will be recorded during module execution.

This function subscribes to nodes so that whenever the value of the node changes while the module is executing the new value will be accumulated and then read using [ziAPIModRead](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one leaf can be subscribed to with one function call.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module in which the nodes should be subscribed to.

[in] path

Path specifying the nodes to subscribe to, may contain wildcards.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or a general error occurred, enable ziAPI's log for detailed information, see [ziAPISetDebugLevel](#).
- ZI_ERROR_LENGTH If the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW When a FIFO overflow occurred.
- ZI_ERROR_COMMAND On an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL If an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND If the given path could not be resolved or no node given by path is able to hold values.
- ZI_ERROR_TIMEOUT When communication timed out.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModUnSubscribe](#), [ziAPIModRead](#)

ziAPIModUnSubscribe

ZIResult_enum ziAPIModUnSubscribe (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path)

Unsubscribes to the nodes specified by path.

This function is the complement to [ziAPIModSubscribe](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one node can be unsubscribed with one function call.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module in which the nodes should be unsubscribed from.

[in] path

Path specifying the nodes to unsubscribe from, may contain wildcards.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH If the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW When a FIFO overflow occurred.
- ZI_ERROR_COMMAND On an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL If an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND If the given path could not be resolved or no node given by path is able to hold values.
- ZI_ERROR_TIMEOUT When communication timed out.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSubscribe](#), [ziAPIModRead](#)

ziAPIModExecute

ZIResult_enum ziAPIModExecute ([ZIConnection](#) conn, [ZIModuleHandle](#) handle)

Starts the module's thread and its associated measurement task.

Once the module's parameters has been configured as required via, e.g. [ziAPIModSetDoubleData](#), this function starts the module's thread. This starts the module's main measurement task which will run asynchronously. The thread will run until either the module has completed its task or until [ziAPIModFinish](#) is called. Subscription or unsubscription is not possible while the module is executing. The status of the module can be obtained with either [ziAPIModFinished](#) or [ziAPIModProgress](#).

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModCreate](#), [ziAPIModProgress](#), [ziAPIModFinish](#)

ziAPIModTrigger

ZIResult_enum ziAPIModTrigger ([ZIConnection](#) conn, [ZIModuleHandle](#) handle)

Manually issue a trigger forcing data recording (SW Trigger Module only).

This function is used with the Software Trigger Module in order to manually issue a trigger in order to force recording of data. A burst of subscribed data will be recorded as configured via the SW Trigger's parameters as would a regular trigger event.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIModProgress

ZIResult_enum ziAPIModProgress (**ZIConnection** conn, **ZIModuleHandle** handle, **ZIDoubleData*** progress)

Queries the current state of progress of the module's measurement task.

This function can be used to query the module's progress in performing its current measurement task, the progress is returned as a double in [0, 1], where 1 indicates task completion.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to execute.

[out] progress

A pointer to **ZIDoubleData** indicating the current progress of the module.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_GENERAL** If a general error occurred, use **ziAPIGetLastError** for a detailed error message.
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

ziAPIModExecute, **ziAPIModFinish**, **ziAPIModFinished**

ziAPIModFinished

ZIResult_enum ziAPIModFinished (**ZIConnection** conn, **ZIModuleHandle** handle, **ZIIntegerData*** finished)

Queries whether the module has finished its measurement task.

This function can be used to query whether the module has finished its task or not.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to execute.

[out] finished

A pointer to **ZIIntegerData**, upon return this will be 0 if the module is still executing or 1 if has finished executing.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_GENERAL** If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModExecute](#), [ziAPIModFinish](#), [ziAPIModProgress](#)

ziAPIModFinish

ZIResult_enum ziAPIModFinish (**ZIConnection** conn, **ZIModuleHandle** handle)

Stops the module performing its measurement task.

This function stops the module performing its associated measurement task and stops recording any data. The task and data recording may be restarted by calling [ziAPIModExecute](#) again.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModProgress](#), [ziAPIModFinished](#)

ziAPIModSave

ZIResult_enum ziAPIModSave ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, const char* fileName)

Saves the currently accumulated data to file.

This function saves the currently accumulated data to a file. The path of the file to save data to is specified via the module's directory parameter.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

[in] fileName

The basename of the file to save the data in.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModExecute](#), [ziAPIModFinish](#), [ziAPIModFinished](#)

ziAPIModRead

ZIResult_enum ziAPIModRead (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path)

Make the currently accumulated data available for use in the C program.

This function can be used to either read (get) module parameters, in this case a path that addresses the module must be specified, or it can be used to read out the currently accumulated data from subscribed nodes in the module. In either case the actual data must then be accessed by the user using [ziAPIModNextNode](#) and [ziAPIModGetChunk](#).

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

[in] path

The path specifying the module parameter(s) to get, specify NULL to obtain all subscribed data.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModGetChunk](#), [ziAPIModNextNode](#)

ziAPIModNextNode

ZIResult_enum ziAPIModNextNode (**ZIConnection** conn, **ZIModuleHandle** handle, char* path, uint32_t bufferSize, **ZIValueType_enum*** valueType, uint64_t* chunks)

Make the data for the next node available for reading with [ziAPIModGetChunk](#).

After calling [ziAPIModRead](#), subscribed data (or module parameters) may now be read out on a node-by-node and chunk-by-chunk basis. All nodes with data available in the module can be iterated over by using [ziAPIModNextNode](#), then for each node the chunks of data available are read out using [ziAPIModGetChunk](#). Calling this function makes the data from the next node available for read.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

[out] path

A string specifying the node's path whose data chunk points to.

[in] bufferSize

The length of the buffer specified as the path output parameter.

[out] valueType

The [ZIValueType_enum](#) of the node's data.

[out] chunks

The number of chunks of data available for the node.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModRead](#), [ziAPIModGetChunk](#), [ziAPIModEventDeallocate](#)

ziAPIModGetChunk

ZIResult_enum ziAPIModGetChunk (**ZIConnection** conn, **ZIModuleHandle** handle, **uint64_t** chunkIndex, **ZIModuleEventPtr*** ev)

Get the specified data chunk from the current node.

Data is read out node-by-node and then chunk-by-chunk. This function can be used to obtain specific data chunks from the current node that data is being read from. More precisely, it preallocates space for an event structure big enough to hold the node's data at the specified chunk index, updates **ZIModuleEventPtr** to point to this space and then copies the chunk data to this space.

Note, before the very first call to ziAPIModGetChunk, the **ZIModuleEventPtr** should be initialized to NULL and then left untouched for all subsequent calls (even after calling **ziAPIModNextNode** to get data from the next node). This is because ziAPIModGetChunk internally manages the required space allocation for the event and then in subsequent calls only reallocates space when it is required. It is optimized to reduce the number of required space reallocations for the event.

The **ZIModuleEventPtr** should be deallocated using **ziAPIModEventDeallocate**, otherwise the lifetime of the **ZIModuleEventPtr** is the same as the lifetime of the module. Indeed, the same **ZIModuleEventPtr** can be used, even for subsequent reads. It is also possible to work with multiple **ZIModuleEventPtr** so that some pointers can be kept for later processing.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to execute.

[out] chunkIndex

The index of the data chunk to update the pointer to.

[out] ev

The module's **ZIModuleEventPtr** that points to the currently available data chunk.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_GENERAL** If a general error occurred, use **ziAPIGetLastError** for a detailed error message.
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

[ziAPIModRead](#), [ziAPIModNextNode](#), [ziAPIModEventDeallocate](#)

ziAPIModEventDeallocate

ZIResult_enum ziAPIModEventDeallocate ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, [ZIModuleEventPtr](#) ev)

Deallocate the [ZIModuleEventPtr](#) being used by the module.

This function deallocates the [ZIModuleEventPtr](#). Since a module event's allocated space is managed internally by [ziAPIModGetChunk](#), when the user no longer requires the event (all data has been read out) it must be deallocated by the user with this function.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

[in] ev

The [ZIModuleEventPtr](#) to deallocate.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModGetChunk](#), [ziAPIModRead](#)

ziAPIModClear

ZIResult_enum ziAPIModClear ([ZIConnection](#) conn, [ZIModuleHandle](#) handle)

Terminates the module's thread and destroys the module.

This function terminates the module's thread. After calling ziAPIModClear the module's handle may not be used any more. A new instance of the module must be initialized if required.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModExecute](#), [ziAPIModFinish](#)

6.2.8. Vector Write

Functions for working with vector data writing.

Enumerations

- enum [ZIVectorWriteStatus_enum](#)
{ ZI_VECTOR_WRITE_STATUS_IDLE,
ZI_VECTOR_WRITE_STATUS_PENDING }

Functions

- [ZIResult_enum](#) [ziAPIVectorWriteBlock](#) ([ZIConnection](#) conn,
const char* path, [ZIVectorData*](#) vectorBlock)
- [ZIResult_enum](#) [ziAPIVectorWriteGetStatus](#) ([ZIConnection](#)
conn, const char* path, uint8_t* status)
status - see [ZIVectorWriteStatus_enum](#)
- [ZIResult_enum](#) [ziAPIVectorWrite](#) ([ZIConnection](#) conn, const
char* path, const void* vectorPtr, uint8_t vectorElementType,
uint64_t vectorSizeElements)
vectorElementType - see [ZIVectorElementType_enum](#)

Enumeration Type Documentation

enum ZIVectorWriteStatus_enum

Enumerator:

- ZI_VECTOR_WRITE_STATUS_IDLE
- ZI_VECTOR_WRITE_STATUS_PENDING

Function Documentation

ziAPIVectorWriteBlock

[ZIResult_enum](#) ziAPIVectorWriteBlock ([ZIConnection](#) conn, const char* path, [ZIVectorData*](#) vectorBlock)

ziAPIVectorWriteGetStatus

ZIResult_enum ziAPIVectorWriteGetStatus (**ZIConnection** conn, const char* path, uint8_t* status)

status - see ZIVectorWriteStatus_enum

ziAPIVectorWrite

ZIResult_enum ziAPIVectorWrite (**ZIConnection** conn, const char* path, const void* vectorPtr, uint8_t vectorElementType, uint64_t vectorSizeElements)

vectorElementType - see ZIVectorElementType_enum

6.2.9. Device discovery

Functions for working with device discovery.

Functions

- `ZIResult_enum` `ziAPIDiscoveryFind (ZIConnection conn, const char* deviceAddress, const char** deviceId)`
Returns the device id for a given device address.
Attention! Invalidates all pointers previously returned by `ziAPIDiscovery*` calls.
- `ZIResult_enum` `ziAPIDiscoveryGet (ZIConnection conn, const char* deviceId, const char** propsJSON)`
Returns the device properties for a given device id in JSON format.
- `ZIResult_enum` `ziAPIDiscoveryGetValueI (ZIConnection conn, const char* deviceId, const char* propName, ZIIntegerData* value)`
Returns given integer property value for a given device id.
- `ZIResult_enum` `ziAPIDiscoveryGetValueS (ZIConnection conn, const char* deviceId, const char* propName, const char** value)`
Returns given string property value for a given device id.

Function Documentation

ziAPIDiscoveryFind

ZIResult_enum ziAPIDiscoveryFind (**ZIConnection** conn, const char* deviceAddress, const char** deviceId)

Returns the device id for a given device address. Attention! Invalidates all pointers previously returned by ziAPIDiscovery* calls.

ziAPIDiscoveryGet

ZIResult_enum ziAPIDiscoveryGet (**ZIConnection** conn, const char* deviceId, const char** propsJSON)

Returns the device properties for a given device id in JSON format.

ziAPIDiscoveryGetValueI

ZIResult_enum ziAPIDiscoveryGetValueI (**ZIConnection** conn, const char* deviceId, const char* propName, ZIIntegerData* value)

Returns given integer property value for a given device id.

ziAPIDiscoveryGetValueS

ZIResult_enum ziAPIDiscoveryGetValueS (**ZIConnection** conn, const char* deviceId, const char* propName, const char** value)

Returns given string property value for a given device id.

6.3. Data Structure Documentation

6.3.1. struct AuxInSample

The [AuxInSample](#) struct holds data for the ZI_DATA_AUXINSAMPLE data type. Deprecated: See [ZIAuxInSample](#).

```
#include "ziAPI.h"

typedef struct AuxInSample {
    ziTimeStampType TimeStamp;
    double Ch0;
    double Ch1;
} AuxInSample;
```

Data Fields

- `ziTimeStampType TimeStamp`
- `double Ch0`
- `double Ch1`

6.3.2. struct ByteArrayData

The [ByteArrayData](#) struct holds data for the ZI_DATA_BYTEARRAY data type. Deprecated: See [ZIByteArray](#).

```
#include "ziAPI.h"

typedef struct ByteArrayData {
    unsigned int Len;
    unsigned char Bytes[0];
} ByteArrayData;
```

Data Fields

- unsigned int Len
- unsigned char Bytes

6.3.3. struct DemodSample

The [DemodSample](#) struct holds data for the ZI_DATA_DEMODSAMPLE data type. Deprecated: See [ZIDemodSample](#).

```
#include "ziAPI.h"

typedef struct DemodSample {
    ziTimeStampType TimeStamp;
    double X;
    double Y;
    double Frequency;
    double Phase;
    unsigned int DIOBits;
    unsigned int Reserved;
    double AuxIn0;
    double AuxIn1;
} DemodSample;
```

Data Fields

- `ziTimeStampType TimeStamp`
- `double X`
- `double Y`
- `double Frequency`
- `double Phase`
- `unsigned int DIOBits`
- `unsigned int Reserved`
- `double AuxIn0`
- `double AuxIn1`

6.3.4. struct DIOSample

The [DIOSample](#) struct holds data for the ZI_DATA_DIOSAMPLE data type. Deprecated: See [ZIDIOSample](#).

```
#include "ziAPI.h"

typedef struct DIOSample {
    ziTimeStampType TimeStamp;
    unsigned int Bits;
    unsigned int Reserved;
} DIOSample;
```

Data Fields

- `ziTimeStampType TimeStamp`
- `unsigned int Bits`
- `unsigned int Reserved`

6.3.5. struct ScopeWave

The structure used to hold a single scope shot (API Level 1). If the client is connected to the Data Server using API Level 4 (recommended if supported by your device class) please see [ZIScopeWave](#) instead ([ZIScopeWaveEx](#) for API Level 5 and above).

```
#include "ziAPI.h"

typedef struct ScopeWave {
    double dt;
    uint32_t ScopeChannel;
    uint32_t TriggerChannel;
    uint32_t BWLimit;
    uint32_t Count;
    int16_t Data[0];
} ScopeWave;
```

Data Fields

- double dt
Time difference between samples.
- uint32_t ScopeChannel
Scope channel of the represented data.
- uint32_t TriggerChannel
Trigger channel of the represented data.
- uint32_t BWLimit
Bandwidth-limit flag.
- uint32_t Count
Count of samples.
- int16_t Data
First wave data.

6.3.6. struct TreeChange

The structure used to hold info about added or removed nodes. This is the version without timestamp used in API v1 compatibility mode.

```
#include "ziAPI.h"

typedef struct TreeChange {
    uint32_t Action;
    char Name[32];
} TreeChange;
```

Data Fields

- uint32_t Action
field indicating which action occurred on the tree. A value of the ZITreeAction_enum (TREE_ACTION) enum.
- char Name
Name of the Path that has been added, removed or changed.

6.3.7. union ziEvent::Val

```
typedef union ziEvent::Val {  
    void* Void;  
    DemodSample* SampleDemod;  
    AuxInSample* SampleAuxIn;  
    DIOSample* SampleDIO;  
    ziDoubleType* Double;  
    ziIntegerType* Integer;  
    TreeChange* Tree;  
    ByteArrayData* ByteArray;  
    ScopeWave* Wave;  
    uint64_t alignment;  
} ziEvent::Val;
```

Data Fields

- void* Void
- DemodSample* SampleDemod
- AuxInSample* SampleAuxIn
- DIOSample* SampleDIO
- ziDoubleType* Double
- ziIntegerType* Integer
- TreeChange* Tree
- ByteArrayData* ByteArray
- ScopeWave* Wave
- uint64_t alignment

6.3.8. struct ZIAdvisorHeader

```
typedef struct ZIAdvisorHeader {  
    uint64_t sampleCount;  
    uint8_t flags;  
    uint8_t sampleFormat;  
    uint8_t reserved0[6];  
    uint8_t reserved1[8];  
} ZIAdvisorHeader;
```

Data Fields

- `uint64_t sampleCount`
Total sample count considered for advisor.
- `uint8_t flags`
Flags.
- `uint8_t sampleFormat`
Sample format Bode = 0, Step = 1, Impulse = 2.
- `uint8_t reserved0`
Reserved space for future use.
- `uint8_t reserved1`
Reserved space for future use.

6.3.9. struct ZIAdvisorSample

```
typedef struct ZIAdvisorSample {  
    double grid;  
    double x;  
    double y;  
} ZIAdvisorSample;
```

Data Fields

- double grid
Grid.
- double x
X.
- double y
Y.

6.3.10. struct ZIAdvisorWave

```
typedef struct ZIAdvisorWave {  
    ZITimeStamp timeStamp;  
    ZIAdvisorHeader header;  
    ZIAdvisorSample data[0];  
    union ZIAdvisorWave::@4 data;  
} ZIAdvisorWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- ZIAdvisorHeader header
- ZIAdvisorSample data
- union ZIAdvisorWave::@4 data
Advisor data vector.

6.3.11. struct ZIAsyncReply

```
typedef struct ZIAsyncReply {  
    ZITimeStamp timeStamp;  
    ZITimeStamp sampleTimeStamp;  
    uint16_t command;  
    uint16_t resultCode;  
    ZIAsyncTag tag;  
} ZIAsyncReply;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp of the reply (server clock)
- ZITimeStamp sampleTimeStamp
Time stamp of the target node sample, to which the reply belongs.
- uint16_t command
Command: 1 - ziAPIAsyncSetDoubleData 2 - ziAPIAsyncSetIntegerData 3 - ziAPIAsyncSetByteArray 4 - ziAPIAsyncSubscribe 5 - ziAPIAsyncUnSubscribe 6 - ziAPIAsyncGetValueAsPollData.
- uint16_t resultCode
Command result code (cast to ZIResult_enum)
- ZIAsyncTag tag
Tag sent along with the async command.

6.3.12. struct ZIAuxInSample

The structure used to hold data for a single auxiliary inputs sample.

```
#include "ziAPI.h"

typedef struct ZIAuxInSample {
    ZITimeStamp timeStamp;
    double ch0;
    double ch1;
} ZIAuxInSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the values have been measured.
- double ch0
Channel 0 voltage.
- double ch1
Channel 1 voltage.

6.3.13. struct ZByteArray

The structure used to hold an arbitrary array of bytes. This is the version without timestamp used in API Level 1 compatibility mode.

```
#include "ziAPI.h"

typedef struct ZByteArray {
    uint32_t length;
    uint8_t bytes[0];
} ZByteArray;
```

Data Fields

- `uint32_t length`
Length of the data readable from the Bytes field.
- `uint8_t bytes`
The data itself. The array has the size given in length.

6.3.14. struct ZIByteArrayTS

The structure used to hold an arbitrary array of bytes. This is the same as [ZIByteArray](#), but with timestamp.

```
#include "ziAPI.h"

typedef struct ZIByteArrayTS {
    ZITimeStamp timeStamp;
    uint32_t length;
    uint8_t bytes[0];
} ZIByteArrayTS;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- uint32_t length
length of the data readable from the bytes field
- uint8_t bytes
the data itself. The array has the size given in length

6.3.15. struct ZICntSample

The structure used to hold data for a single counter sample.

```
#include "ziAPI.h"

typedef struct ZICntSample {
    ZITimeStamp timeStamp;
    uint16_t counter;
    uint16_t reserved0;
    uint16_t id;
    uint16_t reserved1;
} ZICntSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the values have been measured.
- uint16_t counter
Counter value.
- uint16_t reserved0
Reserved.
- uint16_t id
Trigger id.
- uint16_t reserved1
Reserved.

6.3.16. struct ZIDemodSample

The structure used to hold data for a single demodulator sample.

```
#include "ziAPI.h"

typedef struct ZIDemodSample {
    ZITimeStamp timeStamp;
    double x;
    double y;
    double frequency;
    double phase;
    uint32_t dioBits;
    uint32_t trigger;
    double auxIn0;
    double auxIn1;
} ZIDemodSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the sample has been measured.
- double x
X part of the sample.
- double y
Y part of the sample.
- double frequency
Frequency at that sample.
- double phase
Phase at that sample.
- uint32_t dioBits
the current bits of the DIO.
- uint32_t trigger
trigger bits
- double auxIn0
value of Aux input 0.
- double auxIn1
value of Aux input 1.

6.3.17. struct ZIDIOSample

The structure used to hold data for a single digital I/O sample.

```
#include "ziAPI.h"

typedef struct ZIDIOSample {
    ZITimeStamp timeStamp;
    uint32_t bits;
    uint32_t reserved;
} ZIDIOSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the values have been measured.
- uint32_t bits
The digital I/O values.
- uint32_t reserved
Filler to keep 8 bytes alignment in the array of [ZIDIOSample](#) structures.

6.3.18. struct ZIDoubleDataTS

The structure used to hold a single IEEE double value. Same as ZIDoubleData, but with timestamp.

```
#include "ziAPI.h"

typedef struct ZIDoubleDataTS {
    ZITimeStamp timeStamp;
    ZIDoubleData value;
} ZIDoubleDataTS;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the value has changed.
- ZIDoubleData value

6.3.19. struct ZIEvent

This struct holds event data forwarded by the Data Server.

```
#include "ziAPI.h"

typedef struct ZIEvent {
    uint32_t valueType;
    uint32_t count;
    uint8_t path[256];
    void* untyped;
    ZIDoubleData* doubleData;
    ZIDoubleDataTS* doubleDataTS;
    ZIIntegerData* integerData;
    ZIIntegerDataTS* integerDataTS;
    ZIByteArray* byteArray;
    ZIByteArrayTS* byteArrayTS;
    ZICntSample* cntSample;
    ZITreeChangeData* treeChangeData;
    TreeChange* treeChangeDataOld;
    ZIDemodSample* demodSample;
    ZIAuxInSample* auxInSample;
    ZIDIOSample* dioSample;
    ZIScopeWave* scopeWave;
    ZIScopeWaveEx* scopeWaveEx;
    ScopeWave* scopeWaveOld;
    ZIPWAWave* pwaWave;
    ZISweeperWave* sweeperWave;
    ZISpectrumWave* spectrumWave;
    ZIAdvisorWave* advisorWave;
    ZIAsyncReply* asyncReply;
    ZIVectorData* vectorData;
    ZIImpedanceSample* impedanceSample;
    uint64_t alignment;
    union ZIEvent::@6 value;
    uint8_t data[0x400000];
} ZIEvent;
```

Data Fields

- `uint32_t valueType`
Specifies the type of the data held by the [ZIEvent](#), see [ZIValueType_enum](#).
- `uint32_t count`
Number of values available in this event.
- `uint8_t path`
The path to the node from which the event originates.
- `void* untyped`
For convenience. The void field doesn't have a corresponding data type.
- `ZIDoubleData* doubleData`
when `valueType == ZI_VALUE_TYPE_DOUBLE_DATA`
- `ZIDoubleDataTS* doubleDataTS`
when `valueType == ZI_VALUE_TYPE_DOUBLE_DATA_TS`
- `ZIIntegerData* integerData`

- when valueType == ZI_VALUE_TYPE_INTEGER_DATA
- [ZIIntegerDataTS](#)* integerDataTS
when valueType == ZI_VALUE_TYPE_INTEGER_DATA_TS
- [ZIByteArray](#)* byteArray
when valueType == ZI_VALUE_TYPE_BYTE_ARRAY
- [ZIByteArrayTS](#)* byteArrayTS
when valueType == ZI_VALUE_TYPE_BYTE_ARRAY_TS
- [ZICntSample](#)* cntSample
when valueType == ZI_VALUE_TYPE_CNT_SAMPLE
- [ZITreeChangeData](#)* treeChangeData
when valueType == ZI_VALUE_TYPE_TREE_CHANGE_DATA
- [TreeChange](#)* treeChangeDataOld
when valueType ==
ZI_VALUE_TYPE_TREE_CHANGE_DATA_OLD
- [ZIDemodSample](#)* demodSample
when valueType == ZI_VALUE_TYPE_DEMOD_SAMPLE
- [ZIAuxInSample](#)* auxInSample
when valueType == ZI_VALUE_TYPE_AUXIN_SAMPLE
- [ZIDIOSample](#)* dioSample
when valueType == ZI_VALUE_TYPE_DIO_SAMPLE
- [ZIScopeWave](#)* scopeWave
when valueType == ZI_VALUE_TYPE_SCOPE_WAVE
- [ZIScopeWaveEx](#)* scopeWaveEx
when valueType == ZI_VALUE_TYPE_SCOPE_WAVE_EX
- [ScopeWave](#)* scopeWaveOld
when valueType == ZI_VALUE_TYPE_SCOPE_WAVE_OLD
- [ZIPWAWave](#)* pwaWave
when valueType == ZI_VALUE_TYPE_PWA_WAVE
- [ZISweeperWave](#)* sweeperWave
when valueType == ZI_VALUE_TYPE_SWEEPER_WAVE
- [ZISpectrumWave](#)* spectrumWave
when valueType == ZI_VALUE_TYPE_SPECTRUM_WAVE
- [ZIAdvisorWave](#)* advisorWave
when valueType == ZI_VALUE_TYPE_ADVISOR_WAVE
- [ZIAsyncReply](#)* asyncReply
when valueType == ZI_VALUE_TYPE_ASYNC_REPLY
- [ZIVectorData](#)* vectorData

- when valueType == ZI_VALUE_TYPE_VECTOR_DATA
- [ZIImpedanceSample](#)* impedanceSample
when valueType == ZI_VALUE_TYPE_IMPEDANCE_SAMPLE
- uint64_t alignment
ensure union size is 8 bytes
- union ZIEvent::@6 value
Convenience pointer to allow for access to the first entry in Data using the correct type according to [ZIEvent.valueType](#) field.
- uint8_t data
The raw value data.

Detailed Description

[ZIEvent](#) is used to give out events like value changes or errors to the user. Event handling functionality is provided by [ziAPISubscribe](#) and [ziAPIUnSubscribe](#) as well as [ziAPIPollDataEx](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

void ProcessEvent(ZIEvent* Event) {
    unsigned int j;

    switch (Event->valueType) {
    case ZI_VALUE_TYPE_DOUBLE_DATA:

        printf("%u elements of double data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", Event->value.doubleData[j]);

        break;

    case ZI_VALUE_TYPE_INTEGER_DATA:

        printf("%u elements of integer data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", (float)Event->value.integerData[j]);

        break;

    case ZI_VALUE_TYPE_DEMOD_SAMPLE:

        printf("%u elements of sample data %s\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("TS=%f, X=%f, Y=%f.\n",
                (float)Event->value.demodSample[j].timeStamp,
                Event->value.demodSample[j].x,
                Event->value.demodSample[j].y);
```

```
        break;

    case ZI_VALUE_TYPE_TREE_CHANGE_DATA:

        printf("%u elements of tree-changed data, %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++) {
            switch (Event->value.treeChangeDataOld[j].Action) {
                case ZI_TREE_ACTION_REMOVE:
                    printf("Tree removed: %s\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_ADD:
                    printf("treeChangeDataOld added: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_CHANGE:
                    printf("treeChangeDataOld changed: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;
            }
        }

        break;

    default:

        printf("Unexpected event value type: %d.\n", Event->valueType);
        break;
    }
}
```

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

6.3.20. struct ziEvent

This struct holds event data forwarded by the Data Server. Deprecated: See [ZiEvent](#).

```
#include "ziAPI.h"

typedef struct ziEvent {
    uint32_t Type;
    uint32_t Count;
    unsigned char Path[256];
    union ziEvent::Val Val;
    unsigned char Data[0x400000];
} ziEvent;
```

Data Structures

- union [ziEvent::Val](#)

Data Fields

- uint32_t Type
- uint32_t Count
- unsigned char Path
- union ziEvent::Val Val
- unsigned char Data

Detailed Description

[ziEvent](#) is used to give out events like value changes or errors to the user. Event handling functionality is provided by [ziAPISubscribe](#) and [ziAPIUnSubscribe](#) as well as [ziAPIPollDataEx](#).

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

void ProcessEvent(ZiEvent* Event) {
    unsigned int j;

    switch (Event->valueType) {
    case ZI_VALUE_TYPE_DOUBLE_DATA:

        printf("%u elements of double data: %s.\n",
            Event->count,
            Event->path);
```

```
    for (j = 0; j < Event->count; j++)
        printf("%f\n", Event->value.doubleData[j]);

    break;

case ZI_VALUE_TYPE_INTEGER_DATA:

    printf("%u elements of integer data: %s.\n",
        Event->count,
        Event->path);

    for (j = 0; j < Event->count; j++)
        printf("%f\n", (float)Event->value.integerData[j]);

    break;

case ZI_VALUE_TYPE_DEMOD_SAMPLE:

    printf("%u elements of sample data %s\n",
        Event->count,
        Event->path);

    for (j = 0; j < Event->count; j++)
        printf("TS=%f, X=%f, Y=%f.\n",
            (float)Event->value.demodSample[j].timeStamp,
            Event->value.demodSample[j].x,
            Event->value.demodSample[j].y);

    break;

case ZI_VALUE_TYPE_TREE_CHANGE_DATA:

    printf("%u elements of tree-changed data, %s.\n",
        Event->count,
        Event->path);

    for (j = 0; j < Event->count; j++) {
        switch (Event->value.treeChangeDataOld[j].Action) {
            case ZI_TREE_ACTION_REMOVE:
                printf("Tree removed: %s\n",
                    Event->value.treeChangeDataOld[j].Name);
                break;

            case ZI_TREE_ACTION_ADD:
                printf("treeChangeDataOld added: %s.\n",
                    Event->value.treeChangeDataOld[j].Name);
                break;

            case ZI_TREE_ACTION_CHANGE:
                printf("treeChangeDataOld changed: %s.\n",
                    Event->value.treeChangeDataOld[j].Name);
                break;
        }
    }

    break;

default:

    printf("Unexpected event value type: %d.\n", Event->valueType);
    break;
}
```

Data Structure Documentation

union ziEvent::Val

```
typedef union ziEvent::Val {  
    void* Void;  
    DemodSample* SampleDemod;  
    AuxInSample* SampleAuxIn;  
    DIOSample* SampleDIO;  
    ziDoubleType* Double;  
    ziIntegerType* Integer;  
    TreeChange* Tree;  
    ByteArrayData* ByteArray;  
    ScopeWave* Wave;  
    uint64_t alignment;  
} ziEvent::Val;
```

Data Fields

- void* Void
- [DemodSample](#)* SampleDemod
- [AuxInSample](#)* SampleAuxIn
- [DIOSample](#)* SampleDIO
- ziDoubleType* Double
- ziIntegerType* Integer
- [TreeChange](#)* Tree
- [ByteArrayData](#)* ByteArray
- [ScopeWave](#)* Wave
- uint64_t alignment

6.3.21. struct ZIImpedanceSample

The structure used to hold data for a single impedance sample.

```
#include "ziAPI.h"

typedef struct ZIImpedanceSample {
    ZITimeStamp timeStamp;
    double realz;
    double imagz;
    double frequency;
    double phase;
    uint32_t flags;
    uint32_t trigger;
    double param0;
    double param1;
    double drive;
    double bias;
} ZIImpedanceSample;
```

Data Fields

- ZITimeStamp timeStamp
Timestamp at which the sample has been measured.
- double realz
Real part of the impedance sample.
- double imagz
Imaginary part of the impedance sample.
- double frequency
Frequency at that sample.
- double phase
Phase at that sample.
- uint32_t flags
Flags (see ZIImpFlags_enum)
- uint32_t trigger
Trigger bits.
- double param0
Value of model parameter 0.
- double param1
Value of model parameter 1.
- double drive
Drive amplitude.
- double bias
Bias voltage.

6.3.22. struct ZIIntegerDataTS

The structure used to hold a single 64bit signed integer value. Same as ZIIntegerData, but with timestamp.

```
#include "ziAPI.h"

typedef struct ZIIntegerDataTS {
    ZITimeStamp timeStamp;
    ZIIntegerData value;
} ZIIntegerDataTS;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the value has changed.
- ZIIntegerData value

6.3.23. struct ZIModuleEvent

This struct holds data of a single chunk from module lookup.

```
#include "ziAPI.h"

typedef struct ZIModuleEvent {
    uint64_t allocatedSize;
    ZIModuleHeader
        header;
    ZIEvent
        value[0];
} ZIModuleEvent;
```

Data Fields

- `uint64_t allocatedSize`
For internal use - never modify!
- `ZIModuleHeader` header
Module-specific event header.
- `ZIEvent` value
Defines location of stored `ZIEvent`.

6.3.24. struct ZIModuleHeader

Module-specific event header.

```
#include "ziAPI.h"

typedef struct ZIModuleHeader {
    ZIModuleHeaderType_enum
        type;
    void* untyped;
    ZISWTriggerHeader* swTrigger;
    ZISweeperHeader* sweeper;
    union ZIModuleHeader::@7 ptr;
} ZIModuleHeader;
```

Data Fields

- `ZIModuleHeaderType_enum` type
- `void* untyped`
- `ZISWTriggerHeader* swTrigger`
- `ZISweeperHeader* sweeper`
- `union ZIModuleHeader::@7 ptr`

6.3.25. struct ZIModuleHeaderSweeper

Structure to hold information about data returned from the Sweep Module.

```
#include "ziAPI.h"

typedef struct ZIModuleHeaderSweeper {
    char traceName[256];
} ZIModuleHeaderSweeper;
```

Data Fields

- char traceName

6.3.26. struct ZIPWASample

Single PWA sample value.

```
#include "ziAPI.h"

typedef struct ZIPWASample {
    double binPhase;
    double x;
    double y;
    uint32_t countBin;
    uint32_t reserved;
} ZIPWASample;
```

Data Fields

- double binPhase
Phase position of each bin.
- double x
Real PWA result or X component of a demod PWA.
- double y
Y component of the demod PWA.
- uint32_t countBin
Number of events per bin.
- uint32_t reserved
Reserved.

6.3.27. struct ZIPWAWave

PWA Wave.

```
#include "ziAPI.h"

typedef struct ZIPWAWave {
    ZITimeStamp timeStamp;
    uint64_t sampleCount;
    uint32_t inputSelect;
    uint32_t oscSelect;
    uint32_t harmonic;
    uint32_t binCount;
    double frequency;
    uint8_t pwaType;
    uint8_t mode;
    uint8_t overflow;
    uint8_t commensurable;
    uint32_t reservedUInt;
    ZIPWASample
        data[0];
} ZIPWAWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- uint64_t sampleCount
Total sample count considered for PWA.
- uint32_t inputSelect
Input selection used for the PWA.
- uint32_t oscSelect
Oscillator used for the PWA.
- uint32_t harmonic
Harmonic setting.
- uint32_t binCount
Bin count of the PWA.
- double frequency
Frequency during PWA accumulation.
- uint8_t pwaType
Type of the PWA.
- uint8_t mode
PWA Mode [0: zoom PWA, 1: harmonic PWA].
- uint8_t overflow
Overflow indicators. overflow[0]: Data accumulator overflow, overflow[1]: Counter at limit, overflow[6..2]: Reserved, overflow[7]: Invalid (missing frames).
- uint8_t commensurable

Commensurability of the data.

- `uint32_t reservedUInt`
Reserved 32bit.
- `ZIPWASample` data
PWA data vector.

6.3.28. struct ZIScopeWave

The structure used to hold scope data. The data may be formatted differently, depending on settings. See the description of the structure members for details.

```
#include "ziAPI.h"

typedef struct ZIScopeWave {
    ZITimeStamp timeStamp;
    ZITimeStamp triggerTimeStamp;
    double dt;
    uint8_t channelEnable[4];
    uint8_t channelInput[4];
    uint8_t triggerEnable;
    uint8_t triggerInput;
    uint8_t reserved0[2];
    uint8_t channelBWLimit[4];
    uint8_t channelMath[4];
    float channelScaling[4];
    uint32_t sequenceNumber;
    uint32_t segmentNumber;
    uint32_t blockNumber;
    uint64_t totalSamples;
    uint8_t dataTransferMode;
    uint8_t blockMarker;
    uint8_t flags;
    uint8_t sampleFormat;
    uint32_t sampleCount;
    int16_t dataInt16[0];
    int32_t dataInt32[0];
    float dataFloat[0];
    union ZIScopeWave::@0 data;
} ZIScopeWave;
```

Data Fields

- **ZITimeStamp timeStamp**
Time stamp of the last sample in this data block.
- **ZITimeStamp triggerTimeStamp**
Time stamp of the trigger (may also fall between samples and in another block)
- **double dt**
Time difference between samples in seconds.
- **uint8_t channelEnable**
Up to four channels: if channel is enabled, corresponding element is non-zero.
- **uint8_t channelInput**
Specifies the input source for each of the scope four channels: 0 = Signal Input 1, 1 = Signal Input 2, 2 = Trigger Input 1, 3 = Trigger Input 2, 4 = Aux Output 1, 5 = Aux Output 2, 6 = Aux Output 3, 7 = Aux Output 4, 8 = Aux Input 1, 9 = Aux Input 2.
- **uint8_t triggerEnable**
Non-zero if trigger is enabled: Bit(0): rising edge trigger off = 0, on = 1. Bit(1): falling edge trigger off = 0, on = 1.

- `uint8_t triggerInput`
Trigger source (same values as for channel input)
- `uint8_t reserved0`
- `uint8_t channelBWLimit`
Bandwidth-limit flag, per channel. Bit(0): off = 0, on = 1
Bit(7...1): Reserved.
- `uint8_t channelMath`
Math Operation (e.g averaging) Bit (7..0): Reserved.
- `float channelScaling`
Data scaling factors for up to 4 channels.
- `uint32_t sequenceNumber`
Current scope shot sequence number. Identifies a scope shot.
- `uint32_t segmentNumber`
Current segment number.
- `uint32_t blockNumber`
Current block number from the beginning of a scope shot.
Large scope shots are split into blocks, which need to be concatenated to obtain the complete scope shot.
- `uint64_t totalSamples`
Total number of samples in one channel in the current scope shot, same for all channels.
- `uint8_t dataTransferMode`
Data transfer mode SingleTransfer = 0, BlockTransfer = 1, ContinuousTransfer = 3, FFTSingleTransfer = 4. Other values are reserved.
- `uint8_t blockMarker`
Block marker: Bit (0): 1 = End marker for continuous or multi-block transfer Bit (7..0): Reserved.
- `uint8_t flags`
Indicator Flags. Bit (0): 1 = Data loss detected (samples are 0), Bit (1): 1 = Missed trigger, Bit (2): 1 = Transfer failure (corrupted data).
- `uint8_t sampleFormat`
Data format of samples: Int16 = 0, Int32 = 1, Float = 2, Int16Interleaved = 4, Int32Interleaved = 5, FloatInterleaved = 6.
- `uint32_t sampleCount`
Number of samples in one channel in the current block, same for all channels.
- `int16_t dataInt16`

- `int32_t dataInt32`
- `float dataFloat`
- `union ZIScopeWave::@0 data`
First wave data.

6.3.29. struct ZIScopeWaveEx

The structure used to hold scope data (extended). The data may be formatted differently, depending on settings. See the description of the structure members for details.

```
#include "ziAPI.h"

typedef struct ZIScopeWaveEx {
    ZITimeStamp timeStamp;
    ZITimeStamp triggerTimeStamp;
    double dt;
    uint8_t channelEnable[4];
    uint8_t channelInput[4];
    uint8_t triggerEnable;
    uint8_t triggerInput;
    uint8_t reserved0[2];
    uint8_t channelBWLimit[4];
    uint8_t channelMath[4];
    float channelScaling[4];
    uint32_t sequenceNumber;
    uint32_t segmentNumber;
    uint32_t blockNumber;
    uint64_t totalSamples;
    uint8_t dataTransferMode;
    uint8_t blockMarker;
    uint8_t flags;
    uint8_t sampleFormat;
    uint32_t sampleCount;
    double channelOffset[4];
    uint64_t reserved1[32];
    int16_t dataInt16[0];
    int32_t dataInt32[0];
    float dataFloat[0];
    union ZIScopeWaveEx::@1 data;
} ZIScopeWaveEx;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp of the last sample in this data block.
- ZITimeStamp triggerTimeStamp
Time stamp of the trigger (may also fall between samples and in another block)
- double dt
Time difference between samples in seconds.
- uint8_t channelEnable
Up to four channels: if channel is enabled, corresponding element is non-zero.
- uint8_t channelInput
Specifies the input source for each of the scope four channels: 0 = Signal Input 1, 1 = Signal Input 2, 2 = Trigger Input 1, 3 = Trigger Input 2, 4 = Aux Output 1, 5 = Aux Output 2, 6 = Aux Output 3, 7 = Aux Output 4, 8 = Aux Input 1, 9 = Aux Input 2.
- uint8_t triggerEnable

Non-zero if trigger is enabled: Bit(0): rising edge trigger off = 0, on = 1. Bit(1): falling edge trigger off = 0, on = 1.

- `uint8_t triggerInput`
Trigger source (same values as for channel input)
- `uint8_t reserved0`
- `uint8_t channelBWLimit`
Bandwidth-limit flag, per channel. Bit(0): off = 0, on = 1
Bit(7...1): Reserved.
- `uint8_t channelMath`
Math Operation (e.g averaging) Bit (7..0): Reserved.
- `float channelScaling`
Data scaling factors for up to 4 channels.
- `uint32_t sequenceNumber`
Current scope shot sequence number. Identifies a scope shot.
- `uint32_t segmentNumber`
Current segment number.
- `uint32_t blockNumber`
Current block number from the beginning of a scope shot.
Large scope shots are split into blocks, which need to be concatenated to obtain the complete scope shot.
- `uint64_t totalSamples`
Total number of samples in one channel in the current scope shot, same for all channels.
- `uint8_t dataTransferMode`
Data transfer mode SingleTransfer = 0, BlockTransfer = 1, ContinuousTransfer = 3, FFTSingleTransfer = 4. Other values are reserved.
- `uint8_t blockMarker`
Block marker: Bit (0): 1 = End marker for continuous or multi-block transfer Bit (7..0): Reserved.
- `uint8_t flags`
Indicator Flags. Bit (0): 1 = Data loss detected (samples are 0), Bit (1): 1 = Missed trigger, Bit (2): 1 = Transfer failure (corrupted data).
- `uint8_t sampleFormat`
Data format of samples: Int16 = 0, Int32 = 1, Float = 2, Int16Interleaved = 4, Int32Interleaved = 5, FloatInterleaved = 6.
- `uint32_t sampleCount`
Number of samples in one channel in the current block, same for all channels.

- `double channelOffset`
Data offset (scaled) for up to 4 channels.
- `uint64_t reserved1`
- `int16_t dataInt16`
- `int32_t dataInt32`
- `float dataFloat`
- `union ZIScopeWaveEx::@1 data`
First wave data.

6.3.30. struct ZISpectrumDemodSample

```
typedef struct ZISpectrumDemodSample {  
    double grid;  
    double filter;  
    double x;  
    double y;  
    double r;  
} ZISpectrumDemodSample;
```

Data Fields

- double grid
Grid.
- double filter
Filter strength at the specific grid point.
- double x
X.
- double y
Y.
- double r
R.

6.3.31. struct ZISpectrumHeader

```
typedef struct ZISpectrumHeader {  
    uint64_t sampleCount;  
    uint8_t flags;  
    uint8_t sampleFormat;  
    uint8_t spectrumMode;  
    uint8_t window;  
    uint8_t reserved0[4];  
    uint8_t reserved1[8];  
    double bandwidth;  
    double rate;  
    double center;  
    double resolution;  
    double aliasingReject;  
    double nenbw;  
    double overlap;  
} ZISpectrumHeader;
```

Data Fields

- `uint64_t sampleCount`
Total sample count considered for spectrum.
- `uint8_t flags`
Flags Bit 0: Power Bit 1: Spectral density Bit 2: Absolute frequency Bit 3: Full span.
- `uint8_t sampleFormat`
Sample format Demodulator = 0.
- `uint8_t spectrumMode`
Spectrum mode FFT(x+iy) = 0, FFT(r) = 1, FFT(theta) = 2, FFT(freq) = 3, FFT(dtheta/dt)/2pi = 4.
- `uint8_t window`
Window Rectangular = 0, Hann = 1, Hamming = 2, Blackman Harris = 3.
- `uint8_t reserved0`
Reserved space for future use.
- `uint8_t reserved1`
Reserved space for future use.
- `double bandwidth`
Filter bandwidth.
- `double rate`
Rate of the sampled data.
- `double center`
FFT center value.
- `double resolution`
FFT bin resolution.
- `double aliasingReject`

Aliasing reject (dB)

- double nenbw
Correction factor for the used window when calculating spectral density.
- double overlap
FFT overlap [0 .. 1[.

6.3.32. struct ZISpectrumWave

```
typedef struct ZISpectrumWave {  
    ZITimeStamp timeStamp;  
    ZISpectrumHeader header;  
    ZISpectrumDemodSample dataDemod[0];  
    union ZISpectrumWave::@3 data;  
} ZISpectrumWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- ZISpectrumHeader header
- ZISpectrumDemodSample dataDemod
- union ZISpectrumWave::@3 data
Spectrum data vector.

6.3.33. struct ZIStatisticSample

```
typedef struct ZIStatisticSample {  
    double avg;  
    double stddev;  
    double pwr;  
} ZIStatisticSample;
```

Data Fields

- double avg
Average value or single value.
- double stddev
Standard deviation.
- double pwr
Power value.

6.3.34. struct ZISweeperDemodSample

```
typedef struct ZISweeperDemodSample {  
    double grid;  
    double bandwidth;  
    uint64_t count;  
    double tc;  
    double tcMeas;  
    double settling;  
    ZITimeStamp setTimeStamp;  
    ZITimeStamp nextTimeStamp;  
    ZIStatisticSample x;  
    ZIStatisticSample y;  
    ZIStatisticSample r;  
    ZIStatisticSample phase;  
    ZIStatisticSample frequency;  
    ZIStatisticSample auxin0;  
    ZIStatisticSample auxin1;  
} ZISweeperDemodSample;
```

Data Fields

- double grid
Grid value (x-axis)
- double bandwidth
Demodulator bandwidth used for the specific sweep point.
- uint64_t count
Sample count used for statistic calculation.
- double tc
Time constant calculated for the specific sweep point.
- double tcMeas
Time constant used by the device.
- double settling
Settling time (s) used to wait until averaging operation is started.
- ZITimeStamp setTimeStamp
Time stamp when the grid value was set on the device.
- ZITimeStamp nextTimeStamp
Time stamp when the first statistic value was recorded.
- ZIStatisticSample x
Sweep point statistic result of X.
- ZIStatisticSample y
Sweep point statistic result of Y.
- ZIStatisticSample r
Sweep point statistic result of R.
- ZIStatisticSample phase
Sweep point statistic result of phase.

- ZIStatisticSample frequency
Sweep point statistic result of frequency.
- ZIStatisticSample auxin0
Sweep point statistic result of auxin0.
- ZIStatisticSample auxin1
Sweep point statistic result of auxin1.

6.3.35. struct ZISweeperDoubleSample

```
typedef struct ZISweeperDoubleSample {  
    double grid;  
    double bandwidth;  
    uint64_t count;  
    ZIStatisticSample value;  
} ZISweeperDoubleSample;
```

Data Fields

- double grid
Grid value (x-axis)
- double bandwidth
Bandwidth.
- uint64_t count
Sample count used for statistic calculation.
- ZIStatisticSample value
Result value (y-axis)

6.3.36. struct ZISweeperHeader

```
typedef struct ZISweeperHeader {  
    uint64_t sampleCount;  
    uint8_t flags;  
    uint8_t sampleFormat;  
    uint8_t sweepMode;  
    uint8_t bandwidthMode;  
    uint8_t reserved0[4];  
    uint8_t reserved1[8];  
} ZISweeperHeader;
```

Data Fields

- `uint64_t sampleCount`
Total sample count considered for sweeper.
- `uint8_t flags`
Flags Bit 0: Phase unwrap Bit 1: Sinc filter.
- `uint8_t sampleFormat`
Sample format Double = 0, Demodulator = 1.
- `uint8_t sweepMode`
Sweep mode Sequential = 0, Binary = 1, Bidirectional = 2, Reverse = 3.
- `uint8_t bandwidthMode`
Bandwidth mode Manual = 0, Fixed = 1, Auto = 2.
- `uint8_t reserved0`
Reserved space for future use.
- `uint8_t reserved1`
Reserved space for future use.

6.3.37. struct ZISweeperImpedanceSample

```
typedef struct ZISweeperImpedanceSample {  
    double grid;  
    double bandwidth;  
    uint64_t count;  
    double tc;  
    double tcMeas;  
    double settling;  
    ZITimeStamp setTimeStamp;  
    ZITimeStamp nextTimeStamp;  
    ZIStatisticSample realz;  
    ZIStatisticSample imagz;  
    ZIStatisticSample absz;  
    ZIStatisticSample phasez;  
    ZIStatisticSample frequency;  
    ZIStatisticSample param0;  
    ZIStatisticSample param1;  
    ZIStatisticSample drive;  
    ZIStatisticSample bias;  
} ZISweeperImpedanceSample;
```

Data Fields

- double grid
Grid value (x-axis)
- double bandwidth
Demodulator bandwidth used for the specific sweep point.
- uint64_t count
Sample count used for statistic calculation.
- double tc
Time constant calculated for the specific sweep point.
- double tcMeas
Time constant used by the device.
- double settling
Settling time (s) used to wait until averaging operation is started.
- ZITimeStamp setTimeStamp
Time stamp when the grid value was set on the device.
- ZITimeStamp nextTimeStamp
Time stamp when the first statistic value was recorded.
- ZIStatisticSample realz
Sweep point statistic result of X.
- ZIStatisticSample imagz
Sweep point statistic result of Y.
- ZIStatisticSample absz
Sweep point statistic result of R.
- ZIStatisticSample phasez

Sweep point statistic result of phase.

- ZIStatisticSample frequency
Sweep point statistic result of frequency.

- ZIStatisticSample param0
Sweep point statistic result of param0.

- ZIStatisticSample param1
Sweep point statistic result of param1.

- ZIStatisticSample drive
Sweep point statistic result of drive amplitude.

- ZIStatisticSample bias
Sweep point statistic result of bias.

6.3.38. struct ZISweeperWave

```
typedef struct ZISweeperWave {  
    ZITimeStamp timeStamp;  
    ZISweeperHeader header;  
    ZISweeperDoubleSample dataDouble[0];  
    ZISweeperDemodSample dataDemod[0];  
    ZISweeperImpedanceSample dataImpedance[0];  
    union ZISweeperWave::@2 data;  
} ZISweeperWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- ZISweeperHeader header
- ZISweeperDoubleSample dataDouble
- ZISweeperDemodSample dataDemod
- ZISweeperImpedanceSample dataImpedance
- union ZISweeperWave::@2 data
Sweeper data vector.

6.3.39. struct ZISWTriggerHeader

Structure to hold information about data returned from the SW Trigger Module.

```
#include "ziAPI.h"

typedef struct ZISWTriggerHeader {
    ZITimeStamp triggerStart;
    uint64_t triggerNumber;
    uint32_t cols;
    uint32_t rows;
    uint32_t repetitions;
    uint32_t flags;
    uint8_t reserved0[32];
} ZISWTriggerHeader;
```

Data Fields

- ZITimeStamp triggerStart
Trigger timestamp.
- uint64_t triggerNumber
Trigger counter since execution start.
- uint32_t cols
Number of columns.
- uint32_t rows
Number of rows.
- uint32_t repetitions
Number of repetitions in grid mode.
- uint32_t flags
Flags Bit (0): Finished (all repetitions recorded)
- uint8_t reserved0
Reserved space for future use.

6.3.40. struct ZITreeChangeData

The struct is holding info about added or removed nodes.

```
#include "ziAPI.h"

typedef struct ZITreeChangeData {
    ZITimeStamp timeStamp;
    uint32_t action;
    char name[32];
} ZITreeChangeData;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- uint32_t action
field indicating which action occurred on the tree. A value of the ZITreeAction_enum.
- char name
Name of the Path that has been added, removed or changed.

6.3.41. struct ZIVectorData

The structure used to hold vector data block. See the description of the structure members for details.

```
#include "ziAPI.h"

typedef struct ZIVectorData {
    ZITimeStamp timeStamp;
    uint32_t sequenceNumber;
    uint32_t blockNumber;
    uint64_t totalElements;
    uint64_t blockOffset;
    uint32_t blockElements;
    uint8_t flags;
    uint8_t elementType;
    uint8_t reserved0[2];
    uint64_t reserved1[32];
    uint8_t dataUInt8[0];
    uint16_t dataUInt16[0];
    uint32_t dataUInt32[0];
    uint64_t dataUInt64[0];
    int8_t dataInt8[0];
    int16_t dataInt16[0];
    int32_t dataInt32[0];
    int64_t dataInt64[0];
    double dataDouble[0];
    float dataFloat[0];
    union ZIVectorData::@5 data;
} ZIVectorData;
```

Data Fields

- **ZITimeStamp timeStamp**
Time stamp of this array data block.
- **uint32_t sequenceNumber**
Current array transfer sequence number. Incremented for each new transfer. Stays same for all blocks of a single array transfer.
- **uint32_t blockNumber**
Current block number from the beginning of an array transfer. Large array transfers are split into blocks, which need to be concatenated to obtain the complete array.
- **uint64_t totalElements**
Total number of elements in the array.
- **uint64_t blockOffset**
Offset of the current block first element from the beginning of the array.
- **uint32_t blockElements**
Number of elements in the current block.
- **uint8_t flags**
Block marker: Bit (0): 1 = End marker for multi-block transfer
Bit (1): 1 = Transfer failure Bit (7..2): Reserved.

- `uint8_t elementType`
Vector element type, see `ZIVectorElementType_enum`.
- `uint8_t reserved0`
- `uint64_t reserved1`
- `uint8_t dataUInt8`
- `uint16_t dataUInt16`
- `uint32_t dataUInt32`
- `uint64_t dataUInt64`
- `int8_t dataInt8`
- `int16_t dataInt16`
- `int32_t dataInt32`
- `int64_t dataInt64`
- `double dataDouble`
- `float dataFloat`
- `union ZIVectorData::@5 data`
First data element of the current block.

6.4. File Documentation

6.4.1. File ziAPI.h

Header File for the LabOne C/C++ API.

Data Structures

- struct [ZIDoubleDataTS](#)
The structure used to hold a single IEEE double value. Same as [ZIDoubleData](#), but with timestamp.
- struct [ZIIntegerDataTS](#)
The structure used to hold a single 64bit signed integer value. Same as [ZIIntegerData](#), but with timestamp.
- struct [ZITreeChangeData](#)
The struct is holding info about added or removed nodes.
- struct [TreeChange](#)
The structure used to hold info about added or removed nodes. This is the version without timestamp used in API v1 compatibility mode.
- struct [ZIDemodSample](#)
The structure used to hold data for a single demodulator sample.
- struct [ZIAuxInSample](#)
The structure used to hold data for a single auxiliary inputs sample.
- struct [ZIDIOSample](#)
The structure used to hold data for a single digital I/O sample.
- struct [ZIByteArray](#)
The structure used to hold an arbitrary array of bytes. This is the version without timestamp used in API Level 1 compatibility mode.
- struct [ZIByteArrayTS](#)
The structure used to hold an arbitrary array of bytes. This is the same as [ZIByteArray](#), but with timestamp.
- struct [ZICntSample](#)
The structure used to hold data for a single counter sample.
- struct [ScopeWave](#)
The structure used to hold a single scope shot (API Level 1). If the client is connected to the Data Server using API Level 4 (recommended if supported by your device class) please see [ZIScopeWave](#) instead ([ZIScopeWaveEx](#) for API Level 5 and above).

- struct [ZIScopeWave](#)
The structure used to hold scope data. The data may be formatted differently, depending on settings. See the description of the structure members for details.
- struct [ZIScopeWaveEx](#)
The structure used to hold scope data (extended). The data may be formatted differently, depending on settings. See the description of the structure members for details.
- struct [ZIPWASample](#)
Single PWA sample value.
- struct [ZIPWAWave](#)
PWA Wave.
- struct [ZIImpedanceSample](#)
The structure used to hold data for a single impedance sample.
- struct [ZIStatisticSample](#)
- struct [ZISweeperDoubleSample](#)
- struct [ZISweeperDemodSample](#)
- struct [ZISweeperImpedanceSample](#)
- struct [ZISweeperHeader](#)
- struct [ZISweeperWave](#)
- struct [ZISpectrumDemodSample](#)
- struct [ZISpectrumHeader](#)
- struct [ZISpectrumWave](#)
- struct [ZIAdvisorSample](#)
- struct [ZIAdvisorHeader](#)
- struct [ZIAdvisorWave](#)

- struct [ZIVectorData](#)
The structure used to hold vector data block. See the description of the structure members for details.
- struct [ZIAsyncReply](#)
- struct [ZIEvent](#)
This struct holds event data forwarded by the Data Server.
- struct [ZISWTriggerHeader](#)
Structure to hold information about data returned from the SW Trigger Module.
- struct [ZIModuleHeaderSweeper](#)
Structure to hold information about data returned from the Sweep Module.
- struct [ZIModuleHeader](#)
Module-specific event header.
- struct [ZIModuleEvent](#)
This struct holds data of a single chunk from module lookup.
- struct [DemodSample](#)
The [DemodSample](#) struct holds data for the `ZI_DATA_DEMODSAMPLE` data type. Deprecated: See [ZIDemodSample](#).
- struct [AuxInSample](#)
The [AuxInSample](#) struct holds data for the `ZI_DATA_AUXINSAMPLE` data type. Deprecated: See [ZIAuxInSample](#).
- struct [DIOSample](#)
The [DIOSample](#) struct holds data for the `ZI_DATA_DIOSAMPLE` data type. Deprecated: See [ZIDIOSample](#).
- struct [ByteArrayData](#)
The [ByteArrayData](#) struct holds data for the `ZI_DATA_BYTEARRAY` data type. Deprecated: See [ZIByteArray](#).
- struct [ziEvent](#)
This struct holds event data forwarded by the Data Server. Deprecated: See [ZIEvent](#).
- union [ziEvent::Val](#)

Defines

- `#define MAX_PATH_LEN 256`

The maximum length that has to be used for passing paths to functions (including terminating zero)

- `#define MAX_EVENT_SIZE 0x400000`
The maximum size of an event's data block.
- `#define MAX_NAME_LEN 32`
The maximum length of the node name (in tree change event)

Typedefs

- `typedef ZIModuleHandle`
A handle with which to reference an instance of a `ziCore` Module created with [ziAPIModCreate](#).
- `typedef ZIConnection`
The `ZIConnection` is a connection reference; it holds information and helper variables about a connection to the Data Server. There is nothing in this reference which the user user may use, so it is hidden and instead a dummy pointer is used. See [ziAPIInit](#) for how to create a `ZIConnection`.
- `typedef ZIModuleEventPtr`
The pointer to a Module's data chunk to read out, updated via [ziAPIModGetChunk](#).

Enumerations

- `enum ZIResult_enum { ZI_INFO_BASE,`
`ZI_INFO_SUCCESS, ZI_INFO_MAX, ZI_WARNING_BASE,`
`ZI_WARNING_GENERAL, ZI_WARNING_UNDERRUN,`
`ZI_WARNING_OVERFLOW, ZI_WARNING_NOTFOUND,`
`ZI_WARNING_NO_ASYNC, ZI_WARNING_MAX,`
`ZI_ERROR_BASE, ZI_ERROR_GENERAL, ZI_ERROR_USB,`
`ZI_ERROR_MALLOC, ZI_ERROR_MUTEX_INIT,`
`ZI_ERROR_MUTEX_DESTROY, ZI_ERROR_MUTEX_LOCK,`
`ZI_ERROR_MUTEX_UNLOCK, ZI_ERROR_THREAD_START,`
`ZI_ERROR_THREAD_JOIN, ZI_ERROR_SOCKET_INIT,`
`ZI_ERROR_SOCKET_CONNECT, ZI_ERROR_HOSTNAME,`
`ZI_ERROR_CONNECTION, ZI_ERROR_TIMEOUT,`
`ZI_ERROR_COMMAND, ZI_ERROR_SERVER_INTERNAL,`
`ZI_ERROR_LENGTH, ZI_ERROR_FILE, ZI_ERROR_DUPLICATE,`
`ZI_ERROR_READONLY, ZI_ERROR_DEVICE_NOT_VISIBLE,`
`ZI_ERROR_DEVICE_IN_USE, ZI_ERROR_DEVICE_INTERFACE,`
`ZI_ERROR_DEVICE_CONNECTION_TIMEOUT,`
`ZI_ERROR_DEVICE_DIFFERENT_INTERFACE,`
`ZI_ERROR_DEVICE_NEEDS_FW_UPGRADE,`
`ZI_ERROR_ZIEVENT_DATATYPE_MISMATCH,`
`ZI_ERROR_DEVICE_NOT_FOUND,`
`ZI_ERROR_NOT_SUPPORTED,`
`ZI_ERROR_TOO_MANY_CONNECTIONS, ZI_ERROR_MAX,`
`ZI_SUCCESS, ZI_MAX_INFO, ZI_WARNING, ZI_UNDERRUN,`
`ZI_OVERFLOW, ZI_NOTFOUND, ZI_MAX_WARNING,`
`ZI_ERROR, ZI_USB, ZI_MALLOC, ZI_MUTEX_INIT,`
`ZI_MUTEX_DESTROY, ZI_MUTEX_LOCK, ZI_MUTEX_UNLOCK,`
`ZI_THREAD_START, ZI_THREAD_JOIN, ZI_SOCKET_INIT,`

```
ZI_SOCKET_CONNECT, ZI_HOSTNAME, ZI_CONNECTION,  
ZI_TIMEOUT, ZI_COMMAND, ZI_SERVER_INTERNAL,  
ZI_LENGTH, ZI_FILE, ZI_DUPLICATE, ZI_READONLY,  
ZI_MAX_ERROR }
```

Defines return value for all ziAPI functions. Divided into 3 regions: info, warning and error.

- enum [ZIValueType_enum](#) { ZI_VALUE_TYPE_NONE,
ZI_VALUE_TYPE_DOUBLE_DATA,
ZI_VALUE_TYPE_INTEGER_DATA,
ZI_VALUE_TYPE_DEMOD_SAMPLE,
ZI_VALUE_TYPE_SCOPE_WAVE_OLD,
ZI_VALUE_TYPE_AUXIN_SAMPLE,
ZI_VALUE_TYPE_DIO_SAMPLE,
ZI_VALUE_TYPE_BYTE_ARRAY, ZI_VALUE_TYPE_PWA_WAVE,
ZI_VALUE_TYPE_TREE_CHANGE_DATA_OLD,
ZI_VALUE_TYPE_DOUBLE_DATA_TS,
ZI_VALUE_TYPE_INTEGER_DATA_TS,
ZI_VALUE_TYPE_SCOPE_WAVE,
ZI_VALUE_TYPE_SCOPE_WAVE_EX,
ZI_VALUE_TYPE_BYTE_ARRAY_TS,
ZI_VALUE_TYPE_CNT_SAMPLE,
ZI_VALUE_TYPE_TREE_CHANGE_DATA,
ZI_VALUE_TYPE_ASYNC_REPLY,
ZI_VALUE_TYPE_SWEEPER_WAVE,
ZI_VALUE_TYPE_SPECTRUM_WAVE,
ZI_VALUE_TYPE_ADVISOR_WAVE,
ZI_VALUE_TYPE_VECTOR_DATA,
ZI_VALUE_TYPE_IMPEDANCE_SAMPLE,
ZI_DATA_NONE, ZI_DATA_DOUBLE, ZI_DATA_INTEGER,
ZI_DATA_DEMODSAMPLE, ZI_DATA_SCOPEWAVE,
ZI_DATA_AUXINSAMPLE, ZI_DATA_DIOSAMPLE,
ZI_DATA_BYTEARRAY, ZI_DATA_TREE_CHANGED }

Enumerates all types that data in a [ZIEvent](#) may have.

- enum [ZITreeAction_enum](#) { ZI_TREE_ACTION_REMOVE,
ZI_TREE_ACTION_ADD, ZI_TREE_ACTION_CHANGE }

Defines the actions that are performed on a tree,
as returned in the [ZITreeChangeData::action](#) or
[ZITreeChangeDataOld::action](#).

- enum [ZIImpFlags_enum](#) { ZI_IMP_FLAGS_NONE,
ZI_IMP_FLAGS_VALID_INTERNAL,
ZI_IMP_FLAGS_VALID_USER,
ZI_IMP_FLAGS_AUTORANGE_GATING,
ZI_IMP_FLAGS_OVERFLOW_VOLTAGE,
ZI_IMP_FLAGS_OVERFLOW_CURRENT,
ZI_IMP_FLAGS_UNDERFLOW_VOLTAGE,
ZI_IMP_FLAGS_UNDERFLOW_CURRENT,
ZI_IMP_FLAGS_FREQ_EXACT,
ZI_IMP_FLAGS_FREQ_INTERPOLATION,
ZI_IMP_FLAGS_FREQ_EXTRAPOLATION,
ZI_IMP_FLAGS_SUPPRESSION_PARAM0,
ZI_IMP_FLAGS_SUPPRESSION_PARAM1,
ZI_IMP_FLAGS_STRONGCOMPENSATION_PARAM0,
ZI_IMP_FLAGS_STRONGCOMPENSATION_PARAM1,
ZI_IMP_FLAGS_BWC_BIT0, ZI_IMP_FLAGS_BWC_BIT1,

```
ZI_IMP_FLAGS_BWC_BIT2, ZI_IMP_FLAGS_BWC_BIT3,  
ZI_IMP_FLAGS_BWC_MASK,  
ZI_IMP_FLAGS_OPEN_DETECTION }
```

Enumerates the bits set in an [ZIImpedanceSample](#)'s flags.

- `enum ZIVectorElementType_enum`
{ `ZI_VECTOR_ELEMENT_TYPE_UINT8`,
`ZI_VECTOR_ELEMENT_TYPE_UINT16`,
`ZI_VECTOR_ELEMENT_TYPE_UINT32`,
`ZI_VECTOR_ELEMENT_TYPE_UINT64`,
`ZI_VECTOR_ELEMENT_TYPE_FLOAT`,
`ZI_VECTOR_ELEMENT_TYPE_DOUBLE`,
`ZI_VECTOR_ELEMENT_TYPE_ASCII` }
Enumerates all the types that a `::elementType` may have.
- `enum ZIAPIVersion_enum` { `ZI_API_VERSION_0`,
`ZI_API_VERSION_1`, `ZI_API_VERSION_4`, `ZI_API_VERSION_5` }
- `enum ZIListNodes_enum` { `ZI_LIST_NODES_NONE`,
`ZI_LIST_NODES_RECURSIVE`, `ZI_LIST_NODES_ABSOLUTE`,
`ZI_LIST_NODES_LEAFONLY`,
`ZI_LIST_NODES_SETTINGSONLY`, `ZI_LIST_NONE`,
`ZI_LIST_RECURSIVE`, `ZI_LIST_ABSOLUTE`,
`ZI_LIST_LEAFONLY`, `ZI_LIST_SETTINGSONLY` }
Defines the values of the flags used in [ziAPIListNodes](#).
- `enum ZIModuleHeaderType_enum`
{ `ZI_MODULE_HEADER_TYPE_NONE`,
`ZI_MODULE_HEADER_TYPE_SWTRIGGER`,
`ZI_MODULE_HEADER_TYPE_SWEEPER` }
Enumerates all module header types.
- `enum ZIVectorWriteStatus_enum`
{ `ZI_VECTOR_WRITE_STATUS_IDLE`,
`ZI_VECTOR_WRITE_STATUS_PENDING` }
- `enum TREE_ACTION` { `TREE_ACTION_REMOVE`,
`TREE_ACTION_ADD`, `TREE_ACTION_CHANGE` }
`TREE_ACTION` defines the values for the [TreeChange::Action](#) Variable.

Functions

- `ZIResult_enum ziAPIInit (ZIConnection* conn)`
Initializes a [ZIConnection](#) structure.
- `ZIResult_enum ziAPIDestroy (ZIConnection conn)`
Destroys a [ZIConnection](#) structure.
- `ZIResult_enum ziAPIConnect (ZIConnection conn, const
char* hostname, uint16_t port)`
Connects the [ZIConnection](#) to Data Server.
- `ZIResult_enum ziAPIDisconnect (ZIConnection conn)`

Disconnects an established connection.

- `ZIResult_enum` `ziAPIListImplementations (char* implementations, uint32_t bufferSize)`
Returns the list of supported implementations.
- `ZIResult_enum` `ziAPIConnectEx (ZIConnection conn, const char* hostname, uint16_t port, ZIAPIVersion_enum apiLevel, const char* implementation)`
Connects to Data Server and enables extended ziAPI.
- `ZIResult_enum` `ziAPIGetConnectionAPILevel (ZIConnection conn, ZIAPIVersion_enum* apiLevel)`
Returns ziAPI level used for the connection conn.
- `ZIResult_enum` `ziAPIGetRevision (unsigned int* revision)`
Retrieves the revision of ziAPI.
- `ZIResult_enum` `ziAPIListNodes (ZIConnection conn, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)`
Returns all child nodes found at the specified path.
- `ZIResult_enum` `ziAPIUpdateDevices (ZIConnection conn)`
Search for the newly connected devices and update the tree.
- `ZIResult_enum` `ziAPIConnectDevice (ZIConnection conn, const char* deviceSerial, const char* deviceInterface, const char* interfaceParams)`
Connect a device to the server.
- `ZIResult_enum` `ziAPIDisconnectDevice (ZIConnection conn, const char* deviceSerial)`
Disconnect a device from the server.
- `ZIResult_enum` `ziAPIGetValueD (ZIConnection conn, const char* path, ZIDoubleData* value)`
gets the double-type value of the specified node
- `ZIResult_enum` `ziAPIGetValueI (ZIConnection conn, const char* path, ZIIntegerData* value)`
gets the integer-type value of the specified node
- `ZIResult_enum` `ziAPIGetDemodSample (ZIConnection conn, const char* path, ZIDemodSample* value)`
Gets the demodulator sample value of the specified node.
- `ZIResult_enum` `ziAPIGetDIOSample (ZIConnection conn, const char* path, ZIDIOSample* value)`
Gets the Digital I/O sample of the specified node.
- `ZIResult_enum` `ziAPIGetAuxInSample (ZIConnection conn, const char* path, ZIAuxInSample* value)`
gets the AuxIn sample of the specified node
- `ZIResult_enum` `ziAPIGetValueB (ZIConnection conn, const char* path, unsigned char* buffer, unsigned int* length, unsigned int bufferSize)`

gets the ByteArray value of the specified node

- **ZIResult_enum** `ziAPISetValueD (ZIConnection conn, const char* path, ZIDoubleData value)`
asynchronously sets a double-type value to one or more nodes specified in the path
- **ZIResult_enum** `ziAPISetValueI (ZIConnection conn, const char* path, ZIIntegerData value)`
asynchronously sets an integer-type value to one or more nodes specified in a path
- **ZIResult_enum** `ziAPISetValueB (ZIConnection conn, const char* path, unsigned char* buffer, unsigned int length)`
asynchronously sets the binary-type value of one ore more nodes specified in the path
- **ZIResult_enum** `ziAPISyncSetValueD (ZIConnection conn, const char* path, ZIDoubleData* value)`
synchronously sets a double-type value to one or more nodes specified in the path
- **ZIResult_enum** `ziAPISyncSetValueI (ZIConnection conn, const char* path, ZIIntegerData* value)`
synchronously sets an integer-type value to one or more nodes specified in a path
- **ZIResult_enum** `ziAPISyncSetValueB (ZIConnection conn, const char* path, uint8_t* buffer, uint32_t* length, uint32_t bufferSize)`
Synchronously sets the binary-type value of one ore more nodes specified in the path.
- **ZIResult_enum** `ziAPISync (ZIConnection conn)`
Synchronizes the session by dropping all pending data.
- **ZIResult_enum** `ziAPIEchoDevice (ZIConnection conn, const char* deviceSerial)`
Sends an echo command to a device and blocks until answer is received.
- **ZIEvent*** `ziAPIAllocateEventEx ()`
Allocates **ZIEvent** structure and returns the pointer to it. Attention!!! It is the client code responsibility to deallocate the structure by calling `ziAPIDeallocateEventEx!`
- `void ziAPIDeallocateEventEx (ZIEvent* ev)`
Deallocates **ZIEvent** structure created with `ziAPIAllocateEventEx()`.
- **ZIResult_enum** `ziAPISubscribe (ZIConnection conn, const char* path)`
subscribes the nodes given by path for `ziAPIPollDataEx`
- **ZIResult_enum** `ziAPIUnSubscribe (ZIConnection conn, const char* path)`

unsubscribes to the nodes given by path

- `ZIResult_enum` `ziAPIPollDataEx (ZIConnection conn, ZIEvent* ev, uint32_t timeOutMilliseconds)`
checks if an event is available to read
- `ZIResult_enum` `ziAPIGetValueAsPollData (ZIConnection conn, const char* path)`
triggers a value request, which will be given back on the poll event queue
- `ZIResult_enum` `ziAPIAsyncSetDoubleData (ZIConnection conn, const char* path, ZIDoubleData value)`
- `ZIResult_enum` `ziAPIAsyncSetIntegerData (ZIConnection conn, const char* path, ZIntegerData value)`
- `ZIResult_enum` `ziAPIAsyncSetByteArray (ZIConnection conn, const char* path, uint8_t* buffer, uint32_t length)`
- `ZIResult_enum` `ziAPIAsyncSubscribe (ZIConnection conn, const char* path, ZIAsyncTag tag)`
- `ZIResult_enum` `ziAPIAsyncUnSubscribe (ZIConnection conn, const char* path, ZIAsyncTag tag)`
- `ZIResult_enum` `ziAPIAsyncGetValueAsPollData (ZIConnection conn, const char* path, ZIAsyncTag tag)`
- `ZIResult_enum` `ziAPIGetError (ZIResult_enum result, char** buffer, int* base)`
Returns a description and the severity for a `ZIResult_enum`.
- `ZIResult_enum` `ziAPIGetLastError (ZIConnection conn, char* buffer, uint32_t bufferSize)`
Returns the message from the last error that occurred.
- `void` `ziAPISetDebugLevel (int32_t debugLevel)`
Enable ziAPI's log and set the severity level of entries to be included in the log.
- `void` `ziAPIWriteDebugLog (int32_t debugLevel, const char* message)`
Write a message to ziAPI's log with the specified severity.
- `ZIResult_enum` `ReadMEMFile (const char* filename, char* buffer, int32_t bufferSize, int32_t* bytesUsed)`
- `ZIResult_enum` `ziAPIModCreate (ZIConnection conn, ZIModuleHandle* handle, const char* moduleId)`

Create a [ZIModuleHandle](#) that can be used for asynchronous measurement tasks.

- [ZIResult_enum](#) `ziAPIModSetDoubleData (ZIConnection conn, ZIModuleHandle handle, const char* path, ZIDoubleData value)`
Sets a module parameter to the specified double type.
- [ZIResult_enum](#) `ziAPIModSetIntegerData (ZIConnection conn, ZIModuleHandle handle, const char* path, ZIIntegerData value)`
Sets a module parameter to the specified integer type.
- [ZIResult_enum](#) `ziAPIModSetByteArray (ZIConnection conn, ZIModuleHandle handle, const char* path, uint8_t* buffer, uint32_t length)`
Sets a module parameter to the specified byte array.
- [ZIResult_enum](#) `ziAPIModListNodes (ZIConnection conn, ZIModuleHandle handle, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)`
Returns all child parameter node paths found under the specified parent module parameter path.
- [ZIResult_enum](#) `ziAPIModSubscribe (ZIConnection conn, ZIModuleHandle handle, const char* path)`
Subscribes to the nodes specified by path, these nodes will be recorded during module execution.
- [ZIResult_enum](#) `ziAPIModUnSubscribe (ZIConnection conn, ZIModuleHandle handle, const char* path)`
Unsubscribes to the nodes specified by path.
- [ZIResult_enum](#) `ziAPIModExecute (ZIConnection conn, ZIModuleHandle handle)`
Starts the module's thread and its associated measurement task.
- [ZIResult_enum](#) `ziAPIModTrigger (ZIConnection conn, ZIModuleHandle handle)`
Manually issue a trigger forcing data recording (SW Trigger Module only).
- [ZIResult_enum](#) `ziAPIModProgress (ZIConnection conn, ZIModuleHandle handle, ZIDoubleData* progress)`
Queries the current state of progress of the module's measurement task.
- [ZIResult_enum](#) `ziAPIModFinished (ZIConnection conn, ZIModuleHandle handle, ZIIntegerData* finished)`
Queries whether the module has finished its measurement task.
- [ZIResult_enum](#) `ziAPIModFinish (ZIConnection conn, ZIModuleHandle handle)`
Stops the module performing its measurement task.

- `ZIResult_enum` `ziAPIModSave (ZIConnection conn, ZIModuleHandle handle, const char* fileName)`
Saves the currently accumulated data to file.
- `ZIResult_enum` `ziAPIModRead (ZIConnection conn, ZIModuleHandle handle, const char* path)`
Make the currently accumulated data available for use in the C program.
- `ZIResult_enum` `ziAPIModNextNode (ZIConnection conn, ZIModuleHandle handle, char* path, uint32_t bufferSize, ZIValueType_enum* valueType, uint64_t* chunks)`
Make the data for the next node available for reading with `ziAPIModGetChunk`.
- `ZIResult_enum` `ziAPIModGetChunk (ZIConnection conn, ZIModuleHandle handle, uint64_t chunkIndex, ZIModuleEventPtr* ev)`
Get the specified data chunk from the current node.
- `ZIResult_enum` `ziAPIModEventDeallocate (ZIConnection conn, ZIModuleHandle handle, ZIModuleEventPtr ev)`
Deallocate the `ZIModuleEventPtr` being used by the module.
- `ZIResult_enum` `ziAPIModClear (ZIConnection conn, ZIModuleHandle handle)`
Terminates the module's thread and destroys the module.
- `ZIResult_enum` `ziAPIVectorWriteBlock (ZIConnection conn, const char* path, ZIVectorData* vectorBlock)`
- `ZIResult_enum` `ziAPIVectorWriteGetStatus (ZIConnection conn, const char* path, uint8_t* status)`
status - see `ZIVectorWriteStatus_enum`
- `ZIResult_enum` `ziAPIVectorWrite (ZIConnection conn, const char* path, const void* vectorPtr, uint8_t vectorElementType, uint64_t vectorSizeElements)`
vectorElementType - see `ZIVectorElementType_enum`
- `ZIResult_enum` `ziAPIDiscoveryFind (ZIConnection conn, const char* deviceAddress, const char** deviceId)`
Returns the device id for a given device address.
Attention! Invalidates all pointers previously returned by `ziAPIDiscovery*` calls.
- `ZIResult_enum` `ziAPIDiscoveryGet (ZIConnection conn, const char* deviceId, const char** propsJSON)`
Returns the device properties for a given device id in JSON format.
- `ZIResult_enum` `ziAPIDiscoveryGetValue1 (ZIConnection conn, const char* deviceId, const char* propName, ZIIntegerData* value)`
Returns given integer property value for a given device id.

- `ZIResult_enum` `ziAPIDiscoveryGetValueS (ZIConnection conn, const char* deviceId, const char* propName, const char** value)`
Returns given string property value for a given device id.
- `__inline ziEvent* ziAPIAllocateEvent ()`
Deprecated: See `ziAPIAllocateEventEx()`.
- `__inline void ziAPIDeallocateEvent (ziEvent* ev)`
Deprecated: See `ziAPIDeallocateEventEx()`.
- `__inline ZIResult_enum` `ziAPIPollData (ZIConnection conn, ziEvent* ev, int timeout)`
Checks if an event is available to read. Deprecated: See `ziAPIPollDataEx()`.
- `__inline ZIResult_enum` `ziAPIGetValueS (ZIConnection conn, char* path, DemodSample* value)`
- `__inline ZIResult_enum` `ziAPIGetValueDIO (ZIConnection conn, char* path, DIOSample* value)`
- `__inline ZIResult_enum` `ziAPIGetValueAuxIn (ZIConnection conn, char* path, AuxInSample* value)`
- `double` `ziAPISecondsTimeStamp (ziTimeStampType TS)`

Detailed Description

ziAPI provides all functionality to establish a connection with the Data Server and to communicate with it. It has functions for setting and getting parameters in a single call as well as an event-framework with which the user may subscribe the parameter tree and receive the events which occur when values change.

- All functions do not check passed pointers if they're NULL pointers. In that case a segmentation fault will occur.
- The ZIConnection is not thread-safe. One connection can only be used in one thread. If you want to use the ziAPI in a multi-threaded program you will have to use one ZIConnection for each thread that is communicating or implement a mutual exclusion.
- The Data Server is able to handle connections from threads simultaneously. The Data Server takes over the synchronization.

Data Structure Documentation

struct ZIDoubleDataTS

The structure used to hold a single IEEE double value. Same as ZIDoubleData, but with timestamp.

```
#include "ziAPI.h"

typedef struct ZIDoubleDataTS {
    ZITimeStamp timeStamp;
    ZIDoubleData value;
} ZIDoubleDataTS;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the value has changed.
- ZIDoubleData value

struct ZIIntegerDataTS

The structure used to hold a single 64bit signed integer value. Same as ZIIntegerData, but with timestamp.

```
#include "ziAPI.h"

typedef struct ZIIntegerDataTS {
    ZITimeStamp timeStamps;
    ZIIntegerData value;
} ZIIntegerDataTS;
```

Data Fields

- ZITimeStamp timeStamps
Time stamp at which the value has changed.
- ZIIntegerData value

struct ZITreeChangeData

The struct is holding info about added or removed nodes.

```
#include "ziAPI.h"

typedef struct ZITreeChangeData {
    ZITimeStamp timeStamp;
    uint32_t action;
    char name[32];
} ZITreeChangeData;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- uint32_t action
field indicating which action occurred on the tree. A value of the ZITreeAction_enum.
- char name
Name of the Path that has been added, removed or changed.

struct TreeChange

The structure used to hold info about added or removed nodes. This is the version without timestamp used in API v1 compatibility mode.

```
#include "ziAPI.h"

typedef struct TreeChange {
    uint32_t Action;
    char Name[32];
} TreeChange;
```

Data Fields

- uint32_t Action
field indicating which action occurred on the tree. A value of the ZITreeAction_enum (TREE_ACTION) enum.
- char Name
Name of the Path that has been added, removed or changed.

struct ZIDemodSample

The structure used to hold data for a single demodulator sample.

```
#include "ziAPI.h"

typedef struct ZIDemodSample {
    ZITimeStamp timeStamp;
    double x;
    double y;
    double frequency;
    double phase;
    uint32_t dioBits;
    uint32_t trigger;
    double auxIn0;
    double auxIn1;
} ZIDemodSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the sample has been measured.
- double x
X part of the sample.
- double y
Y part of the sample.
- double frequency
Frequency at that sample.
- double phase
Phase at that sample.
- uint32_t dioBits
the current bits of the DIO.
- uint32_t trigger
trigger bits
- double auxIn0
value of Aux input 0.
- double auxIn1
value of Aux input 1.

struct ZIAuxInSample

The structure used to hold data for a single auxiliary inputs sample.

```
#include "ziAPI.h"

typedef struct ZIAuxInSample {
    ZITimeStamp timeStamp;
    double ch0;
    double ch1;
} ZIAuxInSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the values have been measured.
- double ch0
Channel 0 voltage.
- double ch1
Channel 1 voltage.

struct ZIDIOSample

The structure used to hold data for a single digital I/O sample.

```
#include "ziAPI.h"

typedef struct ZIDIOSample {
    ZITimeStamp timeStamp;
    uint32_t bits;
    uint32_t reserved;
} ZIDIOSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the values have been measured.
- uint32_t bits
The digital I/O values.
- uint32_t reserved
Filler to keep 8 bytes alignment in the array of [ZIDIOSample](#) structures.

struct ZIByteArray

The structure used to hold an arbitrary array of bytes. This is the version without timestamp used in API Level 1 compatibility mode.

```
#include "ziAPI.h"

typedef struct ZIByteArray {
    uint32_t length;
    uint8_t bytes[0];
} ZIByteArray;
```

Data Fields

- `uint32_t length`
Length of the data readable from the Bytes field.
- `uint8_t bytes`
The data itself. The array has the size given in length.

struct ZIByteArrayTS

The structure used to hold an arbitrary array of bytes. This is the same as [ZIByteArray](#), but with timestamp.

```
#include "ziAPI.h"

typedef struct ZIByteArrayTS {
    ZITimeStamp timeStamp;
    uint32_t length;
    uint8_t bytes[0];
} ZIByteArrayTS;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- uint32_t length
length of the data readable from the bytes field
- uint8_t bytes
the data itself. The array has the size given in length

struct ZICntSample

The structure used to hold data for a single counter sample.

```
#include "ziAPI.h"

typedef struct ZICntSample {
    ZITimeStamp timeStamp;
    uint16_t counter;
    uint16_t reserved0;
    uint16_t id;
    uint16_t reserved1;
} ZICntSample;
```

Data Fields

- ZITimeStamp timeStamp
The timestamp at which the values have been measured.
- uint16_t counter
Counter value.
- uint16_t reserved0
Reserved.
- uint16_t id
Trigger id.
- uint16_t reserved1
Reserved.

struct ScopeWave

The structure used to hold a single scope shot (API Level 1). If the client is connected to the Data Server using API Level 4 (recommended if supported by your device class) please see [ZIScopeWave](#) instead ([ZIScopeWaveEx](#) for API Level 5 and above).

```
#include "ziAPI.h"

typedef struct ScopeWave {
    double dt;
    uint32_t ScopeChannel;
    uint32_t TriggerChannel;
    uint32_t BWLimit;
    uint32_t Count;
    int16_t Data[0];
} ScopeWave;
```

Data Fields

- double dt
Time difference between samples.
- uint32_t ScopeChannel
Scope channel of the represented data.
- uint32_t TriggerChannel
Trigger channel of the represented data.
- uint32_t BWLimit
Bandwidth-limit flag.
- uint32_t Count
Count of samples.
- int16_t Data
First wave data.

struct ZIScopeWave

The structure used to hold scope data. The data may be formatted differently, depending on settings. See the description of the structure members for details.

```
#include "ziAPI.h"

typedef struct ZIScopeWave {
    ZITimeStamp timeStamp;
    ZITimeStamp triggerTimeStamp;
    double dt;
    uint8_t channelEnable[4];
    uint8_t channelInput[4];
    uint8_t triggerEnable;
    uint8_t triggerInput;
    uint8_t reserved0[2];
    uint8_t channelBWLimit[4];
    uint8_t channelMath[4];
    float channelScaling[4];
    uint32_t sequenceNumber;
    uint32_t segmentNumber;
    uint32_t blockNumber;
    uint64_t totalSamples;
    uint8_t dataTransferMode;
    uint8_t blockMarker;
    uint8_t flags;
    uint8_t sampleFormat;
    uint32_t sampleCount;
    int16_t dataInt16[0];
    int32_t dataInt32[0];
    float dataFloat[0];
    union ZIScopeWave::@0 data;
} ZIScopeWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp of the last sample in this data block.
- ZITimeStamp triggerTimeStamp
Time stamp of the trigger (may also fall between samples and in another block)
- double dt
Time difference between samples in seconds.
- uint8_t channelEnable
Up to four channels: if channel is enabled, corresponding element is non-zero.
- uint8_t channelInput
Specifies the input source for each of the scope four channels: 0 = Signal Input 1, 1 = Signal Input 2, 2 = Trigger Input 1, 3 = Trigger Input 2, 4 = Aux Output 1, 5 = Aux Output 2, 6 = Aux Output 3, 7 = Aux Output 4, 8 = Aux Input 1, 9 = Aux Input 2.
- uint8_t triggerEnable
Non-zero if trigger is enabled: Bit(0): rising edge trigger off = 0, on = 1. Bit(1): falling edge trigger off = 0, on = 1.

- `uint8_t triggerInput`
Trigger source (same values as for channel input)
- `uint8_t reserved0`
- `uint8_t channelBWLimit`
Bandwidth-limit flag, per channel. Bit(0): off = 0, on = 1
Bit(7...1): Reserved.
- `uint8_t channelMath`
Math Operation (e.g averaging) Bit (7..0): Reserved.
- `float channelScaling`
Data scaling factors for up to 4 channels.
- `uint32_t sequenceNumber`
Current scope shot sequence number. Identifies a scope shot.
- `uint32_t segmentNumber`
Current segment number.
- `uint32_t blockNumber`
Current block number from the beginning of a scope shot.
Large scope shots are split into blocks, which need to be concatenated to obtain the complete scope shot.
- `uint64_t totalSamples`
Total number of samples in one channel in the current scope shot, same for all channels.
- `uint8_t dataTransferMode`
Data transfer mode SingleTransfer = 0, BlockTransfer = 1, ContinuousTransfer = 3, FFTSingleTransfer = 4. Other values are reserved.
- `uint8_t blockMarker`
Block marker: Bit (0): 1 = End marker for continuous or multi-block transfer Bit (7..0): Reserved.
- `uint8_t flags`
Indicator Flags. Bit (0): 1 = Data loss detected (samples are 0), Bit (1): 1 = Missed trigger, Bit (2): 1 = Transfer failure (corrupted data).
- `uint8_t sampleFormat`
Data format of samples: Int16 = 0, Int32 = 1, Float = 2, Int16Interleaved = 4, Int32Interleaved = 5, FloatInterleaved = 6.
- `uint32_t sampleCount`
Number of samples in one channel in the current block, same for all channels.
- `int16_t dataInt16`

- `int32_t dataInt32`
- `float dataFloat`
- `union ZIScopeWave::@0 data`
First wave data.

struct ZIScopeWaveEx

The structure used to hold scope data (extended). The data may be formatted differently, depending on settings. See the description of the structure members for details.

```
#include "ziAPI.h"

typedef struct ZIScopeWaveEx {
    ZITimeStamp timeStamp;
    ZITimeStamp triggerTimeStamp;
    double dt;
    uint8_t channelEnable[4];
    uint8_t channelInput[4];
    uint8_t triggerEnable;
    uint8_t triggerInput;
    uint8_t reserved0[2];
    uint8_t channelBWLimit[4];
    uint8_t channelMath[4];
    float channelScaling[4];
    uint32_t sequenceNumber;
    uint32_t segmentNumber;
    uint32_t blockNumber;
    uint64_t totalSamples;
    uint8_t dataTransferMode;
    uint8_t blockMarker;
    uint8_t flags;
    uint8_t sampleFormat;
    uint32_t sampleCount;
    double channelOffset[4];
    uint64_t reserved1[32];
    int16_t dataInt16[0];
    int32_t dataInt32[0];
    float dataFloat[0];
    union ZIScopeWaveEx::@1 data;
} ZIScopeWaveEx;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp of the last sample in this data block.
- ZITimeStamp triggerTimeStamp
Time stamp of the trigger (may also fall between samples and in another block)
- double dt
Time difference between samples in seconds.
- uint8_t channelEnable
Up to four channels: if channel is enabled, corresponding element is non-zero.
- uint8_t channelInput
Specifies the input source for each of the scope four channels: 0 = Signal Input 1, 1 = Signal Input 2, 2 = Trigger Input 1, 3 = Trigger Input 2, 4 = Aux Output 1, 5 = Aux Output 2, 6 = Aux Output 3, 7 = Aux Output 4, 8 = Aux Input 1, 9 = Aux Input 2.
- uint8_t triggerEnable
Non-zero if trigger is enabled: Bit(0): rising edge trigger off = 0, on = 1. Bit(1): falling edge trigger off = 0, on = 1.

- `uint8_t triggerInput`
Trigger source (same values as for channel input)
- `uint8_t reserved0`
- `uint8_t channelBWLimit`
Bandwidth-limit flag, per channel. Bit(0): off = 0, on = 1
Bit(7...1): Reserved.
- `uint8_t channelMath`
Math Operation (e.g averaging) Bit (7..0): Reserved.
- `float channelScaling`
Data scaling factors for up to 4 channels.
- `uint32_t sequenceNumber`
Current scope shot sequence number. Identifies a scope shot.
- `uint32_t segmentNumber`
Current segment number.
- `uint32_t blockNumber`
Current block number from the beginning of a scope shot.
Large scope shots are split into blocks, which need to be concatenated to obtain the complete scope shot.
- `uint64_t totalSamples`
Total number of samples in one channel in the current scope shot, same for all channels.
- `uint8_t dataTransferMode`
Data transfer mode SingleTransfer = 0, BlockTransfer = 1, ContinuousTransfer = 3, FFTSingleTransfer = 4. Other values are reserved.
- `uint8_t blockMarker`
Block marker: Bit (0): 1 = End marker for continuous or multi-block transfer Bit (7..0): Reserved.
- `uint8_t flags`
Indicator Flags. Bit (0): 1 = Data loss detected (samples are 0), Bit (1): 1 = Missed trigger, Bit (2): 1 = Transfer failure (corrupted data).
- `uint8_t sampleFormat`
Data format of samples: Int16 = 0, Int32 = 1, Float = 2, Int16Interleaved = 4, Int32Interleaved = 5, FloatInterleaved = 6.
- `uint32_t sampleCount`
Number of samples in one channel in the current block, same for all channels.
- `double channelOffset`

Data offset (scaled) for up to 4 channels.

- `uint64_t reserved1`

- `int16_t dataInt16`

- `int32_t dataInt32`

- `float dataFloat`

- `union ZIScopeWaveEx::@1 data`
First wave data.

struct ZIPWASample

Single PWA sample value.

```
#include "ziAPI.h"

typedef struct ZIPWASample {
    double binPhase;
    double x;
    double y;
    uint32_t countBin;
    uint32_t reserved;
} ZIPWASample;
```

Data Fields

- double binPhase
Phase position of each bin.
- double x
Real PWA result or X component of a demod PWA.
- double y
Y component of the demod PWA.
- uint32_t countBin
Number of events per bin.
- uint32_t reserved
Reserved.

struct ZIPWAWave

PWA Wave.

```
#include "ziAPI.h"

typedef struct ZIPWAWave {
    ZITimeStamp timeStamp;
    uint64_t sampleCount;
    uint32_t inputSelect;
    uint32_t oscSelect;
    uint32_t harmonic;
    uint32_t binCount;
    double frequency;
    uint8_t pwaType;
    uint8_t mode;
    uint8_t overflow;
    uint8_t commensurable;
    uint32_t reservedUInt;
    ZIPWASample
        data[0];
} ZIPWAWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- uint64_t sampleCount
Total sample count considered for PWA.
- uint32_t inputSelect
Input selection used for the PWA.
- uint32_t oscSelect
Oscillator used for the PWA.
- uint32_t harmonic
Harmonic setting.
- uint32_t binCount
Bin count of the PWA.
- double frequency
Frequency during PWA accumulation.
- uint8_t pwaType
Type of the PWA.
- uint8_t mode
PWA Mode [0: zoom PWA, 1: harmonic PWA].
- uint8_t overflow
Overflow indicators. overflow[0]: Data accumulator overflow, overflow[1]: Counter at limit, overflow[6..2]: Reserved, overflow[7]: Invalid (missing frames).
- uint8_t commensurable
Commensurability of the data.

- `uint32_t reservedUInt`
Reserved 32bit.
- `ZIPWASample` data
PWA data vector.

struct ZIImpedanceSample

The structure used to hold data for a single impedance sample.

```
#include "ziAPI.h"

typedef struct ZIImpedanceSample {
    ZITimeStamp timeStamp;
    double realz;
    double imagz;
    double frequency;
    double phase;
    uint32_t flags;
    uint32_t trigger;
    double param0;
    double param1;
    double drive;
    double bias;
} ZIImpedanceSample;
```

Data Fields

- ZITimeStamp timeStamp
Timestamp at which the sample has been measured.
- double realz
Real part of the impedance sample.
- double imagz
Imaginary part of the impedance sample.
- double frequency
Frequency at that sample.
- double phase
Phase at that sample.
- uint32_t flags
Flags (see ZIImpFlags_enum)
- uint32_t trigger
Trigger bits.
- double param0
Value of model parameter 0.
- double param1
Value of model parameter 1.
- double drive
Drive amplitude.
- double bias
Bias voltage.

struct ZIStatisticSample

```
typedef struct ZIStatisticSample {  
    double avg;  
    double stddev;  
    double pwr;  
} ZIStatisticSample;
```

Data Fields

- double avg
Average value or single value.
- double stddev
Standard deviation.
- double pwr
Power value.

struct ZISweeperDoubleSample

```
typedef struct ZISweeperDoubleSample {  
    double grid;  
    double bandwidth;  
    uint64_t count;  
    ZIStatisticSample value;  
} ZISweeperDoubleSample;
```

Data Fields

- double grid
Grid value (x-axis)
- double bandwidth
Bandwidth.
- uint64_t count
Sample count used for statistic calculation.
- ZIStatisticSample value
Result value (y-axis)

struct ZISweeperDemodSample

```
typedef struct ZISweeperDemodSample {  
    double grid;  
    double bandwidth;  
    uint64_t count;  
    double tc;  
    double tcMeas;  
    double settling;  
    ZITimeStamp setTimeStamp;  
    ZITimeStamp nextTimeStamp;  
    ZIStatisticSample x;  
    ZIStatisticSample y;  
    ZIStatisticSample r;  
    ZIStatisticSample phase;  
    ZIStatisticSample frequency;  
    ZIStatisticSample auxin0;  
    ZIStatisticSample auxin1;  
} ZISweeperDemodSample;
```

Data Fields

- double grid
Grid value (x-axis)
- double bandwidth
Demodulator bandwidth used for the specific sweep point.
- uint64_t count
Sample count used for statistic calculation.
- double tc
Time constant calculated for the specific sweep point.
- double tcMeas
Time constant used by the device.
- double settling
Settling time (s) used to wait until averaging operation is started.
- ZITimeStamp setTimeStamp
Time stamp when the grid value was set on the device.
- ZITimeStamp nextTimeStamp
Time stamp when the first statistic value was recorded.
- ZIStatisticSample x
Sweep point statistic result of X.
- ZIStatisticSample y
Sweep point statistic result of Y.
- ZIStatisticSample r
Sweep point statistic result of R.
- ZIStatisticSample phase
Sweep point statistic result of phase.

- ZIStatisticSample frequency
Sweep point statistic result of frequency.
- ZIStatisticSample auxin0
Sweep point statistic result of auxin0.
- ZIStatisticSample auxin1
Sweep point statistic result of auxin1.

struct ZISweeperImpedanceSample

```
typedef struct ZISweeperImpedanceSample {  
    double grid;  
    double bandwidth;  
    uint64_t count;  
    double tc;  
    double tcMeas;  
    double settling;  
    ZITimeStamp setTimeStamp;  
    ZITimeStamp nextTimeStamp;  
    ZIStatisticSample realz;  
    ZIStatisticSample imagz;  
    ZIStatisticSample absz;  
    ZIStatisticSample phasez;  
    ZIStatisticSample frequency;  
    ZIStatisticSample param0;  
    ZIStatisticSample param1;  
    ZIStatisticSample drive;  
    ZIStatisticSample bias;  
} ZISweeperImpedanceSample;
```

Data Fields

- double grid
Grid value (x-axis)
- double bandwidth
Demodulator bandwidth used for the specific sweep point.
- uint64_t count
Sample count used for statistic calculation.
- double tc
Time constant calculated for the specific sweep point.
- double tcMeas
Time constant used by the device.
- double settling
Settling time (s) used to wait until averaging operation is started.
- ZITimeStamp setTimeStamp
Time stamp when the grid value was set on the device.
- ZITimeStamp nextTimeStamp
Time stamp when the first statistic value was recorded.
- ZIStatisticSample realz
Sweep point statistic result of X.
- ZIStatisticSample imagz
Sweep point statistic result of Y.
- ZIStatisticSample absz
Sweep point statistic result of R.
- ZIStatisticSample phasez

Sweep point statistic result of phase.

- ZIStatisticSample frequency
Sweep point statistic result of frequency.

- ZIStatisticSample param0
Sweep point statistic result of param0.

- ZIStatisticSample param1
Sweep point statistic result of param1.

- ZIStatisticSample drive
Sweep point statistic result of drive amplitude.

- ZIStatisticSample bias
Sweep point statistic result of bias.

struct ZISweeperHeader

```
typedef struct ZISweeperHeader {  
    uint64_t sampleCount;  
    uint8_t flags;  
    uint8_t sampleFormat;  
    uint8_t sweepMode;  
    uint8_t bandwidthMode;  
    uint8_t reserved0[4];  
    uint8_t reserved1[8];  
} ZISweeperHeader;
```

Data Fields

- `uint64_t sampleCount`
Total sample count considered for sweeper.
- `uint8_t flags`
Flags Bit 0: Phase unwrap Bit 1: Sinc filter.
- `uint8_t sampleFormat`
Sample format Double = 0, Demodulator = 1.
- `uint8_t sweepMode`
Sweep mode Sequential = 0, Binary = 1, Bidirectional = 2, Reverse = 3.
- `uint8_t bandwidthMode`
Bandwidth mode Manual = 0, Fixed = 1, Auto = 2.
- `uint8_t reserved0`
Reserved space for future use.
- `uint8_t reserved1`
Reserved space for future use.

struct ZISweeperWave

```
typedef struct ZISweeperWave {  
    ZITimeStamp timeStamp;  
    ZISweeperHeader header;  
    ZISweeperDoubleSample dataDouble[0];  
    ZISweeperDemodSample dataDemod[0];  
    ZISweeperImpedanceSample dataImpedance[0];  
    union ZISweeperWave::@2 data;  
} ZISweeperWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- ZISweeperHeader header
- ZISweeperDoubleSample dataDouble
- ZISweeperDemodSample dataDemod
- ZISweeperImpedanceSample dataImpedance
- union ZISweeperWave::@2 data
Sweeper data vector.

struct ZISpectrumDemodSample

```
typedef struct ZISpectrumDemodSample {  
    double grid;  
    double filter;  
    double x;  
    double y;  
    double r;  
} ZISpectrumDemodSample;
```

Data Fields

- double grid
Grid.
- double filter
Filter strength at the specific grid point.
- double x
X.
- double y
Y.
- double r
R.

struct ZISpectrumHeader

```
typedef struct ZISpectrumHeader {
    uint64_t sampleCount;
    uint8_t flags;
    uint8_t sampleFormat;
    uint8_t spectrumMode;
    uint8_t window;
    uint8_t reserved0[4];
    uint8_t reserved1[8];
    double bandwidth;
    double rate;
    double center;
    double resolution;
    double aliasingReject;
    double nenbw;
    double overlap;
} ZISpectrumHeader;
```

Data Fields

- `uint64_t sampleCount`
Total sample count considered for spectrum.
- `uint8_t flags`
Flags Bit 0: Power Bit 1: Spectral density Bit 2: Absolute frequency Bit 3: Full span.
- `uint8_t sampleFormat`
Sample format Demodulator = 0.
- `uint8_t spectrumMode`
Spectrum mode FFT(x+iy) = 0, FFT(r) = 1, FFT(theta) = 2, FFT(freq) = 3, FFT(dtheta/dt)/2pi = 4.
- `uint8_t window`
Window Rectangular = 0, Hann = 1, Hamming = 2, Blackman Harris = 3.
- `uint8_t reserved0`
Reserved space for future use.
- `uint8_t reserved1`
Reserved space for future use.
- `double bandwidth`
Filter bandwidth.
- `double rate`
Rate of the sampled data.
- `double center`
FFT center value.
- `double resolution`
FFT bin resolution.
- `double aliasingReject`

Aliasing reject (dB)

- double nenbw
Correction factor for the used window when calculating spectral density.
- double overlap
FFT overlap [0 .. 1[.

struct ZISpectrumWave

```
typedef struct ZISpectrumWave {  
    ZITimeStamp timeStamp;  
    ZISpectrumHeader header;  
    ZISpectrumDemodSample dataDemod[0];  
    union ZISpectrumWave::@3 data;  
} ZISpectrumWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- ZISpectrumHeader header
- ZISpectrumDemodSample dataDemod
- union ZISpectrumWave::@3 data
Spectrum data vector.

struct ZIAdvisorSample

```
typedef struct ZIAdvisorSample {  
    double grid;  
    double x;  
    double y;  
} ZIAdvisorSample;
```

Data Fields

- double grid
Grid.
- double x
X.
- double y
Y.

struct ZIAdvisorHeader

```
typedef struct ZIAdvisorHeader {  
    uint64_t sampleCount;  
    uint8_t flags;  
    uint8_t sampleFormat;  
    uint8_t reserved0[6];  
    uint8_t reserved1[8];  
} ZIAdvisorHeader;
```

Data Fields

- `uint64_t sampleCount`
Total sample count considered for advisor.
- `uint8_t flags`
Flags.
- `uint8_t sampleFormat`
Sample format Bode = 0, Step = 1, Impulse = 2.
- `uint8_t reserved0`
Reserved space for future use.
- `uint8_t reserved1`
Reserved space for future use.

struct ZIAdvisorWave

```
typedef struct ZIAdvisorWave {  
    ZITimeStamp timeStamp;  
    ZIAdvisorHeader header;  
    ZIAdvisorSample data[0];  
    union ZIAdvisorWave::@4 data;  
} ZIAdvisorWave;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp at which the data was updated.
- ZIAdvisorHeader header
- ZIAdvisorSample data
- union ZIAdvisorWave::@4 data
Advisor data vector.

struct ZIVectorData

The structure used to hold vector data block. See the description of the structure members for details.

```
#include "ziAPI.h"

typedef struct ZIVectorData {
    ZITimeStamp timeStamp;
    uint32_t sequenceNumber;
    uint32_t blockNumber;
    uint64_t totalElements;
    uint64_t blockOffset;
    uint32_t blockElements;
    uint8_t flags;
    uint8_t elementType;
    uint8_t reserved0[2];
    uint64_t reserved1[32];
    uint8_t dataUInt8[0];
    uint16_t dataUInt16[0];
    uint32_t dataUInt32[0];
    uint64_t dataUInt64[0];
    int8_t dataInt8[0];
    int16_t dataInt16[0];
    int32_t dataInt32[0];
    int64_t dataInt64[0];
    double dataDouble[0];
    float dataFloat[0];
    union ZIVectorData::@5 data;
} ZIVectorData;
```

Data Fields

- **ZITimeStamp timeStamp**
Time stamp of this array data block.
- **uint32_t sequenceNumber**
Current array transfer sequence number. Incremented for each new transfer. Stays same for all blocks of a single array transfer.
- **uint32_t blockNumber**
Current block number from the beginning of an array transfer. Large array transfers are split into blocks, which need to be concatenated to obtain the complete array.
- **uint64_t totalElements**
Total number of elements in the array.
- **uint64_t blockOffset**
Offset of the current block first element from the beginning of the array.
- **uint32_t blockElements**
Number of elements in the current block.
- **uint8_t flags**
Block marker: Bit (0): 1 = End marker for multi-block transfer
Bit (1): 1 = Transfer failure Bit (7..2): Reserved.
- **uint8_t elementType**

Vector element type, see `ZIVectorElementType_enum`.

- `uint8_t reserved0`
- `uint64_t reserved1`
- `uint8_t dataUInt8`
- `uint16_t dataUInt16`
- `uint32_t dataUInt32`
- `uint64_t dataUInt64`
- `int8_t dataInt8`
- `int16_t dataInt16`
- `int32_t dataInt32`
- `int64_t dataInt64`
- `double dataDouble`
- `float dataFloat`
- `union ZIVectorData::@5 data`
First data element of the current block.

struct ZIAsyncReply

```
typedef struct ZIAsyncReply {  
    ZITimeStamp timeStamp;  
    ZITimeStamp sampleTimeStamp;  
    uint16_t command;  
    uint16_t resultCode;  
    ZIAsyncTag tag;  
} ZIAsyncReply;
```

Data Fields

- ZITimeStamp timeStamp
Time stamp of the reply (server clock)
- ZITimeStamp sampleTimeStamp
Time stamp of the target node sample, to which the reply belongs.
- uint16_t command
Command: 1 - ziAPIAsyncSetDoubleData 2 - ziAPIAsyncSetIntegerData 3 - ziAPIAsyncSetByteArray 4 - ziAPIAsyncSubscribe 5 - ziAPIAsyncUnSubscribe 6 - ziAPIAsyncGetValueAsPollData.
- uint16_t resultCode
Command result code (cast to ZIResult_enum)
- ZIAsyncTag tag
Tag sent along with the async command.

struct ZIEvent

This struct holds event data forwarded by the Data Server.

```
#include "ziAPI.h"

typedef struct ZIEvent {
    uint32_t valueType;
    uint32_t count;
    uint8_t path[256];
    void* untyped;
    ZIDoubleData* doubleData;
    ZIDoubleDataTS* doubleDataTS;
    ZIIntegerData* integerData;
    ZIIntegerDataTS* integerDataTS;
    ZIByteArray* byteArray;
    ZIByteArrayTS* byteArrayTS;
    ZICntSample* cntSample;
    ZITreeChangeData* treeChangeData;
    TreeChange* treeChangeDataOld;
    ZIDemodSample* demodSample;
    ZIAuxInSample* auxInSample;
    ZIDIOSample* dioSample;
    ZIScopeWave* scopeWave;
    ZIScopeWaveEx* scopeWaveEx;
    ScopeWave* scopeWaveOld;
    ZIPWAWave* pwaWave;
    ZISweeperWave* sweeperWave;
    ZISpectrumWave* spectrumWave;
    ZIAdvisorWave* advisorWave;
    ZIAsyncReply* asyncReply;
    ZIVectorData* vectorData;
    ZIImpedanceSample* impedanceSample;
    uint64_t alignment;
    union ZIEvent::@6 value;
    uint8_t data[0x400000];
} ZIEvent;
```

Data Fields

- `uint32_t valueType`
Specifies the type of the data held by the [ZIEvent](#), see [ZIValueType_enum](#).
- `uint32_t count`
Number of values available in this event.
- `uint8_t path`
The path to the node from which the event originates.
- `void* untyped`
For convenience. The void field doesn't have a corresponding data type.
- `ZIDoubleData* doubleData`
when `valueType == ZI_VALUE_TYPE_DOUBLE_DATA`
- `ZIDoubleDataTS* doubleDataTS`
when `valueType == ZI_VALUE_TYPE_DOUBLE_DATA_TS`
- `ZIIntegerData* integerData`

- when valueType == ZI_VALUE_TYPE_INTEGER_DATA
- [ZIIntegerDataTS](#)* integerDataTS
when valueType == ZI_VALUE_TYPE_INTEGER_DATA_TS
- [ZIByteArray](#)* byteArray
when valueType == ZI_VALUE_TYPE_BYTE_ARRAY
- [ZIByteArrayTS](#)* byteArrayTS
when valueType == ZI_VALUE_TYPE_BYTE_ARRAY_TS
- [ZICntSample](#)* cntSample
when valueType == ZI_VALUE_TYPE_CNT_SAMPLE
- [ZITreeChangeData](#)* treeChangeData
when valueType == ZI_VALUE_TYPE_TREE_CHANGE_DATA
- [TreeChange](#)* treeChangeDataOld
when valueType ==
ZI_VALUE_TYPE_TREE_CHANGE_DATA_OLD
- [ZIDemodSample](#)* demodSample
when valueType == ZI_VALUE_TYPE_DEMOD_SAMPLE
- [ZIAuxInSample](#)* auxInSample
when valueType == ZI_VALUE_TYPE_AUXIN_SAMPLE
- [ZIDIOSample](#)* dioSample
when valueType == ZI_VALUE_TYPE_DIO_SAMPLE
- [ZIScopeWave](#)* scopeWave
when valueType == ZI_VALUE_TYPE_SCOPE_WAVE
- [ZIScopeWaveEx](#)* scopeWaveEx
when valueType == ZI_VALUE_TYPE_SCOPE_WAVE_EX
- [ScopeWave](#)* scopeWaveOld
when valueType == ZI_VALUE_TYPE_SCOPE_WAVE_OLD
- [ZIPWAWave](#)* pwaWave
when valueType == ZI_VALUE_TYPE_PWA_WAVE
- [ZISweeperWave](#)* sweeperWave
when valueType == ZI_VALUE_TYPE_SWEEPER_WAVE
- [ZISpectrumWave](#)* spectrumWave
when valueType == ZI_VALUE_TYPE_SPECTRUM_WAVE
- [ZIAdvisorWave](#)* advisorWave
when valueType == ZI_VALUE_TYPE_ADVISOR_WAVE
- [ZIAsyncReply](#)* asyncReply
when valueType == ZI_VALUE_TYPE_ASYNC_REPLY
- [ZIVectorData](#)* vectorData

- when valueType == ZI_VALUE_TYPE_VECTOR_DATA
- [ZIImpedanceSample](#)* impedanceSample
when valueType == ZI_VALUE_TYPE_IMPEDANCE_SAMPLE
- uint64_t alignment
ensure union size is 8 bytes
- union ZIEvent::@6 value
Convenience pointer to allow for access to the first entry in Data using the correct type according to [ZIEvent.valueType](#) field.
- uint8_t data
The raw value data.

Detailed Description

[ZIEvent](#) is used to give out events like value changes or errors to the user. Event handling functionality is provided by [ziAPISubscribe](#) and [ziAPIUnSubscribe](#) as well as [ziAPIPollDataEx](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

void ProcessEvent(ZIEvent* Event) {
    unsigned int j;

    switch (Event->valueType) {
    case ZI_VALUE_TYPE_DOUBLE_DATA:

        printf("%u elements of double data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", Event->value.doubleData[j]);

        break;

    case ZI_VALUE_TYPE_INTEGER_DATA:

        printf("%u elements of integer data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", (float)Event->value.integerData[j]);

        break;

    case ZI_VALUE_TYPE_DEMOD_SAMPLE:

        printf("%u elements of sample data %s\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("TS=%f, X=%f, Y=%f.\n",
                (float)Event->value.demodSample[j].timeStamp,
                Event->value.demodSample[j].x,
                Event->value.demodSample[j].y);
```

```
        break;

    case ZI_VALUE_TYPE_TREE_CHANGE_DATA:

        printf("%u elements of tree-changed data, %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++) {
            switch (Event->value.treeChangeDataOld[j].Action) {
                case ZI_TREE_ACTION_REMOVE:
                    printf("Tree removed: %s\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_ADD:
                    printf("treeChangeDataOld added: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_CHANGE:
                    printf("treeChangeDataOld changed: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;
            }
        }

        break;

    default:

        printf("Unexpected event value type: %d.\n", Event->valueType);
        break;
    }
}
```

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

struct ZISWTriggerHeader

Structure to hold information about data returned from the SW Trigger Module.

```
#include "ziAPI.h"

typedef struct ZISWTriggerHeader {
    ZITimeStamp triggerStart;
    uint64_t triggerNumber;
    uint32_t cols;
    uint32_t rows;
    uint32_t repetitions;
    uint32_t flags;
    uint8_t reserved0[32];
} ZISWTriggerHeader;
```

Data Fields

- ZITimeStamp triggerStart
Trigger timestamp.
- uint64_t triggerNumber
Trigger counter since execution start.
- uint32_t cols
Number of columns.
- uint32_t rows
Number of rows.
- uint32_t repetitions
Number of repetitions in grid mode.
- uint32_t flags
Flags Bit (0): Finished (all repetitions recorded)
- uint8_t reserved0
Reserved space for future use.

struct ZIModuleHeaderSweeper

Structure to hold information about data returned from the Sweep Module.

```
#include "ziAPI.h"

typedef struct ZIModuleHeaderSweeper {
    char traceName[256];
} ZIModuleHeaderSweeper;
```

Data Fields

- char traceName

struct ZIModuleHeader

Module-specific event header.

```
#include "ziAPI.h"

typedef struct ZIModuleHeader {
    ZIModuleHeaderType_enum
        type;
    void* untyped;
    ZISWTriggerHeader* swTrigger;
    ZISweeperHeader* sweeper;
    union ZIModuleHeader::@7 ptr;
} ZIModuleHeader;
```

Data Fields

- `ZIModuleHeaderType_enum` type
- `void*` `untyped`
- `ZISWTriggerHeader*` `swTrigger`
- `ZISweeperHeader*` `sweeper`
- `union ZIModuleHeader::@7` `ptr`

struct ZIModuleEvent

This struct holds data of a single chunk from module lookup.

```
#include "ziAPI.h"

typedef struct ZIModuleEvent {
    uint64_t allocatedSize;
    ZIModuleHeader header;

    ZIEvent value[0];
} ZIModuleEvent;
```

Data Fields

- `uint64_t allocatedSize`
For internal use - never modify!
- `ZIModuleHeader header`
Module-specific event header.
- `ZIEvent value`
Defines location of stored `ZIEvent`.

struct DemodSample

The [DemodSample](#) struct holds data for the ZI_DATA_DEMODSAMPLE data type. Deprecated: See [ZIDemodSample](#).

```
#include "ziAPI.h"

typedef struct DemodSample {
    ziTimeStampType TimeStamp;
    double X;
    double Y;
    double Frequency;
    double Phase;
    unsigned int DIOBits;
    unsigned int Reserved;
    double AuxIn0;
    double AuxIn1;
} DemodSample;
```

Data Fields

- `ziTimeStampType TimeStamp`
- `double X`
- `double Y`
- `double Frequency`
- `double Phase`
- `unsigned int DIOBits`
- `unsigned int Reserved`
- `double AuxIn0`
- `double AuxIn1`

struct AuxInSample

The [AuxInSample](#) struct holds data for the ZI_DATA_AUXINSAMPLE data type. Deprecated: See [ZIAuxInSample](#).

```
#include "ziAPI.h"

typedef struct AuxInSample {
    ziTimeStampType TimeStamp;
    double Ch0;
    double Ch1;
} AuxInSample;
```

Data Fields

- `ziTimeStampType TimeStamp`
- `double Ch0`
- `double Ch1`

struct DIOSample

The [DIOSample](#) struct holds data for the ZI_DATA_DIOSAMPLE data type. Deprecated: See [ZIDIOSample](#).

```
#include "ziAPI.h"

typedef struct DIOSample {
    ziTimeStampType TimeStamp;
    unsigned int Bits;
    unsigned int Reserved;
} DIOSample;
```

Data Fields

- `ziTimeStampType TimeStamp`
- `unsigned int Bits`
- `unsigned int Reserved`

struct ByteArrayData

The [ByteArrayData](#) struct holds data for the ZI_DATA_BYTEARRAY data type. Deprecated: See [ZIByteArray](#).

```
#include "ziAPI.h"

typedef struct ByteArrayData {
    unsigned int Len;
    unsigned char Bytes[0];
} ByteArrayData;
```

Data Fields

- unsigned int Len
- unsigned char Bytes

struct `ziEvent`

This struct holds event data forwarded by the Data Server. Deprecated: See [ZIEvent](#).

```
#include "ziAPI.h"

typedef struct ziEvent {
    uint32_t Type;
    uint32_t Count;
    unsigned char Path[256];
    union ziEvent::Val Val;
    unsigned char Data[0x400000];
} ziEvent;
```

Data Structures

- union `ziEvent::Val`

Data Fields

- `uint32_t` Type
- `uint32_t` Count
- unsigned char Path
- union `ziEvent::Val` Val
- unsigned char Data

Detailed Description

[ziEvent](#) is used to give out events like value changes or errors to the user. Event handling functionality is provided by `ziAPISubscribe` and `ziAPIUnSubscribe` as well as `ziAPIPollDataEx`.

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

```
// Copyright [2016] Zurich Instruments AG
#include <stdio.h>

#include "ziAPI.h"

void ProcessEvent(ZIEvent* Event) {
    unsigned int j;

    switch (Event->valueType) {
    case ZI_VALUE_TYPE_DOUBLE_DATA:

        printf("%u elements of double data: %s.\n",
            Event->count,
            Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", Event->value.doubleData[j]);
    }
```

```
        break;

    case ZI_VALUE_TYPE_INTEGER_DATA:

        printf("%u elements of integer data: %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++)
            printf("%f\n", (float)Event->value.integerData[j]);

        break;

    case ZI_VALUE_TYPE_DEMOD_SAMPLE:

        printf("%u elements of sample data %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++)
            printf("TS=%f, X=%f, Y=%f.\n",
                   (float)Event->value.demodSample[j].timeStamp,
                   Event->value.demodSample[j].x,
                   Event->value.demodSample[j].y);

        break;

    case ZI_VALUE_TYPE_TREE_CHANGE_DATA:

        printf("%u elements of tree-changed data, %s.\n",
               Event->count,
               Event->path);

        for (j = 0; j < Event->count; j++) {
            switch (Event->value.treeChangeDataOld[j].Action) {
                case ZI_TREE_ACTION_REMOVE:
                    printf("Tree removed: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_ADD:
                    printf("treeChangeDataOld added: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;

                case ZI_TREE_ACTION_CHANGE:
                    printf("treeChangeDataOld changed: %s.\n",
                           Event->value.treeChangeDataOld[j].Name);
                    break;
            }
        }

        break;

    default:

        printf("Unexpected event value type: %d.\n", Event->valueType);
        break;
    }
}
```

Data Structure Documentation

union ziEvent::Val

```
typedef union ziEvent::Val {  
    void* Void;  
    DemodSample* SampleDemod;  
    AuxInSample* SampleAuxIn;  
    DIOSample* SampleDIO;  
    ziDoubleType* Double;  
    ziIntegerType* Integer;  
    TreeChange* Tree;  
    ByteArrayData* ByteArray;  
    ScopeWave* Wave;  
    uint64_t alignment;  
} ziEvent::Val;
```

Data Fields

- void* Void
- DemodSample* SampleDemod
- AuxInSample* SampleAuxIn
- DIOSample* SampleDIO
- ziDoubleType* Double
- ziIntegerType* Integer
- TreeChange* Tree
- ByteArrayData* ByteArray
- ScopeWave* Wave
- uint64_t alignment

union ziEvent::Val

```
typedef union ziEvent::Val {  
    void* Void;  
    DemodSample* SampleDemod;  
    AuxInSample* SampleAuxIn;  
    DIOSample* SampleDIO;  
    ziDoubleType* Double;  
    ziIntegerType* Integer;  
    TreeChange* Tree;  
    ByteArrayData* ByteArray;  
    ScopeWave* Wave;  
    uint64_t alignment;  
} ziEvent::Val;
```

Data Fields

- void* Void
- DemodSample* SampleDemod
- AuxInSample* SampleAuxIn
- DIOSample* SampleDIO
- ziDoubleType* Double
- ziIntegerType* Integer
- TreeChange* Tree
- ByteArrayData* ByteArray
- ScopeWave* Wave
- uint64_t alignment

Enumeration Type Documentation

enum ZIResult_enum

Defines return value for all ziAPI functions. Divided into 3 regions: info, warning and error.

Enumerator:

- ZI_INFO_BASE
- ZI_INFO_SUCCESS
Success (no error)
- ZI_INFO_MAX
- ZI_WARNING_BASE
- ZI_WARNING_GENERAL
Warning (general);.
- ZI_WARNING_UNDERRUN
FIFO Underrun.
- ZI_WARNING_OVERFLOW
FIFO Overflow.
- ZI_WARNING_NOTFOUND
Value or Node not found.
- ZI_WARNING_NO_ASYNC
Async command executed in sync mode (will be no async reply)
- ZI_WARNING_MAX
- ZI_ERROR_BASE
- ZI_ERROR_GENERAL
Error (general)
- ZI_ERROR_USB
USB Communication failed.
- ZI_ERROR_MALLOC
Memory allocation failed.
- ZI_ERROR_MUTEX_INIT
Unable to initialize mutex.
- ZI_ERROR_MUTEX_DESTROY
Unable to destroy mutex.
- ZI_ERROR_MUTEX_LOCK
Unable to lock mutex.
- ZI_ERROR_MUTEX_UNLOCK
Unable to unlock mutex.

- `ZI_ERROR_THREAD_START`
Unable to start thread.
- `ZI_ERROR_THREAD_JOIN`
Unable to join thread.
- `ZI_ERROR_SOCKET_INIT`
Can't initialize socket.
- `ZI_ERROR_SOCKET_CONNECT`
Unable to connect socket.
- `ZI_ERROR_HOSTNAME`
Hostname not found.
- `ZI_ERROR_CONNECTION`
Connection invalid.
- `ZI_ERROR_TIMEOUT`
Command timed out.
- `ZI_ERROR_COMMAND`
Command internally failed.
- `ZI_ERROR_SERVER_INTERNAL`
Command failed in server.
- `ZI_ERROR_LENGTH`
Provided Buffer length is too small.
- `ZI_ERROR_FILE`
Can't open file or read from it.
- `ZI_ERROR_DUPLICATE`
There is already a similar entry.
- `ZI_ERROR_READONLY`
Attempt to set a read-only node.
- `ZI_ERROR_DEVICE_NOT_VISIBLE`
Device is not visible to the server.
- `ZI_ERROR_DEVICE_IN_USE`
Device is already connected by a different server.
- `ZI_ERROR_DEVICE_INTERFACE`
Device does currently not support the specified interface.
- `ZI_ERROR_DEVICE_CONNECTION_TIMEOUT`
Device connection timeout.
- `ZI_ERROR_DEVICE_DIFFERENT_INTERFACE`
Device already connected over a different Interface.
- `ZI_ERROR_DEVICE_NEEDS_FW_UPGRADE`

Device needs FW upgrade.

- ZI_ERROR_ZIEVENT_DATATYPE_MISMATCH
Trying to get data from a poll event with wrong target data type.
- ZI_ERROR_DEVICE_NOT_FOUND
Device not found.
- ZI_ERROR_NOT_SUPPORTED
Provided arguments are not supported for the command.
- ZI_ERROR_TOO_MANY_CONNECTIONS
Connection invalid.
- ZI_ERROR_MAX
- ZI_SUCCESS
Success (no error)
- ZI_MAX_INFO
- ZI_WARNING
Warning (general);.
- ZI_UNDERRUN
FIFO Underrun.
- ZI_OVERFLOW
FIFO Overflow.
- ZI_NOTFOUND
Value or Node not found.
- ZI_MAX_WARNING
- ZI_ERROR
Error (general)
- ZI_USB
USB Communication failed.
- ZI_MALLOC
Memory allocation failed.
- ZI_MUTEX_INIT
Unable to initialize mutex.
- ZI_MUTEX_DESTROY
Unable to destroy mutex.
- ZI_MUTEX_LOCK
Unable to lock mutex.
- ZI_MUTEX_UNLOCK
Unable to unlock mutex.

- `ZI_THREAD_START`
Unable to start thread.
- `ZI_THREAD_JOIN`
Unable to join thread.
- `ZI_SOCKET_INIT`
Can't initialize socket.
- `ZI_SOCKET_CONNECT`
Unable to connect socket.
- `ZI_HOSTNAME`
Hostname not found.
- `ZI_CONNECTION`
Connection invalid.
- `ZI_TIMEOUT`
Command timed out.
- `ZI_COMMAND`
Command internally failed.
- `ZI_SERVER_INTERNAL`
Command failed in server.
- `ZI_LENGTH`
Provided Buffer length doesn't reach.
- `ZI_FILE`
Can't open file or read from it.
- `ZI_DUPLICATE`
There is already a similar entry.
- `ZI_READONLY`
Attempt to set a read-only node.
- `ZI_MAX_ERROR`

enum ZIValueType_enum

Enumerates all types that data in a [ZIEvent](#) may have.

Enumerator:

- `ZI_VALUE_TYPE_NONE`
No data type, event is invalid.
- `ZI_VALUE_TYPE_DOUBLE_DATA`
[ZIDoubleData](#) type. Use the `ZIEvent.value.doubleData` pointer to read the data of the event.
- `ZI_VALUE_TYPE_INTEGER_DATA`
[ZIIntegerData](#) type. Use the `ZIEvent.value.integerData` pointer to read the data of the event.
- `ZI_VALUE_TYPE_DEMOD_SAMPLE`
[ZIDemodSample](#) type. Use the `ZIEvent.value.demodSample` pointer to read the data of the event.
- `ZI_VALUE_TYPE_SCOPE_WAVE_OLD`
[ScopeWave](#) type, used in v1 compatibility mode. use the `ZIEvent.value.scopeWaveOld` pointer to read the data of the event.
- `ZI_VALUE_TYPE_AUXIN_SAMPLE`
[ZIAuxInSample](#) type. Use the `ZIEvent.value.auxInSample` pointer to read the data of the event.
- `ZI_VALUE_TYPE_DIO_SAMPLE`
[ZIDIOSample](#) type. Use the `ZIEvent.value.dioSample` pointer to read the data of the event.
- `ZI_VALUE_TYPE_BYTE_ARRAY`
[ZIByteArray](#) type. Use the `ZIEvent.value.byteArray` pointer to read the data of the event.
- `ZI_VALUE_TYPE_PWA_WAVE`
[ZIPWAWave](#) type. Use the `ZIEvent.value.pwaWave` pointer to read the data of the event.
- `ZI_VALUE_TYPE_TREE_CHANGE_DATA_OLD`
[TreeChange](#) type - a list of added or removed nodes, used in v1 compatibility mode. Use the `ZIEvent.value.treeChangeDataOld` pointer to read the data of the event.
- `ZI_VALUE_TYPE_DOUBLE_DATA_TS`
[ZIDoubleDataTS](#) type. Use the `ZIEvent.value.doubleDataTS` pointer to read the data of the event.
- `ZI_VALUE_TYPE_INTEGER_DATA_TS`
[ZIIntegerDataTS](#) type. Use the `ZIEvent.value.integerDataTS` pointer to read the data of the event.

- `ZI_VALUE_TYPE_SCOPE_WAVE`
`ZIScopeWave` type. Use the `ZIEvent.value.scopeWave` pointer to read the data of the event.
- `ZI_VALUE_TYPE_SCOPE_WAVE_EX`
`ZIScopeWaveEx` type. Use the `ZIEvent.value.scopeWaveEx` pointer to read the data of the event.
- `ZI_VALUE_TYPE_BYTE_ARRAY_TS`
`ZIByteArrayTS` type. Use the `ZIEvent.value.byteArrayTS` pointer to read the data of the event.
- `ZI_VALUE_TYPE_CNT_SAMPLE`
`ZICntSample` type. Use the `ZIEvent.value.cntSample` pointer to read the data of the event.
- `ZI_VALUE_TYPE_TREE_CHANGE_DATA`
`ZITreeChangeData` type - a list of added or removed nodes. Use the `ZIEvent.value.treeChangeData` pointer to read the data of the event.
- `ZI_VALUE_TYPE_ASYNC_REPLY`
`ZIAsyncReply` type. Use the `ZIEvent.value.asyncReply` pointer to read the data of the event.
- `ZI_VALUE_TYPE_SWEEPER_WAVE`
`ZISweeperWave` type. Use the `ZIEvent.value.sweeperWave` pointer to read the data of the event.
- `ZI_VALUE_TYPE_SPECTRUM_WAVE`
`ZISpectrumWave` type. Use the `ZIEvent.value.spectrumWave` pointer to read the data of the event.
- `ZI_VALUE_TYPE_ADVISOR_WAVE`
`ZIAdvisorWave` type. Use the `ZIEvent.value.advisorWave` pointer to read the data of the event.
- `ZI_VALUE_TYPE_VECTOR_DATA`
`ZIVectorData` type. Use the `ZIEvent.value.vectorData` pointer to access the data of the event.
- `ZI_VALUE_TYPE_IMPEDANCE_SAMPLE`
`ZIImpedanceSample` type. Use the `ZIEvent.value.impedanceSample` pointer to access the data of the event.
- `ZI_DATA_NONE`
no data type. the `ziEvent` is invalid.
- `ZI_DATA_DOUBLE`
double data type. use the `ziEvent::Val.Double` Pointer to read the data of the event.
- `ZI_DATA_INTEGER`
integer data type. use the `ziEvent::Val.Integer` Pointer to read the data of the event.

- `ZI_DATA_DEMODSAMPLE`
`DemodSample` data type. use the `ziEvent::Val.Sample` Pointer to read the data of the event.
- `ZI_DATA_SCOPEWAVE`
`ScopeWave` data type. use the `ziEvent::Val.Wave` Pointer to read the data of the event.
- `ZI_DATA_AUXINSAMPLE`
`MiscADValue` data type. use the `ziEvent::Val.ADValue` Pointer to read the data of the event.
- `ZI_DATA_DIOSAMPLE`
`DIOValue` data type. use the `ziEvent::Val.DIOValue` Pointer to read the data of the event.
- `ZI_DATA_BYTEARRAY`
`ByteArray` data type. use the `ziEvent::Val.ByteArray` Pointer to read the data of the event.
- `ZI_DATA_TREE_CHANGED`
a list of added or removed trees. use the `ziEvent::Val.Tree` Pointer to read the data of the event.

enum ZITreeAction_enum

Defines the actions that are performed on a tree, as returned in the [ZITreeChangeData::action](#) or [ZITreeChangeDataOld::action](#).

Enumerator:

- `ZI_TREE_ACTION_REMOVE`
A node has been removed.
- `ZI_TREE_ACTION_ADD`
A node has been added.
- `ZI_TREE_ACTION_CHANGE`
A node has been changed.

enum ZImpFlags_enum

Enumerates the bits set in an [ZImpedanceSample](#)'s flags.

Enumerator:

- ZI_IMP_FLAGS_NONE
- ZI_IMP_FLAGS_VALID_INTERNAL
- ZI_IMP_FLAGS_VALID_USER
- ZI_IMP_FLAGS_AUTORANGE_GATING
- ZI_IMP_FLAGS_OVERFLOW_VOLTAGE
- ZI_IMP_FLAGS_OVERFLOW_CURRENT
- ZI_IMP_FLAGS_UNDERFLOW_VOLTAGE
- ZI_IMP_FLAGS_UNDERFLOW_CURRENT
- ZI_IMP_FLAGS_FREQ_EXACT
- ZI_IMP_FLAGS_FREQ_INTERPOLATION
- ZI_IMP_FLAGS_FREQ_EXTRAPOLATION
- ZI_IMP_FLAGS_SUPPRESSION_PARAM0
- ZI_IMP_FLAGS_SUPPRESSION_PARAM1
- ZI_IMP_FLAGS_STRONGCOMPENSATION_PARAM0
- ZI_IMP_FLAGS_STRONGCOMPENSATION_PARAM1
- ZI_IMP_FLAGS_BWC_BIT0
- ZI_IMP_FLAGS_BWC_BIT1
- ZI_IMP_FLAGS_BWC_BIT2
- ZI_IMP_FLAGS_BWC_BIT3
- ZI_IMP_FLAGS_BWC_MASK
- ZI_IMP_FLAGS_OPEN_DETECTION

enum ZIVectorElementType_enum

Enumerates all the types that a ::elementType may have.

Enumerator:

- ZI_VECTOR_ELEMENT_TYPE_UINT8
- ZI_VECTOR_ELEMENT_TYPE_UINT16
- ZI_VECTOR_ELEMENT_TYPE_UINT32
- ZI_VECTOR_ELEMENT_TYPE_UINT64
- ZI_VECTOR_ELEMENT_TYPE_FLOAT
- ZI_VECTOR_ELEMENT_TYPE_DOUBLE
- ZI_VECTOR_ELEMENT_TYPE_ASCII_Z
NULL-terminated string.

enum ZIAPIVersion_enum

Enumerator:

- ZI_API_VERSION_0
- ZI_API_VERSION_1
- ZI_API_VERSION_4
- ZI_API_VERSION_5

enum ZIListNodes_enum

Defines the values of the flags used in [ziAPIListNodes](#).

Enumerator:

- `ZI_LIST_NODES_NONE`
Default, return a simple listing of the given node immediate descendants.
- `ZI_LIST_NODES_RECURSIVE`
List the nodes recursively.
- `ZI_LIST_NODES_ABSOLUTE`
Return absolute paths.
- `ZI_LIST_NODES_LEAFONLY`
Return only leaf nodes, which means the nodes at the outermost level of the tree.
- `ZI_LIST_NODES_SETTINGSONLY`
Return only nodes which are marked as setting.
- `ZI_LIST_NONE`
Default, return a simple listing of the given node immediate descendants.
- `ZI_LIST_RECURSIVE`
List the nodes recursively.
- `ZI_LIST_ABSOLUTE`
Return absolute paths.
- `ZI_LIST_LEAFONLY`
Return only leaf nodes, which means the nodes at the outermost level of the tree.
- `ZI_LIST_SETTINGSONLY`
Return only nodes which are marked as setting.

enum ZIModuleHeaderType_enum

Enumerates all module header types.

Enumerator:

- ZI_MODULE_HEADER_TYPE_NONE
- ZI_MODULE_HEADER_TYPE_SWTRIGGER
- ZI_MODULE_HEADER_TYPE_SWEEPER

enum ZIVectorWriteStatus_enum

Enumerator:

- ZI_VECTOR_WRITE_STATUS_IDLE
- ZI_VECTOR_WRITE_STATUS_PENDING

enum TREE_ACTION

TREE_ACTION defines the values for the [TreeChange::Action](#) Variable.

Enumerator:

- TREE_ACTION_REMOVE
a tree has been removed
- TREE_ACTION_ADD
a tree has been added
- TREE_ACTION_CHANGE
a tree has changed

Function Documentation

ziAPIInit

ZIResult_enum ziAPIInit (**ZIConnection*** conn)

Initializes a **ZIConnection** structure.

This function initializes the structure so that it is ready to connect to Data Server. It allocates memory and sets up the infrastructure needed.

Parameters:

[out] conn

Pointer to **ZIConnection** that is to be initialized

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_MALLOC on memory allocation failure
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIDestroy](#), [ziAPIConnect](#), [ziAPIDisconnect](#)

See [Connection](#) for an example

ziAPIDestroy

ZIResult_enum ziAPIDestroy (ZIConnection conn)

Destroys a [ZIConnection](#) structure.

This function frees all memory that has been allocated by [ziAPIInit](#). If it is called with an uninitialized [ZIConnection](#) struct it may result in segmentation faults as well when it is called with a struct for which [ZIAPIDestroy](#) already has been called.

Parameters:

[in] conn

Pointer to [ZIConnection](#) struct that has to be destroyed

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIInit](#), [ziAPIConnect](#), [ziAPIDisconnect](#)

See [Connection](#) for an example

ziAPIConnect

ZIResult_enum ziAPIConnect ([ZIConnection](#) conn, const char* hostname, uint16_t port)

Connects the [ZIConnection](#) to Data Server.

Connects to Data Server using a [ZIConnection](#) and prepares for data exchange. For most cases it is enough to just give a reference to the connection and give NULL for hostname and 0 for the port, so it connects to localhost on the default port.

Parameters:

[in] conn

Pointer to [ZIConnection](#) with which the connection should be established

[in] hostname

Name of the Host to which it should be connected, if NULL "localhost" will be used as default

[in] port

The Number of the port to connect to. If 0, default port of the local Data Server will be used (8005)

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_HOSTNAME if the given host name could not be found
- ZI_ERROR_SOCKET_CONNECT if no connection could be established
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_SOCKET_INIT if initialization of the socket failed
- ZI_ERROR_CONNECTION when the Data Server didn't return the correct answer
- ZI_ERROR_TIMEOUT when initial communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIDisconnect](#), [ziAPIInit](#), [ziAPIDestroy](#)

See [Connection](#) for an example

ziAPIDisconnect

ZIResult_enum `ziAPIDisconnect (ZIConnection conn)`

Disconnects an established connection.

Disconnects from Data Server. If the connection has not been established and the function is called it returns without doing anything.

Parameters:

[in] conn
Pointer to ZIConnection to be disconnected

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIConnect](#), [ziAPIInit](#), [ziAPIDestroy](#)

See [Connection](#) for an example

ziAPIListImplementations

ZIResult_enum **ziAPIListImplementations** (char* implementations, uint32_t bufferSize)

Returns the list of supported implementations.

Returned names are defined by implementations in the linked library and may change depending on software version.

Parameters:

[out] implementations

Pointer to a buffer receiving a newline-delimited list of the names of all the supported ziAPI implementations. The string is zero-terminated.

[in] bufferSize

The size of the buffer assigned to the implementations parameter

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_LENGTH if the length of the char-buffer given by MaxLen is too small for all elements
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIConnectEx](#)

ziAPIConnectEx

ZIResult_enum ziAPIConnectEx (**ZIConnection** conn, const char* hostname, uint16_t port, **ZIAPIVersion_enum** apiLevel, const char* implementation)

Connects to Data Server and enables extended ziAPI.

With apiLevel=ZI_API_VERSION_1 and implementation=NULL, this call is equivalent to plain [ziAPIConnect](#). With other version and implementation values enables corresponding ziAPI extension and connection using different implementation.

Parameters:

[in] conn

Pointer to the ZIConnection with which the connection should be established

[in] hostname

Name of the host to which it should be connected, if NULL "localhost" will be used as default

[in] port

The number of the port to connect to. If 0 the port of the local Data Server will be used

[in] apiLevel

Specifies the ziAPI compatibility level to use for this connection (1 or 4).

[in] implementation

Specifies implementation to use for a connection, must be one of the returned by [ziAPIListImplementations](#) or NULL to select default implementation

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_HOSTNAME if the given host name could not be found
- ZI_ERROR_SOCKET_CONNECT if no connection could be established
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_SOCKET_INIT if initialization of the socket failed
- ZI_ERROR_CONNECTION when the Data Server didn't return the correct answer or requested implementation is not found or doesn't support requested ziAPI level
- ZI_ERROR_TIMEOUT when initial communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIListImplementations](#), [ziAPIConnect](#), [ziAPIDisconnect](#), [ziAPIInit](#), [ziAPIDestroy](#), [ziAPIGetConnectionVersion](#)

See [Connection](#) for an example

ziAPIGetConnectionAPILevel

ZIResult_enum ziAPIGetConnectionAPILevel (**ZIConnection** conn, ZIAPIVersion_enum* apiLevel)

Returns ziAPI level used for the connection conn.

Parameters:

[in] conn

Pointer to ZIConnection

[out] apiLevel

Pointer to preallocated ZIAPIVersion_enum, receiving the ziAPI level

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION if level can not be determined due to conn is not connected

See Also:

[ziAPIConnectEx](#), [ziAPIGetVersion](#)

ziAPIGetRevision

ZIResult_enum ziAPIGetRevision (unsigned int* revision)

Retrieves the revision of ziAPI.

Sets an unsigned int with the revision (build number) of the ziAPI you are using.

Parameters:

[in] revision

Pointer to an unsigned int to fill up with the revision.

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIListNodes

ZIResult_enum ziAPIListNodes (**ZIConnection** conn, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)

Returns all child nodes found at the specified path.

This function returns a list of node names found at the specified path. The path may contain wildcards so that the returned nodes do not necessarily have to have the same parents. The list is returned in a null-terminated char-buffer, each element delimited by a newline. If the maximum length of the buffer (bufferSize) is not sufficient for all elements, nothing will be returned and the return value will be [ZI_LENGTH](#).

Parameters:

[in] conn

Pointer to the ZIConnection for which the node names should be retrieved.

[in] path

Path for which all children will be returned. The path may contain wildcard characters.

[out] nodes

Upon call filled with newline-delimited list of the names of all the children found. The string is zero-terminated.

[in] bufferSize

The length of the buffer used for the nodes output parameter.

[in] flags

A combination of flags (applied bitwise) as defined in [ZIListNodes_enum](#).

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the path's length exceeds [MAX_PATH_LEN](#) or the length of the char-buffer for the nodes given by bufferSize is too small for all elements
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Tree Listing](#) for an example

See Also:

[ziAPIUpdate](#)

ziAPIUpdateDevices

[ZIResult_enum](#) ziAPIUpdateDevices ([ZIConnection](#) conn)

Search for the newly connected devices and update the tree.

This function forces the Data Server to search for newly connected devices and to connect to run them

Parameters:

[in] conn
Pointer to ZIConnection

Returns:

- ZI_INFO_SUCCESS
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIListNodes](#)

ziAPIConnectDevice

ZIResult_enum ziAPIConnectDevice (**ZIConnection** conn, const char* deviceSerial, const char* deviceInterface, const char* interfaceParams)

Connect a device to the server.

This function connects a device with deviceSerial via the specified deviceInterface for use with the server.

Parameters:

[in] conn

Pointer to the ZIConnection with which the connection should be established

[in] deviceSerial

The serial of the device to connect to, e.g., dev2100

[in] deviceInterface

The interface to use for the connection, e.g., USB|1GbE

[in] interfaceParams

Parameters for interface configuration

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIDisconnectDevice](#), [ziAPIConnect](#), [ziAPIDisconnect](#), [ziAPIInit](#)

ziAPIDisconnectDevice

ZIResult_enum `ziAPIDisconnectDevice (ZIConnection conn, const char* deviceSerial)`

Disconnect a device from the server.

This function disconnects a device specified by deviceSerial from the server.

Parameters:

[in] conn

Pointer to the [ZIConnection](#) with which the connection should be established

[in] deviceSerial

The serial of the device to connect to, e.g., dev2100

Returns:

- [ZI_INFO_SUCCESS](#) on success
- [ZI_ERROR_TIMEOUT](#) when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIConnectDevice](#), [ziAPIConnect](#), [ziAPIDisconnect](#), [ziAPIInit](#)

ziAPIGetValueD

ZIResult_enum ziAPIGetValueD (**ZIConnection** conn, const char* path, ZIDoubleData* value)

gets the double-type value of the specified node

This function retrieves the numerical value of the specified node as an double-type value. The value first found is returned if more than one value is available (a wildcard is used in the path).

Parameters:

[in] conn

Pointer to ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to a double in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the path's length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIDoubleData ValueD;

    if ((RetVal = ziAPISetValueI(Conn,
                                "/dev1046/demods/*/rate",
                                100)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }

    if ((RetVal = ziAPIGetValueD(Conn,
                                "/dev1046/demods/0/rate",
```

```
                                &ValueD)) != ZI_INFO_SUCCESS) {  
    ziAPIGetError(RetVal, &ErrBuffer, NULL);  
    fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);  
} else {  
    printf("Value = %f\n", ValueD);  
}  
}
```

See Also:

[ziAPISetValueD](#), [ziAPIGetValueAsPollData](#)

ziAPIGetValueI

ZIResult_enum ziAPIGetValueI (**ZIConnection** conn, const char* path, ZIIntegerData* value)

gets the integer-type value of the specified node

This function retrieves the numerical value of the specified node as an integer-type value. The value first found is returned if more than one value is available (a wildcard is used in the path).

Parameters:

[in] conn

Pointer to ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to an 64bit integer in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the path's length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIIntegerData ValueI;

    if ((RetVal = ziAPISetValueD(Conn,
                                "/dev1046/demods/*/rate",
                                5.53)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }

    if ((RetVal = ziAPIGetValueI(Conn,
                                "/dev1046/demods/0/rate",
```

```
                                &ValueI)) != ZI_INFO_SUCCESS) {  
    ziAPIGetError(RetVal, &ErrBuffer, NULL);  
    fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);  
} else {  
    printf("Value = %f\n", (float)ValueI);  
}  
}
```

See Also:

[ziAPISetValueI](#), [ziAPIGetValueAsPollData](#)

ziAPIGetDemodSample

ZIResult_enum ziAPIGetDemodSample (**ZIConnection** conn, const char* path, **ZIDemodSample*** value)

Gets the demodulator sample value of the specified node.

This function retrieves the value of the specified node as an **DemodSample** struct. The value first found is returned if more than one value is available (a wildcard is used in the path). This function is only applicable to paths matching DEMODS/[0-9]+/SAMPLE.

Parameters:

[in] conn

Pointer to ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to a **ZIDemodSample** struct in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the path's length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void GetSample(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    ZIDemodSample DemodSample;

    if ((RetVal = ziAPIGetDemodSample(Conn,
                                     "/dev1046/demods/0/sample",
                                     &DemodSample)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
```

```
        printf("TS = %f, X=%f, Y=%f\n",  
              (float)DemodSample.timeStamp,  
              DemodSample.x,  
              DemodSample.y);  
    }  
}
```

See Also:

[ziAPIGetValueAsPollData](#)

ziAPIGetDIOSample

ZIResult_enum ziAPIGetDIOSample (**ZIConnection** conn, const char* path, **ZIDIOSample*** value)

Gets the Digital I/O sample of the specified node.

This function retrieves the newest available DIO sample from the specified node. The value first found is returned if more than one value is available (a wildcard is used in the path). This function is only applicable to nodes ending in "/DIOS/[0-9]+/INPUT".

Parameters:

[in] conn

Pointer to the ZIConnection with which the value should be retrieved

[in] path

Path to the node holding the value

[out] value

Pointer to a **ZIDIOSample** struct in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void GetSample(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    ZIDIOSample DIOSample;

    if ((RetVal = ziAPIGetDIOSample(Conn,
                                    "/dev1046/dios/0/output",
                                    &DIOSample)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
```

```
        printf("TS = %f, bits=%08x\n",  
              (float)DIOSample.timeStamp,  
              DIOSample.bits);  
    }  
}
```

See Also:

[ziAPIGetValueAsPollData](#)

ziAPIGetAuxInSample

ZIResult_enum ziAPIGetAuxInSample (**ZIConnection** conn, const char* path, **ZIAuxInSample*** value)

gets the AuxIn sample of the specified node

This function retrieves the newest available AuxIn sample from the specified node. The value first found is returned if more than one value is available (a wildcard is used in the path). This function is only applicable to nodes ending in "/AUXINS/[0-9]+/SAMPLE".

Parameters:

[in] conn

Pointer to the ziConnection with which the Value should be retrieved

[in] path

Path to the Node holding the value

[out] value

Pointer to an **ZIAuxInSample** struct in which the value should be written

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void GetSample(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    ZIAuxInSample AuxInSample;

    if ((RetVal = ziAPIGetAuxInSample(Conn,
                                     "/dev1046/auxins/0/sample",
                                     &AuxInSample)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s\n", ErrBuffer);
    }
}
```

```
    } else {  
        printf("TS = %f, ch0=%f, ch1=%f\n",  
              (float)AuxInSample.timeStamp,  
              AuxInSample.ch0,  
              AuxInSample.ch1);  
    }  
}
```

See Also:

[ziAPIGetValueAsPollData](#)

ziAPIGetValueB

ZIResult_enum ziAPIGetValueB (**ZIConnection** conn, const char* path, unsigned char* buffer, unsigned int* length, unsigned int bufferSize)

gets the Bytearray value of the specified node

This function retrieves the newest available DIO sample from the specified node. The value first found is returned if more than one value is available (a wildcard is used in the path).

Parameters:

[in] conn

Pointer to the ziConnection with which the value should be retrieved

[in] path

Path to the Node holding the value

[out] buffer

Pointer to a buffer to store the retrieved data in

[out] length

Pointer to an unsigned int to store the length of data in. if an error occurred or the length of the passed buffer doesn't reach a zero will be returned

[in] bufferSize

The length of the passed buffer

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void PrintVersion(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
```

```
const char* Path = "ZI/ABOUT/VERSION";
unsigned char Buffer[0xff];
unsigned int Length;

if ((RetVal = ziAPIGetValueB(Conn,
                             Path,
                             Buffer,
                             &Length,
                             sizeof(Buffer) - 1)) != ZI_INFO_SUCCESS) {
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't get value: %s.\n", ErrBuffer);
} else {
    Buffer[Length] = 0;
    printf("%s=\"%s\"\n", Path, Buffer);
}
}
```

See Also:

[ziAPISetValueB](#), [ziAPIGetValueAsPollData](#)

ziAPISetValueD

ZIResult_enum ziAPISetValueD (**ZIConnection** conn, const char* path, ZIDoubleData value)

asynchronously sets a double-type value to one or more nodes specified in the path

This function sets the values of the nodes specified in path to Value. More than one value can be set if a wildcard is used. The function sets the value asynchronously which means that after the function returns you have no security to which value it is finally set nor at what point in time it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set.

[in] path

Path to the Node(s) for which the value(s) will be set to Value.

[in] value

The double-type value that will be written to the node(s).

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred.
- ZI_ERROR_READONLY on attempt to set a read-only node.
- ZI_ERROR_COMMAND on an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIIntegerData ValueI;

    if ((RetVal = ziAPISetValueD(Conn,
                                "/dev1046/demods/*/rate",
                                5.53)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }
}
```

```
    }

    if ((RetVal = ziAPIGetValueI(Conn,
                                "/dev1046/demods/0/rate",
                                &ValueI)) != ZI\_INFO\_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
        printf("Value = %f\n", (float)ValueI);
    }
}
```

See Also:

[ziAPIGetValueD](#). [ziAPISyncSetValueD](#)

ziAPISetValueI

ZIResult_enum ziAPISetValueI (**ZIConnection** conn, const char* path, ZIIntegerData value)

asynchronously sets an integer-type value to one or more nodes specified in a path

This function sets the values of the nodes specified in path to Value. More than one value can be set if a wildcard is used. The function sets the value asynchronously which means that after the function returns you have no security to which value it is finally set nor at what point in time it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set

[in] value

The int-type value that will be written to the node(s)

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred.
- ZI_ERROR_READONLY on attempt to set a read-only node.
- ZI_ERROR_COMMAND on an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void UpdateValue(ZIConnection Conn) {
    ZIResult_enum RetVal;
    char* ErrBuffer;
    ZIDoubleData ValueD;

    if ((RetVal = ziAPISetValueI(Conn,
                                "/dev1046/demods/*/rate",
                                100)) != ZI_INFO_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
    }
}
```

```
    }

    if ((RetVal = ziAPIGetValueD(Conn,
                                "/dev1046/demods/0/rate",
                                &ValueD)) != ZI\_INFO\_SUCCESS) {
        ziAPIGetError(RetVal, &ErrBuffer, NULL);
        fprintf(stderr, "Error, can't get Parameter: %s.\n", ErrBuffer);
    } else {
        printf("Value = %f\n", ValueD);
    }
}
```

See Also:

[ziAPIGetValueI](#). [ziAPISyncSetValueI](#)

ziAPISetValueB

ZIResult_enum ziAPISetValueB (**ZIConnection** conn, const char* path, unsigned char* buffer, unsigned int length)

asynchronously sets the binary-type value of one ore more nodes specified in the path

This function sets the values at the nodes specified in a path. More than one value can be set if a wildcard is used. The function sets the value asynchronously which means that after the function returns you have no security to which value it is finally set nor at what point in time it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set

[in] buffer

Pointer to the byte array with the data

[in] length

Length of the data in the buffer

Returns:

- ZI_INFO_SUCCESS on success.
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred.
- ZI_ERROR_READONLY on attempt to set a read-only node.
- ZI_ERROR_COMMAND on an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values.
- ZI_ERROR_TIMEOUT when communication timed out.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

```
// Copyright [2016] Zurich Instruments AG
#include <stdlib.h>
#include <stdio.h>

#include "ziAPI.h"

void ProgramCPU(ZIConnection Conn,
               unsigned char* Buffer,
               int Len) {
    ZIResult_enum RetVal;
    char* ErrBuffer;

    if ((RetVal = ziAPISetValueB(Conn,
```

```
        "/dev1046/cpus/0/program",
        Buffer,
        Len)) != ZI_INFO_SUCCESS) {
    ziAPIGetError(RetVal, &ErrBuffer, NULL);
    fprintf(stderr, "Error, can't set Parameter: %s.\n", ErrBuffer);
}
}
```

See Also:

[ziAPIGetValueB](#). [ziAPISyncSetValueB](#)

ziAPISyncSetValueD

ZIResult_enum ziAPISyncSetValueD (**ZIConnection** conn, const char* path, ZIDoubleData* value)

synchronously sets a double-type value to one or more nodes specified in the path

This function sets the values of the nodes specified in path to Value. More than one value can be set if a wildcard is used. The function sets the value synchronously. After returning you know that it is set and to which value it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set to value

[in] value

Pointer to a double-type containing the value to be written. When the function returns value holds the effectively written value.

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_READONLY on attempt to set a read-only node
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIGetValueD](#), [ziAPISetValueD](#)

ziAPISyncSetValueI

ZIResult_enum ziAPISyncSetValueI (**ZIConnection** conn, const char* path, **ZIIntegerData*** value)

synchronously sets an integer-type value to one or more nodes specified in a path

This function sets the values of the nodes specified in path to value. More than one value can be set if a wildcard is used. The function sets the value synchronously. After returning you know that it is set and to which value it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the node(s) for which the value(s) will be set

[in] value

Pointer to a int-type containing then value to be written. when the function returns value holds the effectively written value.

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_READONLY on attempt to set a read-only node
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIGetValueI](#), [ziAPISetValueI](#)

ziAPISyncSetValueB

ZIResult_enum ziAPISyncSetValueB (**ZIConnection** conn, const char* path, uint8_t* buffer, uint32_t* length, uint32_t bufferSize)

Synchronously sets the binary-type value of one ore more nodes specified in the path.

This function sets the values at the nodes specified in a path. More than one value can be set if a wildcard is used. This function sets the value synchronously. After returning you know that it is set and to which value it is set.

Parameters:

[in] conn

Pointer to the ziConnection for which the value(s) will be set

[in] path

Path to the Node(s) for which the value(s) will be set

[in] buffer

Pointer to the byte array with the data

[in] length

Length of the data in the buffer

[in] bufferSize

Length of the data in the buffer

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_READONLY on attempt to set a read-only node
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIGetValueB](#), [ziAPISetValueB](#)

ziAPISync

ZIResult_enum ziAPISync (ZIConnection conn)

Synchronizes the session by dropping all pending data.

This function drops any data that is pending for transfer. Any data (including poll data) retrieved afterwards is guaranteed to be produced not earlier than the call to ziAPISync. This ensures in particular that any settings made prior to the call to ziAPISync have been propagated to the device, and the data retrieved afterwards is produced with the new settings already set to the hardware. Note, however, that this does not include any required settling time.

Parameters:

[in] conn

Pointer to the ZIConnection that is to be synchronized

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIEchoDevice

ZIResult_enum ziAPIEchoDevice (**ZIConnection** conn, const char* deviceSerial)

Sends an echo command to a device and blocks until answer is received.

This is useful to flush all buffers between API and device to enforce that further code is only executed after the device executed a previous command. Per device echo is only implemented for HF2. For other device types it is a synonym to ziAPISync, and deviceSerial parameter is ignored.

Parameters:

[in] conn

Pointer to the ZIConnection that is to be synchronized

[in] deviceSerial

The serial of the device to get the echo from, e.g., dev2100

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIAllocateEventEx

ZIEvent* ziAPIAllocateEventEx ()

Allocates [ZIEvent](#) structure and returns the pointer to it. Attention!!! It is the client code responsibility to deallocate the structure by calling [ziAPIDeallocateEventEx](#)!

This function allocates a [ZIEvent](#) structure and returns the pointer to it. Free the memory using [ziAPIDeallocateEventEx](#).

See Also:

[ziAPIDeallocateEventEx](#)

ziAPIDeallocateEventEx

void ziAPIDeallocateEventEx ([ZIEvent*](#) ev)

Deallocates [ZIEvent](#) structure created with [ziAPIAllocateEventEx\(\)](#).

Parameters:

[in] ev

Pointer to [ZIEvent](#) structure to be deallocated..

See Also:

[ziAPIAllocateEventEx](#)

This function is the compliment to [ziAPIAllocateEventEx\(\)](#)

ziAPISubscribe

ZIResult_enum ziAPISubscribe (**ZIConnection** conn, const char* path)

subscribes the nodes given by path for [ziAPIPollDataEx](#)

This function subscribes to nodes so that whenever the value of the node changes the new value can be polled using [ziAPIPollDataEx](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one leaf can be subscribed to with one function call.

Parameters:

[in] conn

Pointer to the ziConnection for which to subscribe for

[in] path

Path to the nodes to subscribe

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPIUnSubscribe](#), [ziAPIPollDataEx](#), [ziAPIGetValueAsPollData](#)

ziAPIUnSubscribe

ZIResult_enum ziAPIUnSubscribe (**ZIConnection** conn, const char* path)

unsubscribes to the nodes given by path

This function is the complement to [ziAPISubscribe](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one node can be unsubscribed with one function call.

Parameters:

[in] conn

Pointer to the ziConnection for which to unsubscribe for

[in] path

Path to the Nodes to unsubscribe

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no node given by path is able to hold values
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIPollDataEx](#), [ziAPIGetValueAsPollData](#)

ziAPIPollDataEx

ZIResult_enum ziAPIPollDataEx (**ZIConnection** conn, **ZIEvent*** ev, uint32_t timeOutMilliseconds)

checks if an event is available to read

This function returns immediately if an event is pending. Otherwise it waits for an event for up to timeOutMilliseconds. All value changes that occur in nodes that have been subscribed to or in children of nodes that have been subscribed to are sent from the Data Server to the ziAPI session. For a description of how the data are available in the struct, refer to the documentation of struct [ziEvent](#). When no event was available within timeOutMilliseconds, the ziEvent::Type field will be ZI_DATA_NONE and the ziEvent::Count field will be zero. Otherwise these fields hold the values corresponding to the event that occurred.

Parameters:

[in] conn

Pointer to the [ZIConnection](#) for which events should be received

[out] ev

Pointer to a [ZIEvent](#) struct in which the received event will be written

[in] timeOutMilliseconds

Time to wait for an event in milliseconds. If -1 it will wait forever, if 0 the function returns immediately.

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIGetValueAsPollData](#), [ziEvent](#)

ziAPIGetValueAsPollData

ZIResult_enum ziAPIGetValueAsPollData ([ZIConnection](#) conn, const char* path)

triggers a value request, which will be given back on the poll event queue

Use this function to receive the value of one or more nodes as one or more events using [ziAPIPollDataEx](#), even when the node is not subscribed or no value change has occurred.

Parameters:

[in] conn

Pointer to the [ZIConnection](#) with which the value should be retrieved

[in] path

Path to the Node holding the value

Returns:

- ZI_INFO_SUCCESS on success
- ZI_ERROR_CONNECTION when the connection is invalid (not connected) or when a communication error occurred
- ZI_ERROR_LENGTH if the Path's Length exceeds MAX_PATH_LEN or the length of the char-buffer for the nodes given by MaxLen is too small for all elements
- ZI_WARNING_OVERFLOW when a FIFO overflow occurred
- ZI_ERROR_COMMAND on an incorrect answer of the server
- ZI_ERROR_SERVER_INTERNAL if an internal error occurred in the Data Server
- ZI_WARNING_NOTFOUND if the given path could not be resolved or no value is attached to the node
- ZI_ERROR_TIMEOUT when communication timed out
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIPollDataEx](#)

ziAPIAsyncSetDoubleData

ZIResult_enum ziAPIAsyncSetDoubleData (**ZIConnection** conn, const char* path, ZIDoubleData value)

ziAPIAsyncSetIntegerData

ZIResult_enum ziAPIAsyncSetIntegerData (**ZIConnection** conn, const char* path, ZIIntegerData value)

ziAPIAsyncSetByteArray

ZIResult_enum ziAPIAsyncSetByteArray (**ZIConnection** conn, const char* path, uint8_t* buffer, uint32_t length)

ziAPIAsyncSubscribe

[ZIResult_enum](#) ziAPIAsyncSubscribe ([ZIConnection](#) conn, const char* path, ZIAsyncTag tag)

ziAPIAsyncUnSubscribe

[ZIResult_enum](#) ziAPIAsyncUnSubscribe ([ZIConnection](#) conn, const char* path, ZIAsyncTag tag)

ziAPIAsyncGetValueAsPollData

[ZIResult_enum](#) ziAPIAsyncGetValueAsPollData ([ZIConnection](#) conn, const char* path, ZIAsyncTag tag)

ziAPIGetError

ZIResult_enum ziAPIGetError (**ZIResult_enum** result, char** buffer, int* base)

Returns a description and the severity for a **ZIResult_enum**.

This function returns a static char pointer to a description string for the given **ZIResult_enum** error code. It also provides a parameter returning the severity (info, warning, error). If the given error code does not exist a description for an unknown error and the base for an error will be returned. If a description or the base is not needed NULL may be passed. In general, it's recommended to use **ziAPIGetLastError** instead to get detailed error messages.

Parameters:

[in] result

A **ZIResult_enum** for which the description or base will be returned

[out] buffer

A pointer to a char array to return the description. May be NULL if no description is needed.

[out] base

The severity for the provided Status parameter:

- ZI_INFO_BASE For infos.
- ZI_WARNING_BASE For warnings.
- ZI_ERROR_BASE For errors.

Returns:

- ZI_INFO_SUCCESS Upon success.

ziAPIGetLastError

ZIResult_enum ziAPIGetLastError (**ZIConnection** conn, char* buffer, uint32_t bufferSize)

Returns the message from the last error that occurred.

This function can be used to obtain the error message from the last error that occurred associated with the provided **ZIConnection**. If the last ziAPI call is successful, then the last error message returned by ziAPIGetError is empty. Only ziAPI function calls that take **ZIConnection** as an input argument influence the message returned by ziAPIGetLastError, if they do not take **ZIConnection** as an input argument the last error message will neither be reset to be empty or set to an error message (in the case of the error). There are some exceptions to this rule, ziAPIGetLastError can also not be used with **ziAPIInit**, **ziAPIConnect**, **ziAPIConnectEx** and **ziAPIDestroy**. Note, a call to ziAPIGetLastError will also reset the last error message to empty if its call was successful. Since the buffer is left unchanged in the case of an error occurring in the call to ziAPIGetLastError it is safest to initialize the buffer with a known value, for example, "ziAPIGetLastError was not successful".

Parameters:

[in] conn

The **ZIConnection** from which to get the error message.

[out] buffer

A pointer to a char array to return the message.

[in] bufferSize

The length of the provided buffer.

Returns:

- **ZI_INFO_SUCCESS** Upon success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred. In this case the provided buffer is left unchanged.
- **ZI_ERROR_LENGTH** If the message's length exceeds the provided bufferSize, the message is truncated and written to buffer.

ziAPISetDebugLevel

void ziAPISetDebugLevel (int32_t debugLevel)

Enable ziAPI's log and set the severity level of entries to be included in the log.

Calling this function enables ziAPI's log at the specified severity level. On Windows the logs can be found by navigating to the Zurich Instruments "Logs" folder entry in the Windows Start Menu: Programs -> Zurich Instruments -> LabOne Servers -> Logs. This will open an Explorer window displaying folders containing log files from various LabOne components, in particular, the `ziAPILog` folder contains logs from ziAPI. On Linux, the logs can be found at `/tmp/ziAPILog_USERNAME`, where `USERNAME` is the same as the output of the `whoami` command.

Parameters:

[in] debugLevel

An integer specifying the log's severity level:

- trace: 0,
- info: 1,
- debug: 2,
- warning: 3,
- error: 4,
- fatal: 5,
- status: 6.

See Also:

[ziAPIWriteDebugLog](#)

ziAPIWriteDebugLog

void ziAPIWriteDebugLog (int32_t debugLevel, const char* message)

Write a message to ziAPI's log with the specified severity.

This function may be used to write a message to ziAPI's log from client code to assist with debugging. Note, this function is only available if the implementation used in [ziAPIConnectEx](#) is "ziAPI_Core" (the default implementation). Also logging must be first enabled using [ziAPISetDebugLevel](#).

Parameters:

[in] debugLevel

An integer specifying the severity of the message to write in the log:

- trace: 0,
- info: 1,
- debug: 2,
- warning: 3,
- error: 4,
- fatal: 5,
- status: 6.

[in] message

A character array comprising of the message to be written.

See Also:

[ziAPISetDebugLevel](#)

ReadMEMFile

[ZIResult_enum](#) ReadMEMFile (const char* filename, char* buffer, int32_t bufferSize, int32_t* bytesUsed)

ziAPIModCreate

ZIResult_enum ziAPIModCreate (**ZIConnection** conn, **ZIModuleHandle*** handle, const char* moduleId)

Create a **ZIModuleHandle** that can be used for asynchronous measurement tasks.

This function initializes a ziCore module and provides a pointer (handle) with which to access and work with it. Note that this function does not start the module's thread. Before the thread can be started (with **ziAPIModExecute**):

- the device serial (e.g., "dev100") to be used with module must be specified via **ziAPIModSetByteArray**.
- the desired data (node paths) to record during the measurement must be specified via **ziAPIModSubscribe**. The module's thread is stopped with **ziAPIModClear**.

Parameters:

[in] conn

The **ZIConnection** which should be used to initialize the module.

[out] handle

Pointer to the initialized **ZIModuleHandle**, which from then on can be used to reference the module.

[in] moduleId

The name specifying the type the module to create (only the following ziCore Modules are currently supported in ziAPI):

- "sweep" to initialize an instance of the Sweeper Module.
- "record" to initialize an instance of the Software Trigger (Recorder) Module.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** when the connection is invalid (not connected) or when a communication error occurred.
- **ZI_WARNING_NOTFOUND** if the provided moduleId was invalid.
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

ziAPIModExecute, **ziAPIModClear**

ziAPIModSetDoubleData

ZIResult_enum ziAPIModSetDoubleData (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, **ZIDoubleData** value)

Sets a module parameter to the specified double type.

This function is used to configure (set) module parameters which have double types.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to set data on.

[in] path

Path to the module parameter path.

[in] value

The double data to write to the path.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSetIntegerData](#), [ziAPIModSetByteArray](#)

ziAPIModSetIntegerData

ZIResult_enum ziAPIModSetIntegerData (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, **ZIntegerData** value)

Sets a module parameter to the specified integer type.

This function is used to configure (set) module parameters which have integer types.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to set data on.

[in] path

Path to the module parameter path.

[in] value

The integer data to write to the path.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSetDoubleData](#), [ziAPIModSetByteArray](#)

ziAPIModSetByteArray

ZIResult_enum ziAPIModSetByteArray (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, uint8_t* buffer, uint32_t length)

Sets a module parameter to the specified byte array.

This function is used to configure (set) module parameters which have byte array types.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to set data on.

[in] path

Path to the module parameter path.

[in] buffer

Pointer to the byte array with the data.

[in] length

Length of the data in the buffer.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSetDoubleData](#), [ziAPIModSetIntegerData](#)

ziAPIModListNodes

ZIResult_enum ziAPIModListNodes (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path, char* nodes, uint32_t bufferSize, uint32_t flags)

Returns all child parameter node paths found under the specified parent module parameter path.

This function returns a list of parameter names found at the specified path. The path may contain wildcards. The list is returned in a null-terminated char-buffer, each element delimited by a newline. If the maximum length of the buffer (bufferSize) is not sufficient for all elements, nothing will be returned and the return value will be **ZI_ERROR_LENGTH**. Note, the provided path must match the module being addressed, i.e., path must exactly start with "sweep/" for the Sweeper Module.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** from which the parameter names should be retrieved.

[in] path

Path for which all children will be returned. The path may contain wildcard characters.

[out] nodes

Upon call filled with newline-delimited list of the names of all the children found. The string is zero-terminated.

[in] bufferSize

The length of the buffer specified as the nodes output parameter.

[in] flags

A combination of flags (applied bitwise) as defined in **ZIListNodes_enum**.

Returns:

- **ZI_INFO_SUCCESS** On success
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_LENGTH** If the path's length exceeds **MAX_PATH_LEN** or the length of the char-buffer for the nodes given by bufferSize is too small for all elements.
- **ZI_WARNING_OVERFLOW** When a FIFO overflow occurred.
- **ZI_ERROR_COMMAND** On an incorrect answer of the server.
- **ZI_ERROR_SERVER_INTERNAL** If an internal error occurred in Data Server.
- **ZI_WARNING_NOTFOUND** If the given path could not be resolved.
- **ZI_ERROR_TIMEOUT** When communication timed out.

- `ZI_ERROR_GENERAL` If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIModSubscribe

ZIResult_enum ziAPIModSubscribe (**ZIConnection** conn, **ZIModuleHandle** handle, const char* path)

Subscribes to the nodes specified by path, these nodes will be recorded during module execution.

This function subscribes to nodes so that whenever the value of the node changes while the module is executing the new value will be accumulated and then read using [ziAPIModRead](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one leaf can be subscribed to with one function call.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module in which the nodes should be subscribed to.

[in] path

Path specifying the nodes to subscribe to, may contain wildcards.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or a general error occurred, enable ziAPI's log for detailed information, see [ziAPISetDebugLevel](#).
- ZI_ERROR_LENGTH If the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW When a FIFO overflow occurred.
- ZI_ERROR_COMMAND On an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL If an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND If the given path could not be resolved or no node given by path is able to hold values.
- ZI_ERROR_TIMEOUT When communication timed out.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModUnSubscribe](#), [ziAPIModRead](#)

ziAPIModUnSubscribe

ZIResult_enum ziAPIModUnSubscribe ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, const char* path)

Unsubscribes to the nodes specified by path.

This function is the complement to [ziAPIModSubscribe](#). By using wildcards or by using a path that is not a leaf node but contains sub nodes, more than one node can be unsubscribed with one function call.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module in which the nodes should be unsubscribed from.

[in] path

Path specifying the nodes to unsubscribe from, may contain wildcards.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_LENGTH If the Path's Length exceeds MAX_PATH_LEN.
- ZI_WARNING_OVERFLOW When a FIFO overflow occurred.
- ZI_ERROR_COMMAND On an incorrect answer of the server.
- ZI_ERROR_SERVER_INTERNAL If an internal error occurred in the Data Server.
- ZI_WARNING_NOTFOUND If the given path could not be resolved or no node given by path is able to hold values.
- ZI_ERROR_TIMEOUT When communication timed out.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModSubscribe](#), [ziAPIModRead](#)

ziAPIModExecute

ZIResult_enum ziAPIModExecute ([ZIConnection](#) conn, [ZIModuleHandle](#) handle)

Starts the module's thread and its associated measurement task.

Once the module's parameters has been configured as required via, e.g. [ziAPIModSetDoubleData](#), this function starts the module's thread. This starts the module's main measurement task which will run asynchronously. The thread will run until either the module has completed its task or until [ziAPIModFinish](#) is called. Subscription or unsubscription is not possible while the module is executing. The status of the module can be obtained with either [ziAPIModFinished](#) or [ziAPIModProgress](#).

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModCreate](#), [ziAPIModProgress](#), [ziAPIModFinish](#)

ziAPIModTrigger

ZIResult_enum ziAPIModTrigger ([ZIConnection](#) conn, [ZIModuleHandle](#) handle)

Manually issue a trigger forcing data recording (SW Trigger Module only).

This function is used with the Software Trigger Module in order to manually issue a trigger in order to force recording of data. A burst of subscribed data will be recorded as configured via the SW Trigger's parameters as would a regular trigger event.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

ziAPIModProgress

ZIResult_enum ziAPIModProgress (**ZIConnection** conn, **ZIModuleHandle** handle, **ZIDoubleData*** progress)

Queries the current state of progress of the module's measurement task.

This function can be used to query the module's progress in performing its current measurement task, the progress is returned as a double in [0, 1], where 1 indicates task completion.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to execute.

[out] progress

A pointer to **ZIDoubleData** indicating the current progress of the module.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_GENERAL** If a general error occurred, use **ziAPIGetLastError** for a detailed error message.
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

ziAPIModExecute, **ziAPIModFinish**, **ziAPIModFinished**

ziAPIModFinished

ZIResult_enum ziAPIModFinished (**ZIConnection** conn, **ZIModuleHandle** handle, **ZIIntegerData*** finished)

Queries whether the module has finished its measurement task.

This function can be used to query whether the module has finished its task or not.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to execute.

[out] finished

A pointer to **ZIIntegerData**, upon return this will be 0 if the module is still executing or 1 if has finished executing.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModExecute](#), [ziAPIModFinish](#), [ziAPIModProgress](#)

ziAPIModFinish

ZIResult_enum ziAPIModFinish (**ZIConnection** conn, **ZIModuleHandle** handle)

Stops the module performing its measurement task.

This function stops the module performing its associated measurement task and stops recording any data. The task and data recording may be restarted by calling [ziAPIModExecute](#) again.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModProgress](#), [ziAPIModFinished](#)

ziAPIModSave

ZIResult_enum ziAPIModSave ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, const char* fileName)

Saves the currently accumulated data to file.

This function saves the currently accumulated data to a file. The path of the file to save data to is specified via the module's directory parameter.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

[in] fileName

The basename of the file to save the data in.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModExecute](#), [ziAPIModFinish](#), [ziAPIModFinished](#)

ziAPIModRead

ZIResult_enum ziAPIModRead ([ZIConnection](#) conn, [ZIModuleHandle](#) handle, const char* path)

Make the currently accumulated data available for use in the C program.

This function can be used to either read (get) module parameters, in this case a path that addresses the module must be specified, or it can be used to read out the currently accumulated data from subscribed nodes in the module. In either case the actual data must then be accessed by the user using [ziAPIModNextNode](#) and [ziAPIModGetChunk](#).

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

[in] path

The path specifying the module parameter(s) to get, specify NULL to obtain all subscribed data.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModGetChunk](#), [ziAPIModNextNode](#)

ziAPIModNextNode

ZIResult_enum ziAPIModNextNode (**ZIConnection** conn, **ZIModuleHandle** handle, char* path, uint32_t bufferSize, **ZIValueType_enum*** valueType, uint64_t* chunks)

Make the data for the next node available for reading with [ziAPIModGetChunk](#).

After calling [ziAPIModRead](#), subscribed data (or module parameters) may now be read out on a node-by-node and chunk-by-chunk basis. All nodes with data available in the module can be iterated over by using [ziAPIModNextNode](#), then for each node the chunks of data available are read out using [ziAPIModGetChunk](#). Calling this function makes the data from the next node available for read.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

[out] path

A string specifying the node's path whose data chunk points to.

[in] bufferSize

The length of the buffer specified as the path output parameter.

[out] valueType

The [ZIValueType_enum](#) of the node's data.

[out] chunks

The number of chunks of data available for the node.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModRead](#), [ziAPIModGetChunk](#), [ziAPIModEventDeallocate](#)

ziAPIModGetChunk

ZIResult_enum ziAPIModGetChunk (**ZIConnection** conn, **ZIModuleHandle** handle, **uint64_t** chunkIndex, **ZIModuleEventPtr*** ev)

Get the specified data chunk from the current node.

Data is read out node-by-node and then chunk-by-chunk. This function can be used to obtain specific data chunks from the current node that data is being read from. More precisely, it preallocates space for an event structure big enough to hold the node's data at the specified chunk index, updates **ZIModuleEventPtr** to point to this space and then copies the chunk data to this space.

Note, before the very first call to ziAPIModGetChunk, the **ZIModuleEventPtr** should be initialized to NULL and then left untouched for all subsequent calls (even after calling **ziAPIModNextNode** to get data from the next node). This is because ziAPIModGetChunk internally manages the required space allocation for the event and then in subsequent calls only reallocates space when it is required. It is optimized to reduce the number of required space reallocations for the event.

The **ZIModuleEventPtr** should be deallocated using **ziAPIModEventDeallocate**, otherwise the lifetime of the **ZIModuleEventPtr** is the same as the lifetime of the module. Indeed, the same **ZIModuleEventPtr** can be used, even for subsequent reads. It is also possible to work with multiple **ZIModuleEventPtr** so that some pointers can be kept for later processing.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to execute.

[out] chunkIndex

The index of the data chunk to update the pointer to.

[out] ev

The module's **ZIModuleEventPtr** that points to the currently available data chunk.

Returns:

- **ZI_INFO_SUCCESS** On success.
- **ZI_ERROR_CONNECTION** When the connection is invalid (not connected) or when a communication error occurred.
- **ZI_ERROR_GENERAL** If a general error occurred, use **ziAPIGetLastError** for a detailed error message.
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

[ziAPIModRead](#), [ziAPIModNextNode](#), [ziAPIModEventDeallocate](#)

ziAPIModEventDeallocate

ZIResult_enum ziAPIModEventDeallocate (**ZIConnection** conn, **ZIModuleHandle** handle, **ZIModuleEventPtr** ev)

Deallocate the **ZIModuleEventPtr** being used by the module.

This function deallocates the **ZIModuleEventPtr**. Since a module event's allocated space is managed internally by **ziAPIModGetChunk**, when the user no longer requires the event (all data has been read out) it must be deallocated by the user with this function.

Parameters:

[in] conn

The **ZIConnection** from which the module was created.

[in] handle

The **ZIModuleHandle** specifying the module to execute.

[in] ev

The **ZIModuleEventPtr** to deallocate.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use **ziAPIGetLastError** for a detailed error message.
- Other return codes may also be returned, for a detailed error message use **ziAPIGetLastError**.

See Also:

ziAPIModGetChunk, **ziAPIModRead**

ziAPIModClear

ZIResult_enum ziAPIModClear ([ZIConnection](#) conn, [ZIModuleHandle](#) handle)

Terminates the module's thread and destroys the module.

This function terminates the module's thread. After calling ziAPIModClear the module's handle may not be used any more. A new instance of the module must be initialized if required.

Parameters:

[in] conn

The [ZIConnection](#) from which the module was created.

[in] handle

The [ZIModuleHandle](#) specifying the module to execute.

Returns:

- ZI_INFO_SUCCESS On success.
- ZI_ERROR_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_ERROR_GENERAL If a general error occurred, use [ziAPIGetLastError](#) for a detailed error message.
- Other return codes may also be returned, for a detailed error message use [ziAPIGetLastError](#).

See Also:

[ziAPIModExecute](#), [ziAPIModFinish](#)

ziAPIVectorWriteBlock

ZIResult_enum ziAPIVectorWriteBlock (**ZIConnection** conn, const char* path, **ZIVectorData*** vectorBlock)

ziAPIVectorWriteGetStatus

ZIResult_enum ziAPIVectorWriteGetStatus (**ZIConnection** conn, const char* path, uint8_t* status)

status - see ZIVectorWriteStatus_enum

ziAPIVectorWrite

ZIResult_enum ziAPIVectorWrite (**ZIConnection** conn, const char* path, const void* vectorPtr, uint8_t vectorElementType, uint64_t vectorSizeElements)

vectorElementType - see ZIVectorElementType_enum

ziAPIDiscoveryFind

ZIResult_enum ziAPIDiscoveryFind (**ZIConnection** conn, const char* deviceAddress, const char** deviceId)

Returns the device id for a given device address. Attention! Invalidates all pointers previously returned by ziAPIDiscovery* calls.

ziAPIDiscoveryGet

ZIResult_enum ziAPIDiscoveryGet (**ZIConnection** conn, const char* deviceId, const char** propsJSON)

Returns the device properties for a given device id in JSON format.

ziAPIDiscoveryGetValueI

ZIResult_enum ziAPIDiscoveryGetValueI (**ZIConnection** conn, const char* deviceId, const char* propName, ZIIntegerData* value)

Returns given integer property value for a given device id.

ziAPIDiscoveryGetValueS

ZIResult_enum ziAPIDiscoveryGetValueS (**ZIConnection** conn, const char* deviceId, const char* propName, const char** value)

Returns given string property value for a given device id.

ziAPIAllocateEvent

`__inline ziEvent* ziAPIAllocateEvent ()`

Deprecated: See [ziAPIAllocateEventEx\(\)](#).

ziAPIDeallocateEvent

```
__inline void ziAPIDeallocateEvent ( ziEvent* ev )
```

Deprecated: See [ziAPIDeallocateEventEx\(\)](#).

ziAPIPollData

`__inline ZIResult_enum ziAPIPollData (ZIConnection conn, ziEvent* ev, int timeOut)`

Checks if an event is available to read. Deprecated: See [ziAPIPollDataEx\(\)](#).

Parameters:

[in] conn

Pointer to the [ZIConnection](#) for which events should be received

[out] ev

Pointer to a [ziEvent](#) struct in which the received event will be written

[in] timeOut

Time to wait for an event in milliseconds. If -1 it will wait forever, if 0 the function returns immediately.

Returns:

- ZI_SUCCESS On success.
- ZI_CONNECTION When the connection is invalid (not connected) or when a communication error occurred.
- ZI_OVERFLOW When a FIFO overflow occurred.

See [Data Handling](#) for an example

See Also:

[ziAPISubscribe](#), [ziAPIUnSubscribe](#), [ziAPIGetValueAsPollData](#), [ziEvent](#)

ziAPIGetValueS

```
__inline ZIResult_enum ziAPIGetValueS ( ZIConnection conn, char* path,  
DemodSample* value )
```

ziAPIGetValueDIO

```
__inline ZIResult_enum ziAPIGetValueDIO ( ZIConnection conn, char* path,  
DIOSample* value )
```

ziAPIGetValueAuxIn

```
__inline ZIResult_enum ziAPIGetValueAuxIn ( ZIConnection conn, char* path,  
AuxInSample* value )
```

ziAPISecondsTimeStamp

double ziAPISecondsTimeStamp (ziTimeStampType TS)

Deprecated: timestamps should instead be converted to seconds by dividing by the instrument's "clockbase". This is available as an leaf under the instrument's root "device" branch in the node hierarchy, e.g., /dev2001/clockbase.

Parameters:

[in] TS
the timestamp to convert to seconds

Returns:

The timestamp in seconds as a double

Glossary

This glossary provides easy to understand descriptions for many terms related to measurement instrumentation including the abbreviations used inside this user manual.

A

A/D	Analog to Digital See Also ADC .
AC	Alternate Current
ADC	Analog to Digital Converter
AM	Amplitude Modulation
Amplitude Modulated AFM (AM-AFM)	AFM mode where the amplitude change between drive and measured signal encodes the topography or the measured AFM variable. See Also Atomic Force Microscope .
API	Application Programming Interface
ASCII	American Standard Code for Information Interchange
Atomic Force Microscope (AFM)	Microscope that scans surfaces by means an oscillating mechanical structure (e.g. cantilever, tuning fork) whose oscillating tip gets so close to the surface to enter in interaction because of electrostatic, chemical, magnetic or other forces. With an AFM it is possible to produce images with atomic resolution. See Also Amplitude Modulated AFM , Frequency Modulated AFM , Phase modulation AFM .
AVAR	Allen Variance

B

Bandwidth (BW)	<p>The signal bandwidth represents the highest frequency components of interest in a signal. For filters the signal bandwidth is the cut-off point, where the transfer function of a system shows 3 dB attenuation versus DC. In this context the bandwidth is a synonym of cut-off frequency $f_{\text{cut-off}}$ or 3dB frequency $f_{-3\text{dB}}$. The concept of bandwidth is used when the dynamic behavior of a signal is important or separation of different signals is required.</p> <p>In the context of a open-loop or closed-loop system, the bandwidth can be used to indicate the fastest speed of the system, or the highest signal update change rate that is possible with the system.</p> <p>Sometimes the term bandwidth is erroneously used as synonym of frequency range. See Also Noise Equivalent Power Bandwidth.</p>
BNC	Bayonet Neill-Concelman Connector

C

CF	Clock Fail (internal processor clock missing)
----	---

Common Mode Rejection Ratio (CMRR)	Specification of a differential amplifier (or other device) indicating the ability of an amplifier to obtain the difference between two inputs while rejecting the components that do not differ from the signal (common mode). A high CMRR is important in applications where the signal of interest is represented by a small voltage fluctuation superimposed on a (possibly large) voltage offset, or when relevant information is contained in the voltage difference between two signals. The simplest mathematical definition of common-mode rejection ratio is: $CMRR = 20 * \log(\text{differential gain} / \text{common mode gain})$.
------------------------------------	--

CSV	Comma Separated Values
-----	------------------------

D

D/A	Digital to Analog
-----	-------------------

DAC	Digital to Analog Converter
-----	-----------------------------

DC	Direct Current
----	----------------

DDS	Direct Digital Synthesis
-----	--------------------------

DHCP	Dynamic Host Configuration Protocol
------	-------------------------------------

DIO	Digital Input/Output
-----	----------------------

DNS	Domain Name Server
-----	--------------------

DSP	Digital Signal Processor
-----	--------------------------

DUT	Device Under Test
-----	-------------------

Dynamic Reserve (DR)	The measure of a lock-in amplifier's capability to withstand the disturbing signals and noise at non-reference frequencies, while maintaining the specified measurement accuracy within the signal bandwidth.
----------------------	---

E

XML	Extensible Markup Language. See Also XML .
-----	---

F

FFT	Fast Fourier Transform
-----	------------------------

FIFO	First In First Out
------	--------------------

FM	Frequency Modulation
----	----------------------

Frequency Accuracy (FA)	Measure of an instrument's ability to faithfully indicate the correct frequency versus a traceable standard.
-------------------------	--

Frequency Modulated AFM (FM-AFM)	AFM mode where the frequency change between drive and measured signal encodes the topography or the measured AFM variable. See Also Atomic Force Microscope .
----------------------------------	--

Frequency Response Analyzer (FRA)	Instrument capable to stimulate a device under test and plot the frequency response over a selectable frequency range with a fine granularity.
-----------------------------------	--

Frequency Sweeper

See Also [Frequency Response Analyzer](#).

G

Gain Phase Meter

See Also [Vector Network Analyzer](#).

GPIO

General Purpose Interface Bus

GUI

Graphical User Interface

I

I/O

Input / Output

Impedance Spectroscope
(IS)

Instrument suited to stimulate a device under test and to measure the impedance (by means of a current measurement) at a selectable frequency and its amplitude and phase change over time. The output is both amplitude and phase information referred to the stimulus signal.

Input Amplitude Accuracy
(IAA)

Measure of instrument's capability to faithfully indicate the signal amplitude at the input channel versus a traceable standard.

Input voltage noise (IVN)

Total noise generated by the instrument and referred to the signal input, thus expressed as additional source of noise for the measured signal.

IP

Internet Protocol

L

LAN

Local Area Network

LED

Light Emitting Diode

Lock-in Amplifier (LI, LIA)

Instrument suited for the acquisition of small signals in noisy environments, or quickly changing signal with good signal to noise ratio - lock-in amplifiers recover the signal of interest knowing the frequency of the signal by demodulation with the suited reference frequency - the result of the demodulation are amplitude and phase of the signal compared to the reference: these are value pairs in the complex plane (X, Y) , (R, Θ) .

M

Media Access Control
address (MAC address)

Refers to the unique identifier assigned to network adapters for physical network communication.

Multi-frequency (MF)

Refers to the simultaneous measurement of signals modulated at arbitrary frequencies. The objective of multi-frequency is to increase the information that can be derived from a measurement which is particularly important for one-time, non-repeating events, and to increase the speed of a measurement since different frequencies do not have to be applied one after the other.

See Also [Multi-harmonic](#).

Multi-harmonic (MH)

Refers to the simultaneous measurement of modulated signals at various harmonic frequencies. The objective of multi-frequency is to increase the

information that can be derived from a measurement which is particularly important for one-time, non-repeating events, and to increase the speed of a measurement since different frequencies do not have to be applied one after the other.

See Also [Multi-frequency](#).

N

Noise Equivalent Power Bandwidth (NEPBW)

Effective bandwidth considering the area below the transfer function of a low-pass filter in the frequency spectrum. NEPBW is used when the amount of power within a certain bandwidth is important, such as noise measurements. This unit corresponds to a perfect filter with infinite steepness at the equivalent frequency.

See Also [Bandwidth](#).

Nyquist Frequency (NF)

For sampled analog signals, the Nyquist frequency corresponds to two times the highest frequency component that is being correctly represented after the signal conversion.

O

Output Amplitude Accuracy (OAA)

Measure of an instrument's ability to faithfully output a set voltage at a given frequency versus a traceable standard.

OV

Over Volt (signal input saturation and clipping of signal)

P

PC

Personal Computer

PD

Phase Detector

Phase-locked Loop (PLL)

Electronic circuit that serves to track and control a defined frequency. For this purpose a copy of the external signal is generated such that it is in phase with the original signal, but with usually better spectral characteristics. It can act as frequency stabilization, frequency multiplication, or as frequency recovery. In both analog and digital implementations it consists of a phase detector, a loop filter, a controller, and an oscillator.

Phase modulation AFM (PM-AFM)

AFM mode where the phase between drive and measured signal encodes the topography or the measured AFM variable.

See Also [Atomic Force Microscope](#).

PID

Proportional-Integral-Derivative

PL

Packet Loss (loss of packets of data between the instruments and the host computer)

R

RISC

Reduced Instruction Set Computer

Root Mean Square (RMS)

Statistical measure of the magnitude of a varying quantity. It is especially useful when variates are positive and negative, e.g., sinusoids, sawtooth, square waves. For a sine wave the following relation holds between the

amplitude and the RMS value: $U_{\text{RMS}} = U_{\text{PK}} / \sqrt{2} = U_{\text{PK}} / 1.41$. The RMS is also called quadratic mean.

RT Real-time

S

Scalar Network Analyzer (SNA) Instrument that measures the voltage of an analog input signal providing just the amplitude (gain) information.
See Also [Spectrum Analyzer](#), [Vector Network Analyzer](#).

SL Sample Loss (loss of samples between the instrument and the host computer)

Spectrum Analyzer (SA) Instrument that measures the voltage of an analog input signal providing just the amplitude (gain) information over a defined spectrum.
See Also [Scalar Network Analyzer](#).

SSH Secure Shell

T

TC Time Constant

TCP/IP Transmission Control Protocol / Internet Protocol

Thread An independent sequence of instructions to be executed by a processor.

Total Harmonic Distortion (THD) Measure of the non-linearity of signal channels (input and output)

TTL Transistor to Transistor Logic level

U

UHF Ultra-High Frequency

UHS Ultra-High Stability

USB Universal Serial Bus

V

VCO Voltage Controlled Oscillator

Vector Network Analyzer (VNA) Instrument that measures the network parameters of electrical networks, commonly expressed as s-parameters. For this purpose it measures the voltage of an input signal providing both amplitude (gain) and phase information. For this characteristic an older name was gain phase meter.
See Also [Gain Phase Meter](#), [Scalar Network Analyzer](#).

X

XML Extensible Markup Language: Markup language that defines a set of rules for encoding documents in a format that is both human-readable and machine-readable.

Z

ZCtrl	Zurich Instruments Control bus
ZoomFFT	This technique performs FFT processing on demodulated samples, for instance after a lock-in amplifier. Since the resolution of an FFT depends on the number of point acquired and the spanned time (not the sample rate), it is possible to obtain very highly resolution spectral analysis.
ZSync	Zurich Instruments Synchronization bus

Index

A

API

- Compatibility, 11
- Levels, 11
- Versions, 11

Asynchronous commands, 19

C

C API (see ziAPI)

C Programming Language (see ziAPI)

Comparison of LabOne APIs, 8

D

Data Server, 6

Device Settings Module, 40-40

L

LabOne

- API overview, 8
- Comparison of APIs, 8

LabVIEW, 123-131

- Comparison to other interfaces, 8
- Concepts, 126
- Examples, finding, 128
- Examples, running, 128
- Finding examples, 128
- Finding help, 127
- Getting started, 126
- Installing the API, 124
 - Linux, 124
 - Windows, 124
- LabOne VI Palette, 126
- Modules, 127
- Palette, LabOne, 126
- Running examples, 128
- Tips and tricks, 131
- VI Palette, 126

M

Matlab, 45-75

- Built-in help, 50
- Command reference, 58
- Comparison to other interfaces, 8
- Contents of the API package, 49
- Examples, running, 51
- Getting started, 49
- Help, accessing, 50
- Installing the API, 46
- List of Examples, 49
- List of Utility functions, 49
- Logging, 53
- Modules, 52
- Modules, configuring, 52

Reference, 58

Requirements, 46

Running examples, 51

Tips and tricks, 54

Troubleshooting, 56

Verifying correct configuration, 47

N

Node

- Concept, 13
- Leaf, 13
- Node hierarchy, 13

P

PLL Advisor Module, 41

Polling Data

- Concept, 13

Python, 76-122

- Built-in help, 81
- Command reference, 86
 - zhinst package, 86
 - zhinst's utility functions, 86
 - ziDAQRecorder class, 111
 - ziDAQServer class, 94
 - ziDAQSweeper class, 103
 - ziDAQZoomFFT class, 108
 - ziDeviceSettings class, 101
 - ziPidAdvisor class, 118
 - ziPLLAdvisor class, 115
- Comparison to other interfaces, 8
- Contents of the API package, 81
- Examples, running, 81
- Exploring the available examples, 82
- Getting started, 81
- Help, accessing, 81
- Installing the API, 77
 - Linux, 79
 - Windows, 78
- Loading data in Python, 85
- Locating the zhinst installation, 83
- Logging, 84
- Modules, 83
- Modules, configuring, 83
- Recommended python packages for ziPython, 77
- Reference, 86
- Requirements for using ziPython, 77
- Running examples, 81
- Tips and tricks, 85

S

Software Trigger Module, 34-39

Sweeper Module, 21-31

- Bandwidth control, 21
- Measurement data, 23
- Measurement data, averaging, 23
- Scanning mode, 21

Settling time, 22
 Settling time, definition, 22
 Sweep parameter, 21
 Sweep range, 21
 Synchronous commands, 19

U

UHF
 Leaf, 16

Z

ziAPI
 Comparison to other interfaces, 8
 ziAPI, C API functions and data types
 absz (see ZISweeperImpedanceSample)
 Action (see TreeChange)
 action (see ZITreeChangeData)
 advisorWave (see ZIEvent)
 aliasingReject (see ZISpectrumHeader)
 alignment (see ZIEvent) (see ziEvent::Val)
 allocatedSize (see ZIModuleEvent)
 asyncReply (see ZIEvent)
 AuxIn0 (see DemodSample)
 auxIn0 (see ZIDemodSample)
 auxIn1 (see ZISweeperDemodSample)
 AuxIn1 (see DemodSample)
 auxIn1 (see ZIDemodSample)
 auxIn1 (see ZISweeperDemodSample)
 auxInSample (see ZIEvent)
 avg (see ZIStatisticSample)
 bandwidth (see ZISpectrumHeader) (see
 ZISweeperDemodSample) (see
 ZISweeperDoubleSample) (see
 ZISweeperImpedanceSample)
 bandwidthMode (see ZISweeperHeader)
 bias (see ZIImpedanceSample) (see
 ZISweeperImpedanceSample)
 binCount (see ZIPWAWave)
 binPhase (see ZIPWASample)
 Bits (see DIOSample)
 bits (see ZIDIOSample)
 blockElements (see ZIVectorData)
 blockMarker (see ZIScopeWave) (see
 ZIScopeWaveEx)
 blockNumber (see ZIScopeWave) (see
 ZIScopeWaveEx) (see ZIVectorData)
 blockOffset (see ZIVectorData)
 BWLimit (see ScopeWave)
 byteArray (see ZIEvent)
 ByteArray (see ziEvent::Val)
 byteArrayTS (see ZIEvent)
 Bytes (see ByteArrayData)
 bytes (see ZIByteArray) (see ZIByteArrayTS)
 center (see ZISpectrumHeader)
 Ch0 (see AuxInSample)
 ch0 (see ZIAuxInSample)

Ch1 (see AuxInSample)
 ch1 (see ZIAuxInSample)
 channelBWLimit (see ZIScopeWave) (see
 ZIScopeWaveEx)
 channelEnable (see ZIScopeWave) (see
 ZIScopeWaveEx)
 channelInput (see ZIScopeWave) (see
 ZIScopeWaveEx)
 channelMath (see ZIScopeWave) (see
 ZIScopeWaveEx)
 channelOffset (see ZIScopeWaveEx)
 channelScaling (see ZIScopeWave) (see
 ZIScopeWaveEx)
 cntSample (see ZIEvent)
 cols (see ZISWTriggerHeader)
 command (see ZIAsyncReply)
 commensurable (see ZIPWAWave)
 count (see ZIEvent) (see
 ZISweeperDemodSample) (see
 ZISweeperDoubleSample) (see
 ZISweeperImpedanceSample)
 Count (see ScopeWave) (see ziEvent)
 countBin (see ZIPWASample)
 counter (see ZICntSample)
 data (see ZIAdvisorWave) (see ZIEvent)
 (see ZIPWAWave) (see ZIScopeWave) (see
 ZIScopeWaveEx) (see ZISpectrumWave) (see
 ZISweeperWave) (see ZIVectorData)
 Data (see ScopeWave) (see ziEvent)
 dataDemod (see ZISpectrumWave) (see
 ZISweeperWave)
 dataDouble (see ZISweeperWave) (see
 ZIVectorData)
 dataFloat (see ZIScopeWave) (see
 ZIScopeWaveEx) (see ZIVectorData)
 dataImpedance (see ZISweeperWave)
 dataInt16 (see ZIScopeWave) (see
 ZIScopeWaveEx) (see ZIVectorData)
 dataInt32 (see ZIScopeWave) (see
 ZIScopeWaveEx) (see ZIVectorData)
 dataInt64 (see ZIVectorData)
 dataInt8 (see ZIVectorData)
 dataTransferMode (see ZIScopeWave) (see
 ZIScopeWaveEx)
 dataUInt16 (see ZIVectorData)
 dataUInt32 (see ZIVectorData)
 dataUInt64 (see ZIVectorData)
 dataUInt8 (see ZIVectorData)
 demodSample (see ZIEvent)
 DIOBits (see DemodSample)
 dioBits (see ZIDemodSample)
 dioSample (see ZIEvent)
 Double (see ziEvent::Val)
 doubleData (see ZIEvent)
 doubleDataTS (see ZIEvent)
 drive (see ZIImpedanceSample) (see
 ZISweeperImpedanceSample)

dt (see ScopeWave) (see ZIScopeWave) (see ZIScopeWaveEx)
 elementType (see ZIVectorData)
 filter (see ZISpectrumDemodSample)
 flags (see ZIAdvisorHeader) (see ZIImpedanceSample) (see ZIScopeWave) (see ZIScopeWaveEx) (see ZISpectrumHeader) (see ZISweeperHeader) (see ZISWTriggerHeader) (see ZIVectorData)
 Frequency (see DemodSample)
 frequency (see ZIDemodSample) (see ZIImpedanceSample) (see ZIPWAWave) (see ZISweeperDemodSample) (see ZISweeperImpedanceSample)
 grid (see ZIAdvisorSample) (see ZISpectrumDemodSample) (see ZISweeperDemodSample) (see ZISweeperDoubleSample) (see ZISweeperImpedanceSample)
 harmonic (see ZIPWAWave)
 header (see ZIAdvisorWave) (see ZIModuleEvent) (see ZISpectrumWave) (see ZISweeperWave)
 id (see ZICntSample)
 imagz (see ZIImpedanceSample) (see ZISweeperImpedanceSample)
 impedanceSample (see ZIEvent)
 inputSelect (see ZIPWAWave)
 Integer (see ziEvent::Val)
 integerData (see ZIEvent)
 integerDataTS (see ZIEvent)
 Len (see ByteArrayData)
 length (see ZIByteArray) (see ZIByteArrayTS)
 MAX_EVENT_SIZE, 306
 MAX_NAME_LEN, 306
 MAX_PATH_LEN, 305
 mode (see ZIPWAWave)
 Name (see TreeChange)
 name (see ZITreeChangeData)
 nenbw (see ZISpectrumHeader)
 nextTimeStamp (see ZISweeperDemodSample) (see ZISweeperImpedanceSample)
 oscSelect (see ZIPWAWave)
 overflow (see ZIPWAWave)
 overlap (see ZISpectrumHeader)
 param0 (see ZIImpedanceSample) (see ZISweeperImpedanceSample)
 param1 (see ZIImpedanceSample) (see ZISweeperImpedanceSample)
 path (see ZIEvent)
 Path (see ziEvent)
 Phase (see DemodSample)
 phase (see ZIDemodSample) (see ZIImpedanceSample) (see ZISweeperDemodSample)
 phasez (see ZISweeperImpedanceSample)
 ptr (see ZIModuleHeader)
 pwaType (see ZIPWAWave)
 pwaWave (see ZIEvent)
 pwr (see ZIStatisticSample)
 r (see ZISpectrumDemodSample) (see ZISweeperDemodSample)
 rate (see ZISpectrumHeader)
 ReadMEMFile, 436
 realz (see ZIImpedanceSample) (see ZISweeperImpedanceSample)
 repetitions (see ZISWTriggerHeader)
 Reserved (see DemodSample) (see DIOsample)
 reserved (see ZIDIOsample) (see ZIPWASample)
 reserved0 (see ZIAdvisorHeader) (see ZICntSample) (see ZIScopeWave) (see ZIScopeWaveEx) (see ZISpectrumHeader) (see ZISweeperHeader) (see ZISWTriggerHeader) (see ZIVectorData)
 reserved1 (see ZIAdvisorHeader) (see ZICntSample) (see ZIScopeWaveEx) (see ZISpectrumHeader) (see ZISweeperHeader) (see ZIVectorData)
 reservedUInt (see ZIPWAWave)
 resolution (see ZISpectrumHeader)
 resultCode (see ZIAsyncReply)
 rows (see ZISWTriggerHeader)
 SampleAuxIn (see ziEvent::Val)
 sampleCount (see ZIAdvisorHeader) (see ZIPWAWave) (see ZIScopeWave) (see ZIScopeWaveEx) (see ZISpectrumHeader) (see ZISweeperHeader)
 SampleDemod (see ziEvent::Val)
 SampleDIO (see ziEvent::Val)
 sampleFormat (see ZIAdvisorHeader) (see ZIScopeWave) (see ZIScopeWaveEx) (see ZISpectrumHeader) (see ZISweeperHeader)
 sampleTimeStamp (see ZIAsyncReply)
 ScopeChannel (see ScopeWave)
 scopeWave (see ZIEvent)
 scopeWaveEx (see ZIEvent)
 scopeWaveOld (see ZIEvent)
 segmentNumber (see ZIScopeWave) (see ZIScopeWaveEx)
 sequenceNumber (see ZIScopeWave) (see ZIScopeWaveEx) (see ZIVectorData)
 setTimeStamp (see ZISweeperDemodSample) (see ZISweeperImpedanceSample)
 settling (see ZISweeperDemodSample) (see ZISweeperImpedanceSample)
 spectrumMode (see ZISpectrumHeader)
 spectrumWave (see ZIEvent)
 stddev (see ZIStatisticSample)
 sweeper (see ZIModuleHeader)
 sweeperWave (see ZIEvent)
 sweepMode (see ZISweeperHeader)
 swTrigger (see ZIModuleHeader)
 tag (see ZIAsyncReply)
 tc (see ZISweeperDemodSample) (see ZISweeperImpedanceSample)

tcMeas (see ZISweeperDemodSample) (see ZISweeperImpedanceSample)
 TimeStamp (see AuxInSample) (see DemodSample) (see DIOSample)
 timeStamp (see ZIAdvisorWave) (see ZIAsyncReply) (see ZIAuxInSample) (see ZIByteArrayTS) (see ZICntSample) (see ZIDemodSample) (see ZIDIOSample) (see ZIDoubleDataTS) (see ZIImpedanceSample) (see ZIIntegerDataTS) (see ZIPWAWave) (see ZIScopeWave) (see ZIScopeWaveEx) (see ZISpectrumWave) (see ZISweeperWave) (see ZITreeChangeData) (see ZIVectorData)
 totalElements (see ZIVectorData)
 totalSamples (see ZIScopeWave) (see ZIScopeWaveEx)
 traceName (see ZIModuleHeaderSweeper)
 Tree (see ziEvent::Val)
 treeChangeData (see ZIEvent)
 treeChangeDataOld (see ZIEvent)
 TREE_ACTION, 384
 TREE_ACTION_ADD (see TREE_ACTION)
 TREE_ACTION_CHANGE (see TREE_ACTION)
 TREE_ACTION_REMOVE (see TREE_ACTION)
 trigger (see ZIDemodSample) (see ZIImpedanceSample)
 TriggerChannel (see ScopeWave)
 triggerEnable (see ZIScopeWave) (see ZIScopeWaveEx)
 triggerInput (see ZIScopeWave) (see ZIScopeWaveEx)
 triggerNumber (see ZISWTriggerHeader)
 triggerStart (see ZISWTriggerHeader)
 triggerTimeStamp (see ZIScopeWave) (see ZIScopeWaveEx)
 Type (see ziEvent)
 type (see ZIModuleHeader)
 untyped (see ZIEvent) (see ZIModuleHeader)
 Val (see ziEvent)
 value (see ZIDoubleDataTS) (see ZIEvent) (see ZIIntegerDataTS) (see ZIModuleEvent) (see ZISweeperDoubleSample)
 valueType (see ZIEvent)
 vectorData (see ZIEvent)
 Void (see ziEvent::Val)
 Wave (see ziEvent::Val)
 window (see ZISpectrumHeader)
 X (see DemodSample)
 x (see ZIAdvisorSample) (see ZIDemodSample) (see ZIPWASample) (see ZISpectrumDemodSample) (see ZISweeperDemodSample)
 Y (see DemodSample)
 y (see ZIAdvisorSample) (see ZIDemodSample) (see ZIPWASample) (see ZISpectrumDemodSample) (see ZISweeperDemodSample)
 ziAPIAllocateEvent, 463
 ziAPIAllocateEventEx, 192, 420
 ziAPIAsyncGetValueAsPollData, 205, 431
 ziAPIAsyncSetByteArray, 202, 428
 ziAPIAsyncSetDoubleData, 200, 426
 ziAPIAsyncSetIntegerData, 201, 427
 ziAPIAsyncSubscribe, 203, 429
 ziAPIAsyncUnSubscribe, 204, 430
 ziAPIConnect, 140, 387
 ziAPIConnectDevice, 152, 395
 ziAPIConnectEx, 143, 390
 ziAPIDeallocateEvent, 464
 ziAPIDeallocateEventEx, 193, 421
 ziAPIDestroy, 139, 386
 ziAPIDisconnect, 141, 388
 ziAPIDisconnectDevice, 153, 396
 ziAPIDiscoveryFind, 244, 459
 ziAPIDiscoveryGet, 245, 460
 ziAPIDiscoveryGetValueI, 246, 461
 ziAPIDiscoveryGetValueS, 247, 462
 ziAPIEchoDevice, 178, 419
 ziAPIGetAuxInSample, 164, 405
 ziAPIGetConnectionAPILevel, 144, 391
 ziAPIGetDemodSample, 160, 401
 ziAPIGetDIOSample, 162, 403
 ziAPIGetError, 207, 432
 ziAPIGetLastError, 208, 433
 ziAPIGetRevision, 145, 392
 ziAPIGetValueAsPollData, 197, 425
 ziAPIGetValueAuxIn, 181, 468
 ziAPIGetValueB, 166, 407
 ziAPIGetValueD, 156, 397
 ziAPIGetValueDIO, 180, 467
 ziAPIGetValueI, 158, 399
 ziAPIGetValueS, 179, 466
 ziAPIInit, 138, 385
 ziAPIListImplementations, 142, 389
 ziAPIListNodes, 149, 393
 ziAPIModClear, 237, 455
 ziAPIModCreate, 219, 437
 ziAPIModEventDeallocate, 236, 454
 ziAPIModExecute, 227, 445
 ziAPIModFinish, 231, 449
 ziAPIModFinished, 230, 448
 ziAPIModGetChunk, 235, 453
 ziAPIModListNodes, 223, 441
 ziAPIModNextNode, 234, 452
 ziAPIModProgress, 229, 447
 ziAPIModRead, 233, 451
 ziAPIModSave, 232, 450
 ziAPIModSetByteArray, 222, 440
 ziAPIModSetDoubleData, 220, 438
 ziAPIModSetIntegerData, 221, 439
 ziAPIModSubscribe, 225, 443
 ziAPIModTrigger, 228, 446
 ziAPIModUnSubscribe, 226, 444
 ziAPIPollData, 198, 465

ziAPIPollDataEx, 196, 424	ZI_ERROR_DEVICE_IN_USE (see ZIResult_enum)
ziAPISecondsTimeStamp, 469	ZI_ERROR_DEVICE_NEEDS_FW_UPGRADE (see ZIResult_enum)
ziAPISetDebugLevel, 209, 434	ZI_ERROR_DEVICE_NOT_FOUND (see ZIResult_enum)
ziAPISetValueB, 172, 413	ZI_ERROR_DEVICE_NOT_VISIBLE (see ZIResult_enum)
ziAPISetValueD, 168, 409	ZI_ERROR_DUPLICATE (see ZIResult_enum)
ziAPISetValueI, 170, 411	ZI_ERROR_FILE (see ZIResult_enum)
ziAPISubscribe, 194, 422	ZI_ERROR_GENERAL (see ZIResult_enum)
ziAPISync, 177, 418	ZI_ERROR_HOSTNAME (see ZIResult_enum)
ziAPISyncSetValueB, 176, 417	ZI_ERROR_LENGTH (see ZIResult_enum)
ziAPISyncSetValueD, 174, 415	ZI_ERROR_MALLOC (see ZIResult_enum)
ziAPISyncSetValueI, 175, 416	ZI_ERROR_MAX (see ZIResult_enum)
ziAPIUnSubscribe, 195, 423	ZI_ERROR_MUTEX_DESTROY (see ZIResult_enum)
ziAPIUpdateDevices, 151, 394	ZI_ERROR_MUTEX_INIT (see ZIResult_enum)
ziAPIVectorWrite, 242, 458	ZI_ERROR_MUTEX_LOCK (see ZIResult_enum)
ziAPIVectorWriteBlock, 240, 456	ZI_ERROR_MUTEX_UNLOCK (see ZIResult_enum)
ziAPIVectorWriteGetStatus, 241, 457	ZI_ERROR_NOT_SUPPORTED (see ZIResult_enum)
ZIAPIVersion_enum, 137, 380	ZI_ERROR_READONLY (see ZIResult_enum)
ziAPIWriteDebugLog, 210, 435	ZI_ERROR_SERVER_INTERNAL (see ZIResult_enum)
ZIConnection, 135, 306	ZI_ERROR_SOCKET_CONNECT (see ZIResult_enum)
ZIImpFlags_enum, 378	ZI_ERROR_SOCKET_INIT (see ZIResult_enum)
ZIListNode_enum, 148, 381	ZI_ERROR_THREAD_JOIN (see ZIResult_enum)
ZIModuleEventPtr, 211, 306	ZI_ERROR_THREAD_START (see ZIResult_enum)
ZIModuleHandle, 306	ZI_ERROR_TIMEOUT (see ZIResult_enum)
ZIModuleHeaderType_enum, 218, 382	ZI_ERROR_TOO_MANY_CONNECTIONS (see ZIResult_enum)
ZIResult_enum, 370	ZI_ERROR_USB (see ZIResult_enum)
ZITreeAction_enum, 377	ZI_ERROR_ZIEVENT_DATATYPE_MISMATCH (see ZIResult_enum)
ZIValueType_enum, 374	ZI_FILE (see ZIResult_enum)
ZIVectorElementType_enum, 379	ZI_HOSTNAME (see ZIResult_enum)
ZIVectorWriteStatus_enum, 239, 383	ZI_IMP_FLAGS_AUTORANGE_GATING (see ZIImpFlags_enum)
ZI_API_VERSION_0 (see ZIAPIVersion_enum)	ZI_IMP_FLAGS_BWC_BIT0 (see ZIImpFlags_enum)
ZI_API_VERSION_1 (see ZIAPIVersion_enum)	ZI_IMP_FLAGS_BWC_BIT1 (see ZIImpFlags_enum)
ZI_API_VERSION_4 (see ZIAPIVersion_enum)	ZI_IMP_FLAGS_BWC_BIT2 (see ZIImpFlags_enum)
ZI_API_VERSION_5 (see ZIAPIVersion_enum)	ZI_IMP_FLAGS_BWC_BIT3 (see ZIImpFlags_enum)
ZI_COMMAND (see ZIResult_enum)	ZI_IMP_FLAGS_BWC_MASK (see ZIImpFlags_enum)
ZI_CONNECTION (see ZIResult_enum)	ZI_IMP_FLAGS_FREQ_EXACT (see ZIImpFlags_enum)
ZI_DATA_AUXINSAMPLE (see ZIValueType_enum)	ZI_IMP_FLAGS_FREQ_EXTRAPOLATION (see ZIImpFlags_enum)
ZI_DATA_BYTEARRAY (see ZIValueType_enum)	ZI_IMP_FLAGS_FREQ_INTERPOLATION (see ZIImpFlags_enum)
ZI_DATA_DEMODSAMPLE (see ZIValueType_enum)	ZI_IMP_FLAGS_NONE (see ZIImpFlags_enum)
ZI_DATA_DIOSAMPLE (see ZIValueType_enum)	
ZI_DATA_DOUBLE (see ZIValueType_enum)	
ZI_DATA_INTEGER (see ZIValueType_enum)	
ZI_DATA_NONE (see ZIValueType_enum)	
ZI_DATA_SCOPEWAVE (see ZIValueType_enum)	
ZI_DATA_TREE_CHANGED (see ZIValueType_enum)	
ZI_DUPLICATE (see ZIResult_enum)	
ZI_ERROR (see ZIResult_enum)	
ZI_ERROR_BASE (see ZIResult_enum)	
ZI_ERROR_COMMAND (see ZIResult_enum)	
ZI_ERROR_CONNECTION (see ZIResult_enum)	
ZI_ERROR_DEVICE_CONNECTION_TIMEOUT (see ZIResult_enum)	
ZI_ERROR_DEVICE_DIFFERENT_INTERFACE (see ZIResult_enum)	
ZI_ERROR_DEVICE_INTERFACE (see ZIResult_enum)	

ZI_IMP_FLAGS_OPEN_DETECTION	(see	ZI_SERVER_INTERNAL (see ZIResult_enum)
ZIImpFlags_enum)		
ZI_IMP_FLAGS_OVERFLOW_CURRENT	(see	ZI_SOCKET_CONNECT (see ZIResult_enum)
ZIImpFlags_enum)		
ZI_IMP_FLAGS_OVERFLOW_VOLTAGE	(see	ZI_SOCKET_INIT (see ZIResult_enum)
ZIImpFlags_enum)		
ZI_IMP_FLAGS_STRONGCOMPENSATION_PARAM0		ZI_SUCCESS (see ZIResult_enum)
(see ZIImpFlags_enum)		
ZI_IMP_FLAGS_STRONGCOMPENSATION_PARAM1		ZI_THREAD_JOIN (see ZIResult_enum)
(see ZIImpFlags_enum)		
ZI_IMP_FLAGS_SUPPRESSION_PARAM0	(see	ZI_THREAD_START (see ZIResult_enum)
ZIImpFlags_enum)		
ZI_IMP_FLAGS_SUPPRESSION_PARAM1	(see	ZI_TIMEOUT (see ZIResult_enum)
ZIImpFlags_enum)		
ZI_IMP_FLAGS_UNDERFLOW_CURRENT	(see	ZI_TREE_ACTION_ADD (see ZITreeAction_enum)
ZIImpFlags_enum)		
ZI_IMP_FLAGS_UNDERFLOW_VOLTAGE	(see	ZI_TREE_ACTION_CHANGE (see
ZIImpFlags_enum)		ZITreeAction_enum)
ZI_IMP_FLAGS_VALID_INTERNAL	(see	ZI_TREE_ACTION_REMOVE (see
ZIImpFlags_enum)		ZITreeAction_enum)
ZI_IMP_FLAGS_VALID_USER	(see	ZI_UNDERRUN (see ZIResult_enum)
ZIImpFlags_enum)		
ZI_INFO_BASE (see ZIResult_enum)		ZI_USB (see ZIResult_enum)
ZI_INFO_MAX (see ZIResult_enum)		
ZI_INFO_SUCCESS (see ZIResult_enum)		ZI_VALUE_TYPE_ADVISOR_WAVE (see
ZI_LENGTH (see ZIResult_enum)		ZIValueType_enum)
ZI_LIST_ABSOLUTE (see ZIListNodes_enum)		ZI_VALUE_TYPE_ASYNC_REPLY (see
ZI_LIST_LEAFSONLY (see ZIListNodes_enum)		ZIValueType_enum)
ZI_LIST_NODES_ABSOLUTE	(see	ZI_VALUE_TYPE_AUXIN_SAMPLE (see
ZIListNodes_enum)		ZIValueType_enum)
ZI_LIST_NODES_LEAFSONLY	(see	ZI_VALUE_TYPE_BYTE_ARRAY (see
ZIListNodes_enum)		ZIValueType_enum)
ZI_LIST_NODES_NONE (see ZIListNodes_enum)		ZI_VALUE_TYPE_BYTE_ARRAY_TS (see
ZI_LIST_NODES_RECURSIVE	(see	ZIValueType_enum)
ZIListNodes_enum)		
ZI_LIST_NODES_SETTINGSONLY	(see	ZI_VALUE_TYPE_CNT_SAMPLE (see
ZIListNodes_enum)		ZIValueType_enum)
ZI_LIST_NONE (see ZIListNodes_enum)		ZI_VALUE_TYPE_DEMOD_SAMPLE (see
ZI_LIST_RECURSIVE (see ZIListNodes_enum)		ZIValueType_enum)
ZI_LIST_SETTINGSONLY (see ZIListNodes_enum)		ZI_VALUE_TYPE_DIO_SAMPLE (see
ZI_MALLOC (see ZIResult_enum)		ZIValueType_enum)
ZI_MAX_ERROR (see ZIResult_enum)		ZI_VALUE_TYPE_DOUBLE_DATA (see
ZI_MAX_INFO (see ZIResult_enum)		ZIValueType_enum)
ZI_MAX_WARNING (see ZIResult_enum)		ZI_VALUE_TYPE_DOUBLE_DATA_TS (see
ZI_MODULE_HEADER_TYPE_NONE	(see	ZIValueType_enum)
ZIModuleHeaderType_enum)		
ZI_MODULE_HEADER_TYPE_SWEEPER	(see	ZI_VALUE_TYPE_IMPEDANCE_SAMPLE (see
ZIModuleHeaderType_enum)		ZIValueType_enum)
ZI_MODULE_HEADER_TYPE_SWTRIGGER	(see	ZI_VALUE_TYPE_INTEGER_DATA (see
ZIModuleHeaderType_enum)		ZIValueType_enum)
ZI_MUTEX_DESTROY (see ZIResult_enum)		ZI_VALUE_TYPE_INTEGER_DATA_TS (see
ZI_MUTEX_INIT (see ZIResult_enum)		ZIValueType_enum)
ZI_MUTEX_LOCK (see ZIResult_enum)		ZI_VALUE_TYPE_NONE (see ZIValueType_enum)
ZI_MUTEX_UNLOCK (see ZIResult_enum)		
ZI_NOTFOUND (see ZIResult_enum)		ZI_VALUE_TYPE_PWA_WAVE (see
ZI_OVERFLOW (see ZIResult_enum)		ZIValueType_enum)
ZI_READONLY (see ZIResult_enum)		ZI_VALUE_TYPE_SCOPE_WAVE (see
		ZIValueType_enum)
		ZI_VALUE_TYPE_SCOPE_WAVE_EX (see
		ZIValueType_enum)
		ZI_VALUE_TYPE_SCOPE_WAVE_OLD (see
		ZIValueType_enum)
		ZI_VALUE_TYPE_SPECTRUM_WAVE (see
		ZIValueType_enum)
		ZI_VALUE_TYPE_SWEEPER_WAVE (see
		ZIValueType_enum)
		ZI_VALUE_TYPE_TREE_CHANGE_DATA (see
		ZIValueType_enum)
		ZI_VALUE_TYPE_TREE_CHANGE_DATA_OLD (see
		ZIValueType_enum)

ZI_VALUE_TYPE_VECTOR_DATA (see
ZIValueType_enum)
ZI_VECTOR_ELEMENT_TYPE_ASCII_Z (see
ZIVectorElementType_enum)
ZI_VECTOR_ELEMENT_TYPE_DOUBLE (see
ZIVectorElementType_enum)
ZI_VECTOR_ELEMENT_TYPE_FLOAT (see
ZIVectorElementType_enum)
ZI_VECTOR_ELEMENT_TYPE_UINT16 (see
ZIVectorElementType_enum)
ZI_VECTOR_ELEMENT_TYPE_UINT32 (see
ZIVectorElementType_enum)
ZI_VECTOR_ELEMENT_TYPE_UINT64 (see
ZIVectorElementType_enum)
ZI_VECTOR_ELEMENT_TYPE_UINT8 (see
ZIVectorElementType_enum)
ZI_VECTOR_WRITE_STATUS_IDLE (see
ZIVectorWriteStatus_enum)
ZI_VECTOR_WRITE_STATUS_PENDING (see
ZIVectorWriteStatus_enum)
ZI_WARNING (see ZIResult_enum)
ZI_WARNING_BASE (see ZIResult_enum)
ZI_WARNING_GENERAL (see ZIResult_enum)
ZI_WARNING_MAX (see ZIResult_enum)
ZI_WARNING_NOTFOUND (see ZIResult_enum)
ZI_WARNING_NO_ASYNC (see ZIResult_enum)
ZI_WARNING_OVERFLOW (see ZIResult_enum)
ZI_WARNING_UNDERRUN (see ZIResult_enum)
ziCore
 Tips and Tricks, 43-43
zoomFFT Module, 32-33